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Nº EX-2021-86572011 -APN-DGA#ANMAT

AM

Digitally signed by Gestion Documental Electronica  
Date: 2022.11.03 18:25:06 -03:00

Digitally signed by Gestion Documental  
Electronica  
Date: 2022.11.03 18:25:07 -03:00

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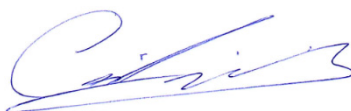
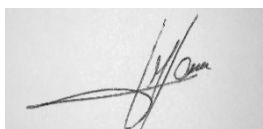
The image shows a rectangular label for the CHEMLIS® Chagas R-iELISA kit. The label is white with a blue border and contains the following information:

- Top Right:** CHEMTEST logo (three blue curved lines).
- Top Left:** CHEMLIS® Chagas R-iELISA in large black font.
- Center:** Kit para la detección de anticuerpos anti *Trypanosoma cruzi* en muestras de suero y plasma de humanos. Below it, the English translation: *Kit for the detection of antibodies against Trypanosoma cruzi in human serum and plasma*.
- Below Center:** Autorizado por ANMAT PM 2360-02. Below that, a blue box with the text: IVD Diagnóstico in vitro / In vitro diagnostics.
- Bottom Section (Blue Background):**
  - Left:** ELISA logo (four white circles) and Humanos logo (white person icon).
  - Center:** REF 15-CL01-2P and LOT 15-CL01/2003-1.
  - Right:** 192 (with a downward arrow icon), 03/2021 (with a calendar icon), and 20°C to 8°C (with a thermometer icon).
  - Bottom Right:** Chemtest Argentina S.A. Director Técnico: Andrés E. Ciocchini. Av. 25 de Mayo 1021, San Martín (C.P. 1650), Buenos Aires, Argentina. info@chemtest.net - chemtest.net.

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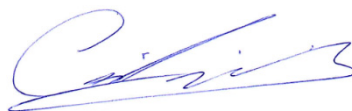
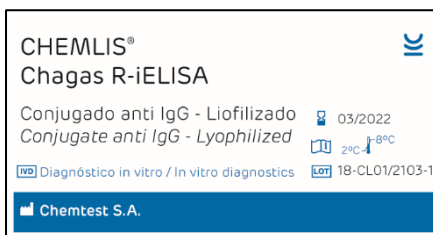
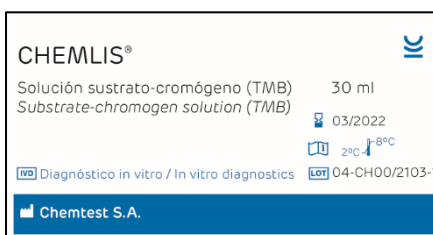
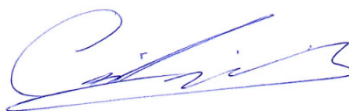
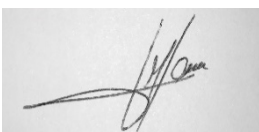
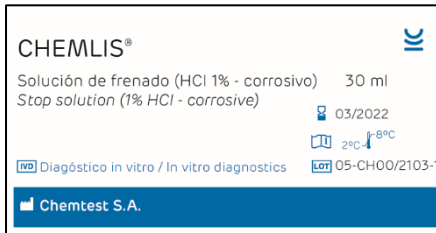
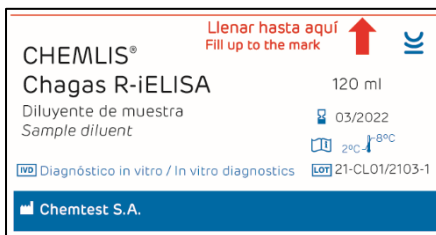
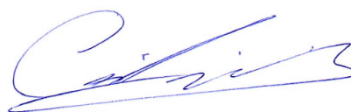
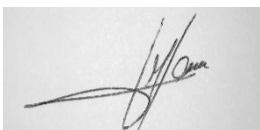
*Kit for the detection of antibodies against Trypanosoma cruzi in human serum and plasma*

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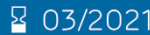
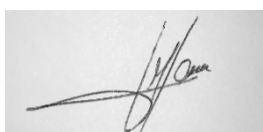
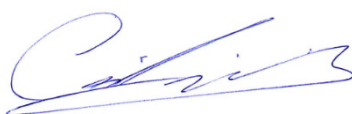
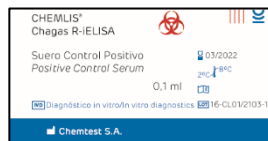
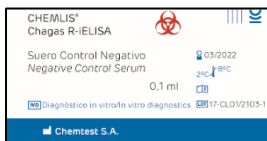
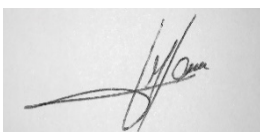
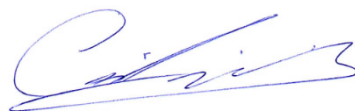
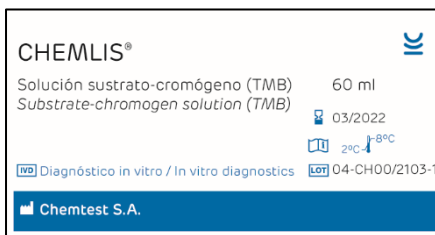
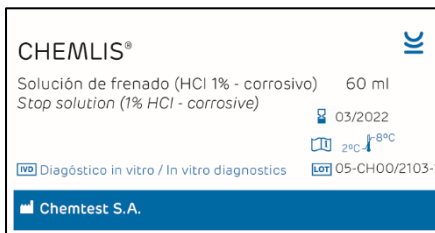
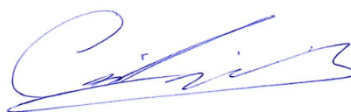
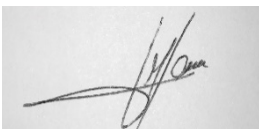
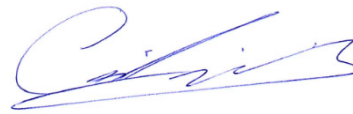
*Indirect ELISA kit for the detection of anti *Trypanosoma cruzi* antibodies in human serum and plasma.*

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Electrochemical magnetic microbeads-based biosensor for point-of-care serodiagnosis of infectious diseases. *Biosensors & Bioelectronics*; 80; 6-2016; 24-33. <https://doi.org/10.1016/j.bios.2016.01.021>.

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**CHEMLIS®**  
**Chagas R-iELISA**

Kit components / Presentations

Components	Description	REF: 15-CL01-2P Presentation 2 microplates	REF: 15-CL01-4P Presentation 4 microplates
Microplate	96-well microtitre plate (12 x 8 strip wells with strip holder) coated with the antigens, sealed and stored dry	2 microplates (192 tests)	4 microplates (384 tests)
Conjugate	Lyophilized (horseradish peroxidase conjugated anti-bovine IgG antibodies)	2 units	4 units
Substrate-chromogen solution	Ready to use [3,3',5,5'-Tetramethylbenzidine in substrate buffer containing hydrogen peroxide (H <sub>2</sub> O <sub>2</sub> )]- STORE IN THE DARK	1 x 30 ml	1 x 60 ml
Stop solution	Ready to use (contains HCl 1%) - CORROSIVE	1 x 30 ml	1 x 60 ml
Wash solution	10X concentrate	1 x 120 ml, 10X (sufficient for 1.2 liters)	1 x 120 ml, 10X (sufficient for 1.2 liters)
Sample diluent	To reconstitute - CONTAINS PRESERVATIVE [Thimerosal (C <sub>9</sub> H <sub>9</sub> HgNaO <sub>2</sub> S) 0.01%].	1 x 120 ml	2 x 120 ml
Positive control	Inactivated serum – CONTAINS PRESERVATIVE [Thimerosal (C <sub>9</sub> H <sub>9</sub> HgNaO <sub>2</sub> S) 0.01%]	1 x 0.1 ml	1 x 0.1 ml
Negative control	Inactivated serum – CONTAINS PRESERVATIVE [Thimerosal (C <sub>9</sub> H <sub>9</sub> HgNaO <sub>2</sub> S) 0.01%]	1 x 0.1 ml	1 x 0.1 ml
Polyethylene bag	With hermetic closure. Reusable.	1 unit	1 unit
Instruction manual		1 unit	1 unit

## Name and application

### CHEMLIS® Chagas R-iELISA

Indirect ELISA kit for the detection of antibodies against *Trypanosoma cruzi* in human serum and plasma.

CHEMLIS® Chagas R-iELISA is an indirect ELISA kit (Enzyme Linked Immunosorbent Assay) for the detection of IgG specific antibodies against *Trypanosoma cruzi*, causative agent of Chagas disease, in human serum samples or plasma obtained with heparin, sodium citrate or EDTA. Thanks to its exclusive combination of four recombinant antigens, the kit has an excellent diagnostic performance. CHEMLIS® Chagas R-iELISA can be used as a screening and/or a confirmatory test.

### Introduction

American Trypanosomiasis or Chagas disease is a major health and economic problem in Latin America caused by the protozoan parasite *Trypanosoma cruzi*. This parasite is found throughout the American continent in a variety of wild and domestic mammalian reservoirs and transmitted by the bite of infected, hematophagous triatomine insects. Apart from vectorial transmission, humans can become infected with *T. cruzi* by ingestion of contaminated food and fluids, from mother to child during pregnancy and through contaminated blood transfusion or organ transplantation. It is estimated that 10-12 million people are currently infected and that up to 120 million individuals living in endemic areas are at risk of infection, and thus at risk of developing cardiac or gut pathology normally associated with chronic Chagas Disease. Increasing travel and immigration has also brought the risk of *T. cruzi* infection into non-endemic countries such as Spain, Australia and the U.S., where some autochthonous cases of transmis-

sion have already been reported. Diagnosis of Chagas disease is challenging because it is often asymptomatic during the acute phase, and evolves into a chronic stage that has different clinical forms. In addition, and due to a major decline in parasitemia during the chronic phase of the disease, *T. cruzi* detection in blood samples by direct examination, blood-culture or xenodiagnosis is difficult as well as labor- and time-consuming. Therefore, detection of anti-*T. cruzi* antibodies in body fluids (typically blood samples) is still the most effective method for demonstrating direct exposure to the parasite. At present, the most widely used serologic methods are indirect hemagglutination assay (IHA), indirect immunofluorescence (IIF), and enzyme-linked immunosorbent assay (ELISA) using total parasite homogenates or semi-purified antigenic fractions. In spite of their simplicity and low cost, these tests show variations in their reproducibility and reliability that can be attributed to the poor standardization of the antigens. The advent of recombinant DNA and peptide synthesis technologies allowed the production and one-step purification of large amounts of highly pure *T. cruzi* immunodominant antigens. An additional advantage of both recombinant antigens and synthetic peptides is that they minimize specificity problems, one of the major drawbacks of immunodiagnosis of Chagas disease.

### Principles of the technique

CHEMLIS® Chagas R-iELISA kit is a solid phase enzyme immunoassay for the detection IgG antibodies against *Trypanosoma cruzi* in human serum or plasma samples. In this assay the antibodies present in the sample react with the recombinant antigens (exclusive combination of four recombinant proteins) that coat the wells. If the sample contains antibodies against these *Trypanosoma cruzi* proteins, they bind to the well. Upon addition of the horse-

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radish peroxidase (HRP) conjugated antibody, a complex is formed with the IgG antibodies attached to the antigen. Free antibodies are removed by washing before adding the substrate-chromogen solution. The blue color that appears is the product of the reaction of the substrate-chromogen with the conjugate and color intensity is directly proportional to the amount of antibodies present in the sample reacting with the antigens. The higher the antibody titer present in the sample, the more intense color develops in the test wells. Adding the stop solution, upon which a color change to yellow is observed, stops the reaction. The results are measured in a microplate spectrophotometer determining the absorbance (Abs) at 450 nm.

### Storage and expiration

Storage: between 2 and 8°C.

Transport temperature: 4 to 15°C for 72 hours.

Proper storage of the kit and all reagents will allow the use until the expiration date.

Shelf life: 12 months. Do not use reagents from different lots or after the expiration date.

### Materials needed but not provided

- High precision multichannel pipette and micropipettes (10, 100 and 1000 µl).
- Containers for multichannel pipette.
- Disposable pipette tips.
- Microplate reader for absorbance determination at 450 nm.
- Orbital rotator.
- Timer.
- Tubes for the dilutions of the controls and samples.
- Container (graduated cylinder) for the preparation of the wash solution.
- Distilled water (Type II) or ultrapure water (Type I).
- Microplate covers (lid, aluminium foil or adhesive).


### Precautions

1. Carefully read and strictly follow all instructions detailed in the leaflet included in the kit.
2. Store the kit and all reagents at 2-8 °C.
3. All reagents should equilibrate at room temperature (20-25°C) before use, and return to 2-8 °C following use.
4. Handle all materials and reagents according to the Good Laboratory Practices.
5. The samples, as well as positive and negative controls, should be considered potentially infectious and, therefore, the kit components that contact them should be handled with gloves and disposed as pathological waste according to legal regulations. The stop solution contains hydrochloric acid (corrosive). Skin contact generates irritation and should be handled with gloves. The rest of the kit components does not represent any potential risk for health and the environment.
6. Do not use the kit beyond expiration date or if the kit components were not kept in the conditions indicated above.
7. Do not mix components or instruction manuals from different test kits batches.
8. Care should be taken to avoid contamination of kit components.
9. If multichannel pipette containers are reused, they should be washed thoroughly with high quality water for reuse and identified to be used with the same solution.
10. Use different pipettes tips for each reagent and for each sample.
11. Include positive and negative controls, and the target sample in duplicate in each test.
12. For reconstitution of the reagents use only distilled water (Type II) or ultrapure water (Type I).

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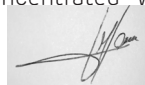should be brought to room temperature (20-25°C) and mixed to ensure dissolution of any precipitated salts.

The 10X concentrated Wash solution must be diluted 1/10 with distilled (Type II) or ultrapure (Type I) water; for example, adding 30 ml of 10X Wash solution to 270 ml of water (sufficient quantity of Wash solution 1X to process one microplate: 13 ml to reconstitute the IgG conjugate, 120 ml to reconstitute the Sample diluent and the rest for the washes). Mix very well before using. Once prepared, the 1X Wash solution can be stored at 2-8 °C for up to 7 days

It is recommended to prepare sufficient quantity to process the plates of the day.

#### b. Conjugate

Reconstitute the conjugate by adding 13 ml (per flask) of 1X Wash solution previously prepared (see previous item). Add the Washing solution carefully. Allow to stand for a minute and mix well. Prepare immediately before use.

For greater accuracy it is recommended to pipet the 13 ml using two graduation marks of the pipette (for example: pipetting twice from zero to 6.5 ml).

The remaining conjugate can be aliquoted in eppendorf type tubes and stored at -20°C (cannot be stored in the refrigerator) for up to 3 months. During this period of time each aliquot can be thaw and re-frozen one time.

#### c. Sample diluent

Reconstitute the Sample diluent by adding the 1X Wash solution previously prepared (see previous item) up to the level indicated in the corresponding bottle (120 ml) and mix very well. Prepare immediately before use.

The remaining Sample diluent can be stored in the refrigerator (2 to 8°C) up to 3 months or at -20°C for a year. Thaw it completely, bring to room temperature

and mix very well until complete dissolution before using.

### Preparation of controls

#### a. Positive (PC) and negative (NC) controls

Controls (PC and NC) should be diluted 1:200 using the Sample diluent.

Add 2 µl of the corresponding control to 400 µl of Sample diluent (dilution 1:200) and mix well. Use 100 µl of each diluted control per well and in duplicate.

#### b. Sample blank (SB)

On each plate a sample blank consisting of Sample diluent should be added. Dispense directly and in duplicate 100 µl per well of Sample diluent.

### Sample preparation

Samples (serum or plasma) should be diluted 1:200 using the Sample diluent.

Add 2 µl of the sample (serum or plasma) to 400 µl of Sample diluent (dilution 1:200) and mix well. Use 100 µl of the diluted sample per well.

#### Notes:

-It is recommended to analyze serum or plasma (obtained using heparin, sodium citrate or EDTA as anticoagulant) samples recently obtained (fresh) and free of turbidity.

-If the samples contain particles or a high content of triglycerides, they should be clarified by centrifugation.


-If the samples are not analyzed immediately, they can be kept refrigerated (2 to 8 °C) for up to 72 hs. For longer storage, store at -20°C (-15 to -20°C), taking into account that during storage and as a consequence of the freezing and thawing process, the antibody titer may change. If samples are frozen, they must be completely thawed, homogenized and centrifuged before testing.

-It is recommended to analyze the samples at least in duplicate.

-Samples contaminated and/or not well conserved may produce erroneous results.

### Test procedure

1. Allow the reagents equilibrate to room temperature (20-25°C). Do not remove the

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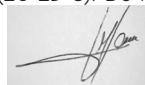plate from its pouch until it has reached room temperature. If the entire plate is not used, separate only the 8-well strips needed to analyze the controls and samples. Store the rest of the strips, together with the desiccants, in the polyethylene bag with hermetic seal included in the kit and store it again in the refrigerator (2-8°C).

Note: once the original packaging of the plate has been opened, the strips can be stored in the hermetically sealed polyethylene bag with the desiccants, and in the refrigerator (2-8°C) for up to 4 months, as long as the indicated expiration date is not exceeded.

2. Add in duplicate 100 µl per well of the positive control, negative control and sample blank. See Item Preparation of controls.

3. Add 100 µl per well of each sample. See Item Sample preparation.

4. Cover the plate and incubate at room temperature (20-25°C) with gentle orbital shaking for 1 hour.

5. Wash/rinse the plates 4 times with 200 µl of 1X Wash Solution per well. See item Reagents preparation.

Note: remove the liquid completely after each washing by dump and tapping onto absorbent material or aspiration with pipette. Avoid plate drying between plate washings and before adding the next reagent.

6. Add 100 µl of Conjugate to each well. See item Preparation of the conjugate. Cover the plate and incubate at room temperature (20-25°C) with gentle orbital shaking for 1 hour.

7. Repeat step 5.

8. Add 100 µl of Substrate-chromogen solution to each well, and incubated with gentle orbital shaking for 10 minutes (± 1 min) at room temperature (20-25°C) and in the dark (for example, cover the plate with aluminum foil). Start counting the time after filling the first well.

9. Stop the reaction by adding 100 µl of Stop solution into each well and mix thoroughly by lightly tapping the edges of the plate. This solution must be added in the same order and speed as the Substrate-chromogen solution (see Step 8).

10. Measure the absorbance at 450 nm using a plate spectrophotometer. There is no need to perform a blank reading. Perform the reading within 15 minutes after adding the stop solution.

Note: samples with high reactivity values can lead to the appearance of a black precipitate due to the precipitation of the substrate.

## Calculations

### Expression of results

Results are expressed as the percentage of reactivity ( $PR_{PC} = 100\%$ ) with respect to the positive control included in each test. The absorbance values measured at 450 nm ( $Abs_{450}$ ) of each sample are related to the  $Abs_{450}$  value of the positive control (PC) as follows:

$$PR_{Sample} = \frac{Abs_{450} \text{ Sample}^*}{Average Abs_{450} \text{ PC}} \times 100$$


\* If the sample was run in duplicate, consider the average of the  $Abs_{450}$  value.

### Validity criteria

To confirm the validity of the test the following criteria must be met:

- The average value of the PC ( $Abs_{450}$ ) must be higher than 1.6 ( $Abs_{450} \text{ PC} > 1.6$ ).
- The  $Abs_{450}$  values of the controls (PC and NC) and samples duplicates must not differ by more than 20% of the corresponding average.
- The  $Abs_{450}$  values of the blank sample must be less than 0.1 ( $Abs_{450} \text{ Blank sample} < 0.1$ ).
- The average  $Abs_{450}$  value of the NC must be less than 0.15 (Average  $Abs_{450} \text{ NC} < 0.15$ ).
- The Average  $Abs_{450} \text{ PC} / \text{Average } Abs_{450} \text{ NC}$  ratio must be higher than 10.5 (Average  $Abs_{450} \text{ PC} / \text{Average } Abs_{450} \text{ NC} > 10.5$ ).

If any of these criteria is not fulfilled, the test is considered invalid. In this case, the test should be repeated strictly fol-

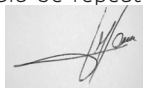lowing the instructions detailed in the Instruction manual included in the kit.

### Interpretation of results

The test sample results should be interpreted as follows:

Sample (dilution 1:200)

PR	Interpretation
≤ 21 %	Non-Reactive
> 21 ≤ 25,5 %	Indeterminate
> 25,5 %	Reactive

In the case of obtaining an "Indeterminate" result, the sample should be re-tested. If the result is still indeterminate, it is recommended to test a second sample obtained after a period of at least three weeks. If the new sample is indeterminate again, the epidemiological situation should be considered and the sample should be analyzed with another technique if available.

### Performance characteristics

#### *Diagnostic sensitivity and specificity*

More than 390 serum samples obtained from patients with acute or chronic Chagas disease and from healthy donors were analyzed to determine the diagnostic sensitivity and specificity of the assay. Based on these results, a ROC (receiver-operating-analysis) curve analysis was performed and it was determined that taking a cut-off value of 23%, CHEMLIS® Chagas R-iELISA has a diagnostic sensitivity of 100% and a specificity of 100%.

#### *Hook effect*

Based on the analysis of samples with high titers of specific anti-T. cruzi antibodies, it was determined that the CHEMLIS® Chagas R-iELISA kit does not have a hook effect in the detection of IgG antibodies.

#### *Cross reactivity*

Cross-reactivity of CHEMLIS® Chagas R-iELISA was evaluated by testing clinical

samples obtained from patients with other conditions. Samples with positive serology obtained from patients with the following diseases or bacterial, parasitic, and viral infections were analyzed:

-Bacterial: Escherichia coli O157, O145 and O121 infections, Human Brucellosis, Tuberculosis, Syphilis, Streptococcus pneumoniae, Klebsiella pneumoniae and Acinetobacter baumannii infections.

-Viral: Dengue, VHB, VHC, HTLV I/II, HIV, Hantavirus.

-Parasitic: Leishmaniasis y Toxoplasmosis.

-Autoimmune diseases: rheumatoid arthritis, Systemic Lupus Erythematosus and autoimmune vasculitis.

CHEMLIS® Chagas R-iELISA has a high level of specificity since no cross-reactivity was observed with the exception of a sample from a patient with tuberculosis and another patient with Dengue.

#### *Interfering substances*

The following potentially interfering substances have no impact on the performance of CHEMLIS® Chagas R-iELISA (final test concentrations of interfering substances are indicated in parentheses):

-Oxalic acid (0.6 mg/ml)

-Bovine albumin (20 mg/ml)

-Human Immunoglobulin G (hIgG) (2 mg/ml)

-Human Immunoglobulin M (hIgG) (2 mg/ml)

-Human Immunoglobulin A (hIgG) (2 mg/ml)

-Cholesterol (5 mg/ml)

-Hemoglobin (10 mg/ml)

-Bilirubin (2 mg/ml)

#### *Repeatability and Reproducibility*

Repeatability and reproducibility of CHEMLIS® Chagas R-iELISA was established using in-house reference panels containing negative specimens and a range of positive specimens. There were no differences observed within-run, between-run, between-lots, between-sites, and between-days.

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Electrochemical magnetic microbeads-based biosensor for point-of-care serodiagnosis of infectious diseases.

Biosensors & Bioelectronics; 80; 6-2016; 24-33. <https://doi.org/10.1016/j.bios.2016.01.021>.

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Digitally signed by Gestion Documental Electronica  
Date: 2022.10.25 10:00:37 -03:00

Digitally signed by Gestion Documental  
Electronica  
Date: 2022.10.25 10:00:38 -03:00