



República Argentina - Poder Ejecutivo Nacional
2019 - Año de la Exportación

Disposición

Número:

Referencia: 1-47-3110-1051-19-0

VISTO el Expediente N° 1-47-3110-1051-19-0 del Registro de esta Administración Nacional de Medicamentos, Alimentos y Tecnología Médica (ANMAT), y

CONSIDERANDO:

Que por las presentes actuaciones BIOIMPLANT S.R.L., solicita se autorice la inscripción en el Registro Productores y Productos de Tecnología Médica (RPPTM) de esta Administración Nacional, de un nuevo producto médico.

Que las actividades de elaboración y comercialización de productos médicos se encuentran contempladas por la Ley 16463, el Decreto 9763/64, y MERCOSUR/GMC/RES. N° 40/00, incorporada al ordenamiento jurídico nacional por Disposición ANMAT N° 2318/02 (TO 2004), y normas complementarias.

Que consta la evaluación técnica producida por la Dirección Nacional de Productos Médicos, en la que informa que el producto estudiado reúne los requisitos técnicos que contempla la norma legal vigente, y que los establecimientos declarados demuestran aptitud para la elaboración y el control de calidad del producto cuya inscripción en el Registro se solicita.

Que corresponde autorizar la inscripción en el RPPTM del producto médico objeto de la solicitud.

Que se actúa en virtud de las facultades conferidas por el Decreto N° 1490/92 y sus modificatorios.

Por ello;

EL ADMINISTRADOR NACIONAL DE LA ADMINISTRACIÓN NACIONAL DE
MEDICAMENTOS, ALIMENTOS Y TECNOLOGÍA MÉDICA

DISPONE:

ARTÍCULO 1º.- Autorízase la inscripción en el Registro Nacional de Productores y Productos de Tecnología Médica (RPPTM) de la Administración Nacional de Medicamentos, Alimentos y Tecnología Médica (ANMAT) del producto médico marca DOUBLE MEDICAL nombre descriptivo PLACAS BLOQUEADAS DE TITANIO y nombre técnico Placas, para Huesos, de acuerdo con lo solicitado por BIOIMPLANT S.R.L., con los Datos Identificatorios Característicos que figuran al pie de la presente.

ARTÍCULO 2º.- Autorízase los textos de los proyectos de rótulo/s y de instrucciones de uso que obran en documento N° IF-2019-65396542-APN-DNPM#ANMAT.

ARTÍCULO 3º.- En los rótulos e instrucciones de uso autorizados deberá figurar la leyenda “Autorizado por la ANMAT PM-2342-01”, con exclusión de toda otra leyenda no contemplada en la normativa vigente.

ARTICULO 4º.- Extiéndase el Certificado de Autorización e Inscripción en el RPPTM con los datos característicos mencionados en esta disposición.

ARTÍCULO 5º.- La vigencia del Certificado de Autorización será de cinco (5) años, a partir de la fecha de la presente disposición.

ARTÍCULO 6º.- Regístrese. Inscríbase en el Registro Nacional de Productores y Productos de Tecnología Médica al nuevo producto. Por el Departamento de Mesa de Entrada, notifíquese al interesado, haciéndole entrega de la presente Disposición, conjuntamente con rótulos e instrucciones de uso autorizados y el Certificado mencionado en el artículo 4º. Gírese a la Dirección de Gestión de Información Técnica a los fines de confeccionar el legajo correspondiente. Cumplido, archívese.

DATOS IDENTIFICATORIOS CARACTERISTICOS

Nombre descriptivo: PLACAS BLOQUEADAS DE TITANIO

Código de identificación y nombre técnico UMDNS: 13-050- Placas, para Huesos

Marca(s) de (los) producto(s) médico(s): DOUBLE MEDICAL

Clase de Riesgo: III

Indicación/es autorizada/s: para tratar fracturas de huesos grandes (clavícula, escapula, pelvis, húmero, cúbito, radio, fémur, tibia y peroné) y pequeños (metacarpo, metatarso y falanges).

Modelo/s:

Nº Referencia	Descripción Perforac.
020140003	T-Locking Plate, Combi-holes 3 Holes
020140004	T-Locking Plate, Combi-holes 4 Holes
020140005	T-Locking Plate, Combi-holes 5 Holes
020140006	T-Locking Plate, Combi-holes 6 Holes
020140007	T-Locking Plate, Combi-holes 7 Holes

020140008 T-Locking Plate, Combi-holes 8 Holes
020150004 T-Buttress Locking Plate, Combi-holes 4 Holes
020150005 T-Buttress Locking Plate, Combi-holes 5 Holes
020150006 T-Buttress Locking Plate, Combi-holes 6 Holes
020150003 T-Buttress Locking Plate, Combi-holes 3 Holes
020150007 T-Buttress Locking Plate, Combi-holes 7 Holes
020150008 T-Buttress Locking Plate, Combi-holes 8 Holes
020161003 L-Buttress Locking Plate left 3 Holes
020161004 L-Buttress Locking Plate left 4 Holes
020161005 L-Buttress Locking Plate left 5 Holes
020161006 L-Buttress Locking Plate left 6 Holes
020161007 L-Buttress Locking Plate left 7 Holes
020161008 L-Buttress Locking Plate left 8 Holes
020162003 L-Buttress Locking Plate, right 3 Holes
020162004 L-Buttress Locking Plate, right 4 Holes
020162005 L-Buttress Locking Plate, right 5 Holes
020162006 L-Buttress Locking Plate, right 6 Holes
020162007 L-Buttress Locking Plate, right 7 Holes
020162008 L-Buttress Locking Plate, right 8 Holes
020171004 Distal Femoral Sub-Condylar Locking Plate, Left 4 Holes
020171005 Distal Femoral Sub-Condylar Locking Plate, Left 5 Holes
020171006 Distal Femoral Sub-Condylar Locking Plate, Left 6 Holes
020171007 Distal Femoral Sub-Condylar Locking Plate, Left 7 Holes
020171008 Distal Femoral Sub-Condylar Locking Plate, Left 8 Holes
020171009 Distal Femoral Sub-Condylar Locking Plate, Left 9 Holes
020171010 Distal Femoral Sub-Condylar Locking Plate, Left 10 Holes
020171011 Distal Femoral Sub-Condylar Locking Plate, Left 11 Holes
020171012 Distal Femoral Sub-Condylar Locking Plate, Left 12 Holes

020172004 Distal Femoral Sub-Condylar Locking Plate, Right 4 Holes
020172005 Distal Femoral Sub-Condylar Locking Plate, Right 5 Holes
020172006 Distal Femoral Sub-Condylar Locking Plate, Right 6 Holes
020172007 Distal Femoral Sub-Condylar Locking Plate, Right 7 Holes
020172008 Distal Femoral Sub-Condylar Locking Plate, Right 8 Holes
020172009 Distal Femoral Sub-Condylar Locking Plate, Right 9 Holes
020172010 Distal Femoral Sub-Condylar Locking Plate, Right 10 Holes
020172011 Distal Femoral Sub-Condylar Locking Plate, Right 11 Holes
020172012 Distal Femoral Sub-Condylar Locking Plate, Right 12 Holes
020211060 Calcaneal Locking Plate , Left 8 Holes
020211075 Calcaneal Locking Plate , Left 10 Holes
020212060 Calcaneal Locking Plate , Right 8 Holes
020212075 Calcaneal Locking Plate , Right 10 Holes
020221015 Calcaneal Locking Plate , Left 15 Holes
020222015 Calcaneal Locking Plate , Right 15 Holes
020231015 Calcaneal Locking Plate , Left 15 Holes
020232015 Calcaneal Locking Plate , Right 15 Holes
020251003 Proximal Femoral Locking Plate, Multi-holes, Left 3 Holes
020251004 Proximal Femoral Locking Plate, Multi-holes, Left 4 Holes
020251005 Proximal Femoral Locking Plate, Multi-holes, Left 5 Holes
020251006 Proximal Femoral Locking Plate, Multi-holes, Left 6 Holes
020251007 Proximal Femoral Locking Plate, Multi-holes, Left 7 Holes
020251008 Proximal Femoral Locking Plate, Multi-holes, Left 8 Holes
020251009 Proximal Femoral Locking Plate, Multi-holes, Left 9 Holes
020251010 Proximal Femoral Locking Plate, Multi-holes, Left 10 Holes
020251011 Proximal Femoral Locking Plate, Multi-holes, Left 11 Holes
020251012 Proximal Femoral Locking Plate, Multi-holes, Left 12 Holes
020252003 Proximal Femoral Locking Plate, Multi-holes, Right 3 Holes
020252004 Proximal Femoral Locking Plate, Multi-holes, Right 4 Holes

020252005 Proximal Femoral Locking Plate, Multi-holes, Right 5 Holes

020252006 Proximal Femoral Locking Plate, Multi-holes, Right 6 Holes

020252007 Proximal Femoral Locking Plate, Multi-holes, Right 7 Holes

020252008 Proximal Femoral Locking Plate, Multi-holes, Right 8 Holes

020252009 Proximal Femoral Locking Plate, Multi-holes, Right 9 Holes

020252010 Proximal Femoral Locking Plate, Multi-holes, Right 10 Holes

020252011 Proximal Femoral Locking Plate, Multi-holes, Right 11 Holes

020252012 Proximal Femoral Locking Plate, Multi-holes, Right 12 Holes

020280004 Pelvic Reconstruction Locking Plate 4 Holes

020280005 Pelvic Reconstruction Locking Plate 5 Holes

020280006 Pelvic Reconstruction Locking Plate 6 Holes

020280007 Pelvic Reconstruction Locking Plate 7 Holes

020280008 Pelvic Reconstruction Locking Plate 8 Holes

020280009 Pelvic Reconstruction Locking Plate 9 Holes

020280010 Pelvic Reconstruction Locking Plate 10 Holes

020280011 Pelvic Reconstruction Locking Plate 11 Holes

020280012 Pelvic Reconstruction Locking Plate 12 Holes

020280013 Pelvic Reconstruction Locking Plate 13 Holes

020280014 Pelvic Reconstruction Locking Plate 14 Holes

020320000 3.5mm Locking Spacer

020330000 5.0mm Locking Spacer

020360000 3.5mm Cable Spacer

020370000 5.0mm Cable Spacer

020380003 Proximal Humeral Locking Plate , Combi-holes 3 Holes

020380004 Proximal Humeral Locking Plate , Combi-holes 4 Holes

020380005 Proximal Humeral Locking Plate , Combi-holes 5 Holes

020380006 Proximal Humeral Locking Plate , Combi-holes 6 Holes

020380007 Proximal Humeral Locking Plate , Combi-holes 7 Holes

020380008 Proximal Humeral Locking Plate , Combi-holes	8 Holes
020380009 Proximal Humeral Locking Plate , Combi-holes	9 Holes
020380010 Proximal Humeral Locking Plate , Combi-holes	10 Holes
020380011 Proximal Humeral Locking Plate , Combi-holes	11 Holes
020380012 Proximal Humeral Locking Plate , Combi-holes	12 Holes
020380013 Proximal Humeral Locking Plate , Combi-holes	13 Holes
020400006 Broad Limited Contact Locking Plate, Combi-holes	6 Holes
020400007 Broad Limited Contact Locking Plate, Combi-holes	7 Holes
020400008 Broad Limited Contact Locking Plate, Combi-holes	8 Holes
020400009 Broad Limited Contact Locking Plate, Combi-holes	9 Holes
020400010 Broad Limited Contact Locking Plate, Combi-holes	10 Holes
020400011 Broad Limited Contact Locking Plate, Combi-holes	11 Holes
020400012 Broad Limited Contact Locking Plate, Combi-holes	12 Holes
020400013 Broad Limited Contact Locking Plate, Combi-holes	13 Holes
020400014 Broad Limited Contact Locking Plate, Combi-holes	14 Holes
020400015 Broad Limited Contact Locking Plate, Combi-holes	15 Holes
020400016 Broad Limited Contact Locking Plate, Combi-holes	16 Holes
020410004 Narrow Limited Contact Locking Plate, Combi-holes	4 Holes
020410005 Narrow Limited Contact Locking Plate, Combi-holes	5 Holes
020410006 Narrow Limited Contact Locking Plate, Combi-holes	6 Holes
020410007 Narrow Limited Contact Locking Plate, Combi-holes	7 Holes
020410008 Narrow Limited Contact Locking Plate, Combi-holes	8 Holes
020410009 Narrow Limited Contact Locking Plate, Combi-holes	9 Holes
020410010 Narrow Limited Contact Locking Plate, Combi-holes	10 Holes
020410011 Narrow Limited Contact Locking Plate, Combi-holes	11 Holes
020410012 Narrow Limited Contact Locking Plate, Combi-holes	12 Holes
020410013 Narrow Limited Contact Locking Plate, Combi-holes	13 Holes
020410014 Narrow Limited Contact Locking Plate, Combi-holes	14 Holes
020410015 Narrow Limited Contact Locking Plate, Combi-holes	15 Holes

020410016	Narrow Limited Contact Locking Plate, Combi-holes	16 Holes
020421003	Proximal Femoral Locking Plate, Combi-holes, Left	3 Holes
020421004	Proximal Femoral Locking Plate, Combi-holes, Left	4 Holes
020421005	Proximal Femoral Locking Plate, Combi-holes, Left	5 Holes
020421006	Proximal Femoral Locking Plate, Combi-holes, Left	6 Holes
020421007	Proximal Femoral Locking Plate, Combi-holes, Left	7 Holes
020421008	Proximal Femoral Locking Plate, Combi-holes, Left	8 Holes
020421009	Proximal Femoral Locking Plate, Combi-holes, Left	9 Holes
020421010	Proximal Femoral Locking Plate, Combi-holes, Left	10 Holes
020421011	Proximal Femoral Locking Plate, Combi-holes, Left	11 Holes
020421012	Proximal Femoral Locking Plate, Combi-holes, Left	12 Holes
020421013	Proximal Femoral Locking Plate, Combi-holes, Left	13 Holes
020422003	Proximal Femoral Locking Plate, Combi-holes, Right	3 Holes
020422004	Proximal Femoral Locking Plate, Combi-holes, Right	4 Holes
020422005	Proximal Femoral Locking Plate, Combi-holes, Right	5 Holes
020422006	Proximal Femoral Locking Plate, Combi-holes, Right	6 Holes
020422007	Proximal Femoral Locking Plate, Combi-holes, Right	7 Holes
020422008	Proximal Femoral Locking Plate, Combi-holes, Right	8 Holes
020422009	Proximal Femoral Locking Plate, Combi-holes, Right	9 Holes
020422010	Proximal Femoral Locking Plate, Combi-holes, Right	10 Holes
020422011	Proximal Femoral Locking Plate, Combi-holes, Right	11 Holes
020422012	Proximal Femoral Locking Plate, Combi-holes, Right	12 Holes
020422013	Proximal Femoral Locking Plate, Combi-holes, Right	13 Holes
020451003	Distal Medial Tibial Locking Plate, Combi-holes, Left	3 Holes
020451004	Distal Medial Tibial Locking Plate, Combi-holes, Left	4 Holes
020451005	Distal Medial Tibial Locking Plate, Combi-holes, Left	5 Holes
020451006	Distal Medial Tibial Locking Plate, Combi-holes, Left	6 Holes
020451007	Distal Medial Tibial Locking Plate, Combi-holes, Left	7 Holes

020451008	Distal Medial Tibial Locking Plate, Combi-holes, Left	8 Holes
020451009	Distal Medial Tibial Locking Plate, Combi-holes, Left	9 Holes
020451010	Distal Medial Tibial Locking Plate, Combi-holes, Left	10 Holes
020451011	Distal Medial Tibial Locking Plate, Combi-holes, Left	11 Holes
020451012	Distal Medial Tibial Locking Plate, Combi-holes, Left	12 Holes
020451013	Distal Medial Tibial Locking Plate, Combi-holes, Left	13 Holes
020452003	Distal Medial Tibial Locking Plate, Combi-holes, Right	3 Holes
020452004	Distal Medial Tibial Locking Plate, Combi-holes, Right	4 Holes
020452005	Distal Medial Tibial Locking Plate, Combi-holes, Right	5 Holes
020452006	Distal Medial Tibial Locking Plate, Combi-holes, Right	6 Holes
020452007	Distal Medial Tibial Locking Plate, Combi-holes, Right	7 Holes
020452008	Distal Medial Tibial Locking Plate, Combi-holes, Right	8 Holes
020452009	Distal Medial Tibial Locking Plate, Combi-holes, Right	9 Holes
020452010	Distal Medial Tibial Locking Plate, Combi-holes, Right	10 Holes
020452011	Distal Medial Tibial Locking Plate, Combi-holes, Right	11 Holes
020452012	Distal Medial Tibial Locking Plate, Combi-holes, Right	12 Holes
020452013	Distal Medial Tibial Locking Plate, Combi-holes, Right	13 Holes
020470004	Small Limited Contact Locking Plate, Combi-holes	4 Holes
020470005	Small Limited Contact Locking Plate, Combi-holes	5 Holes
020470006	Small Limited Contact Locking Plate, Combi-holes	6 Holes
020470007	Small Limited Contact Locking Plate, Combi-holes	7 Holes
020470008	Small Limited Contact Locking Plate, Combi-holes	8 Holes
020470009	Small Limited Contact Locking Plate, Combi-holes	9 Holes
020470010	Small Limited Contact Locking Plate, Combi-holes	10 Holes
020470011	Small Limited Contact Locking Plate, Combi-holes	11 Holes
020470012	Small Limited Contact Locking Plate, Combi-holes	12 Holes
020470013	Small Limited Contact Locking Plate, Combi-holes	13 Holes
020470014	Small Limited Contact Locking Plate, Combi-holes	14 Holes
020480004	Reconstruction Locking Plate, combi-holes	4 Holes

020480005 Reconstruction Locking Plate, combi-holes	5 Holes
020480006 Reconstruction Locking Plate, combi-holes	6 Holes
020480007 Reconstruction Locking Plate, combi-holes	7 Holes
020480008 Reconstruction Locking Plate, combi-holes	8 Holes
020480009 Reconstruction Locking Plate, combi-holes	9 Holes
020480010 Reconstruction Locking Plate, combi-holes	10 Holes
020480011 Reconstruction Locking Plate, combi-holes	11 Holes
020480012 Reconstruction Locking Plate, combi-holes	12 Holes
020480013 Reconstruction Locking Plate, combi-holes	13 Holes
020480014 Reconstruction Locking Plate, combi-holes	14 Holes
020490004 One-third Tubular Locking Plate, Combi-holes	4 Holes
020490005 One-third Tubular Locking Plate, Combi-holes	5 Holes
020490006 One-third Tubular Locking Plate, Combi-holes	6 Holes
020490007 One-third Tubular Locking Plate, Combi-holes	7 Holes
020490008 One-third Tubular Locking Plate, Combi-holes	8 Holes
020490009 One-third Tubular Locking Plate, Combi-holes	9 Holes
020490010 One-third Tubular Locking Plate, Combi-holes	10 Holes
020490011 One-third Tubular Locking Plate, Combi-holes	11 Holes
020490012 One-third Tubular Locking Plate, Combi-holes	12 Holes
020501004 Clavicular Locking Plate, Combi-holes, Left	4 Holes
020501005 Clavicular Locking Plate, Combi-holes, Left	5 Holes
020501006 Clavicular Locking Plate, Combi-holes, Left	6 Holes
020501007 Clavicular Locking Plate, Combi-holes, Left	7 Holes
020501008 Clavicular Locking Plate, Combi-holes, Left	8 Holes
020501009 Clavicular Locking Plate, Combi-holes, Left	9 Holes
020501010 Clavicular Locking Plate, Combi-holes, Left	10 Holes
020501011 Clavicular Locking Plate, Combi-holes, Left	11 Holes
020501012 Clavicular Locking Plate, Combi-holes, Left	12 Holes

020502004 Clavicular Locking Plate, Combi-holes, Right 4 Holes
020502005 Clavicular Locking Plate, Combi-holes, Right 5 Holes
020502006 Clavicular Locking Plate, Combi-holes, Right 6 Holes
020502007 Clavicular Locking Plate, Combi-holes, Right 7 Holes
020502008 Clavicular Locking Plate, Combi-holes, Right 8 Holes
020502009 Clavicular Locking Plate, Combi-holes, Right 9 Holes
020502010 Clavicular Locking Plate, Combi-holes, Right 10 Holes
020502011 Clavicular Locking Plate, Combi-holes, Right 11 Holes
020502012 Clavicular Locking Plate, Combi-holes, Right 12 Holes
020511002 Distal Fibular Locking Plate, Combi-holes, Left 2 Holes
020511003 Distal Fibular Locking Plate, Combi-holes, Left 3 Holes
020511004 Distal Fibular Locking Plate, Combi-holes, Left 4 Holes
020511005 Distal Fibular Locking Plate, Combi-holes, Left 5 Holes
020511006 Distal Fibular Locking Plate, Combi-holes, Left 6 Holes
020511007 Distal Fibular Locking Plate, Combi-holes, Left 7 Holes
020511008 Distal Fibular Locking Plate, Combi-holes, Left 8 Holes
020511009 Distal Fibular Locking Plate, Combi-holes, Left 9 Holes
020512002 Distal Fibular Locking Plate, Combi-holes, Right 2 Holes
020512003 Distal Fibular Locking Plate, Combi-holes, Right 3 Holes
020512004 Distal Fibular Locking Plate, Combi-holes, Right 4 Holes
020512005 Distal Fibular Locking Plate, Combi-holes, Right 5 Holes
020512006 Distal Fibular Locking Plate, Combi-holes, Right 6 Holes
020512007 Distal Fibular Locking Plate, Combi-holes, Right 7 Holes
020512008 Distal Fibular Locking Plate, Combi-holes, Right 8 Holes
020512009 Distal Fibular Locking Plate, Combi-holes, Right 9 Holes
020520003 Olecranon Locking Plate , Combi-holes 3 Holes
020520004 Olecranon Locking Plate , Combi-holes 4 Holes
020520005 Olecranon Locking Plate , Combi-holes 5 Holes
020520006 Olecranon Locking Plate , Combi-holes 6 Holes

020520007 Olecranon Locking Plate , Combi-holes 7 Holes
020520008 Olecranon Locking Plate , Combi-holes 8 Holes
020541000 Proximal Femoral Locking Plate for 7.3mm Cannulated Locking Screw, Left 51mm
020542000 Proximal Femoral Locking Plate for 7.3mm Cannulated Locking Screw, Right 51mm
020621004 Distal Humeral Sub-Condylar Locking Plate, Left 4 Holes
020621007 Distal Humeral Sub-Condylar Locking Plate, Left 7 Holes
020621010 Distal Humeral Sub-Condylar Locking Plate, Left 10 Holes
020622004 Distal Humeral Sub-Condylar Locking Plate, Right 4 Holes
020622007 Distal Humeral Sub-Condylar Locking Plate, Right 7 Holes
020622010 Distal Humeral Sub-Condylar Locking Plate, Right 10 Holes
020631002 Distal Radius Locking Plate , Small, Combi-holes, Left 2 Holes
020631003 Distal Radius Locking Plate , Small, Combi-holes, Left 3 Holes
020631004 Distal Radius Locking Plate , Small, Combi-holes, Left 4 Holes
020631005 Distal Radius Locking Plate , Small, Combi-holes, Left 5 Holes
020632002 Distal Radius Locking Plate , Small, Combi-holes, Right 2 Holes
020632003 Distal Radius Locking Plate , Small, Combi-holes, Right 3 Holes
020632004 Distal Radius Locking Plate , Small, Combi-holes, Right 4 Holes
020632005 Distal Radius Locking Plate , Small, Combi-holes, Right 5 Holes
020641002 Distal Radius Locking Plate , Large, Combi-holes, Left 2 Holes
020641003 Distal Radius Locking Plate , Large, Combi-holes, Left 3 Holes
020641004 Distal Radius Locking Plate , Large, Combi-holes, Left 4 Holes
020641005 Distal Radius Locking Plate , Large, Combi-holes, Left 5 Holes
020642002 Distal Radius Locking Plate , Large, Combi-holes, Right 2 Holes
020642003 Distal Radius Locking Plate , Large, Combi-holes, Right 3 Holes
020642004 Distal Radius Locking Plate , Large, Combi-holes, Right 4 Holes
020642005 Distal Radius Locking Plate , Large, Combi-holes, Right 5 Holes
020651002 Proximal Radius Locking Plate, left 2 Holes
020651003 Proximal Radius Locking Plate, left 3 Holes

020651004 Proximal Radius Locking Plate, left 4 Holes
020652002 Proximal Radius Locking Plate, right 2 Holes
020652003 Proximal Radius Locking Plate, right 3 Holes
020652004 Proximal Radius Locking Plate, right 4 Holes
020661003 Distal Clavicular Locking Plate, Combi-holes, Left 3 Holes
020661004 Distal Clavicular Locking Plate, Combi-holes, Left 4 Holes
020661005 Distal Clavicular Locking Plate, Combi-holes, Left 5 Holes
020661006 Distal Clavicular Locking Plate, Combi-holes, Left 6 Holes
020661007 Distal Clavicular Locking Plate, Combi-holes, Left 7 Holes
020661008 Distal Clavicular Locking Plate, Combi-holes, Left 8 Holes
020662003 Distal Clavicular Locking Plate, Combi-holes, Right 3 Holes
020662004 Distal Clavicular Locking Plate, Combi-holes, Right 4 Holes
020662005 Distal Clavicular Locking Plate, Combi-holes, Right 5 Holes
020662006 Distal Clavicular Locking Plate, Combi-holes, Right 6 Holes
020662007 Distal Clavicular Locking Plate, Combi-holes, Right 7 Holes
020662008 Distal Clavicular Locking Plate, Combi-holes, Right 8 Holes
020811003 Clavicular Hook Locking Plate , Height 15mm, Combi-holes, Left 3 Holes
020811004 Clavicular Hook Locking Plate , Height 15mm, Combi-holes, Left 4 Holes
020811005 Clavicular Hook Locking Plate , Height 15mm, Combi-holes, Left 5 Holes
020811006 Clavicular Hook Locking Plate , Height 15mm, Combi-holes, Left 6 Holes
020811007 Clavicular Hook Locking Plate , Height 15mm, Combi-holes, Left 7 Holes
020811008 Clavicular Hook Locking Plate , Height 15mm, Combi-holes, Left 8 Holes
020812003 Clavicular Hook Locking Plate , Height 15mm, Combi-holes, Right 3 Holes
020812004 Clavicular Hook Locking Plate , Height 15mm, Combi-holes, Right 4 Holes
020812005 Clavicular Hook Locking Plate , Height 15mm, Combi-holes, Right 5 Holes
020812006 Clavicular Hook Locking Plate , Height 15mm, Combi-holes, Right 6 Holes
020812007 Clavicular Hook Locking Plate , Height 15mm, Combi-holes, Right 7 Holes
020812008 Clavicular Hook Locking Plate , Height 15mm, Combi-holes, Right 8 Holes
020821003 Clavicular Hook Locking Plate , Height 18mm, Combi-holes, Left 3 Holes

020821004 Clavicular Hook Locking Plate , Height 18mm, Combi-holes, Left 4 Holes

020821005 Clavicular Hook Locking Plate , Height 18mm, Combi-holes, Left 5 Holes

020821006 Clavicular Hook Locking Plate , Height 18mm, Combi-holes, Left 6 Holes

020821007 Clavicular Hook Locking Plate , Height 18mm, Combi-holes, Left 7 Holes

020821008 Clavicular Hook Locking Plate , Height 18mm, Combi-holes, Left 8 Holes

020822003 Clavicular Hook Locking Plate , Height 18mm, Combi-holes, Right 3 Holes

020822004 Clavicular Hook Locking Plate , Height 18mm, Combi-holes, Right 4 Holes

020822005 Clavicular Hook Locking Plate , Height 18mm, Combi-holes, Right 5 Holes

020822006 Clavicular Hook Locking Plate , Height 18mm, Combi-holes, Right 6 Holes

020822007 Clavicular Hook Locking Plate , Height 18mm, Combi-holes, Right 7 Holes

020822008 Clavicular Hook Locking Plate , Height 18mm, Combi-holes, Right 8 Holes

020961005 Proximal Humeral Greater Tuberosity Locking Compression Plate, left 5 Holes

020962005 Proximal Humeral Greater Tuberosity Locking Compression Plate, right 5 Holes

020971000 Proximal Femoral Locking Plate, for 7.3mm Cannulated Locking Screw, small, left 43mm

020972000 Proximal Femoral Locking Plate, for 7.3mm Cannulated Locking Screw, small, right 43mm

021301003 Distal Lateral Femoral Condylar Locking Plate, left 3 Holes

021302003 Distal Lateral Femoral Condylar Locking Plate, right 3 Holes

021311003 Distal Medial Femoral Condylar Locking Plate, left 3 Holes

021312003 Distal Medial Femoral Condylar Locking Plate, right 3 Holes

021321004 Distal Medial Femoral Sub-condylar Locking Plate, large, left 4 Holes

021321007 Distal Medial Femoral Sub-condylar Locking Plate, large, left 7 Holes

021322004 Distal Medial Femoral Sub-condylar Locking Plate, large, right 4 Holes

021322007 Distal Medial Femoral Sub-condylar Locking Plate, large, right 7 Holes

021331004 Distal Medial Femoral Sub-condylar Locking Plate, small, left 4 Holes

021331007 Distal Medial Femoral Sub-condylar Locking Plate, small, left 7 Holes

021332004 Distal Medial Femoral Sub-condylar Locking Plate, small, right 4 Holes

021332007 Distal Medial Femoral Sub-condylar Locking Plate, small, right 7 Holes

021371003 Distal Lateral Humeral Locking Plate, Combi-holes, Left 3 Holes

021371004 Distal Lateral Humeral Locking Plate, Combi-holes, Left 4 Holes
021371005 Distal Lateral Humeral Locking Plate, Combi-holes, Left 5 Holes
021371006 Distal Lateral Humeral Locking Plate, Combi-holes, Left 6 Holes
021371007 Distal Lateral Humeral Locking Plate, Combi-holes, Left 7 Holes
021371008 Distal Lateral Humeral Locking Plate, Combi-holes, Left 8 Holes
021371009 Distal Lateral Humeral Locking Plate, Combi-holes, Left 9 Holes
021372003 Distal Lateral Humeral Locking Plate, Combi-holes, Right 3 Holes
021372004 Distal Lateral Humeral Locking Plate, Combi-holes, Right 4 Holes
021372005 Distal Lateral Humeral Locking Plate, Combi-holes, Right 5 Holes
021372006 Distal Lateral Humeral Locking Plate, Combi-holes, Right 6 Holes
021372007 Distal Lateral Humeral Locking Plate, Combi-holes, Right 7 Holes
021372008 Distal Lateral Humeral Locking Plate, Combi-holes, Right 8 Holes
021372009 Distal Lateral Humeral Locking Plate, Combi-holes, Right 9 Holes
021381003 Distal Medial Humeral Locking Plate, Combi-holes, Left 3 Holes
021381004 Distal Medial Humeral Locking Plate, Combi-holes, Left 4 Holes
021381005 Distal Medial Humeral Locking Plate, Combi-holes, Left 5 Holes
021381006 Distal Medial Humeral Locking Plate, Combi-holes, Left 6 Holes
021381007 Distal Medial Humeral Locking Plate, Combi-holes, Left 7 Holes
021381008 Distal Medial Humeral Locking Plate, Combi-holes, Left 8 Holes
021381009 Distal Medial Humeral Locking Plate, Combi-holes, Left 9 Holes
021382003 Distal Medial Humeral Locking Plate, Combi-holes, Right 3 Holes
021382004 Distal Medial Humeral Locking Plate, Combi-holes, Right 4 Holes
021382005 Distal Medial Humeral Locking Plate, Combi-holes, Right 5 Holes
021382006 Distal Medial Humeral Locking Plate, Combi-holes, Right 6 Holes
021382007 Distal Medial Humeral Locking Plate, Combi-holes, Right 7 Holes
021382008 Distal Medial Humeral Locking Plate, Combi-holes, Right 8 Holes
021382009 Distal Medial Humeral Locking Plate, Combi-holes, Right 9 Holes
021390003 T-Locking Plate, Right Angled, Head 3 Holes, Combi-holes 3 Holes

021400004 T-Locking Plate, Right Angled, Head 4 Holes, Combi-holes	4 Holes
021390005 T-Locking Plate, Right Angled, Head 3 Holes, Combi-holes	5 Holes
021400006 T-Locking Plate, Right Angled, Head 4 Holes, Combi-holes	6 Holes
021411003 T-Locking Plate, Oblique Angled, Combi-holes, Left 3 Holes	
021411004 T-Locking Plate, Oblique Angled, Combi-holes, Left 4 Holes	
021411005 T-Locking Plate, Oblique Angled, Combi-holes, Left 5 Holes	
021411006 T-Locking Plate, Oblique Angled, Combi-holes, Left 6 Holes	
021412003 T-Locking Plate, Oblique Angled, Combi-holes, Right	3 Holes
021412004 T-Locking Plate, Oblique Angled, Combi-holes, Right	4 Holes
021412005 T-Locking Plate, Oblique Angled, Combi-holes, Right	5 Holes
021412006 T-Locking Plate, Oblique Angled, Combi-holes, Right	6 Holes
021421003 Olecranon Locking Plate , Combi-holes, Left	3 Holes
021421004 Olecranon Locking Plate , Combi-holes, Left	4 Holes
021421005 Olecranon Locking Plate , Combi-holes, Left	5 Holes
021421006 Olecranon Locking Plate , Combi-holes, Left	6 Holes
021421007 Olecranon Locking Plate , Combi-holes, Left	7 Holes
021421008 Olecranon Locking Plate , Combi-holes, Left	8 Holes
021422003 Olecranon Locking Plate , Combi-holes, Right	3 Holes
021422004 Olecranon Locking Plate , Combi-holes, Right	4 Holes
021422005 Olecranon Locking Plate , Combi-holes, Right	5 Holes
021422006 Olecranon Locking Plate , Combi-holes, Right	6 Holes
021422007 Olecranon Locking Plate , Combi-holes, Right	7 Holes
021422008 Olecranon Locking Plate , Combi-holes, Right	8 Holes
021431005 Distal Femoral Locking Plate , Combi-holes, Left	5 Holes
021431006 Distal Femoral Locking Plate , Combi-holes, Left	6 Holes
021431007 Distal Femoral Locking Plate , Combi-holes, Left	7 Holes
021431008 Distal Femoral Locking Plate , Combi-holes, Left	8 Holes
021431009 Distal Femoral Locking Plate , Combi-holes, Left	9 Holes
021431010 Distal Femoral Locking Plate , Combi-holes, Left	10 Holes

021431011	Distal Femoral Locking Plate	, Combi-holes, Left	11 Holes
021431012	Distal Femoral Locking Plate	, Combi-holes, Left	12 Holes
021431013	Distal Femoral Locking Plate	, Combi-holes, Left	13 Holes
021432005	Distal Femoral Locking Plate	, Combi-holes, Right	5 Holes
021432006	Distal Femoral Locking Plate	, Combi-holes, Right	6 Holes
021432007	Distal Femoral Locking Plate	, Combi-holes, Right	7 Holes
021432008	Distal Femoral Locking Plate	, Combi-holes, Right	8 Holes
021432009	Distal Femoral Locking Plate	, Combi-holes, Right	9 Holes
021432010	Distal Femoral Locking Plate	, Combi-holes, Right	10 Holes
021432011	Distal Femoral Locking Plate	, Combi-holes, Right	11 Holes
021432012	Distal Femoral Locking Plate	, Combi-holes, Right	12 Holes
021432013	Distal Femoral Locking Plate	, Combi-holes, Right	13 Holes
021441003	Lateral Tibial Head Buttress Locking Plate	, Combi-holes, Left	3 Holes
021441004	Lateral Tibial Head Buttress Locking Plate	, Combi-holes, Left	4 Holes
021441005	Lateral Tibial Head Buttress Locking Plate	, Combi-holes, Left	5 Holes
021441006	Lateral Tibial Head Buttress Locking Plate	, Combi-holes, Left	6 Holes
021441007	Lateral Tibial Head Buttress Locking Plate	, Combi-holes, Left	7 Holes
021441008	Lateral Tibial Head Buttress Locking Plate	, Combi-holes, Left	8 Holes
021441009	Lateral Tibial Head Buttress Locking Plate	, Combi-holes, Left	9 Holes
021441010	Lateral Tibial Head Buttress Locking Plate	, Combi-holes, Left	10 Holes
021441011	Lateral Tibial Head Buttress Locking Plate	, Combi-holes, Left	11 Holes
021441012	Lateral Tibial Head Buttress Locking Plate	, Combi-holes, Left	12 Holes
021441013	Lateral Tibial Head Buttress Locking Plate	, Combi-holes, Left	13 Holes
021442003	Lateral Tibial Head Buttress Locking Plate	, Combi-holes, Right	3 Holes
021442004	Lateral Tibial Head Buttress Locking Plate	, Combi-holes, Right	4 Holes
021442005	Lateral Tibial Head Buttress Locking Plate	, Combi-holes, Right	5 Holes
021442006	Lateral Tibial Head Buttress Locking Plate	, Combi-holes, Right	6 Holes
021442007	Lateral Tibial Head Buttress Locking Plate	, Combi-holes, Right	7 Holes

021442008	Lateral Tibial Head Buttress Locking Plate	, Combi-holes, Right	8 Holes
021442009	Lateral Tibial Head Buttress Locking Plate	, Combi-holes, Right	9 Holes
021442010	Lateral Tibial Head Buttress Locking Plate	, Combi-holes, Right	10 Holes
021442011	Lateral Tibial Head Buttress Locking Plate	, Combi-holes, Right	11 Holes
021442012	Lateral Tibial Head Buttress Locking Plate	, Combi-holes, Right	12 Holes
021442013	Lateral Tibial Head Buttress Locking Plate	, Combi-holes, Right	13 Holes
021450006	1.5S-Mini Wed Locking Plate	34.4mm	
021460008	2.0S-Mini Wed Locking Plate	47.7mm	
021470010	2.4S-Mini Wed Locking Plate	77.2mm	
021501004	Low Profile Distal medial Tibial Locking Plate, Minimally Invasive Technique	left	Holes
021501006	Low Profile Distal medial Tibial Locking Plate, Minimally Invasive Technique	left	Holes
021501008	Low Profile Distal medial Tibial Locking Plate, Minimally Invasive Technique	left	Holes
021501010	Low Profile Distal medial Tibial Locking Plate, Minimally Invasive Technique	left	10 Holes
021501012	Low Profile Distal medial Tibial Locking Plate, Minimally Invasive Technique	left	12 Holes
021501014	Low Profile Distal medial Tibial Locking Plate, Minimally Invasive Technique	left	14 Holes
021502004	Low Profile Distal medial Tibial Locking Plate, Minimally Invasive Technique	right	4 Holes
021502006	Low Profile Distal medial Tibial Locking Plate, Minimally Invasive Technique	right	6 Holes
021502008	Low Profile Distal medial Tibial Locking Plate, Minimally Invasive Technique	right	8 Holes
021502010	Low Profile Distal medial Tibial Locking Plate, Minimally Invasive Technique	right	10 Holes
021502012	Low Profile Distal medial Tibial Locking Plate, Minimally Invasive Technique	right	12 Holes
021502014	Low Profile Distal medial Tibial Locking Plate, Minimally Invasive Technique	right	14 Holes
021510008	Curved Locking Compression Plate, broad		8 Holes

021510009	Curved Locking Compression Plate, broad	9 Holes
021510010	Curved Locking Compression Plate, broad	10 Holes
021510011	Curved Locking Compression Plate, broad	11 Holes
021510012	Curved Locking Compression Plate, broad	12 Holes
021510013	Curved Locking Compression Plate, broad	13 Holes
021510014	Curved Locking Compression Plate, broad	14 Holes
021510015	Curved Locking Compression Plate, broad	15 Holes
021510016	Curved Locking Compression Plate, broad	16 Holes
021510017	Curved Locking Compression Plate, broad	17 Holes
021510018	Curved Locking Compression Plate, broad	18 Holes
021520002	Posterior Medial Proximal Tibial Plate, combi-holes	2 Holes
021520004	Posterior Medial Proximal Tibial Plate, combi-holes	4 Holes
021520006	Posterior Medial Proximal Tibial Plate, combi-holes	6 Holes
021520008	Posterior Medial Proximal Tibial Plate, combi-holes	8 Holes
021531002	Distal Radius Locking Plate , variable angle, small, left	2 Holes
021531003	Distal Radius Locking Plate , variable angle, small, left	3 Holes
021531004	Distal Radius Locking Plate , variable angle, small, left	4 Holes
021531005	Distal Radius Locking Plate , variable angle, small, left	5 Holes
021532002	Distal Radius Locking Plate , variable angle, small, right	2 Holes
021532003	Distal Radius Locking Plate , variable angle, small, right	3 Holes
021532004	Distal Radius Locking Plate , variable angle, small, right	4 Holes
021532005	Distal Radius Locking Plate , variable angle, small, right	5 Holes
021541002	Distal Radius Locking Plate , variable angle, large, left	2 Holes
021541003	Distal Radius Locking Plate , variable angle, large, left	3 Holes
021541004	Distal Radius Locking Plate , variable angle, large, left	4 Holes
021541005	Distal Radius Locking Plate , variable angle, large, left	5 Holes
021542002	Distal Radius Locking Plate , variable angle, large, right	2 Holes
021542003	Distal Radius Locking Plate , variable angle, large, right	3 Holes
021542004	Distal Radius Locking Plate , variable angle, large, right	4 Holes

021542005 Distal Radius Locking Plate , variable angle, large, right 5 Holes

021550006 Cable Locking Plate 6 Holes

021550008 Cable Locking Plate 8 Holes

021550010 Cable Locking Plate 10 Holes

021561002 Cable Grip Locking Plate, Left 2 Holes

021561004 Cable Grip Locking Plate, Left 4 Holes

021561006 Cable Grip Locking Plate, Left 6 Holes

021561008 Cable Grip Locking Plate, Left 8 Holes

021561010 Cable Grip Locking Plate, Left 10 Holes

021561012 Cable Grip Locking Plate, Left 12 Holes

021562002 Cable Grip Locking Plate, Right 2 Holes

021562004 Cable Grip Locking Plate, Right 4 Holes

021562006 Cable Grip Locking Plate, Right 6 Holes

021562008 Cable Grip Locking Plate, Right 8 Holes

021562010 Cable Grip Locking Plate, Right 10 Holes

021562012 Cable Grip Locking Plate, Right 12 Holes

021601003 Superior Anterior Distal Clavicle Locking Plate, with lateral extension, combi-holes, left 3 Holes

021601004 Superior Anterior Distal Clavicle Locking Plate, with lateral extension, combi-holes, left 4 Holes

021601005 Superior Anterior Distal Clavicle Locking Plate, with lateral extension, combi-holes, left 5 Holes

021601006 Superior Anterior Distal Clavicle Locking Plate, with lateral extension, combi-holes, left 6 Holes

021602003 Superior Anterior Distal Clavicle Locking Plate, with lateral extension, combi-holes, right 3 Holes

021602004 Superior Anterior Distal Clavicle Locking Plate, with lateral extension, combi-holes, right 4 Holes

021602005 Superior Anterior Distal Clavicle Locking Plate, with lateral extension, combi-holes, right 5 Holes

021602006 Superior Anterior Distal Clavicle Locking Plate, with lateral extension, combi-holes, right 6 Holes

021611004	Distal Medial Radius Locking Plate, left	4 Holes
021612004	Distal Medial Radius Locking Plate, right	4 Holes
021620005	Distal Lateral Radius Locking Plate	5 Holes
021620006	Distal Lateral Radius Locking Plate	6 Holes
021631003	Proximal Femoral Locking Plate , left	3 Holes
021631004	Proximal Femoral Locking Plate , Left	4 Holes
021631005	Proximal Femoral Locking Plate , Left	5 Holes
021631006	Proximal Femoral Locking Plate , Left	6 Holes
021631007	Proximal Femoral Locking Plate , Left	7 Holes
021631008	Proximal Femoral Locking Plate , Left	8 Holes
021631009	Proximal Femoral Locking Plate , Left	9 Holes
021631010	Proximal Femoral Locking Plate , Left	10 Holes
021631011	Proximal Femoral Locking Plate , Left	11 Holes
021631012	Proximal Femoral Locking Plate , Left	12 Holes
021631013	Proximal Femoral Locking Plate , Left	13 Holes
021632003	Proximal Femoral Locking Plate , Right	3 Holes
021632004	Proximal Femoral Locking Plate , Right	4 Holes
021632005	Proximal Femoral Locking Plate , Right	5 Holes
021632006	Proximal Femoral Locking Plate , Right	6 Holes
021632007	Proximal Femoral Locking Plate , Right	7 Holes
021632008	Proximal Femoral Locking Plate , Right	8 Holes
021632009	Proximal Femoral Locking Plate , Right	9 Holes
021632010	Proximal Femoral Locking Plate , Right	10 Holes
021632011	Proximal Femoral Locking Plate , Right	11 Holes
021632012	Proximal Femoral Locking Plate , Right	12 Holes
021632013	Proximal Femoral Locking Plate , Right	13 Holes
021641004	Ankle Fusion Locking Plate, left	4 Holes
021642004	Ankle Fusion Locking Plate, right	4 Holes
021671005	Proximal Lateral Tibial Locking Plate II left	5 Holes

021671009 Proximal Lateral Tibial Locking Plate II left 9 Holes
021671013 Proximal Lateral Tibial Locking Plate II left 13 Holes
021671017 Proximal Lateral Tibial Locking Plate II left 17 Holes
021672005 Proximal Lateral Tibial Locking Plate II right 5 Holes
021672009 Proximal Lateral Tibial Locking Plate II right 9 Holes
021672013 Proximal Lateral Tibial Locking Plate II right 13 Holes
021672017 Proximal Lateral Tibial Locking Plate II right 17 Holes
021681004 Proximal Medial Tibial Locking Plate II left 4 Holes
021681006 Proximal Medial Tibial Locking Plate II left 6 Holes
021681008 Proximal Medial Tibial Locking Plate II left 8 Holes
021682004 Proximal Medial Tibial Locking Plate II right 4 Holes
021682006 Proximal Medial Tibial Locking Plate II right 6 Holes
021682008 Proximal Medial Tibial Locking Plate II right 8 Holes
021691003 Proximal Posterior Tibial Locking Plate left, II 3 Holes
021691005 Proximal Posterior Tibial Locking Plate left, II 5 Holes
021691007 Proximal Posterior Tibial Locking Plate left, II 7 Holes
021692003 Proximal Posterior Tibial Locking Plate right, II 3 Holes
021692005 Proximal Posterior Tibial Locking Plate right, II 5 Holes
021692007 Proximal Posterior Tibial Locking Plate right, II 7 Holes
021800004 Humeral Shaft Straight Locking Plate 4 Holes
021800005 Humeral Shaft Straight Locking Plate 5 Holes
021800006 Humeral Shaft Straight Locking Plate 6 Holes
021800007 Humeral Shaft Straight Locking Plate 7 Holes
021800008 Humeral Shaft Straight Locking Plate 8 Holes
021800009 Humeral Shaft Straight Locking Plate 9 Holes
021800010 Humeral Shaft Straight Locking Plate 10 Holes
021800011 Humeral Shaft Straight Locking Plate 11 Holes
021800012 Humeral Shaft Straight Locking Plate 12 Holes

021800013 Humeral Shaft Straight Locking Plate 13 Holes
021800014 Humeral Shaft Straight Locking Plate 14 Holes
021800015 Humeral Shaft Straight Locking Plate 15 Holes
021800016 Humeral Shaft Straight Locking Plate 16 Holes
021810127 Sacroiliac Joint Plate 127mm
021810137 Sacroiliac Joint Plate 137mm
021810147 Sacroiliac Joint Plate 147mm
021821003 Proximal Lateral Tibial Locking Plate, left 3 Holes
021821004 Proximal Lateral Tibial Locking Plate, left 4 Holes
021821005 Proximal Lateral Tibial Locking Plate, left 5 Holes
021821006 Proximal Lateral Tibial Locking Plate, left 6 Holes
021821007 Proximal Lateral Tibial Locking Plate, left 7 Holes
021821008 Proximal Lateral Tibial Locking Plate, left 8 Holes
021821009 Proximal Lateral Tibial Locking Plate, left 9 Holes
021821010 Proximal Lateral Tibial Locking Plate, left 10 Holes
021821011 Proximal Lateral Tibial Locking Plate, left 11 Holes
021821012 Proximal Lateral Tibial Locking Plate, left 12 Holes
021821013 Proximal Lateral Tibial Locking Plate, left 13 Holes
021821014 Proximal Lateral Tibial Locking Plate, left 14 Holes
021821015 Proximal Lateral Tibial Locking Plate, left 15 Holes
021822003 Proximal Lateral Tibial Locking Plate, right 3 Holes
021822004 Proximal Lateral Tibial Locking Plate, right 4 Holes
021822005 Proximal Lateral Tibial Locking Plate, right 5 Holes
021822006 Proximal Lateral Tibial Locking Plate, right 6 Holes
021822007 Proximal Lateral Tibial Locking Plate, right 7 Holes
021822008 Proximal Lateral Tibial Locking Plate, right 8 Holes
021822009 Proximal Lateral Tibial Locking Plate, right 9 Holes
021822010 Proximal Lateral Tibial Locking Plate, right 10 Holes
021822011 Proximal Lateral Tibial Locking Plate, right 11 Holes

021822012 Proximal Lateral Tibial Locking Plate, right 12 Holes
021822013 Proximal Lateral Tibial Locking Plate, right 13 Holes
021822014 Proximal Lateral Tibial Locking Plate, right 14 Holes
021822015 Proximal Lateral Tibial Locking Plate, right 15 Holes
021881003 Posterior Lateral Fibular Locking Plate combi holes, left 3 Holes
021881004 Posterior Lateral Fibular Locking Plate combi holes, left 4 Holes
021881005 Posterior Lateral Fibular Locking Plate combi holes, left 5 Holes
021881006 Posterior Lateral Fibular Locking Plate combi holes, left 6 Holes
021881007 Posterior Lateral Fibular Locking Plate combi holes, left 7 Holes
021881008 Posterior Lateral Fibular Locking Plate combi holes, left 8 Holes
021882003 Posterior Lateral Fibular Locking Plate combi holes, right 3 Holes
021882004 Posterior Lateral Fibular Locking Plate combi holes, right 4 Holes
021882005 Posterior Lateral Fibular Locking Plate combi holes, right 5 Holes
021882006 Posterior Lateral Fibular Locking Plate combi holes, right 6 Holes
021882007 Posterior Lateral Fibular Locking Plate combi holes, right 7 Holes
021882008 Posterior Lateral Fibular Locking Plate combi holes, right 8 Holes
021891003 Distal Lateral Humeral Locking Plate, head 2.7 hole, combi-holes, Left 3 Holes
021891004 Distal Lateral Humeral Locking Plate, head 2.7 hole, combi-holes, Left 4 Holes
021891005 Distal Lateral Humeral Locking Plate, head 2.7 hole, combi-holes, Left 5 Holes
021891006 Distal Lateral Humeral Locking Plate, head 2.7 hole, combi-holes, Left 6 Holes
021891007 Distal Lateral Humeral Locking Plate, head 2.7 hole, combi-holes, Left 7 Holes
021891008 Distal Lateral Humeral Locking Plate, head 2.7 hole, combi-holes, Left 8 Holes
021891009 Distal Lateral Humeral Locking Plate, head 2.7 hole, combi-holes, Left 9 Holes
021892003 Distal Lateral Humeral Locking Plate, head 2.7 hole, combi-holes, Right 3 Holes
021892004 Distal Lateral Humeral Locking Plate, head 2.7 hole, combi-holes, Right 4 Holes
021892005 Distal Lateral Humeral Locking Plate, head 2.7 hole, combi-holes, Right 5 Holes
021892006 Distal Lateral Humeral Locking Plate, head 2.7 hole, combi-holes, Right 6 Holes
021892007 Distal Lateral Humeral Locking Plate, head 2.7 hole, combi-holes, Right 7 Holes

021892008 Distal Lateral Humeral Locking Plate, head 2.7 hole, combi-holes, Right 8 Holes
021892009 Distal Lateral Humeral Locking Plate, head 2.7 hole, combi-holes, Right 9 Holes
021901003 Distal Medial Humeral Locking Plate, head 2.7 hole, combi-holes, Left 3 Holes
021901004 Distal Medial Humeral Locking Plate, head 2.7 hole, combi-holes, Left 4 Holes
021901005 Distal Medial Humeral Locking Plate, head 2.7 hole, combi-holes, Left 5 Holes
021901006 Distal Medial Humeral Locking Plate, head 2.7 hole, combi-holes, Left 6 Holes
021901007 Distal Medial Humeral Locking Plate, head 2.7 hole, combi-holes, Left 7 Holes
021901008 Distal Medial Humeral Locking Plate, head 2.7 hole, combi-holes, Left 8 Holes
021901009 Distal Medial Humeral Locking Plate, head 2.7 hole, combi-holes, Left 9 Holes
021902003 Distal Medial Humeral Locking Plate, head 2.7 hole, combi-holes, Right 3 Holes
021902004 Distal Medial Humeral Locking Plate, head 2.7 hole, combi-holes, Right 4 Holes
021902005 Distal Medial Humeral Locking Plate, head 2.7 hole, combi-holes, Right 5 Holes
021902006 Distal Medial Humeral Locking Plate, head 2.7 hole, combi-holes, Right 6 Holes
021902007 Distal Medial Humeral Locking Plate, head 2.7 hole, combi-holes, Right 7 Holes
021902008 Distal Medial Humeral Locking Plate, head 2.7 hole, combi-holes, Right 8 Holes
021902009 Distal Medial Humeral Locking Plate, head 2.7 hole, combi-holes, Right 9 Holes
022100008 Wrist Fusion Locking Plate, length 117mm, standard, combi-holes 8 Holes
022110008 Wrist Fusion Locking Plate, length 117mm, short, combi-holes 8 Holes
022120009 Wrist Fusion Locking Plate, length 112mm, straight, combi-holes 9 Holes
022181002 Distal Radius Locking Plate, head 7 holes variable angle, left 2 Holes
022181003 Distal Radius Locking Plate, head 7 holes variable angle, left 3 Holes
022181004 Distal Radius Locking Plate, head 7 holes variable angle, left 4 Holes
022181005 Distal Radius Locking Plate, head 7 holes variable angle, left 5 Holes
022182002 Distal Radius Locking Plate, head 7 holes variable angle, right 2 Holes
022182003 Distal Radius Locking Plate, head 7 holes variable angle, right 3 Holes
022182004 Distal Radius Locking Plate, head 7 holes variable angle, right 4 Holes
022182005 Distal Radius Locking Plate, head 7 holes variable angle, right 5 Holes
022360004 Proximal Medial Tibial Osteotomy Locking Plate, combi-holes 4 Holes
022981003 Distal Lateral Femoral Osteotomy Locking Plate, combi-holes, left 3 Holes

022981004	Distal Lateral Femoral Osteotomy Locking Plate, combi-holes,left	4 Holes
022981005	Distal Lateral Femoral Osteotomy Locking Plate, combi-holes,left	5 Holes
022981006	Distal Lateral Femoral Osteotomy Locking Plate, combi-holes,left	6 Holes
022981007	Distal Lateral Femoral Osteotomy Locking Plate, combi-holes,left	7 Holes
022981008	Distal Lateral Femoral Osteotomy Locking Plate, combi-holes,left	8 Holes
022981009	Distal Lateral Femoral Osteotomy Locking Plate, combi-holes,left	9 Holes
022981010	Distal Lateral Femoral Osteotomy Locking Plate, combi-holes,left	10 Holes
022981011	Distal Lateral Femoral Osteotomy Locking Plate, combi-holes,left	11 Holes
022981012	Distal Lateral Femoral Osteotomy Locking Plate, combi-holes,left	12 Holes
022981013	Distal Lateral Femoral Osteotomy Locking Plate, combi-holes,left	13 Holes
022982003	Distal Lateral Femoral Osteotomy Locking Plate, combi-holes,right	3 Holes
022982004	Distal Lateral Femoral Osteotomy Locking Plate, combi-holes,right	4 Holes
022982005	Distal Lateral Femoral Osteotomy Locking Plate, combi-holes,right	5 Holes
022982006	Distal Lateral Femoral Osteotomy Locking Plate, combi-holes,right	6 Holes
022982007	Distal Lateral Femoral Osteotomy Locking Plate, combi-holes,right	7 Holes
022982008	Distal Lateral Femoral Osteotomy Locking Plate, combi-holes,right	8 Holes
022982009	Distal Lateral Femoral Osteotomy Locking Plate, combi-holes,right	9 Holes
022982010	Distal Lateral Femoral Osteotomy Locking Plate, combi-holes,right	10 Holes
022982011	Distal Lateral Femoral Osteotomy Locking Plate, combi-holes,right	11 Holes
022982012	Distal Lateral Femoral Osteotomy Locking Plate, combi-holes,right	12 Holes
022982013	Distal Lateral Femoral Osteotomy Locking Plate, combi-holes,right	13 Holes
022991003	Proximal Lateral Tibial Osteotomy Locking Plate,left	3 Holes
022991005	Proximal Lateral Tibial Osteotomy Locking Plate,left	5 Holes
022991007	Proximal Lateral Tibial Osteotomy Locking Plate,left	7 Holes
022991009	Proximal Lateral Tibial Osteotomy Locking Plate,left	9 Holes
022991011	Proximal Lateral Tibial Osteotomy Locking Plate,left	11 Holes
022991013	Proximal Lateral Tibial Osteotomy Locking Plate,left	13 Holes
022992003	Proximal Lateral Tibial Osteotomy Locking Plate,right	3 Holes

022992005 Proximal Lateral Tibial Osteotomy Locking Plate,right	5 Holes
022992007 Proximal Lateral Tibial Osteotomy Locking Plate,right	7 Holes
022992009 Proximal Lateral Tibial Osteotomy Locking Plate,right	9 Holes
022992011 Proximal Lateral Tibial Osteotomy Locking Plate,right	11 Holes
022992013 Proximal Lateral Tibial Osteotomy Locking Plate,right	13 Holes
025331000 Coronoid Locking Plate, Left	left
025332000 Coronoid Locking Plate, Right	right
025341003 Distal Lateral Tibial Locking Plate , combi-holes,left	3 Holes
025341005 Distal Lateral Tibial Locking Plate , combi-holes,Left	5Holes
025341007 Distal Lateral Tibial Locking Plate , combi-holes,left	7 Holes
025341009 Distal Lateral Tibial Locking Plate , combi-holes,Left	9Holes
025341011 Distal Lateral Tibial Locking Plate , combi-holes,left	11 Holes
025341013 Distal Lateral Tibial Locking Plate , combi-holes,Left	13Holes
025341015 Distal Lateral Tibial Locking Plate , combi-holes,left	15 Holes
025341017 Distal Lateral Tibial Locking Plate , combi-holes,Left	17Holes
025342003 Distal Lateral Tibial Locking Plate , combi-holes,right	3 Holes
025342005 Distal Lateral Tibial Locking Plate , combi-holes,Right	5Holes
025342007 Distal Lateral Tibial Locking Plate , combi-holes,right	7 Holes
025342009 Distal Lateral Tibial Locking Plate , combi-holes,Right	9Holes
025342011 Distal Lateral Tibial Locking Plate , combi-holes,right	11 Holes
025342013 Distal Lateral Tibial Locking Plate , combi-holes,Right	13Holes
025342015 Distal Lateral Tibial Locking Plate , combi-holes,right	15 Holes
025342017 Distal Lateral Tibial Locking Plate , combi-holes,Right	17Holes
025351008 Distal Diaphyseal Metacarpal Radius Locking Plate, left	8 Holes
025351010 Distal Diaphyseal Metacarpal Radius Locking Plate, left	10 Holes
025351012 Distal Diaphyseal Metacarpal Radius Locking Plate, left	12 Holes
025352008 Distal Diaphyseal Metacarpal Radius Locking Plate, right	8 Holes
025352010 Distal Diaphyseal Metacarpal Radius Locking Plate, right	10 Holes
025352012 Distal Diaphyseal Metacarpal Radius Locking Plate, right	12 Holes

025551005 Proximal Humeral Greater Tuberosity Locking Compression Plate, suture hole, left	5 Holes
025552005 Proximal Humeral Greater Tuberosity Locking Compression Plate, suture hole, right	5 Holes
025681006 Superior Anterior Clavicle Locking Plate , combi-holes, left	6 Holes
025681007 Superior Anterior Clavicle Locking Plate , combi-holes, left	7 Holes
025681008 Superior Anterior Clavicle Locking Plate , combi-holes, left	8 Holes
025681009 Superior Anterior Clavicle Locking Plate , combi-holes, left	9 Holes
025681010 Superior Anterior Clavicle Locking Plate , combi-holes, left	10 Holes
025681011 Superior Anterior Clavicle Locking Plate , combi-holes, left	11 Holes
025681012 Superior Anterior Clavicle Locking Plate , combi-holes, left	12 Holes
025682006 Superior Anterior Clavicle Locking Plate , combi-holes, right	6 Holes
025682007 Superior Anterior Clavicle Locking Plate , combi-holes, right	7 Holes
025682008 Superior Anterior Clavicle Locking Plate , combi-holes, right	8 Holes
025682009 Superior Anterior Clavicle Locking Plate , combi-holes, right	9 Holes
025682010 Superior Anterior Clavicle Locking Plate , combi-holes, right	10 Holes
025682011 Superior Anterior Clavicle Locking Plate , combi-holes, right	11 Holes
025682012 Superior Anterior Clavicle Locking Plate , combi-holes, right	12 Holes
025850006 Metaphyseal Locking Plate, 3,5	6 Holes
025850007 Metaphyseal Locking Plate, 3,5	7 Holes
025850008 Metaphyseal Locking Plate, 3,5	8 Holes
025850009 Metaphyseal Locking Plate, 3,5	9 Holes
025850010 Metaphyseal Locking Plate, 3,5	10 Holes
025850011 Metaphyseal Locking Plate, 3,5	11 Holes
025850012 Metaphyseal Locking Plate, 3,5	12 Holes
025850014 Metaphyseal Locking Plate, 3,5	14 Holes
025870004 Distal Medial Humeral Reconstruction Locking Plate	4 Holes
025870006 Distal Medial Humeral Reconstruction Locking Plate	6 Holes
025870008 Distal Medial Humeral Reconstruction Locking Plate	8 Holes
025870010 Distal Medial Humeral Reconstruction Locking Plate	10 Holes

025870012 Distal Medial Humeral Reconstruction Locking Plate 12 Holes
026012008 Olecranon Locking Plate, curved, riht 8 holes
026020007 Olecranon Locking Plate, straight 7 holes
026110004 MD Distal Tibial T-Locking Plate 4 holes
026121004 MD Distal Tibial L-Locking Plate, left 4 holes
026122004 MD Distal Tibial L-Locking Plate, right 4 holes
026131006 Distal Diaphyseal Radius Locking Plate II, left 6 Holes
026131008 Distal Diaphyseal Radius Locking Plate II, left 8 Holes
026131010 Distal Diaphyseal Radius Locking Plate II, left 10 Holes
026131012 Distal Diaphyseal Radius Locking Plate II, left 12 Holes
026131014 Distal Diaphyseal Radius Locking Plate II, left 14 Holes
026132006 Distal Diaphyseal Radius Locking Plate II, right 6 Holes
026132008 Distal Diaphyseal Radius Locking Plate II, right 8 Holes
026132010 Distal Diaphyseal Radius Locking Plate II, right 10 Holes
026132012 Distal Diaphyseal Radius Locking Plate II, right 12 Holes
026132014 Distal Diaphyseal Radius Locking Plate II, right 14 Holes
026141000 Acromion Locking Plate, left
026142000 Acromion Locking Plate, right
026151000 Scapula Medial Locking Plate, left
026152000 Scapula Medial Locking Plate, right
026161000 Scapula Lateral Locking Plate, left
026162000 Scapula Lateral Locking Plate, right
026331000 Glenoid Locking Plate, left
026332000 Glenoid Locking Plate, right
026611012 Olecranon U-locking Plate, left 12 Holes
026612012 Olecranon U-locking Plate, right 12 Holes
026620006 MT Calcaneal Step Locking Plate
026620008 MT Calcaneal Step Locking Plate
026620010 MT Calcaneal Step Locking Plate

026651000 MT First Metatarsal Osteotomy Locking Plate, left
026651003 MT First Metatarsal Osteotomy Locking Plate, left 3 Holes
026651004 MT First Metatarsal Osteotomy Locking Plate, left 4 Holes
026651005 MT First Metatarsal Osteotomy Locking Plate, left 5 Holes
026651006 MT First Metatarsal Osteotomy Locking Plate, left 6 Holes
026651007 MT First Metatarsal Osteotomy Locking Plate, left 7 Holes
026652000 MT First Metatarsal Osteotomy Locking Plate, right
026652003 MT First Metatarsal Osteotomy Locking Plate, right 3 Holes
026652004 MT First Metatarsal Osteotomy Locking Plate, right 4 Holes
026652005 MT First Metatarsal Osteotomy Locking Plate, right 5 Holes
026652006 MT First Metatarsal Osteotomy Locking Plate, right 6 Holes
026652007 MT First Metatarsal Osteotomy Locking Plate, right 7 Holes
026660000 MT Lateral Column Lengthening Locking Plate
026660008 MT Lateral Column Lengthening Locking Plate
026660010 MT Lateral Column Lengthening Locking Plate
026660012 MT Lateral Column Lengthening Locking Plate
030300006 1,5 mm Locking Screw, stardrive, self tapping 6 mm
030300007 1,5 mm Locking Screw, stardrive, self tapping 7 mm
030300008 1,5 mm Locking Screw, stardrive, self tapping 8 mm
030300009 1,5 mm Locking Screw, stardrive, self tapping 9 mm
030300010 1,5 mm Locking Screw, stardrive, self tapping 10 mm
030300011 1,5 mm Locking Screw, stardrive, self tapping 11 mm
030300012 1,5 mm Locking Screw, stardrive, self tapping 12 mm
030300013 1,5 mm Locking Screw, stardrive, self tapping 13 mm
030300014 1,5 mm Locking Screw, stardrive, self tapping 14 mm
030300015 1,5 mm Locking Screw, stardrive, self tapping 15 mm
030300016 1,5 mm Locking Screw, stardrive, self tapping 16 mm
030300017 1,5 mm Locking Screw, stardrive, self tapping 17 mm

030300018	1,5 mm Locking Screw, stardrive, self tapping	18 mm
030300019	1,5 mm Locking Screw, stardrive, self tapping	19 mm
030300020	1,5 mm Locking Screw, stardrive, self tapping	20 mm
030320006	2,0mm Locking Screw, stardrive, self tapping	6 mm
030320007	2,0mm Locking Screw, stardrive, self tapping	7 mm
030320008	2,0mm Locking Screw, stardrive, self tapping	8 mm
030320009	2,0mm Locking Screw, stardrive, self tapping	9 mm
030320010	2,0mm Locking Screw, stardrive, self tapping	10 mm
030320011	2,0mm Locking Screw, stardrive, self tapping	11 mm
030320012	2,0mm Locking Screw, stardrive, self tapping	12 mm
030320013	2,0mm Locking Screw, stardrive, self tapping	13 mm
030320014	2,0mm Locking Screw, stardrive, self tapping	14 mm
030320015	2,0mm Locking Screw, stardrive, self tapping	15 mm
030320016	2,0mm Locking Screw, stardrive, self tapping	16 mm
030320017	2,0mm Locking Screw, stardrive, self tapping	17 mm
030320018	2,0mm Locking Screw, stardrive, self tapping	18 mm
030320019	2,0mm Locking Screw, stardrive, self tapping	19 mm
030320020	2,0mm Locking Screw, stardrive, self tapping	20 mm
030320022	2,0mm Locking Screw, stardrive, self tapping	22 mm
030320024	2,0mm Locking Screw, stardrive, self tapping	24 mm
030320026	2,0mm Locking Screw, stardrive, self tapping	26 mm
030320028	2,0mm Locking Screw, stardrive, self tapping	28 mm
030320030	2,0mm Locking Screw, stardrive, self tapping	30 mm
030340006	2,7 mm Locking Screw, stardrive, selft-tapping	6 mm
030340007	2,7 mm Locking Screw, stardrive, selft-tapping	7 mm
030340008	2,7 mm Locking Screw, stardrive, selft-tapping	8 mm
030340009	2,7 mm Locking Screw, stardrive, selft-tapping	9 mm
030340010	2,7 mm Locking Screw, stardrive, selft-tapping	10 mm
030340011	2,7 mm Locking Screw, stardrive, selft-tapping	11 mm

030340012	2,7 mm Locking Screw, stardrive, self-tapping	12 mm
030340013	2,7 mm Locking Screw, stardrive, self-tapping	13 mm
030340014	2,7 mm Locking Screw, stardrive, self-tapping	14 mm
030340015	2,7 mm Locking Screw, stardrive, self-tapping	15 mm
030340016	2,7 mm Locking Screw, stardrive, self-tapping	16 mm
030340017	2,7 mm Locking Screw, stardrive, self-tapping	17 mm
030340018	2,7 mm Locking Screw, stardrive, self-tapping	18 mm
030340019	2,7 mm Locking Screw, stardrive, self-tapping	19 mm
030340020	2,7 mm Locking Screw, stardrive, self-tapping	20 mm
030340022	2,7 mm Locking Screw, stardrive, self-tapping	22 mm
030340024	2,7 mm Locking Screw, stardrive, self-tapping	24 mm
030340026	2,7 mm Locking Screw, stardrive, self-tapping	26 mm
030340028	2,7 mm Locking Screw, stardrive, self-tapping	28 mm
030340030	2,7 mm Locking Screw, stardrive, self-tapping	30 mm
030340035	2.7mm Locking Screw, Stardrive, Self-tapping	35mm
030340040	2.7mm Locking Screw, Stardrive, Self-tapping	40mm
030340045	2.7mm Locking Screw, Stardrive, Self-tapping	45mm
030340050	2.7mm Locking Screw, Stardrive, Self-tapping	50mm
031271205	MD Anterior Process Calcaneal Locking Plate, small, left	5 Holes
031271705	MD Anterior Process Calcaneal Locking Plate, large, left	5 Holes
031272205	MD Anterior Process Calcaneal Locking Plate, small, right	5 Holes
031272705	MD Anterior Process Calcaneal Locking Plate, large, right	5 Holes
031281005	MD Posterior Tuberosity Calcaneal Locking Plate, small, left	5 Holes
031281006	MD Posterior Tuberosity Calcaneal Locking Plate, large, left	6 Holes
031282005	MD Posterior Tuberosity Calcaneal Locking Plate, small, right	5 Holes
031282006	MD Posterior Tuberosity Calcaneal Locking Plate, large, right	6 Holes
031291008	MD AP and PT Calcaneal Locking Plate, medium, left	8 Holes
031291009	MD AP and PT Calcaneal Locking Plate, large, left	9 Holes

031291308 MD AP and PT Calcaneal Locking Plate, small, left 8 Holes
031292008 MD AP and PT Calcaneal Locking Plate, medium, right 8 Holes
031292009 MD AP and PT Calcaneal Locking Plate, large, right 9 Holes
031292308 MD AP and PT Calcaneal Locking Plate, small, right 8 Holes
031301011 MD Calcaneal Locking Plate, small, left 11Holes
031302011 MD Calcaneal Locking Plate, small, right 11Holes
031311011 MD Calcaneal Locking Plate, medium, left 11Holes
031312011 MD Calcaneal Locking Plate, medium, right 11Holes
031321011 MD Calcaneal Locking Plate, large, left 11Holes
031322011 MD Calcaneal Locking Plate, large, right 11Holes
031360002 MT T-Locking Plate for First Metatarsa 2 Holes
031851006 MD Fifth Metatarsal Locking Plate
031852006 MD Fifth Metatarsal Locking Plate
031861009 MT Cuboid Locking Plate, left
031862009 MT Cuboid Locking Plate, right
031870005 MT Talus Locking Plate
031941005 MT First MTC Plantar Fusion Locking Plate, large, left
031942005 MT First MTC Plantar Fusion Locking Plate, large, right
031951005 MT First MTC Plantar Fusion Locking Plate, small, left
031952005 MT First MTC Plantar Fusion Locking Plate, small, right
031961004 Distal Lateral Fibular Locking Plate, left 4 holes
031961005 Distal Lateral Fibular Locking Plate, left 5 holes
031961006 Distal Lateral Fibular Locking Plate, left 6 holes
031961007 Distal Lateral Fibular Locking Plate, left 7 holes
031961009 Distal Lateral Fibular Locking Plate, left 9 holes
031961011 Distal Lateral Fibular Locking Plate, left 11 holes
031961013 Distal Lateral Fibular Locking Plate, left 13 holes
031962004 Distal Lateral Fibular Locking Plate, right 4 holes
031962005 Distal Lateral Fibular Locking Plate, right 5 holes

031962006	Distal Lateral Fibular Locking Plate, right	6 holes
031962007	Distal Lateral Fibular Locking Plate, right	7 holes
031962009	Distal Lateral Fibular Locking Plate, right	9 holes
031962011	Distal Lateral Fibular Locking Plate, right	11 holes
031962013	Distal Lateral Fibular Locking Plate, right	13 holes
031970003	MD Distal Tibial Locking Plate	3 holes
031980004	MD Distal Medial Tibial Locking Plate	4 holes
031990004	MD Distal Anterior Tibial Locking Plate	4 holes
031990006	MD Distal Anterior Tibial Locking Plate	6 holes
032001009	MT Navicular Locking Plate, left	
032002009	MT Navicular Locking Plate, right	
040070010	3.5mm Locking Screw, Self-tapping	10mm
040070012	3.5mm Locking Screw, Self-tapping	12mm
040070014	3.5mm Locking Screw, Self-tapping	14mm
040070016	3.5mm Locking Screw, Self-tapping	16mm
040070018	3.5mm Locking Screw, Self-tapping	18mm
040070020	3.5mm Locking Screw, Self-tapping	20mm
040070022	3.5mm Locking Screw, Self-tapping	22mm
040070024	3.5mm Locking Screw, Self-tapping	24mm
040070026	3.5mm Locking Screw, Self-tapping	26mm
040070028	3.5mm Locking Screw, Self-tapping	28mm
040070030	3.5mm Locking Screw, Self-tapping	30mm
040070032	3.5mm Locking Screw, Self-tapping	32mm
040070034	3.5mm Locking Screw, Self-tapping	34mm
040070035	3.5mm Locking Screw, Self-tapping	35mm
040070036	3.5mm Locking Screw, Self-tapping	36mm
040070038	3.5mm Locking Screw, Self-tapping	38mm
040070040	3.5mm Locking Screw, Self-tapping	40mm

040070042	3.5mm Locking Screw, Self-tapping	42mm
040070044	3.5mm Locking Screw, Self-tapping	44mm
040070045	3.5mm Locking Screw, Self-tapping	45mm
040070046	3.5mm Locking Screw, Self-tapping	46mm
040070048	3.5mm Locking Screw, Self-tapping	48mm
040070050	3.5mm Locking Screw, Self-tapping	50mm
040070052	3.5mm Locking Screw, Self-tapping	52mm
040070054	3.5mm Locking Screw, Self-tapping	54mm
040070055	3.5mm Locking Screw, Self-tapping	55mm
040070056	3.5mm Locking Screw, Self-tapping	56mm
040070058	3.5mm Locking Screw, Self-tapping	58mm
040070060	3.5mm Locking Screw, Self-tapping	60mm
040070065	3.5mm Locking Screw, Self-tapping	65mm
040070070	3.5mm Locking Screw, Self-tapping	70mm
040070075	3.5mm Locking Screw, Self-tapping	75mm
040080012	3.5mm Locking Screw, Self-drilling	12mm
040080014	3.5mm Locking Screw, Self-drilling	14mm
040080016	3.5mm Locking Screw, Self-drilling	16mm
040080018	3.5mm Locking Screw, Self-drilling	18mm
040080020	3.5mm Locking Screw, Self-drilling	20mm
040080022	3.5mm Locking Screw, Self-drilling	22mm
040080024	3.5mm Locking Screw, Self-drilling	24mm
040080026	3.5mm Locking Screw, Self-drilling	26mm
040080028	3.5mm Locking Screw, Self-drilling	28mm
040080030	3.5mm Locking Screw, Self-drilling	30mm
040080032	3.5mm Locking Screw, Self-drilling	32mm
040080034	3.5mm Locking Screw, Self-drilling	34mm
040080035	3.5mm Locking Screw, Self-drilling	35mm
040080036	3.5mm Locking Screw, Self-drilling	36mm

040080038	3.5mm Locking Screw, Self-drilling	38mm
040080040	3.5mm Locking Screw, Self-drilling	40mm
040080042	3.5mm Locking Screw, Self-drilling	42mm
040080044	3.5mm Locking Screw, Self-drilling	44mm
040080045	3.5mm Locking Screw, Self-drilling	45mm
040080046	3.5mm Locking Screw, Self-drilling	46mm
040080048	3.5mm Locking Screw, Self-drilling	48mm
040080050	3.5mm Locking Screw, Self-drilling	50mm
040080052	3.5mm Locking Screw, Self-drilling	52mm
040080054	3.5mm Locking Screw, Self-drilling	54mm
040080055	3.5mm Locking Screw, Self-drilling	55mm
040080056	3.5mm Locking Screw, Self-drilling	56mm
040080058	3.5mm Locking Screw, Self-drilling	58mm
040080060	3.5mm Locking Screw, Self-drilling	60mm
040090016	5.0mm Locking Screw, Self-tapping	16mm
040090018	5.0mm Locking Screw, Self-tapping	18mm
040090020	5.0mm Locking Screw, Self-tapping	20mm
040090022	5.0mm Locking Screw, Self-tapping	22mm
040090024	5.0mm Locking Screw, Self-tapping	24mm
040090026	5.0mm Locking Screw, Self-tapping	26mm
040090028	5.0mm Locking Screw, Self-tapping	28mm
040090030	5.0mm Locking Screw, Self-tapping	30mm
040090032	5.0mm Locking Screw, Self-tapping	32mm
040090034	5.0mm Locking Screw, Self-tapping	34mm
040090036	5.0mm Locking Screw, Self-tapping	36mm
040090038	5.0mm Locking Screw, Self-tapping	38mm
040090040	5.0mm Locking Screw, Self-tapping	40mm
040090042	5.0mm Locking Screw, Self-tapping	42mm

040090044	5.0mm Locking Screw, Self-tapping	44mm
040090045	5.0mm Locking Screw, Self-tapping	45mm
040090046	5.0mm Locking Screw, Self-tapping	46mm
040090048	5.0mm Locking Screw, Self-tapping	48mm
040090050	5.0mm Locking Screw, Self-tapping	50mm
040090055	5.0mm Locking Screw, Self-tapping	55mm
040090060	5.0mm Locking Screw, Self-tapping	60mm
040090065	5.0mm Locking Screw, Self-tapping	65mm
040090070	5.0mm Locking Screw, Self-tapping	70mm
040090075	5.0mm Locking Screw, Self-tapping	75mm
040090080	5.0mm Locking Screw, Self-tapping	80mm
040090085	5.0mm Locking Screw, Self-tapping	85mm
040090090	5.0mm Locking Screw, Self-tapping	90mm
040100020	5.0mm Locking Screw, Self-drilling	20mm
040100022	5.0mm Locking Screw, Self-drilling	22mm
040100024	5.0mm Locking Screw, Self-drilling	24mm
040100026	5.0mm Locking Screw, Self-drilling	26mm
040100028	5.0mm Locking Screw, Self-drilling	28mm
040100030	5.0mm Locking Screw, Self-drilling	30mm
040100032	5.0mm Locking Screw, Self-drilling	32mm
040100034	5.0mm Locking Screw, Self-drilling	34mm
040100036	5.0mm Locking Screw, Self-drilling	36mm
040100038	5.0mm Locking Screw, Self-drilling	38mm
040100040	5.0mm Locking Screw, Self-drilling	40mm
040100042	5.0mm Locking Screw, Self-drilling	42mm
040100044	5.0mm Locking Screw, Self-drilling	44mm
040100046	5.0mm Locking Screw, Self-drilling	46mm
040100048	5.0mm Locking Screw, Self-drilling	48mm
040100050	5.0mm Locking Screw, Self-drilling	50mm

040100055	5.0mm Locking Screw, Self-drilling	55mm
040100060	5.0mm Locking Screw, Self-drilling	60mm
040100065	5.0mm Locking Screw, Self-drilling	65mm
040100070	5.0mm Locking Screw, Self-drilling	70mm
040100075	5.0mm Locking Screw, Self-drilling	75mm
040100080	5.0mm Locking Screw, Self-drilling	80mm
040100085	5.0mm Locking Screw, Self-drilling	85mm
040100090	5.0mm Locking Screw, Self-drilling	90mm
040130060	7.3mm Locking Screw	60mm
040130065	7.3mm Locking Screw	65mm
040130070	7.3mm Locking Screw	70mm
040130075	7.3mm Locking Screw	75mm
040130076	7.3mm Locking Screw, 76.5mm	76.5mm
040130080	7.3mm Locking Screw	80mm
040130083	7.3mm Locking Screw, 83.5mm	83.5mm
040130085	7.3mm Locking Screw	85mm
040130090	7.3mm Locking Screw	90mm
040130095	7.3mm Locking Screw	95mm
040130100	7.3mm Locking Screw	100mm
040130105	7.3mm Locking Screw	105mm
040130110	7.3mm Locking Screw	110mm
040130115	7.3mm Locking Screw	115mm
040140060	6.5mm Cannulated Locking Screw, Fully Threaded	60mm
040140065	6.5mm Cannulated Locking Screw, Fully Threaded	65mm
040140070	6.5mm Cannulated Locking Screw, Fully Threaded	70mm
040140075	6.5mm Cannulated Locking Screw, Fully Threaded	75mm
040140080	6.5mm Cannulated Locking Screw, Fully Threaded	80mm
040140085	6.5mm Cannulated Locking Screw, Fully Threaded	85mm

040140090 6.5mm Cannulated Locking Screw, Fully Threaded 90mm
040140095 6.5mm Cannulated Locking Screw, Fully Threaded 95mm
040140100 6.5mm Cannulated Locking Screw, Fully Threaded 100mm
040140105 6.5mm Cannulated Locking Screw, Fully Threaded 105mm
040140110 6.5mm Cannulated Locking Screw, Fully Threaded 110mm
040140115 6.5mm Cannulated Locking Screw, Fully Threaded 115mm
040140120 6.5mm Cannulated Locking Screw, Fully Threaded 120mm
040240006 2.4mm Locking Screw, Self-tapping, Stardrive 6mm
040240007 2.4mm Locking Screw, Self-tapping, Stardrive 7mm
040240008 2.4mm Locking Screw, Self-tapping, Stardrive 8mm
040240009 2.4mm Locking Screw, Self-tapping, Stardrive 9mm
040240010 2.4mm Locking Screw, Self-tapping, Stardrive 10mm
040240011 2.4mm Locking Screw, Self-tapping, Stardrive 11mm
040240012 2.4mm Locking Screw, Self-tapping, Stardrive 12mm
040240013 2.4mm Locking Screw, Self-tapping, Stardrive 13mm
040240014 2.4mm Locking Screw, Self-tapping, Stardrive 14mm
040240015 2.4mm Locking Screw, Self-tapping, Stardrive 15mm
040240016 2.4mm Locking Screw, Self-tapping, Stardrive 16mm
040240017 2.4mm Locking Screw, Self-tapping, Stardrive 17mm
040240018 2.4mm Locking Screw, Self-tapping, Stardrive 18mm
040240019 2.4mm Locking Screw, Self-tapping, Stardrive 19mm
040240020 2.4mm Locking Screw, Self-tapping, Stardrive 20mm
040240022 2.4mm Locking Screw, Self-tapping, Stardrive 22mm
040240024 2.4mm Locking Screw, Self-tapping, Stardrive 24mm
040240026 2.4mm Locking Screw, Self-tapping, Stardrive 26mm
040240028 2.4mm Locking Screw, Self-tapping, Stardrive 28mm
040240030 2.4mm Locking Screw, Self-tapping, Stardrive 30mm
040260070 7.3mm Cannulated Locking Screw, Thread Length 16mm 70mm
040260075 7.3mm Cannulated Locking Screw, Thread Length 16mm 75mm

040260080 7.3mm Cannulated Locking Screw, Thread Length 16mm 80mm

040260085 7.3mm Cannulated Locking Screw, Thread Length 16mm 85mm

040260090 7.3mm Cannulated Locking Screw, Thread Length 16mm 90mm

040260095 7.3mm Cannulated Locking Screw, Thread Length 16mm 95mm

040260100 7.3mm Cannulated Locking Screw, Thread Length 16mm

100mm

040260105 7.3mm Cannulated Locking Screw, Thread Length 16mm

105mm

040260110 7.3mm Cannulated Locking Screw, Thread Length 16mm

110mm

040270070 7.3mm Cannulated Locking Screw, Thread Length 32mm 70mm

040270075 7.3mm Cannulated Locking Screw, Thread Length 32mm 75mm

040270080 7.3mm Cannulated Locking Screw, Thread Length 32mm 80mm

040270085 7.3mm Cannulated Locking Screw, Thread Length 32mm 85mm

040270090 7.3mm Cannulated Locking Screw, Thread Length 32mm 90mm

040270095 7.3mm Cannulated Locking Screw, Thread Length 32mm 95mm

040270100 7.3mm Cannulated Locking Screw, Thread Length 32mm

100mm

040270105 7.3mm Cannulated Locking Screw, Thread Length 32mm

105mm

040270110 7.3mm Cannulated Locking Screw, Thread Length 32mm

110mm

040270115 7.3mm Cannulated Locking Screw, Thread Length 32mm

115mm

040300010 5.0mm Locking Screw, Self-tapping, Flat Head 10mm

040300014 5.0mm Locking Screw, Self-tapping, Flat Head 14mm

040300018 5.0mm Locking Screw, Self-tapping, Flat Head 18mm

040760008 2.4mm Locking Screw , variable angle, stardrive, self-tapping 8mm

040760010 2.4mm Locking Screw , variable angle, stardrive, self-tapping 10mm
040760012 2.4mm Locking Screw , variable angle, stardrive, self-tapping 12mm
040760014 2.4mm Locking Screw , variable angle, stardrive, self-tapping 14mm
040760016 2.4mm Locking Screw , variable angle, stardrive, self-tapping 16mm
040760018 2.4mm Locking Screw , variable angle, stardrive, self-tapping 18mm
040760020 2.4mm Locking Screw , variable angle, stardrive, self-tapping 20mm
040760022 2.4mm Locking Screw , variable angle, stardrive, self-tapping 22mm
040760024 2.4mm Locking Screw , variable angle, stardrive, self-tapping 24mm
040760026 2.4mm Locking Screw , variable angle, stardrive, self-tapping 26mm
040760028 2.4mm Locking Screw , variable angle, stardrive, self-tapping 28mm
040760030 2.4mm Locking Screw , variable angle, stardrive, self-tapping 30mm
040910060 5.5mm Locking Screw, Self-tapping 60mm
040910065 5.5mm Locking Screw, Self-tapping 65mm
040910070 5.5mm Locking Screw, Self-tapping 70mm
040910075 5.5mm Locking Screw, Self-tapping 75mm
040910080 5.5mm Locking Screw, Self-tapping 80mm
040910085 5.5mm Locking Screw, Self-tapping 85mm
040910090 5.5mm Locking Screw, Self-tapping 90mm
040910095 5.5mm Locking Screw, Self-tapping 95mm
040910100 5.5mm Locking Screw, Self-tapping 100mm
040910105 5.5mm Locking Screw, Self-tapping 105mm
040910110 5.5mm Locking Screw, Self-tapping 110mm
040920060 6.5mm Cannulated Locking Screw, Self-tapping 60mm
040920065 6.5mm Cannulated Locking Screw, Self-tapping 65mm
040920070 6.5mm Cannulated Locking Screw, Self-tapping 70mm
040920075 6.5mm Cannulated Locking Screw, Self-tapping 75mm
040920080 6.5mm Cannulated Locking Screw, Self-tapping 80mm
040920085 6.5mm Cannulated Locking Screw, Self-tapping 85mm
040920090 6.5mm Cannulated Locking Screw, Self-tapping 90mm

040920095	6.5mm Cannulated Locking Screw, Self-tapping	95mm
040920100	6.5mm Cannulated Locking Screw, Self-tapping	100mm
040920105	6.5mm Cannulated Locking Screw, Self-tapping	105mm
040920110	6.5mm Cannulated Locking Screw, Self-tapping	110mm
040940014	4,0 mm Locking Screw, self-tapping	14 mm
040940016	4,0 mm Locking Screw, self-tapping	16 mm
040940018	4,0 mm Locking Screw, self-tapping	18 mm
040940020	4,0 mm Locking Screw, self-tapping	20 mm
040940022	4,0 mm Locking Screw, self-tapping	22 mm
040940024	4,0 mm Locking Screw, self-tapping	24 mm
040940026	4,0 mm Locking Screw, self-tapping	26 mm
040940028	4,0 mm Locking Screw, self-tapping	28 mm
040940030	4,0 mm Locking Screw, self-tapping	30 mm
040940032	4,0 mm Locking Screw, self-tapping	32 mm
040940034	4,0 mm Locking Screw, self-tapping	34 mm
040940036	4,0 mm Locking Screw, self-tapping	36 mm
040940038	4,0 mm Locking Screw, self-tapping	38 mm
040940040	4,0 mm Locking Screw, self-tapping	40 mm
040940045	4,0 mm Locking Screw, self-tapping	45 mm
040940050	4,0 mm Locking Screw, self-tapping	50 mm
040940055	4,0 mm Locking Screw, self-tapping	55 mm
040940060	4,0 mm Locking Screw, self-tapping	60 mm
040940065	4,0 mm Locking Screw, self-tapping	65 mm
040940070	4,0 mm Locking Screw, self-tapping	70 mm
040940075	4,0 mm Locking Screw, self-tapping	75 mm
040940080	4,0 mm Locking Screw, self-tapping	80 mm
041050014	1,8 mm Locking Buttress Screw, variable angle, stardrive	14 mm
041050016	1,8 mm Locking Buttress Screw, variable angle, stardrive	16 mm

041050018 1,8 mm Locking Buttress Screw, variable angle, stardrive 18 mm
041050020 1,8 mm Locking Buttress Screw, variable angle, stardrive 20 mm
041050022 1,8 mm Locking Buttress Screw, variable angle, stardrive 22 mm
041050024 1,8 mm Locking Buttress Screw, variable angle, stardrive 24 mm
041050026 1,8 mm Locking Buttress Screw, variable angle, stardrive 26 mm
041170010 MT 2.7mm Locking Screw , stardrive, self-tapping 10mm
041170012 MT 2.7mm Locking Screw , stardrive, self-tapping 12mm
041170014 MT 2.7mm Locking Screw , stardrive, self-tapping 14mm
041170016 MT 2.7mm Locking Screw , stardrive, self-tapping 16mm
041170018 MT 2.7mm Locking Screw , stardrive, self-tapping 18mm
041170020 MT 2.7mm Locking Screw , stardrive, self-tapping 20mm
041170022 MT 2.7mm Locking Screw , stardrive, self-tapping 22mm
041170024 MT 2.7mm Locking Screw , stardrive, self-tapping 24mm
041170026 MT 2.7mm Locking Screw , stardrive, self-tapping 26mm
041170028 MT 2.7mm Locking Screw , stardrive, self-tapping 28mm
041170030 MT 2.7mm Locking Screw , stardrive, self-tapping 30mm
041170032 MT 2.7mm Locking Screw , stardrive, self-tapping 32mm
041170034 MT 2.7mm Locking Screw , stardrive, self-tapping 34mm
041170036 MT 2.7mm Locking Screw , stardrive, self-tapping 36mm
041170038 MT 2.7mm Locking Screw , stardrive, self-tapping 38mm
041170040 MT 2.7mm Locking Screw , stardrive, self-tapping 40mm
041170042 MT 2.7mm Locking Screw , stardrive, self-tapping 42mm
041170044 MT 2.7mm Locking Screw , stardrive, self-tapping 44mm
041170046 MT 2.7mm Locking Screw , stardrive, self-tapping 46mm
041170048 MT 2.7mm Locking Screw , stardrive, self-tapping 48mm
041170050 MT 2.7mm Locking Screw , stardrive, self-tapping 50mm
041180010 MT 3.0mm Locking Screw, stardrive, self-tapping 10mm
041180012 MT 3.0mm Locking Screw, stardrive, self-tapping 12mm
041180014 MT 3.0mm Locking Screw, stardrive, self-tapping 14mm

041180016 MT 3.0mm Locking Screw, stardrive, self-tapping 16mm
041180018 MT 3.0mm Locking Screw, stardrive, self-tapping 18mm
041180020 MT 3.0mm Locking Screw, stardrive, self-tapping 20mm
041180022 MT 3.0mm Locking Screw, stardrive, self-tapping 22mm
041180024 MT 3.0mm Locking Screw, stardrive, self-tapping 24mm
041180026 MT 3.0mm Locking Screw, stardrive, self-tapping 26mm
041180028 MT 3.0mm Locking Screw, stardrive, self-tapping 28mm
041180030 MT 3.0mm Locking Screw, stardrive, self-tapping 30mm
041180032 MT 3.0mm Locking Screw, stardrive, self-tapping 32mm
041180034 MT 3.0mm Locking Screw, stardrive, self-tapping 34mm
041180036 MT 3.0mm Locking Screw, stardrive, self-tapping 36mm
041180038 MT 3.0mm Locking Screw, stardrive, self-tapping 38mm
041180040 MT 3.0mm Locking Screw, stardrive, self-tapping 40mm
041180042 MT 3.0mm Locking Screw, stardrive, self-tapping 42mm
041180044 MT 3.0mm Locking Screw, stardrive, self-tapping 44mm
041180046 MT 3.0mm Locking Screw, stardrive, self-tapping 46mm
041180048 MT 3.0mm Locking Screw, stardrive, self-tapping 48mm
041180050 MT 3.0mm Locking Screw, stardrive, self-tapping 50mm
041190010 MT 3.5mm Locking Screw, stardrive, self-tapping 10mm
041190012 MT 3.5mm Locking Screw, stardrive, self-tapping 12mm
041190014 MT 3.5mm Locking Screw, stardrive, self-tapping 14mm
041190016 MT 3.5mm Locking Screw, stardrive, self-tapping 16mm
041190018 MT 3.5mm Locking Screw, stardrive, self-tapping 18mm
041190020 MT 3.5mm Locking Screw, stardrive, self-tapping 20mm
041190022 MT 3.5mm Locking Screw, stardrive, self-tapping 22mm
041190024 MT 3.5mm Locking Screw, stardrive, self-tapping 24mm
041190026 MT 3.5mm Locking Screw, stardrive, self-tapping 26mm
041190028 MT 3.5mm Locking Screw, stardrive, self-tapping 28mm

041190030 MT 3.5mm Locking Screw, stardrive, self-tapping 30mm
041190032 MT 3.5mm Locking Screw, stardrive, self-tapping 32mm
041190034 MT 3.5mm Locking Screw, stardrive, self-tapping 34mm
041190036 MT 3.5mm Locking Screw, stardrive, self-tapping 36mm
041190038 MT 3.5mm Locking Screw, stardrive, self-tapping 38mm
041190040 MT 3.5mm Locking Screw, stardrive, self-tapping 40mm
041190042 MT 3.5mm Locking Screw, stardrive, self-tapping 42mm
041190044 MT 3.5mm Locking Screw, stardrive, self-tapping 44mm
041190046 MT 3.5mm Locking Screw, stardrive, self-tapping 46mm
041190048 MT 3.5mm Locking Screw, stardrive, self-tapping 48mm
041190050 MT 3.5mm Locking Screw, stardrive, self-tapping 50mm

Instrumental Asociado

Período de vida útil: NC

Forma de presentación: no estéril, por unidad

Condición de uso: Venta exclusiva a profesionales e instituciones sanitarias

Fabricante:

DOUBLE MEDICAL TECHNOLOGY INC.

Lugar de Elaboración:

No. 18, Shanbianhong East Road, Haicang District, Xiamen, 361026, P.R. China.

Expediente N° 1-47-3110-1051-19-0

Digitally signed by BELLOSO Waldo Horacio
Date: 2019.08.22 16:58:27 ART
Location: Ciudad Autónoma de Buenos Aires

Digitally signed by GESTION DOCUMENTAL ELECTRONICA -
GDE
DN: cn=GESTION DOCUMENTAL ELECTRONICA - GDE, cn=AR,
ou=SECRETARIA DE GOBIERNO DE MODERNIZACION,
ou=SECRETARIA DE MODERNIZACION ADMINISTRATIVA,
serialNumber=CUIT 30715117564
Date: 2019.08.22 16:58:29 -0300'




BIOPROT
IMPLANTES QUIRURGICOS


de bioimplant SRL

ANEXO IIIB

- INFORMACIONES DE LOS ROTULOS E INSTRUCCIONES DE USO DE PRODUCTOS MEDICOS

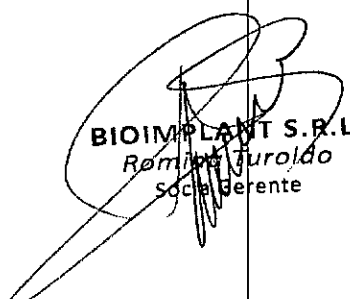

BIOIMPLANT S.R.L.
Román Jurolfo
Socio Gerente



Ma. VICTORIA COZZI
FARMACEUTICA
Mat. 2194






PLACAS BLOQUEADAS DE TITANIO


MODELO DE
ROTULO
ANEXO IIIB

REV.0
FECHA: 01/2019

- MODELO DE ROTULO CON LA INFORMACION DE LA DISPOSICION 2318/02
T.O. 2004 ANEXO IIIB ITEM 2

	FABRICANTE: Double Medical Technology Inc. - N° 18, Shanbianhong East Road Haicang District, Xiamen Fujian, CHINA 361026
NOMBRE GENÉRICO: Placas bloqueadas de titanio MARCA: DOUBLE MEDICAL MODELO: Según corresponda PRESENTACION: una unidad por envase	
REF (Número de serie)	LOT Según corresponda
 PRODUCTO DE UN SOLO USO	 ESTERILIZAR ANTES DE USAR AUTOCLAVE AL VACÍO
	El dispositivo debe ser almacenado en un lugar fresco, bien ventilado y seco a una temperatura por debajo de los 40°C y a una humedad relativa entre 30% y 80%. Proteger de la luz solar directa. Leer las instrucciones de uso antes de usar
VENTA EXCLUSIVA A PROFESIONALES E INSTITUCIONES SANITARIAS	
IMPORTADOR: BIOIMPLANT SRL Montevideo 567 (2000) Rosario - Santa Fe - Argentina // Tel./Fax: (0341) 4485178 / 5049	
Dirección Técnica: Farm. María Victoria Cozzi. Matrícula 2194	
Producto autorizado por la ANMAT PM-2342-1	


BIOIMPLANT S.R.L.
Romina Turoldo
Socia Gerente


Ma. VICTORIA COZZI
FARMACEUTICA
Mat. 2194

• **INSTRUCCIONES DE USO SEGUN DISPOSICION 2318/02 T.O. 2004 ,ANEXO IIIB,
ITEM 3 :**

**3.1. Las indicaciones contempladas en el ítem 2 de éste reglamento (Rótulo), salvo las
que figuran en los ítems 2.4 y 2.5:**

FABRICANTE: Double Medical Technology Inc. - N° 18, Shanbianhong East Road, Haicang District,
Xiamen Fujian, CHINA 361026

NOMBRE GENÉRICO: Placas bloqueadas de titanio

MARCA: DOUBLE MEDICAL

MODELO: Según corresponda

PRESENTACION: una unidad por envase

PRODUCTO DE UN SOLO USO ESTERILIZAR ANTES DE USAR POR AUTOCLAVE AL VACÍO

El dispositivo debe ser almacenado en un lugar fresco, bien ventilado y seco, a una temperatura por debajo de los 40°C y a una humedad relativa entre 30% y 80%. Proteger de la luz solar directa.

Leer las instrucciones de uso antes de usar

VENTA EXCLUSIVA A PROFESIONALES E INSTITUCIONES SANITARIAS

IMPORTADOR: BIOIMPLANT SRL Montevideo 567 (2000) Rosario - Santa Fe - Argentina // Tel./Fax:
(0341) 4485178 / 5049

Dirección Técnica: Farm. María Victoria Cozzi. Matrícula 2194

Producto autorizado por la ANMAT PM-2342-1

La familia de placas bloqueadas Double Medical importadas por Bioprot de Bioimplant SRL se compone de varios modelos que presentan envasados en forma unitaria, no estériles:

NºReferencia	Descripción	Perforac.	Material
020140003	T-Locking Plate, Combi-holes	3 Holes	Aleación de titanio
020140004	T-Locking Plate, Combi-holes	4 Holes	Aleación de titanio
020140005	T-Locking Plate, Combi-holes	5 Holes	Aleación de titanio
020140006	T-Locking Plate, Combi-holes	6 Holes	Aleación de titanio
020140007	T-Locking Plate, Combi-holes	7 Holes	Aleación de titanio
020140008	T-Locking Plate, Combi-holes	8 Holes	Aleación de titanio
020150004	T-Buttress Locking Plate, Combi-holes	4 Holes	Aleación de titanio
020150005	T-Buttress Locking Plate, Combi-holes	5 Holes	Aleación de titanio
020150006	T-Buttress Locking Plate, Combi-holes	6 Holes	Aleación de titanio
020150003	T-Buttress Locking Plate, Combi-holes	3 Holes	Aleación de titanio
020150007	T-Buttress Locking Plate, Combi-holes	7 Holes	Aleación de titanio
020150008	T-Buttress Locking Plate, Combi-holes	8 Holes	Aleación de titanio
020161003	L-Buttress Locking Plate left	3 Holes	Aleación de titanio
020161004	L-Buttress Locking Plate left	4 Holes	Aleación de titanio
020161005	L-Buttress Locking Plate left	5 Holes	Aleación de titanio
020161006	L-Buttress Locking Plate left	6 Holes	Aleación de titanio
020161007	L-Buttress Locking Plate left	7 Holes	Aleación de titanio
020161008	L-Buttress Locking Plate left	8 Holes	Aleación de titanio
020162003	L-Buttress Locking Plate, right	3 Holes	Aleación de titanio
020162004	L-Buttress Locking Plate, right	4 Holes	Aleación de titanio
020162005	L-Buttress Locking Plate, right	5 Holes	Aleación de titanio
020162006	L-Buttress Locking Plate, right	6 Holes	Aleación de titanio

BIOIMPLANT S.R.L.
 Romina Alfredo
 Socia Gerente

Ma. VICTORIA COZZI
 FARMACEUTICA
 Mat. 2194






PLACAS BLOQUEADAS DE TITANIO

INSTRUCCIONES DE
USO
ANEXO IIIB

REV.0
FECHA: 01/2019

020162007	L-Buttress Locking Plate, right		7 Holes	Aleación de titanio
020162008	L-Buttress Locking Plate, right		8 Holes	Aleación de titanio
020171004	Distal Femoral Sub-Condylar Locking Plate, Left		4 Holes	Aleación de titanio
020171005	Distal Femoral Sub-Condylar Locking Plate, Left		5 Holes	Aleación de titanio
020171006	Distal Femoral Sub-Condylar Locking Plate, Left		6 Holes	Aleación de titanio
020171007	Distal Femoral Sub-Condylar Locking Plate, Left		7 Holes	Aleación de titanio
020171008	Distal Femoral Sub-Condylar Locking Plate, Left		8 Holes	Aleación de titanio
020171009	Distal Femoral Sub-Condylar Locking Plate, Left		9 Holes	Aleación de titanio
020171010	Distal Femoral Sub-Condylar Locking Plate, Left		10 Holes	Aleación de titanio
020171011	Distal Femoral Sub-Condylar Locking Plate, Left		11 Holes	Aleación de titanio
020171012	Distal Femoral Sub-Condylar Locking Plate, Left		12 Holes	Aleación de titanio
020172004	Distal Femoral Sub-Condylar Locking Plate, Right		4 Holes	Aleación de titanio
020172005	Distal Femoral Sub-Condylar Locking Plate, Right		5 Holes	Aleación de titanio
020172006	Distal Femoral Sub-Condylar Locking Plate, Right		6 Holes	Aleación de titanio
020172007	Distal Femoral Sub-Condylar Locking Plate, Right		7 Holes	Aleación de titanio
020172008	Distal Femoral Sub-Condylar Locking Plate, Right		8 Holes	Aleación de titanio
020172009	Distal Femoral Sub-Condylar Locking Plate, Right		9 Holes	Aleación de titanio
020172010	Distal Femoral Sub-Condylar Locking Plate, Right		10 Holes	Aleación de titanio
020172011	Distal Femoral Sub-Condylar Locking Plate, Right		11 Holes	Aleación de titanio
020172012	Distal Femoral Sub-Condylar Locking Plate, Right		12 Holes	Aleación de titanio
020211060	Calcaneal Locking Plate I, Left		8 Holes	titanio sin alea
020211075	Calcaneal Locking Plate I, Left		10 Holes	titanio sin alea
020212060	Calcaneal Locking Plate I, Right		8 Holes	titanio sin alea
020212075	Calcaneal Locking Plate I, Right		10 Holes	titanio sin alea
020221015	Calcaneal Locking Plate II, Left		15 Holes	titanio sin alea
020222015	Calcaneal Locking Plate II, Right		15 Holes	titanio sin alea
020231015	Calcaneal Locking Plate II, Left		15 Holes	titanio sin alea
020232015	Calcaneal Locking Plate II, Right		15 Holes	titanio sin alea
020251003	Proximal Femoral Locking Plate, Multi-holes, Left		3 Holes	Aleación de titanio
020251004	Proximal Femoral Locking Plate, Multi-holes, Left		4 Holes	Aleación de titanio
020251005	Proximal Femoral Locking Plate, Multi-holes, Left		5 Holes	Aleación de titanio
020251006	Proximal Femoral Locking Plate, Multi-holes, Left		6 Holes	Aleación de titanio
020251007	Proximal Femoral Locking Plate, Multi-holes, Left		7 Holes	Aleación de titanio
020251008	Proximal Femoral Locking Plate, Multi-holes, Left		8 Holes	Aleación de titanio
020251009	Proximal Femoral Locking Plate, Multi-holes, Left		9 Holes	Aleación de titanio
020251010	Proximal Femoral Locking Plate, Multi-holes, Left		10 Holes	Aleación de titanio
020251011	Proximal Femoral Locking Plate, Multi-holes, Left		11 Holes	Aleación de titanio
020251012	Proximal Femoral Locking Plate, Multi-holes, Left		12 Holes	Aleación de titanio
020252003	Proximal Femoral Locking Plate, Multi-holes, Right		3 Holes	Aleación de titanio
020252004	Proximal Femoral Locking Plate, Multi-holes, Right		4 Holes	Aleación de titanio
020252005	Proximal Femoral Locking Plate, Multi-holes, Right		5 Holes	Aleación de titanio
020252006	Proximal Femoral Locking Plate, Multi-holes, Right		6 Holes	Aleación de titanio
020252007	Proximal Femoral Locking Plate, Multi-holes, Right		7 Holes	Aleación de titanio
020252008	Proximal Femoral Locking Plate, Multi-holes, Right		8 Holes	Aleación de titanio
020252009	Proximal Femoral Locking Plate, Multi-holes, Right		9 Holes	Aleación de titanio
020252010	Proximal Femoral Locking Plate, Multi-holes, Right		10 Holes	Aleación de titanio
020252011	Proximal Femoral Locking Plate, Multi-holes, Right		11 Holes	Aleación de titanio
020252012	Proximal Femoral Locking Plate, Multi-holes, Right		12 Holes	Aleación de titanio
020280004	Pelvic Reconstruction Locking Plate		4 Holes	titanio sin alea
020280005	Pelvic Reconstruction Locking Plate		5 Holes	titanio sin alea


BIOIMPLANT S.R.L.
 Román Buroldo
 Sonda Gerente

Ma. VICTORIA COZZI
 FARMACEUTICA
 Mat. 2194




BIOPROT de bioimplant SRL
 IMPLANTES QUIRURGICOS

PLACAS BLOQUEADAS DE TITANIO

**INSTRUCCIONES DE
 USO
 ANEXO IIIB**

REV.0
 FECHA: 01/2019

020280006	Pelvic Reconstruction Locking Plate	6 Holes	titanio sin alear
020280007	Pelvic Reconstruction Locking Plate	7 Holes	titanio sin alear
020280008	Pelvic Reconstruction Locking Plate	8 Holes	titanio sin alear
020280009	Pelvic Reconstruction Locking Plate	9 Holes	titanio sin alear
020280010	Pelvic Reconstruction Locking Plate	10 Holes	titanio sin alear
020280011	Pelvic Reconstruction Locking Plate	11 Holes	titanio sin alear
020280012	Pelvic Reconstruction Locking Plate	12 Holes	titanio sin alear
020280013	Pelvic Reconstruction Locking Plate	13 Holes	titanio sin alear
020280014	Pelvic Reconstruction Locking Plate	14 Holes	titanio sin alear
020320000	3.5mm Locking Spacer		Aleación de titanio
020330000	5.0mm Locking Spacer		Aleación de titanio
020360000	3.5mm Cable Spacer II		Aleación de titanio
020370000	5.0mm Cable Spacer II		Aleación de titanio
020380003	Proximal Humeral Locking Plate II, Combi-holes	3 Holes	titanio sin alear
020380004	Proximal Humeral Locking Plate II, Combi-holes	4 Holes	titanio sin alear
020380005	Proximal Humeral Locking Plate II, Combi-holes	5 Holes	titanio sin alear
020380006	Proximal Humeral Locking Plate II, Combi-holes	6 Holes	titanio sin alear
020380007	Proximal Humeral Locking Plate II, Combi-holes	7 Holes	titanio sin alear
020380008	Proximal Humeral Locking Plate II, Combi-holes	8 Holes	titanio sin alear
020380009	Proximal Humeral Locking Plate II, Combi-holes	9 Holes	titanio sin alear
020380010	Proximal Humeral Locking Plate II, Combi-holes	10 Holes	titanio sin alear
020380011	Proximal Humeral Locking Plate II, Combi-holes	11 Holes	titanio sin alear
020380012	Proximal Humeral Locking Plate II, Combi-holes	12 Holes	titanio sin alear
020380013	Proximal Humeral Locking Plate II, Combi-holes	13 Holes	titanio sin alear
020400006	Broad Limited Contact Locking Plate, Combi-holes	6 Holes	Aleación de titanio
020400007	Broad Limited Contact Locking Plate, Combi-holes	7 Holes	Aleación de titanio
020400008	Broad Limited Contact Locking Plate, Combi-holes	8 Holes	Aleación de titanio
020400009	Broad Limited Contact Locking Plate, Combi-holes	9 Holes	Aleación de titanio
020400010	Broad Limited Contact Locking Plate, Combi-holes	10 Holes	Aleación de titanio
020400011	Broad Limited Contact Locking Plate, Combi-holes	11 Holes	Aleación de titanio
020400012	Broad Limited Contact Locking Plate, Combi-holes	12 Holes	Aleación de titanio
020400013	Broad Limited Contact Locking Plate, Combi-holes	13 Holes	Aleación de titanio
020400014	Broad Limited Contact Locking Plate, Combi-holes	14 Holes	Aleación de titanio
020400015	Broad Limited Contact Locking Plate, Combi-holes	15 Holes	Aleación de titanio
020400016	Broad Limited Contact Locking Plate, Combi-holes	16 Holes	Aleación de titanio
020410004	Narrow Limited Contact Locking Plate, Combi-holes	4 Holes	Aleación de titanio
020410005	Narrow Limited Contact Locking Plate, Combi-holes	5 Holes	Aleación de titanio
020410006	Narrow Limited Contact Locking Plate, Combi-holes	6 Holes	Aleación de titanio
020410007	Narrow Limited Contact Locking Plate, Combi-holes	7 Holes	Aleación de titanio
020410008	Narrow Limited Contact Locking Plate, Combi-holes	8 Holes	Aleación de titanio
020410009	Narrow Limited Contact Locking Plate, Combi-holes	9 Holes	Aleación de titanio
020410010	Narrow Limited Contact Locking Plate, Combi-holes	10 Holes	Aleación de titanio
020410011	Narrow Limited Contact Locking Plate, Combi-holes	11 Holes	Aleación de titanio
020410012	Narrow Limited Contact Locking Plate, Combi-holes	12 Holes	Aleación de titanio
020410013	Narrow Limited Contact Locking Plate, Combi-holes	13 Holes	Aleación de titanio
020410014	Narrow Limited Contact Locking Plate, Combi-holes	14 Holes	Aleación de titanio
020410015	Narrow Limited Contact Locking Plate, Combi-holes	15 Holes	Aleación de titanio
020410016	Narrow Limited Contact Locking Plate, Combi-holes	16 Holes	Aleación de titanio
020421003	Proximal Femoral Locking Plate, Combi-holes, Left	3 Holes	Aleación de titanio
020421004	Proximal Femoral Locking Plate, Combi-holes, Left	4 Holes	Aleación de titanio

BIOIMPLANT S.R.L.
 Rosanna Turaldo
 Soňa Gerente

Ma. VICTORIA COZZI
 FARMACEUTICA
 Mat. 2194

020421005	Proximal Femoral Locking Plate, Combi-holes, Left	5 Holes	Aleación de titanio
020421006	Proximal Femoral Locking Plate, Combi-holes, Left	6 Holes	Aleación de titanio
020421007	Proximal Femoral Locking Plate, Combi-holes, Left	7 Holes	Aleación de titanio
020421008	Proximal Femoral Locking Plate, Combi-holes, Left	8 Holes	Aleación de titanio
020421009	Proximal Femoral Locking Plate, Combi-holes, Left	9 Holes	Aleación de titanio
020421010	Proximal Femoral Locking Plate, Combi-holes, Left	10 Holes	Aleación de titanio
020421011	Proximal Femoral Locking Plate, Combi-holes, Left	11 Holes	Aleación de titanio
020421012	Proximal Femoral Locking Plate, Combi-holes, Left	12 Holes	Aleación de titanio
020421013	Proximal Femoral Locking Plate, Combi-holes, Left	13 Holes	Aleación de titanio
020422003	Proximal Femoral Locking Plate, Combi-holes, Right	3 Holes	Aleación de titanio
020422004	Proximal Femoral Locking Plate, Combi-holes, Right	4 Holes	Aleación de titanio
020422005	Proximal Femoral Locking Plate, Combi-holes, Right	5 Holes	Aleación de titanio
020422006	Proximal Femoral Locking Plate, Combi-holes, Right	6 Holes	Aleación de titanio
020422007	Proximal Femoral Locking Plate, Combi-holes, Right	7 Holes	Aleación de titanio
020422008	Proximal Femoral Locking Plate, Combi-holes, Right	8 Holes	Aleación de titanio
020422009	Proximal Femoral Locking Plate, Combi-holes, Right	9 Holes	Aleación de titanio
020422010	Proximal Femoral Locking Plate, Combi-holes, Right	10 Holes	Aleación de titanio
020422011	Proximal Femoral Locking Plate, Combi-holes, Right	11 Holes	Aleación de titanio
020422012	Proximal Femoral Locking Plate, Combi-holes, Right	12 Holes	Aleación de titanio
020422013	Proximal Femoral Locking Plate, Combi-holes, Right	13 Holes	Aleación de titanio
020451003	Distal Medial Tibial Locking Plate, Combi-holes, Left	3 Holes	Aleación de titanio
020451004	Distal Medial Tibial Locking Plate, Combi-holes, Left	4 Holes	Aleación de titanio
020451005	Distal Medial Tibial Locking Plate, Combi-holes, Left	5 Holes	Aleación de titanio
020451006	Distal Medial Tibial Locking Plate, Combi-holes, Left	6 Holes	Aleación de titanio
020451007	Distal Medial Tibial Locking Plate, Combi-holes, Left	7 Holes	Aleación de titanio
020451008	Distal Medial Tibial Locking Plate, Combi-holes, Left	8 Holes	Aleación de titanio
020451009	Distal Medial Tibial Locking Plate, Combi-holes, Left	9 Holes	Aleación de titanio
020451010	Distal Medial Tibial Locking Plate, Combi-holes, Left	10 Holes	Aleación de titanio
020451011	Distal Medial Tibial Locking Plate, Combi-holes, Left	11 Holes	Aleación de titanio
020451012	Distal Medial Tibial Locking Plate, Combi-holes, Left	12 Holes	Aleación de titanio
020451013	Distal Medial Tibial Locking Plate, Combi-holes, Left	13 Holes	Aleación de titanio
020452003	Distal Medial Tibial Locking Plate, Combi-holes, Right	3 Holes	Aleación de titanio
020452004	Distal Medial Tibial Locking Plate, Combi-holes, Right	4 Holes	Aleación de titanio
020452005	Distal Medial Tibial Locking Plate, Combi-holes, Right	5 Holes	Aleación de titanio
020452006	Distal Medial Tibial Locking Plate, Combi-holes, Right	6 Holes	Aleación de titanio
020452007	Distal Medial Tibial Locking Plate, Combi-holes, Right	7 Holes	Aleación de titanio
020452008	Distal Medial Tibial Locking Plate, Combi-holes, Right	8 Holes	Aleación de titanio
020452009	Distal Medial Tibial Locking Plate, Combi-holes, Right	9 Holes	Aleación de titanio
020452010	Distal Medial Tibial Locking Plate, Combi-holes, Right	10 Holes	Aleación de titanio
020452011	Distal Medial Tibial Locking Plate, Combi-holes, Right	11 Holes	Aleación de titanio
020452012	Distal Medial Tibial Locking Plate, Combi-holes, Right	12 Holes	Aleación de titanio
020452013	Distal Medial Tibial Locking Plate, Combi-holes, Right	13 Holes	Aleación de titanio
020470004	Small Limited Contact Locking Plate, Combi-holes	4 Holes	Aleación de titanio
020470005	Small Limited Contact Locking Plate, Combi-holes	5 Holes	Aleación de titanio
020470006	Small Limited Contact Locking Plate, Combi-holes	6 Holes	Aleación de titanio
020470007	Small Limited Contact Locking Plate, Combi-holes	7 Holes	Aleación de titanio
020470008	Small Limited Contact Locking Plate, Combi-holes	8 Holes	Aleación de titanio
020470009	Small Limited Contact Locking Plate, Combi-holes	9 Holes	Aleación de titanio
020470010	Small Limited Contact Locking Plate, Combi-holes	10 Holes	Aleación de titanio
020470011	Small Limited Contact Locking Plate, Combi-holes	11 Holes	Aleación de titanio

020470012	Small Limited Contact Locking Plate, Combi-holes	12 Holes	Aleación de titanio
020470013	Small Limited Contact Locking Plate, Combi-holes	13 Holes	Aleación de titanio
020470014	Small Limited Contact Locking Plate, Combi-holes	14 Holes	Aleación de titanio
020480004	Reconstruction Locking Plate, combi-holes	4 Holes	Aleación de titanio
020480005	Reconstruction Locking Plate, combi-holes	5 Holes	Aleación de titanio
020480006	Reconstruction Locking Plate, combi-holes	6 Holes	Aleación de titanio
020480007	Reconstruction Locking Plate, combi-holes	7 Holes	Aleación de titanio
020480008	Reconstruction Locking Plate, combi-holes	8 Holes	Aleación de titanio
020480009	Reconstruction Locking Plate, combi-holes	9 Holes	Aleación de titanio
020480010	Reconstruction Locking Plate, combi-holes	10 Holes	Aleación de titanio
020480011	Reconstruction Locking Plate, combi-holes	11 Holes	Aleación de titanio
020480012	Reconstruction Locking Plate, combi-holes	12 Holes	Aleación de titanio
020480013	Reconstruction Locking Plate, combi-holes	13 Holes	Aleación de titanio
020480014	Reconstruction Locking Plate, combi-holes	14 Holes	Aleación de titanio
020490004	One-third Tubular Locking Plate, Combi-holes	4 Holes	Aleación de titanio
020490005	One-third Tubular Locking Plate, Combi-holes	5 Holes	Aleación de titanio
020490006	One-third Tubular Locking Plate, Combi-holes	6 Holes	Aleación de titanio
020490007	One-third Tubular Locking Plate, Combi-holes	7 Holes	Aleación de titanio
020490008	One-third Tubular Locking Plate, Combi-holes	8 Holes	Aleación de titanio
020490009	One-third Tubular Locking Plate, Combi-holes	9 Holes	Aleación de titanio
020490010	One-third Tubular Locking Plate, Combi-holes	10 Holes	Aleación de titanio
020490011	One-third Tubular Locking Plate, Combi-holes	11 Holes	Aleación de titanio
020490012	One-third Tubular Locking Plate, Combi-holes	12 Holes	Aleación de titanio
020501004	Clavicular Locking Plate, Combi-holes, Left	4 Holes	Aleación de titanio
020501005	Clavicular Locking Plate, Combi-holes, Left	5 Holes	Aleación de titanio
020501006	Clavicular Locking Plate, Combi-holes, Left	6 Holes	Aleación de titanio
020501007	Clavicular Locking Plate, Combi-holes, Left	7 Holes	Aleación de titanio
020501008	Clavicular Locking Plate, Combi-holes, Left	8 Holes	Aleación de titanio
020501009	Clavicular Locking Plate, Combi-holes, Left	9 Holes	Aleación de titanio
020501010	Clavicular Locking Plate, Combi-holes, Left	10 Holes	Aleación de titanio
020501011	Clavicular Locking Plate, Combi-holes, Left	11 Holes	Aleación de titanio
020501012	Clavicular Locking Plate, Combi-holes, Left	12 Holes	Aleación de titanio
020502004	Clavicular Locking Plate, Combi-holes, Right	4 Holes	Aleación de titanio
020502005	Clavicular Locking Plate, Combi-holes, Right	5 Holes	Aleación de titanio
020502006	Clavicular Locking Plate, Combi-holes, Right	6 Holes	Aleación de titanio
020502007	Clavicular Locking Plate, Combi-holes, Right	7 Holes	Aleación de titanio
020502008	Clavicular Locking Plate, Combi-holes, Right	8 Holes	Aleación de titanio
020502009	Clavicular Locking Plate, Combi-holes, Right	9 Holes	Aleación de titanio
020502010	Clavicular Locking Plate, Combi-holes, Right	10 Holes	Aleación de titanio
020502011	Clavicular Locking Plate, Combi-holes, Right	11 Holes	Aleación de titanio
020502012	Clavicular Locking Plate, Combi-holes, Right	12 Holes	Aleación de titanio
020511002	Distal Fibular Locking Plate, Combi-holes, Left	2 Holes	Aleación de titanio
020511003	Distal Fibular Locking Plate, Combi-holes, Left	3 Holes	Aleación de titanio
020511004	Distal Fibular Locking Plate, Combi-holes, Left	4 Holes	Aleación de titanio
020511005	Distal Fibular Locking Plate, Combi-holes, Left	5 Holes	Aleación de titanio
020511006	Distal Fibular Locking Plate, Combi-holes, Left	6 Holes	Aleación de titanio
020511007	Distal Fibular Locking Plate, Combi-holes, Left	7 Holes	Aleación de titanio
020511008	Distal Fibular Locking Plate, Combi-holes, Left	8 Holes	Aleación de titanio
020511009	Distal Fibular Locking Plate, Combi-holes, Left	9 Holes	Aleación de titanio
020512002	Distal Fibular Locking Plate, Combi-holes, Right	2 Holes	Aleación de titanio



BIOPROT de bioimplant SRL
 IMPLANTES Y INSTRUMENTOS

PLACAS BLOQUEADAS DE TITANIO

**INSTRUCCIONES DE
 USO
 ANEXO IIIB**

REV.0
 FECHA: 01/2019

020512003	Distal Fibular Locking Plate, Combi-holes, Right	3 Holes	Aleación de titanio
020512004	Distal Fibular Locking Plate, Combi-holes, Right	4 Holes	Aleación de titanio
020512005	Distal Fibular Locking Plate, Combi-holes, Right	5 Holes	Aleación de titanio
020512006	Distal Fibular Locking Plate, Combi-holes, Right	6 Holes	Aleación de titanio
020512007	Distal Fibular Locking Plate, Combi-holes, Right	7 Holes	Aleación de titanio
020512008	Distal Fibular Locking Plate, Combi-holes, Right	8 Holes	Aleación de titanio
020512009	Distal Fibular Locking Plate, Combi-holes, Right	9 Holes	Aleación de titanio
020520003	Olecranon Locking Plate III, Combi-holes	3 Holes	Aleación de titanio
020520004	Olecranon Locking Plate III, Combi-holes	4 Holes	Aleación de titanio
020520005	Olecranon Locking Plate III, Combi-holes	5 Holes	Aleación de titanio
020520006	Olecranon Locking Plate III, Combi-holes	6 Holes	Aleación de titanio
020520007	Olecranon Locking Plate III, Combi-holes	7 Holes	Aleación de titanio
020520008	Olecranon Locking Plate III, Combi-holes	8 Holes	Aleación de titanio
020541000	Proximal Femoral Locking Plate for 7.3mm Cannulated Locking Screw, Left	51mm	Aleación de titanio
020542000	Proximal Femoral Locking Plate for 7.3mm Cannulated Locking Screw, Right	51mm	Aleación de titanio
020621004	Distal Humeral Sub-Condylar Locking Plate, Left	4 Holes	titanio sin alear
020621007	Distal Humeral Sub-Condylar Locking Plate, Left	7 Holes	titanio sin alear
020621010	Distal Humeral Sub-Condylar Locking Plate, Left	10 Holes	titanio sin alear
020622004	Distal Humeral Sub-Condylar Locking Plate, Right	4 Holes	titanio sin alear
020622007	Distal Humeral Sub-Condylar Locking Plate, Right	7 Holes	titanio sin alear
020622010	Distal Humeral Sub-Condylar Locking Plate, Right	10 Holes	titanio sin alear
020631002	Distal Radius Locking Plate I, Small, Combi-holes, Left	2 Holes	titanio sin alear
020631003	Distal Radius Locking Plate I, Small, Combi-holes, Left	3 Holes	titanio sin alear
020631004	Distal Radius Locking Plate I, Small, Combi-holes, Left	4 Holes	titanio sin alear
020631005	Distal Radius Locking Plate I, Small, Combi-holes, Left	5 Holes	titanio sin alear
020632002	Distal Radius Locking Plate I, Small, Combi-holes, Right	2 Holes	titanio sin alear
020632003	Distal Radius Locking Plate I, Small, Combi-holes, Right	3 Holes	titanio sin alear
020632004	Distal Radius Locking Plate I, Small, Combi-holes, Right	4 Holes	titanio sin alear
020632005	Distal Radius Locking Plate I, Small, Combi-holes, Right	5 Holes	titanio sin alear
020641002	Distal Radius Locking Plate I, Large, Combi-holes, Left	2 Holes	titanio sin alear
020641003	Distal Radius Locking Plate I, Large, Combi-holes, Left	3 Holes	titanio sin alear
020641004	Distal Radius Locking Plate I, Large, Combi-holes, Left	4 Holes	titanio sin alear
020641005	Distal Radius Locking Plate I, Large, Combi-holes, Left	5 Holes	titanio sin alear
020642002	Distal Radius Locking Plate I, Large, Combi-holes, Right	2 Holes	titanio sin alear
020642003	Distal Radius Locking Plate I, Large, Combi-holes, Right	3 Holes	titanio sin alear
020642004	Distal Radius Locking Plate I, Large, Combi-holes, Right	4 Holes	titanio sin alear
020642005	Distal Radius Locking Plate I, Large, Combi-holes, Right	5 Holes	titanio sin alear
020651002	Proximal Radius Locking Plate, left	2 Holes	titanio sin alear
020651003	Proximal Radius Locking Plate, left	3 Holes	titanio sin alear
020651004	Proximal Radius Locking Plate, left	4 Holes	titanio sin alear
020652002	Proximal Radius Locking Plate, right	2 Holes	titanio sin alear
020652003	Proximal Radius Locking Plate, right	3 Holes	titanio sin alear
020652004	Proximal Radius Locking Plate, right	4 Holes	titanio sin alear
020661003	Distal Clavicular Locking Plate, Combi-holes, Left	3 Holes	titanio sin alear
020661004	Distal Clavicular Locking Plate, Combi-holes, Left	4 Holes	titanio sin alear
020661005	Distal Clavicular Locking Plate, Combi-holes, Left	5 Holes	titanio sin alear
020661006	Distal Clavicular Locking Plate, Combi-holes, Left	6 Holes	titanio sin alear
020661007	Distal Clavicular Locking Plate, Combi-holes, Left	7 Holes	titanio sin alear
020661008	Distal Clavicular Locking Plate, Combi-holes, Left	8 Holes	titanio sin alear
020662003	Distal Clavicular Locking Plate, Combi-holes, Right	3 Holes	titanio sin alear

BIOIMPLANT S.R.L.
 Romina Pardo
 Socia Gerente

Ma. VICTORIA COZZI
 FARMACEUTICA
 Mat. 2194

121

020662004	Distal Clavicular Locking Plate, Combi-holes, Right	4 Holes	titanio sin alear
020662005	Distal Clavicular Locking Plate, Combi-holes, Right	5 Holes	titanio sin alear
020662006	Distal Clavicular Locking Plate, Combi-holes, Right	6 Holes	titanio sin alear
020662007	Distal Clavicular Locking Plate, Combi-holes, Right	7 Holes	titanio sin alear
020662008	Distal Clavicular Locking Plate, Combi-holes, Right	8 Holes	titanio sin alear
020811003	Clavicular Hook Locking Plate III, Height 15mm, Combi-holes, Left	3 Holes	titanio sin alear
020811004	Clavicular Hook Locking Plate III, Height 15mm, Combi-holes, Left	4 Holes	titanio sin alear
020811005	Clavicular Hook Locking Plate III, Height 15mm, Combi-holes, Left	5 Holes	titanio sin alear
020811006	Clavicular Hook Locking Plate III, Height 15mm, Combi-holes, Left	6 Holes	titanio sin alear
020811007	Clavicular Hook Locking Plate III, Height 15mm, Combi-holes, Left	7 Holes	titanio sin alear
020811008	Clavicular Hook Locking Plate III, Height 15mm, Combi-holes, Left	8 Holes	titanio sin alear
020812003	Clavicular Hook Locking Plate III, Height 15mm, Combi-holes, Right	3 Holes	titanio sin alear
020812004	Clavicular Hook Locking Plate III, Height 15mm, Combi-holes, Right	4 Holes	titanio sin alear
020812005	Clavicular Hook Locking Plate III, Height 15mm, Combi-holes, Right	5 Holes	titanio sin alear
020812006	Clavicular Hook Locking Plate III, Height 15mm, Combi-holes, Right	6 Holes	titanio sin alear
020812007	Clavicular Hook Locking Plate III, Height 15mm, Combi-holes, Right	7 Holes	titanio sin alear
020812008	Clavicular Hook Locking Plate III, Height 15mm, Combi-holes, Right	8 Holes	titanio sin alear
020821003	Clavicular Hook Locking Plate III, Height 18mm, Combi-holes, Left	3 Holes	titanio sin alear
020821004	Clavicular Hook Locking Plate III, Height 18mm, Combi-holes, Left	4 Holes	titanio sin alear
020821005	Clavicular Hook Locking Plate III, Height 18mm, Combi-holes, Left	5 Holes	titanio sin alear
020821006	Clavicular Hook Locking Plate III, Height 18mm, Combi-holes, Left	6 Holes	titanio sin alear
020821007	Clavicular Hook Locking Plate III, Height 18mm, Combi-holes, Left	7 Holes	titanio sin alear
020821008	Clavicular Hook Locking Plate III, Height 18mm, Combi-holes, Left	8 Holes	titanio sin alear
020822003	Clavicular Hook Locking Plate III, Height 18mm, Combi-holes, Right	3 Holes	titanio sin alear
020822004	Clavicular Hook Locking Plate III, Height 18mm, Combi-holes, Right	4 Holes	titanio sin alear
020822005	Clavicular Hook Locking Plate III, Height 18mm, Combi-holes, Right	5 Holes	titanio sin alear
020822006	Clavicular Hook Locking Plate III, Height 18mm, Combi-holes, Right	6 Holes	titanio sin alear
020822007	Clavicular Hook Locking Plate III, Height 18mm, Combi-holes, Right	7 Holes	titanio sin alear
020822008	Clavicular Hook Locking Plate III, Height 18mm, Combi-holes, Right	8 Holes	titanio sin alear
020961005	Proximal Humeral Greater Tuberosity Locking Compression Plate, left	5 Holes	titanio sin alear
020962005	Proximal Humeral Greater Tuberosity Locking Compression Plate, right	5 Holes	titanio sin alear
020971000	Proximal Femoral Locking Plate, for 7.3mm Cannulated Locking Screw, small, left	43mm	Aleación de titanio
020972000	Proximal Femoral Locking Plate, for 7.3mm Cannulated Locking Screw, small, right	43mm	Aleación de titanio
021301003	Distal Lateral Femoral Condylar Locking Plate, left	3 Holes	Aleación de titanio
021302003	Distal Lateral Femoral Condylar Locking Plate, right	3 Holes	Aleación de titanio
021311003	Distal Medial Femoral Condylar Locking Plate, left	3 Holes	Aleación de titanio
021312003	Distal Medial Femoral Condylar Locking Plate, right	3 Holes	Aleación de titanio
021321004	Distal Medial Femoral Sub-condylar Locking Plate, large, left	4 Holes	Aleación de titanio
021321007	Distal Medial Femoral Sub-condylar Locking Plate, large, left	7 Holes	Aleación de titanio
021322004	Distal Medial Femoral Sub-condylar Locking Plate, large, right	4 Holes	Aleación de titanio
021322007	Distal Medial Femoral Sub-condylar Locking Plate, large, right	7 Holes	Aleación de titanio
021331004	Distal Medial Femoral Sub-condylar Locking Plate, small, left	4 Holes	Aleación de titanio
021331007	Distal Medial Femoral Sub-condylar Locking Plate, small, left	7 Holes	Aleación de titanio
021332004	Distal Medial Femoral Sub-condylar Locking Plate, small, right	4 Holes	Aleación de titanio
021332007	Distal Medial Femoral Sub-condylar Locking Plate, small, right	7 Holes	Aleación de titanio
021371003	Distal Lateral Humeral Locking Plate, Combi-holes, Left	3 Holes	Aleación de titanio
021371004	Distal Lateral Humeral Locking Plate, Combi-holes, Left	4 Holes	Aleación de titanio
021371005	Distal Lateral Humeral Locking Plate, Combi-holes, Left	5 Holes	Aleación de titanio
021371006	Distal Lateral Humeral Locking Plate, Combi-holes, Left	6 Holes	Aleación de titanio
021371007	Distal Lateral Humeral Locking Plate, Combi-holes, Left	7 Holes	Aleación de titanio

021371008	Distal Lateral Humeral Locking Plate, Combi-holes, Left	8 Holes	Aleación de titanio
021371009	Distal Lateral Humeral Locking Plate, Combi-holes, Left	9 Holes	Aleación de titanio
021372003	Distal Lateral Humeral Locking Plate, Combi-holes, Right	3 Holes	Aleación de titanio
021372004	Distal Lateral Humeral Locking Plate, Combi-holes, Right	4 Holes	Aleación de titanio
021372005	Distal Lateral Humeral Locking Plate, Combi-holes, Right	5 Holes	Aleación de titanio
021372006	Distal Lateral Humeral Locking Plate, Combi-holes, Right	6 Holes	Aleación de titanio
021372007	Distal Lateral Humeral Locking Plate, Combi-holes, Right	7 Holes	Aleación de titanio
021372008	Distal Lateral Humeral Locking Plate, Combi-holes, Right	8 Holes	Aleación de titanio
021372009	Distal Lateral Humeral Locking Plate, Combi-holes, Right	9 Holes	Aleación de titanio
021381003	Distal Medial Humeral Locking Plate, Combi-holes, Left	3 Holes	Aleación de titanio
021381004	Distal Medial Humeral Locking Plate, Combi-holes, Left	4 Holes	Aleación de titanio
021381005	Distal Medial Humeral Locking Plate, Combi-holes, Left	5 Holes	Aleación de titanio
021381006	Distal Medial Humeral Locking Plate, Combi-holes, Left	6 Holes	Aleación de titanio
021381007	Distal Medial Humeral Locking Plate, Combi-holes, Left	7 Holes	Aleación de titanio
021381008	Distal Medial Humeral Locking Plate, Combi-holes, Left	8 Holes	Aleación de titanio
021381009	Distal Medial Humeral Locking Plate, Combi-holes, Left	9 Holes	Aleación de titanio
021382003	Distal Medial Humeral Locking Plate, Combi-holes, Right	3 Holes	Aleación de titanio
021382004	Distal Medial Humeral Locking Plate, Combi-holes, Right	4 Holes	Aleación de titanio
021382005	Distal Medial Humeral Locking Plate, Combi-holes, Right	5 Holes	Aleación de titanio
021382006	Distal Medial Humeral Locking Plate, Combi-holes, Right	6 Holes	Aleación de titanio
021382007	Distal Medial Humeral Locking Plate, Combi-holes, Right	7 Holes	Aleación de titanio
021382008	Distal Medial Humeral Locking Plate, Combi-holes, Right	8 Holes	Aleación de titanio
021382009	Distal Medial Humeral Locking Plate, Combi-holes, Right	9 Holes	Aleación de titanio
021390003	T-Locking Plate, Right Angled, Head 3 Holes, Combi-holes	3 Holes	Aleación de titanio
021400004	T-Locking Plate, Right Angled, Head 4 Holes, Combi-holes	4 Holes	Aleación de titanio
021390005	T-Locking Plate, Right Angled, Head 3 Holes, Combi-holes	5 Holes	Aleación de titanio
021400006	T-Locking Plate, Right Angled, Head 4 Holes, Combi-holes	6 Holes	Aleación de titanio
021411003	T-Locking Plate, Oblique Angled, Combi-holes, Left	3 Holes	Aleación de titanio
021411004	T-Locking Plate, Oblique Angled, Combi-holes, Left	4 Holes	Aleación de titanio
021411005	T-Locking Plate, Oblique Angled, Combi-holes, Left	5 Holes	Aleación de titanio
021411006	T-Locking Plate, Oblique Angled, Combi-holes, Left	6 Holes	Aleación de titanio
021412003	T-Locking Plate, Oblique Angled, Combi-holes, Right	3 Holes	Aleación de titanio
021412004	T-Locking Plate, Oblique Angled, Combi-holes, Right	4 Holes	Aleación de titanio
021412005	T-Locking Plate, Oblique Angled, Combi-holes, Right	5 Holes	Aleación de titanio
021412006	T-Locking Plate, Oblique Angled, Combi-holes, Right	6 Holes	Aleación de titanio
021421003	Olecranon Locking Plate II, Combi-holes, Left	3 Holes	Aleación de titanio
021421004	Olecranon Locking Plate II, Combi-holes, Left	4 Holes	Aleación de titanio
021421005	Olecranon Locking Plate II, Combi-holes, Left	5 Holes	Aleación de titanio
021421006	Olecranon Locking Plate II, Combi-holes, Left	6 Holes	Aleación de titanio
021421007	Olecranon Locking Plate II, Combi-holes, Left	7 Holes	Aleación de titanio
021421008	Olecranon Locking Plate II, Combi-holes, Left	8 Holes	Aleación de titanio
021422003	Olecranon Locking Plate II, Combi-holes, Right	3 Holes	Aleación de titanio
021422004	Olecranon Locking Plate II, Combi-holes, Right	4 Holes	Aleación de titanio
021422005	Olecranon Locking Plate II, Combi-holes, Right	5 Holes	Aleación de titanio
021422006	Olecranon Locking Plate II, Combi-holes, Right	6 Holes	Aleación de titanio
021422007	Olecranon Locking Plate II, Combi-holes, Right	7 Holes	Aleación de titanio
021422008	Olecranon Locking Plate II, Combi-holes, Right	8 Holes	Aleación de titanio
021431005	Distal Femoral Locking Plate II, Combi-holes, Left	5 Holes	Aleación de titanio
021431006	Distal Femoral Locking Plate II, Combi-holes, Left	6 Holes	Aleación de titanio
021431007	Distal Femoral Locking Plate II, Combi-holes, Left	7 Holes	Aleación de titanio



BIOPROT de bioimplant S.R.L.
FARMACIA SUBSECCION

PLACAS BLOQUEADAS DE TITANIO

**INSTRUCCIONES DE
USO
ANEXO IIIB**

REV.0
FECHA: 01/2019

021431008	Distal Femoral Locking Plate II, Combi-holes, Left	8 Holes	Aleación de titanio
021431009	Distal Femoral Locking Plate II, Combi-holes, Left	9 Holes	Aleación de titanio
021431010	Distal Femoral Locking Plate II, Combi-holes, Left	10 Holes	Aleación de titanio
021431011	Distal Femoral Locking Plate II, Combi-holes, Left	11 Holes	Aleación de titanio
021431012	Distal Femoral Locking Plate II, Combi-holes, Left	12 Holes	Aleación de titanio
021431013	Distal Femoral Locking Plate II, Combi-holes, Left	13 Holes	Aleación de titanio
021432005	Distal Femoral Locking Plate II, Combi-holes, Right	5 Holes	Aleación de titanio
021432006	Distal Femoral Locking Plate II, Combi-holes, Right	6 Holes	Aleación de titanio
021432007	Distal Femoral Locking Plate II, Combi-holes, Right	7 Holes	Aleación de titanio
021432008	Distal Femoral Locking Plate II, Combi-holes, Right	8 Holes	Aleación de titanio
021432009	Distal Femoral Locking Plate II, Combi-holes, Right	9 Holes	Aleación de titanio
021432010	Distal Femoral Locking Plate II, Combi-holes, Right	10 Holes	Aleación de titanio
021432011	Distal Femoral Locking Plate II, Combi-holes, Right	11 Holes	Aleación de titanio
021432012	Distal Femoral Locking Plate II, Combi-holes, Right	12 Holes	Aleación de titanio
021432013	Distal Femoral Locking Plate II, Combi-holes, Right	13 Holes	Aleación de titanio
021441003	Lateral Tibial Head Buttress Locking Plate II, Combi-holes, Left	3 Holes	Aleación de titanio
021441004	Lateral Tibial Head Buttress Locking Plate II, Combi-holes, Left	4 Holes	Aleación de titanio
021441005	Lateral Tibial Head Buttress Locking Plate II, Combi-holes, Left	5 Holes	Aleación de titanio
021441006	Lateral Tibial Head Buttress Locking Plate II, Combi-holes, Left	6 Holes	Aleación de titanio
021441007	Lateral Tibial Head Buttress Locking Plate II, Combi-holes, Left	7 Holes	Aleación de titanio
021441008	Lateral Tibial Head Buttress Locking Plate II, Combi-holes, Left	8 Holes	Aleación de titanio
021441009	Lateral Tibial Head Buttress Locking Plate II, Combi-holes, Left	9 Holes	Aleación de titanio
021441010	Lateral Tibial Head Buttress Locking Plate II, Combi-holes, Left	10 Holes	Aleación de titanio
021441011	Lateral Tibial Head Buttress Locking Plate II, Combi-holes, Left	11 Holes	Aleación de titanio
021441012	Lateral Tibial Head Buttress Locking Plate II, Combi-holes, Left	12 Holes	Aleación de titanio
021441013	Lateral Tibial Head Buttress Locking Plate II, Combi-holes, Left	13 Holes	Aleación de titanio
021442003	Lateral Tibial Head Buttress Locking Plate II, Combi-holes, Right	3 Holes	Aleación de titanio
021442004	Lateral Tibial Head Buttress Locking Plate II, Combi-holes, Right	4 Holes	Aleación de titanio
021442005	Lateral Tibial Head Buttress Locking Plate II, Combi-holes, Right	5 Holes	Aleación de titanio
021442006	Lateral Tibial Head Buttress Locking Plate II, Combi-holes, Right	6 Holes	Aleación de titanio
021442007	Lateral Tibial Head Buttress Locking Plate II, Combi-holes, Right	7 Holes	Aleación de titanio
021442008	Lateral Tibial Head Buttress Locking Plate II, Combi-holes, Right	8 Holes	Aleación de titanio
021442009	Lateral Tibial Head Buttress Locking Plate II, Combi-holes, Right	9 Holes	Aleación de titanio
021442010	Lateral Tibial Head Buttress Locking Plate II, Combi-holes, Right	10 Holes	Aleación de titanio
021442011	Lateral Tibial Head Buttress Locking Plate II, Combi-holes, Right	11 Holes	Aleación de titanio
021442012	Lateral Tibial Head Buttress Locking Plate II, Combi-holes, Right	12 Holes	Aleación de titanio
021442013	Lateral Tibial Head Buttress Locking Plate II, Combi-holes, Right	13 Holes	Aleación de titanio
021450006	1.5S-Mini Wed Locking Plate	34.4mm	Aleación de titanio
021460008	2.0S-Mini Wed Locking Plate	47.7mm	Aleación de titanio
021470010	2.4S-Mini Wed Locking Plate	77.2mm	Aleación de titanio
021501004	Low Profile Distal medial Tibial Locking Plate, Minimally Invasive Technique left	4 Holes	Aleación de titanio
021501006	Low Profile Distal medial Tibial Locking Plate, Minimally Invasive Technique left	6 Holes	Aleación de titanio
021501008	Low Profile Distal medial Tibial Locking Plate, Minimally Invasive Technique left	8 Holes	Aleación de titanio
021501010	Low Profile Distal medial Tibial Locking Plate, Minimally Invasive Technique left	10 Holes	Aleación de titanio
021501012	Low Profile Distal medial Tibial Locking Plate, Minimally Invasive Technique left	12 Holes	Aleación de titanio
021501014	Low Profile Distal medial Tibial Locking Plate, Minimally Invasive Technique left	14 Holes	Aleación de titanio
021502004	Low Profile Distal medial Tibial Locking Plate, Minimally Invasive Technique right	4 Holes	Aleación de titanio
021502006	Low Profile Distal medial Tibial Locking Plate, Minimally Invasive Technique right	6 Holes	Aleación de titanio
021502008	Low Profile Distal medial Tibial Locking Plate, Minimally Invasive Technique right	8 Holes	Aleación de titanio
021502010	Low Profile Distal medial Tibial Locking Plate, Minimally Invasive Technique right	10 Holes	Aleación de titanio

BIOIMPLANT S.R.L.
Román Urdola
Social Berente

Ma. VICTORIA COZZI
FARMACEUTICA
Mat. 2194

124

021502012	Low Profile Distal medial Tibial Locking Plate, Minimally Invasive Technique right	12 Holes	Aleación de titanio
021502014	Low Profile Distal medial Tibial Locking Plate, Minimally Invasive Technique right	14 Holes	Aleación de titanio
021510008	Curved Locking Compression Plate, broad	8 Holes	Aleación de titanio
021510009	Curved Locking Compression Plate, broad	9 Holes	Aleación de titanio
021510010	Curved Locking Compression Plate, broad	10 Holes	Aleación de titanio
021510011	Curved Locking Compression Plate, broad	11 Holes	Aleación de titanio
021510012	Curved Locking Compression Plate, broad	12 Holes	Aleación de titanio
021510013	Curved Locking Compression Plate, broad	13 Holes	Aleación de titanio
021510014	Curved Locking Compression Plate, broad	14 Holes	Aleación de titanio
021510015	Curved Locking Compression Plate, broad	15 Holes	Aleación de titanio
021510016	Curved Locking Compression Plate, broad	16 Holes	Aleación de titanio
021510017	Curved Locking Compression Plate, broad	17 Holes	Aleación de titanio
021510018	Curved Locking Compression Plate, broad	18 Holes	Aleación de titanio
021520002	Posterior Medial Proximal Tibial Plate, combi-holes	2 Holes	Aleación de titanio
021520004	Posterior Medial Proximal Tibial Plate, combi-holes	4 Holes	Aleación de titanio
021520006	Posterior Medial Proximal Tibial Plate, combi-holes	6 Holes	Aleación de titanio
021520008	Posterior Medial Proximal Tibial Plate, combi-holes	8 Holes	Aleación de titanio
021531002	Distal Radius Locking Plate I , variable angle, small, left	2 Holes	Aleación de titanio
021531003	Distal Radius Locking Plate I , variable angle, small, left	3 Holes	Aleación de titanio
021531004	Distal Radius Locking Plate I , variable angle, small, left	4 Holes	Aleación de titanio
021531005	Distal Radius Locking Plate I , variable angle, small, left	5 Holes	Aleación de titanio
021532002	Distal Radius Locking Plate I , variable angle, small, right	2 Holes	Aleación de titanio
021532003	Distal Radius Locking Plate I , variable angle, small, right	3 Holes	Aleación de titanio
021532004	Distal Radius Locking Plate I , variable angle, small, right	4 Holes	Aleación de titanio
021532005	Distal Radius Locking Plate I , variable angle, small, right	5 Holes	Aleación de titanio
021541002	Distal Radius Locking Plate I , variable angle, large, left	2 Holes	Aleación de titanio
021541003	Distal Radius Locking Plate I , variable angle, large, left	3 Holes	Aleación de titanio
021541004	Distal Radius Locking Plate I , variable angle, large, left	4 Holes	Aleación de titanio
021541005	Distal Radius Locking Plate I , variable angle, large, left	5 Holes	Aleación de titanio
021542002	Distal Radius Locking Plate I , variable angle, large, right	2 Holes	Aleación de titanio
021542003	Distal Radius Locking Plate I , variable angle, large, right	3 Holes	Aleación de titanio
021542004	Distal Radius Locking Plate I , variable angle, large, right	4 Holes	Aleación de titanio
021542005	Distal Radius Locking Plate I , variable angle, large, right	5 Holes	Aleación de titanio
021550006	Cable Locking Plate	6 Holes	Aleación de titanio
021550008	Cable Locking Plate	8 Holes	Aleación de titanio
021550010	Cable Locking Plate	10 Holes	Aleación de titanio
021561002	Cable Grip Locking Plate, Left	2 Holes	Aleación de titanio
021561004	Cable Grip Locking Plate, Left	4 Holes	Aleación de titanio
021561006	Cable Grip Locking Plate, Left	6 Holes	Aleación de titanio
021561008	Cable Grip Locking Plate, Left	8 Holes	Aleación de titanio
021561010	Cable Grip Locking Plate, Left	10 Holes	Aleación de titanio
021561012	Cable Grip Locking Plate, Left	12 Holes	Aleación de titanio
021562002	Cable Grip Locking Plate, Right	2 Holes	Aleación de titanio
021562004	Cable Grip Locking Plate, Right	4 Holes	Aleación de titanio
021562006	Cable Grip Locking Plate, Right	6 Holes	Aleación de titanio
021562008	Cable Grip Locking Plate, Right	8 Holes	Aleación de titanio
021562010	Cable Grip Locking Plate, Right	10 Holes	Aleación de titanio
021562012	Cable Grip Locking Plate, Right	12 Holes	Aleación de titanio
021601003	Superior Anterior Distal Clavicle Locking Plate, with lateral extension, combi-holes, left	3 Holes	Aleación de titanio

021601004	Superior Anterior Distal Clavicle Locking Plate, with lateral extension, combi-holes, left	4 Holes	Aleación de titanio
021601005	Superior Anterior Distal Clavicle Locking Plate, with lateral extension, combi-holes, left	5 Holes	Aleación de titanio
021601006	Superior Anterior Distal Clavicle Locking Plate, with lateral extension, combi-holes, left	6 Holes	Aleación de titanio
021602003	Superior Anterior Distal Clavicle Locking Plate, with lateral extension, combi-holes, right	3 Holes	Aleación de titanio
021602004	Superior Anterior Distal Clavicle Locking Plate, with lateral extension, combi-holes, right	4 Holes	Aleación de titanio
021602005	Superior Anterior Distal Clavicle Locking Plate, with lateral extension, combi-holes, right	5 Holes	Aleación de titanio
021602006	Superior Anterior Distal Clavicle Locking Plate, with lateral extension, combi-holes, right	6 Holes	Aleación de titanio
021611004	Distal Medial Radius Locking Plate, left	4 Holes	Aleación de titanio
021612004	Distal Medial Radius Locking Plate, right	4 Holes	Aleación de titanio
021620005	Distal Lateral Radius Locking Plate	5 Holes	Aleación de titanio
021620006	Distal Lateral Radius Locking Plate	6 Holes	Aleación de titanio
021631003	Proximal Femoral Locking Plate III, left	3 Holes	Aleación de titanio
021631004	Proximal Femoral Locking Plate III, Left	4 Holes	Aleación de titanio
021631005	Proximal Femoral Locking Plate III, Left	5 Holes	Aleación de titanio
021631006	Proximal Femoral Locking Plate III, Left	6 Holes	Aleación de titanio
021631007	Proximal Femoral Locking Plate III, Left	7 Holes	Aleación de titanio
021631008	Proximal Femoral Locking Plate III, Left	8 Holes	Aleación de titanio
021631009	Proximal Femoral Locking Plate III, Left	9 Holes	Aleación de titanio
021631010	Proximal Femoral Locking Plate III, Left	10 Holes	Aleación de titanio
021631011	Proximal Femoral Locking Plate III, Left	11 Holes	Aleación de titanio
021631012	Proximal Femoral Locking Plate III, Left	12 Holes	Aleación de titanio
021631013	Proximal Femoral Locking Plate III, Left	13 Holes	Aleación de titanio
021632003	Proximal Femoral Locking Plate III, Right	3 Holes	Aleación de titanio
021632004	Proximal Femoral Locking Plate III, Right	4 Holes	Aleación de titanio
021632005	Proximal Femoral Locking Plate III, Right	5 Holes	Aleación de titanio
021632006	Proximal Femoral Locking Plate III, Right	6 Holes	Aleación de titanio
021632007	Proximal Femoral Locking Plate III, Right	7 Holes	Aleación de titanio
021632008	Proximal Femoral Locking Plate III, Right	8 Holes	Aleación de titanio
021632009	Proximal Femoral Locking Plate III, Right	9 Holes	Aleación de titanio
021632010	Proximal Femoral Locking Plate III, Right	10 Holes	Aleación de titanio
021632011	Proximal Femoral Locking Plate III, Right	11 Holes	Aleación de titanio
021632012	Proximal Femoral Locking Plate III, Right	12 Holes	Aleación de titanio
021632013	Proximal Femoral Locking Plate III, Right	13 Holes	Aleación de titanio
021641004	Ankle Fusion Locking Plate, left	4 Holes	Aleación de titanio
021642004	Ankle Fusion Locking Plate, right	4 Holes	Aleación de titanio
021671005	Proximal Lateral Tibial Locking Plate II left	5 Holes	Aleación de titanio
021671009	Proximal Lateral Tibial Locking Plate II left	9 Holes	Aleación de titanio
021671013	Proximal Lateral Tibial Locking Plate II left	13 Holes	Aleación de titanio
021671017	Proximal Lateral Tibial Locking Plate II left	17 Holes	Aleación de titanio
021672005	Proximal Lateral Tibial Locking Plate II right	5 Holes	Aleación de titanio
021672009	Proximal Lateral Tibial Locking Plate II right	9 Holes	Aleación de titanio
021672013	Proximal Lateral Tibial Locking Plate II right	13 Holes	Aleación de titanio
021672017	Proximal Lateral Tibial Locking Plate II right	17 Holes	Aleación de titanio
021681004	Proximal Medial Tibial Locking Plate II left	4 Holes	Aleación de titanio
021681006	Proximal Medial Tibial Locking Plate II left	6 Holes	Aleación de titanio
021681008	Proximal Medial Tibial Locking Plate II left	8 Holes	Aleación de titanio

021682004	Proximal Medial Tibial Locking Plate II right	4 Holes	Aleación de titanio
021682006	Proximal Medial Tibial Locking Plate II right	6 Holes	Aleación de titanio
021682008	Proximal Medial Tibial Locking Plate II right	8 Holes	Aleación de titanio
021691003	Proximal Posterior Tibial Locking Plate left, II	3 Holes	Aleación de titanio
021691005	Proximal Posterior Tibial Locking Plate left, II	5 Holes	Aleación de titanio
021691007	Proximal Posterior Tibial Locking Plate left, II	7 Holes	Aleación de titanio
021692003	Proximal Posterior Tibial Locking Plate right, II	3 Holes	Aleación de titanio
021692005	Proximal Posterior Tibial Locking Plate right, II	5 Holes	Aleación de titanio
021692007	Proximal Posterior Tibial Locking Plate right, II	7 Holes	Aleación de titanio
021800004	Humeral Shaft Straight Locking Plate	4 Holes	Aleación de titanio
021800005	Humeral Shaft Straight Locking Plate	5 Holes	Aleación de titanio
021800006	Humeral Shaft Straight Locking Plate	6 Holes	Aleación de titanio
021800007	Humeral Shaft Straight Locking Plate	7 Holes	Aleación de titanio
021800008	Humeral Shaft Straight Locking Plate	8 Holes	Aleación de titanio
021800009	Humeral Shaft Straight Locking Plate	9 Holes	Aleación de titanio
021800010	Humeral Shaft Straight Locking Plate	10 Holes	Aleación de titanio
021800011	Humeral Shaft Straight Locking Plate	11 Holes	Aleación de titanio
021800012	Humeral Shaft Straight Locking Plate	12 Holes	Aleación de titanio
021800013	Humeral Shaft Straight Locking Plate	13 Holes	Aleación de titanio
021800014	Humeral Shaft Straight Locking Plate	14 Holes	Aleación de titanio
021800015	Humeral Shaft Straight Locking Plate	15 Holes	Aleación de titanio
021800016	Humeral Shaft Straight Locking Plate	16 Holes	Aleación de titanio
021810127	Sacroiliac Joint Plate II	127mm	Aleación de titanio
021810137	Sacroiliac Joint Plate II	137mm	Aleación de titanio
021810147	Sacroiliac Joint Plate II	147mm	Aleación de titanio
021821003	Proximal Lateral Tibial Locking Plate, left	3 Holes	Aleación de titanio
021821004	Proximal Lateral Tibial Locking Plate, left	4 Holes	Aleación de titanio
021821005	Proximal Lateral Tibial Locking Plate, left	5 Holes	Aleación de titanio
021821006	Proximal Lateral Tibial Locking Plate, left	6 Holes	Aleación de titanio
021821007	Proximal Lateral Tibial Locking Plate, left	7 Holes	Aleación de titanio
021821008	Proximal Lateral Tibial Locking Plate, left	8 Holes	Aleación de titanio
021821009	Proximal Lateral Tibial Locking Plate, left	9 Holes	Aleación de titanio
021821010	Proximal Lateral Tibial Locking Plate, left	10 Holes	Aleación de titanio
021821011	Proximal Lateral Tibial Locking Plate, left	11 Holes	Aleación de titanio
021821012	Proximal Lateral Tibial Locking Plate, left	12 Holes	Aleación de titanio
021821013	Proximal Lateral Tibial Locking Plate, left	13 Holes	Aleación de titanio
021821014	Proximal Lateral Tibial Locking Plate, left	14 Holes	Aleación de titanio
021821015	Proximal Lateral Tibial Locking Plate, left	15 Holes	Aleación de titanio
021822003	Proximal Lateral Tibial Locking Plate, right	3 Holes	Aleación de titanio
021822004	Proximal Lateral Tibial Locking Plate, right	4 Holes	Aleación de titanio
021822005	Proximal Lateral Tibial Locking Plate, right	5 Holes	Aleación de titanio
021822006	Proximal Lateral Tibial Locking Plate, right	6 Holes	Aleación de titanio
021822007	Proximal Lateral Tibial Locking Plate, right	7 Holes	Aleación de titanio
021822008	Proximal Lateral Tibial Locking Plate, right	8 Holes	Aleación de titanio
021822009	Proximal Lateral Tibial Locking Plate, right	9 Holes	Aleación de titanio
021822010	Proximal Lateral Tibial Locking Plate, right	10 Holes	Aleación de titanio
021822011	Proximal Lateral Tibial Locking Plate, right	11 Holes	Aleación de titanio
021822012	Proximal Lateral Tibial Locking Plate, right	12 Holes	Aleación de titanio
021822013	Proximal Lateral Tibial Locking Plate, right	13 Holes	Aleación de titanio
021822014	Proximal Lateral Tibial Locking Plate, right	14 Holes	Aleación de titanio



BIOPROT de bioimplant SRL
 INFLUENTES QUEJAS

PLACAS BLOQUEADAS DE TITANIO

**INSTRUCCIONES DE
 USO
 ANEXO IIIB**

REV.0
 FECHA: 01/2019

021822015	Proximal Lateral Tibial Locking Plate, right		15 Holes	Aleación de titanio
021881003	Posterior Lateral Fibular Locking Plate combi	holes, left	3 Holes	Aleación de titanio
021881004	Posterior Lateral Fibular Locking Plate combi	holes, left	4 Holes	Aleación de titanio
021881005	Posterior Lateral Fibular Locking Plate combi	holes, left	5 Holes	Aleación de titanio
021881006	Posterior Lateral Fibular Locking Plate combi	holes, left	6 Holes	Aleación de titanio
021881007	Posterior Lateral Fibular Locking Plate combi	holes, left	7 Holes	Aleación de titanio
021881008	Posterior Lateral Fibular Locking Plate combi	holes, left	8 Holes	Aleación de titanio
021882003	Posterior Lateral Fibular Locking Plate combi	holes, right	3 Holes	Aleación de titanio
021882004	Posterior Lateral Fibular Locking Plate combi	holes, right	4 Holes	Aleación de titanio
021882005	Posterior Lateral Fibular Locking Plate combi	holes, right	5 Holes	Aleación de titanio
021882006	Posterior Lateral Fibular Locking Plate combi	holes, right	6 Holes	Aleación de titanio
021882007	Posterior Lateral Fibular Locking Plate combi	holes, right	7 Holes	Aleación de titanio
021882008	Posterior Lateral Fibular Locking Plate combi	holes, right	8 Holes	Aleación de titanio
021891003	Distal Lateral Humeral Locking Plate, head 2	7 hole, combi-holes, Left	3 Holes	Aleación de titanio
021891004	Distal Lateral Humeral Locking Plate, head 2	7 hole, combi-holes, Left	4 Holes	Aleación de titanio
021891005	Distal Lateral Humeral Locking Plate, head 2	7 hole, combi-holes, Left	5 Holes	Aleación de titanio
021891006	Distal Lateral Humeral Locking Plate, head 2	7 hole, combi-holes, Left	6 Holes	Aleación de titanio
021891007	Distal Lateral Humeral Locking Plate, head 2	7 hole, combi-holes, Left	7 Holes	Aleación de titanio
021891008	Distal Lateral Humeral Locking Plate, head 2	7 hole, combi-holes, Left	8 Holes	Aleación de titanio
021891009	Distal Lateral Humeral Locking Plate, head 2	7 hole, combi-holes, Left	9 Holes	Aleación de titanio
021892003	Distal Lateral Humeral Locking Plate, head 2	7 hole, combi-holes, Right	3 Holes	Aleación de titanio
021892004	Distal Lateral Humeral Locking Plate, head 2	7 hole, combi-holes, Right	4 Holes	Aleación de titanio
021892005	Distal Lateral Humeral Locking Plate, head 2	7 hole, combi-holes, Right	5 Holes	Aleación de titanio
021892006	Distal Lateral Humeral Locking Plate, head 2	7 hole, combi-holes, Right	6 Holes	Aleación de titanio
021892007	Distal Lateral Humeral Locking Plate, head 2	7 hole, combi-holes, Right	7 Holes	Aleación de titanio
021892008	Distal Lateral Humeral Locking Plate, head 2	7 hole, combi-holes, Right	8 Holes	Aleación de titanio
021892009	Distal Lateral Humeral Locking Plate, head 2	7 hole, combi-holes, Right	9 Holes	Aleación de titanio
021901003	Distal Medial Humeral Locking Plate, head 2	7 hole, combi-holes, Left	3 Holes	Aleación de titanio
021901004	Distal Medial Humeral Locking Plate, head 2	7 hole, combi-holes, Left	4 Holes	Aleación de titanio
021901005	Distal Medial Humeral Locking Plate, head 2	7 hole, combi-holes, Left	5 Holes	Aleación de titanio
021901006	Distal Medial Humeral Locking Plate, head 2	7 hole, combi-holes, Left	6 Holes	Aleación de titanio
021901007	Distal Medial Humeral Locking Plate, head 2	7 hole, combi-holes, Left	7 Holes	Aleación de titanio
021901008	Distal Medial Humeral Locking Plate, head 2	7 hole, combi-holes, Left	8 Holes	Aleación de titanio
021901009	Distal Medial Humeral Locking Plate, head 2	7 hole, combi-holes, Left	9 Holes	Aleación de titanio
021902003	Distal Medial Humeral Locking Plate, head 2	7 hole, combi-holes, Right	3 Holes	Aleación de titanio
021902004	Distal Medial Humeral Locking Plate, head 2	7 hole, combi-holes, Right	4 Holes	Aleación de titanio
021902005	Distal Medial Humeral Locking Plate, head 2	7 hole, combi-holes, Right	5 Holes	Aleación de titanio
021902006	Distal Medial Humeral Locking Plate, head 2	7 hole, combi-holes, Right	6 Holes	Aleación de titanio
021902007	Distal Medial Humeral Locking Plate, head 2	7 hole, combi-holes, Right	7 Holes	Aleación de titanio
021902008	Distal Medial Humeral Locking Plate, head 2	7 hole, combi-holes, Right	8 Holes	Aleación de titanio
021902009	Distal Medial Humeral Locking Plate, head 2	7 hole, combi-holes, Right	9 Holes	Aleación de titanio
022100008	Wrist Fusion Locking Plate, length 117mm, s	tandard, combi-holes	8 Holes	Aleación de titanio
022110008	Wrist Fusion Locking Plate, length 117mm, s	hort, combi-holes	8 Holes	Aleación de titanio
022120009	Wrist Fusion Locking Plate, length 112mm, s	traight, combi-holes	9 Holes	Aleación de titanio
022181002	Distal Radius Locking Plate, head 7 holes v	ariable angle, left	2 Holes	Aleación de titanio
022181003	Distal Radius Locking Plate, head 7 holes v	ariable angle, left	3 Holes	Aleación de titanio
022181004	Distal Radius Locking Plate, head 7 holes v	ariable angle, left	4 Holes	Aleación de titanio
022181005	Distal Radius Locking Plate, head 7 holes v	ariable angle, left	5 Holes	Aleación de titanio
022182002	Distal Radius Locking Plate, head 7 holes v	ariable angle, right	2 Holes	Aleación de titanio
022182003	Distal Radius Locking Plate, head 7 holes v	ariable angle, right	3 Holes	Aleación de titanio

BIOIMPLANT S.R.L.
 Roxana Turalde
 Jefa Gerente

Ma. VICTORIA COZZI
 FARMACEUTICA
 Mat. 2194

022182004	Distal Radius Locking Plate, head 7 holes variable angle, right	4 Holes	Aleación de titanio
022182005	Distal Radius Locking Plate, head 7 holes variable angle, right	5 Holes	Aleación de titanio
022360004	Proximal Medial Tibial Osteotomy Locking Plate, combi-holes	4 Holes	Aleación de titanio
022981003	Distal Lateral Femoral Osteotomy Locking Plate, combi-holes, left	3 Holes	Aleación de titanio
022981004	Distal Lateral Femoral Osteotomy Locking Plate, combi-holes, left	4 Holes	Aleación de titanio
022981005	Distal Lateral Femoral Osteotomy Locking Plate, combi-holes, left	5 Holes	Aleación de titanio
022981006	Distal Lateral Femoral Osteotomy Locking Plate, combi-holes, left	6 Holes	Aleación de titanio
022981007	Distal Lateral Femoral Osteotomy Locking Plate, combi-holes, left	7 Holes	Aleación de titanio
022981008	Distal Lateral Femoral Osteotomy Locking Plate, combi-holes, left	8 Holes	Aleación de titanio
022981009	Distal Lateral Femoral Osteotomy Locking Plate, combi-holes, left	9 Holes	Aleación de titanio
022981010	Distal Lateral Femoral Osteotomy Locking Plate, combi-holes, left	10 Holes	Aleación de titanio
022981011	Distal Lateral Femoral Osteotomy Locking Plate, combi-holes, left	11 Holes	Aleación de titanio
022981012	Distal Lateral Femoral Osteotomy Locking Plate, combi-holes, left	12 Holes	Aleación de titanio
022981013	Distal Lateral Femoral Osteotomy Locking Plate, combi-holes, left	13 Holes	Aleación de titanio
022982003	Distal Lateral Femoral Osteotomy Locking Plate, combi-holes, right	3 Holes	Aleación de titanio
022982004	Distal Lateral Femoral Osteotomy Locking Plate, combi-holes, right	4 Holes	Aleación de titanio
022982005	Distal Lateral Femoral Osteotomy Locking Plate, combi-holes, right	5 Holes	Aleación de titanio
022982006	Distal Lateral Femoral Osteotomy Locking Plate, combi-holes, right	6 Holes	Aleación de titanio
022982007	Distal Lateral Femoral Osteotomy Locking Plate, combi-holes, right	7 Holes	Aleación de titanio
022982008	Distal Lateral Femoral Osteotomy Locking Plate, combi-holes, right	8 Holes	Aleación de titanio
022982009	Distal Lateral Femoral Osteotomy Locking Plate, combi-holes, right	9 Holes	Aleación de titanio
022982010	Distal Lateral Femoral Osteotomy Locking Plate, combi-holes, right	10 Holes	Aleación de titanio
022982011	Distal Lateral Femoral Osteotomy Locking Plate, combi-holes, right	11 Holes	Aleación de titanio
022982012	Distal Lateral Femoral Osteotomy Locking Plate, combi-holes, right	12 Holes	Aleación de titanio
022982013	Distal Lateral Femoral Osteotomy Locking Plate, combi-holes, right	13 Holes	Aleación de titanio
022991003	Proximal Lateral Tibial Osteotomy Locking Plate, left	3 Holes	Aleación de titanio
022991005	Proximal Lateral Tibial Osteotomy Locking Plate, left	5 Holes	Aleación de titanio
022991007	Proximal Lateral Tibial Osteotomy Locking Plate, left	7 Holes	Aleación de titanio
022991009	Proximal Lateral Tibial Osteotomy Locking Plate, left	9 Holes	Aleación de titanio
022991011	Proximal Lateral Tibial Osteotomy Locking Plate, left	11 Holes	Aleación de titanio
022991013	Proximal Lateral Tibial Osteotomy Locking Plate, left	13 Holes	Aleación de titanio
022992003	Proximal Lateral Tibial Osteotomy Locking Plate, right	3 Holes	Aleación de titanio
022992005	Proximal Lateral Tibial Osteotomy Locking Plate, right	5 Holes	Aleación de titanio
022992007	Proximal Lateral Tibial Osteotomy Locking Plate, right	7 Holes	Aleación de titanio
022992009	Proximal Lateral Tibial Osteotomy Locking Plate, right	9 Holes	Aleación de titanio
022992011	Proximal Lateral Tibial Osteotomy Locking Plate, right	11 Holes	Aleación de titanio
022992013	Proximal Lateral Tibial Osteotomy Locking Plate, right	13 Holes	Aleación de titanio
025331000	Coronoid Locking Plate, Left	left	Aleación de titanio
025332000	Coronoid Locking Plate, Right	right	Aleación de titanio
025341003	Distal Lateral Tibial Locking Plate II, combi-holes, left	3 Holes	Aleación de titanio
025341005	Distal Lateral Tibial Locking Plate II, combi-holes, Left	5Holes	Aleación de titanio
025341007	Distal Lateral Tibial Locking Plate II, combi-holes, left	7 Holes	Aleación de titanio
025341009	Distal Lateral Tibial Locking Plate II, combi-holes, Left	9Holes	Aleación de titanio
025341011	Distal Lateral Tibial Locking Plate II, combi-holes, left	11 Holes	Aleación de titanio
025341013	Distal Lateral Tibial Locking Plate II, combi-holes, Left	13Holes	Aleación de titanio
025341015	Distal Lateral Tibial Locking Plate II, combi-holes, left	15 Holes	Aleación de titanio
025341017	Distal Lateral Tibial Locking Plate II, combi-holes, Left	17Holes	Aleación de titanio
025342003	Distal Lateral Tibial Locking Plate II, combi-holes, right	3 Holes	Aleación de titanio
025342005	Distal Lateral Tibial Locking Plate II, combi-holes, Right	5Holes	Aleación de titanio
025342007	Distal Lateral Tibial Locking Plate II, combi-holes, right	7 Holes	Aleación de titanio



BIOPROT de bioimplant S.R.L.
 IMPLANTES QUIRÚRGICOS

PLACAS BLOQUEADAS DE TITANIO

**INSTRUCCIONES DE
 USO
 ANEXO IIIB**

REV.0
 FECHA: 01/2019

025342009	Distal Lateral Tibial Locking Plate II , combi-holes,Right	9Holes	Aleación de titanio
025342011	Distal Lateral Tibial Locking Plate II , combi-holes,right	11 Holes	Aleación de titanio
025342013	Distal Lateral Tibial Locking Plate II , combi-holes,Right	13Holes	Aleación de titanio
025342015	Distal Lateral Tibial Locking Plate II , combi-holes,right	15 Holes	Aleación de titanio
025342017	Distal Lateral Tibial Locking Plate II , combi-holes,Right	17Holes	Aleación de titanio
025351008	Distal Diaphyseal Metacarpal Radius Locking Plate, left	8 Holes	Aleación de titanio
025351010	Distal Diaphyseal Metacarpal Radius Locking Plate, left	10 Holes	Aleación de titanio
025351012	Distal Diaphyseal Metacarpal Radius Locking Plate, left	12 Holes	Aleación de titanio
025352008	Distal Diaphyseal Metacarpal Radius Locking Plate, right	8 Holes	Aleación de titanio
025352010	Distal Diaphyseal Metacarpal Radius Locking Plate, right	10 Holes	Aleación de titanio
025352012	Distal Diaphyseal Metacarpal Radius Locking Plate, right	12 Holes	Aleación de titanio
025551005	Proximal Humeral Greater Tuberosity Locking Compression Plate, suture hole, left	5 Holes	Aleación de titanio
025552005	Proximal Humeral Greater Tuberosity Locking Compression Plate, suture hole, right	5 Holes	Aleación de titanio
025681006	Superior Anterior Clavicle Locking Plate II , combi-holes, left	6 Holes	Aleación de titanio
025681007	Superior Anterior Clavicle Locking Plate II , combi-holes, left	7 Holes	Aleación de titanio
025681008	Superior Anterior Clavicle Locking Plate II , combi-holes, left	8 Holes	Aleación de titanio
025681009	Superior Anterior Clavicle Locking Plate II , combi-holes, left	9 Holes	Aleación de titanio
025681010	Superior Anterior Clavicle Locking Plate II , combi-holes, left	10 Holes	Aleación de titanio
025681011	Superior Anterior Clavicle Locking Plate II , combi-holes, left	11 Holes	Aleación de titanio
025681012	Superior Anterior Clavicle Locking Plate II , combi-holes, left	12 Holes	Aleación de titanio
025682006	Superior Anterior Clavicle Locking Plate II , combi-holes, right	6 Holes	Aleación de titanio
025682007	Superior Anterior Clavicle Locking Plate II , combi-holes, right	7 Holes	Aleación de titanio
025682008	Superior Anterior Clavicle Locking Plate II , combi-holes, right	8 Holes	Aleación de titanio
025682009	Superior Anterior Clavicle Locking Plate II , combi-holes, right	9 Holes	Aleación de titanio
025682010	Superior Anterior Clavicle Locking Plate II , combi-holes, right	10 Holes	Aleación de titanio
025682011	Superior Anterior Clavicle Locking Plate II , combi-holes, right	11 Holes	Aleación de titanio
025682012	Superior Anterior Clavicle Locking Plate II , combi-holes, right	12 Holes	Aleación de titanio
025850006	Metaphyseal Locking Plate, 3,5	6 Holes	Aleación de titanio
025850007	Metaphyseal Locking Plate, 3,5	7 Holes	Aleación de titanio
025850008	Metaphyseal Locking Plate, 3,5	8 Holes	Aleación de titanio
025850009	Metaphyseal Locking Plate, 3,5	9 Holes	Aleación de titanio
025850010	Metaphyseal Locking Plate, 3,5	10 Holes	Aleación de titanio
025850011	Metaphyseal Locking Plate, 3,5	11 Holes	Aleación de titanio
025850012	Metaphyseal Locking Plate, 3,5	12 Holes	Aleación de titanio
025850014	Metaphyseal Locking Plate, 3,5	14 Holes	Aleación de titanio
025870004	Distal Medial Humeral Reconstruction Locking Plate	4 Holes	Aleación de titanio
025870006	Distal Medial Humeral Reconstruction Locking Plate	6 Holes	Aleación de titanio
025870008	Distal Medial Humeral Reconstruction Locking Plate	8 Holes	Aleación de titanio
025870010	Distal Medial Humeral Reconstruction Locking Plate	10 Holes	Aleación de titanio
025870012	Distal Medial Humeral Reconstruction Locking Plate	12 Holes	Aleación de titanio
026012008	Olecranon Locking Plate, curved, riht	8 holes	Aleación de titanio
026020007	Olecranon Locking Plate, straight	7 holes	Aleación de titanio
026110004	MD Distal Tibial T-Locking Plate	4 holes	Aleación de titanio
026121004	MD Distal Tibial L-Locking Plate, left	4 holes	Aleación de titanio
026122004	MD Distal Tibial L-Locking Plate, right	4 holes	Aleación de titanio
026131006	Distal Diaphyseal Radius Locking Plate II, left	6 Holes	Aleación de titanio
026131008	Distal Diaphyseal Radius Locking Plate II, left	8 Holes	Aleación de titanio
026131010	Distal Diaphyseal Radius Locking Plate II, left	10 Holes	Aleación de titanio
026131012	Distal Diaphyseal Radius Locking Plate II, left	12 Holes	Aleación de titanio
026131014	Distal Diaphyseal Radius Locking Plate II, left	14 Holes	Aleación de titanio

BIOIMPLANT S.R.L.
 Román Ruoldo
 Social Gerente

Ma. VICTORIA COZZI
 FARMACEUTICA
 Mat. 2194



BIOPROT de bioimplant SRL
 IMPLANTES ORTOPÉDICOS

PLACAS BLOQUEADAS DE TITANIO

**INSTRUCCIONES DE
 USO
 ANEXO IIIB**

REV.0
 FECHA: 01/2019

026132006	Distal Diaphyseal Radius Locking Plate II, right	6 Holes	Aleación de titanio
026132008	Distal Diaphyseal Radius Locking Plate II, right	8 Holes	Aleación de titanio
026132010	Distal Diaphyseal Radius Locking Plate II, right	10 Holes	Aleación de titanio
026132012	Distal Diaphyseal Radius Locking Plate II, right	12 Holes	Aleación de titanio
026132014	Distal Diaphyseal Radius Locking Plate II, right	14 Holes	Aleación de titanio
026141000	Acromion Locking Plate, left		Aleación de titanio
026142000	Acromion Locking Plate, right		Aleación de titanio
026151000	Scapula Medial Locking Plate, left		Aleación de titanio
026152000	Scapula Medial Locking Plate, right		Aleación de titanio
026161000	Scapula Lateral Locking Plate, left		Aleación de titanio
026162000	Scapula Lateral Locking Plate, right		Aleación de titanio
026331000	Glenoid Locking Plate, left		Aleación de titanio
026332000	Glenoid Locking Plate, right		Aleación de titanio
026611012	Olecranon U-locking Plate, left	12 Holes	Aleación de titanio
026612012	Olecranon U-locking Plate, right	12 Holes	Aleación de titanio
026620006	MT Calcaneal Step Locking Plate		Aleación de titanio
026620008	MT Calcaneal Step Locking Plate		Aleación de titanio
026620010	MT Calcaneal Step Locking Plate		Aleación de titanio
026651000	MT First Metatarsal Osteotomy Locking Plate, left		Aleación de titanio
026651003	MT First Metatarsal Osteotomy Locking Plate, left	3 Holes	Aleación de titanio
026651004	MT First Metatarsal Osteotomy Locking Plate, left	4 Holes	Aleación de titanio
026651005	MT First Metatarsal Osteotomy Locking Plate, left	5 Holes	Aleación de titanio
026651006	MT First Metatarsal Osteotomy Locking Plate, left	6 Holes	Aleación de titanio
026651007	MT First Metatarsal Osteotomy Locking Plate, left	7 Holes	Aleación de titanio
026652000	MT First Metatarsal Osteotomy Locking Plate, right		Aleación de titanio
026652003	MT First Metatarsal Osteotomy Locking Plate, right	3 Holes	Aleación de titanio
026652004	MT First Metatarsal Osteotomy Locking Plate, right	4 Holes	Aleación de titanio
026652005	MT First Metatarsal Osteotomy Locking Plate, right	5 Holes	Aleación de titanio
026652006	MT First Metatarsal Osteotomy Locking Plate, right	6 Holes	Aleación de titanio
026652007	MT First Metatarsal Osteotomy Locking Plate, right	7 Holes	Aleación de titanio
026660000	MT Lateral Column Lengthening Locking Plate		Aleación de titanio
026660008	MT Lateral Column Lengthening Locking Plate		Aleación de titanio
026660010	MT Lateral Column Lengthening Locking Plate		Aleación de titanio
026660012	MT Lateral Column Lengthening Locking Plate		Aleación de titanio

Tornillos de titanio:

Las placas bloqueadas de titanio Double Medical se fijan exclusivamente con los tornillos de titanio Double Medical según modelo de placa y son los detallados a continuación:

030300006	1,5 mm Locking Screw, stardrive, self tapping	6 mm	Aleación de titanio
030300007	1,5 mm Locking Screw, stardrive, self tapping	7 mm	Aleación de titanio
030300008	1,5 mm Locking Screw, stardrive, self tapping	8 mm	Aleación de titanio
030300009	1,5 mm Locking Screw, stardrive, self tapping	9 mm	Aleación de titanio
030300010	1,5 mm Locking Screw, stardrive, self tapping	10 mm	Aleación de titanio
030300011	1,5 mm Locking Screw, stardrive, self tapping	11 mm	Aleación de titanio
030300012	1,5 mm Locking Screw, stardrive, self tapping	12 mm	Aleación de titanio
030300013	1,5 mm Locking Screw, stardrive, self tapping	13 mm	Aleación de titanio
030300014	1,5 mm Locking Screw, stardrive, self tapping	14 mm	Aleación de titanio
030300015	1,5 mm Locking Screw, stardrive, self tapping	15 mm	Aleación de titanio
030300016	1,5 mm Locking Screw, stardrive, self tapping	16 mm	Aleación de titanio

BIOIMPLANT S.R.L.
 Román Purodo
 Social Gerente

Ma. VICTORIA COZZI
 FARMACEUTICA
 Mat. 2194

030300017	1,5 mm Locking Screw, stardrive, self tapping	17 mm	Aleación de titanio
030300018	1,5 mm Locking Screw, stardrive, self tapping	18 mm	Aleación de titanio
030300019	1,5 mm Locking Screw, stardrive, self tapping	19 mm	Aleación de titanio
030300020	1,5 mm Locking Screw, stardrive, self tapping	20 mm	Aleación de titanio
030320006	2,0mm Locking Screw, stardrive, self tapping	6 mm	Aleación de titanio
030320007	2,0mm Locking Screw, stardrive, self tapping	7 mm	Aleación de titanio
030320008	2,0mm Locking Screw, stardrive, self tapping	8 mm	Aleación de titanio
030320009	2,0mm Locking Screw, stardrive, self tapping	9 mm	Aleación de titanio
030320010	2,0mm Locking Screw, stardrive, self tapping	10 mm	Aleación de titanio
030320011	2,0mm Locking Screw, stardrive, self tapping	11 mm	Aleación de titanio
030320012	2,0mm Locking Screw, stardrive, self tapping	12 mm	Aleación de titanio
030320013	2,0mm Locking Screw, stardrive, self tapping	13 mm	Aleación de titanio
030320014	2,0mm Locking Screw, stardrive, self tapping	14 mm	Aleación de titanio
030320015	2,0mm Locking Screw, stardrive, self tapping	15 mm	Aleación de titanio
030320016	2,0mm Locking Screw, stardrive, self tapping	16 mm	Aleación de titanio
030320017	2,0mm Locking Screw, stardrive, self tapping	17 mm	Aleación de titanio
030320018	2,0mm Locking Screw, stardrive, self tapping	18 mm	Aleación de titanio
030320019	2,0mm Locking Screw, stardrive, self tapping	19 mm	Aleación de titanio
030320020	2,0mm Locking Screw, stardrive, self tapping	20 mm	Aleación de titanio
030320022	2,0mm Locking Screw, stardrive, self tapping	22 mm	Aleación de titanio
030320024	2,0mm Locking Screw, stardrive, self tapping	24 mm	Aleación de titanio
030320026	2,0mm Locking Screw, stardrive, self tapping	26 mm	Aleación de titanio
030320028	2,0mm Locking Screw, stardrive, self tapping	28 mm	Aleación de titanio
030320030	2,0mm Locking Screw, stardrive, self tapping	30 mm	Aleación de titanio
030340006	2,7 mm Locking Screw, stardrive, self-tapping	6 mm	Aleación de titanio
030340007	2,7 mm Locking Screw, stardrive, self-tapping	7 mm	Aleación de titanio
030340008	2,7 mm Locking Screw, stardrive, self-tapping	8 mm	Aleación de titanio
030340009	2,7 mm Locking Screw, stardrive, self-tapping	9 mm	Aleación de titanio
030340010	2,7 mm Locking Screw, stardrive, self-tapping	10 mm	Aleación de titanio
030340011	2,7 mm Locking Screw, stardrive, self-tapping	11 mm	Aleación de titanio
030340012	2,7 mm Locking Screw, stardrive, self-tapping	12 mm	Aleación de titanio
030340013	2,7 mm Locking Screw, stardrive, self-tapping	13 mm	Aleación de titanio
030340014	2,7 mm Locking Screw, stardrive, self-tapping	14 mm	Aleación de titanio
030340015	2,7 mm Locking Screw, stardrive, self-tapping	15 mm	Aleación de titanio
030340016	2,7 mm Locking Screw, stardrive, self-tapping	16 mm	Aleación de titanio
030340017	2,7 mm Locking Screw, stardrive, self-tapping	17 mm	Aleación de titanio
030340018	2,7 mm Locking Screw, stardrive, self-tapping	18 mm	Aleación de titanio
030340019	2,7 mm Locking Screw, stardrive, self-tapping	19 mm	Aleación de titanio
030340020	2,7 mm Locking Screw, stardrive, self-tapping	20 mm	Aleación de titanio
030340022	2,7 mm Locking Screw, stardrive, self-tapping	22 mm	Aleación de titanio
030340024	2,7 mm Locking Screw, stardrive, self-tapping	24 mm	Aleación de titanio
030340026	2,7 mm Locking Screw, stardrive, self-tapping	26 mm	Aleación de titanio
030340028	2,7 mm Locking Screw, stardrive, self-tapping	28 mm	Aleación de titanio
030340030	2,7 mm Locking Screw, stardrive, self-tapping	30 mm	Aleación de titanio
030340035	2,7mm Locking Screw, Stardrive, Self-tapping	35mm	Aleación de titanio
030340040	2,7mm Locking Screw, Stardrive, Self-tapping	40mm	Aleación de titanio
030340045	2,7mm Locking Screw, Stardrive, Self-tapping	45mm	Aleación de titanio
030340050	2,7mm Locking Screw, Stardrive, Self-tapping	50mm	Aleación de titanio
031271205	MD Anterior Process Calcaneal Locking Plate, small, left	5 Holes	Aleación de titanio
031271705	MD Anterior Process Calcaneal Locking Plate, large, left	5 Holes	Aleación de titanio

031272205	MD Anterior Process Calcaneal Locking Plate, small, right	5 Holes	Aleación de titanio
031272705	MD Anterior Process Calcaneal Locking Plate, large, right	5 Holes	Aleación de titanio
031281005	MD Posterior Tuberosity Calcaneal Locking Plate, small, left	5 Holes	Aleación de titanio
031281006	MD Posterior Tuberosity Calcaneal Locking Plate, large, left	6 Holes	Aleación de titanio
031282005	MD Posterior Tuberosity Calcaneal Locking Plate, small, right	5 Holes	Aleación de titanio
031282006	MD Posterior Tuberosity Calcaneal Locking Plate, large, right	6 Holes	Aleación de titanio
031291008	MD AP and PT Calcaneal Locking Plate, medium, left	8 Holes	Aleación de titanio
031291009	MD AP and PT Calcaneal Locking Plate, large, left	9 Holes	Aleación de titanio
031291308	MD AP and PT Calcaneal Locking Plate, small, left	8 Holes	Aleación de titanio
031292008	MD AP and PT Calcaneal Locking Plate, medium, right	8 Holes	Aleación de titanio
031292009	MD AP and PT Calcaneal Locking Plate, large, right	9 Holes	Aleación de titanio
031292308	MD AP and PT Calcaneal Locking Plate, small, right	8 Holes	Aleación de titanio
031301011	MD Calcaneal Locking Plate, small, left	11Holes	Aleación de titanio
031302011	MD Calcaneal Locking Plate, small, right	11Holes	Aleación de titanio
031311011	MD Calcaneal Locking Plate, medium, left	11Holes	Aleación de titanio
031312011	MD Calcaneal Locking Plate, medium, right	11Holes	Aleación de titanio
031321011	MD Calcaneal Locking Plate, large, left	11Holes	Aleación de titanio
031322011	MD Calcaneal Locking Plate, large, right	11Holes	Aleación de titanio
031360002	MT T-Locking Plate for First Metatarsa	2 Holes	Aleación de titanio
031851006	MD Fifth Metatarsal Locking Plate		Aleación de titanio
031852006	MD Fifth Metatarsal Locking Plate		Aleación de titanio
031861009	MT Cuboid Locking Plate, left		Aleación de titanio
031862009	MT Cuboid Locking Plate, right		Aleación de titanio
031870005	MT Talus Locking Plate		Aleación de titanio
031941005	MT First MTC Plantar Fusion Locking Plate, large, left		Aleación de titanio
031942005	MT First MTC Plantar Fusion Locking Plate, large, right		Aleación de titanio
031951005	MT First MTC Plantar Fusion Locking Plate, small, left		Aleación de titanio
031952005	MT First MTC Plantar Fusion Locking Plate, small, right		Aleación de titanio
031961004	Distal Lateral Fibular Locking Plate, left	4 holes	Aleación de titanio
031961005	Distal Lateral Fibular Locking Plate, left	5 holes	Aleación de titanio
031961006	Distal Lateral Fibular Locking Plate, left	6 holes	Aleación de titanio
031961007	Distal Lateral Fibular Locking Plate, left	7 holes	Aleación de titanio
031961009	Distal Lateral Fibular Locking Plate, left	9 holes	Aleación de titanio
031961011	Distal Lateral Fibular Locking Plate, left	11 holes	Aleación de titanio
031961013	Distal Lateral Fibular Locking Plate, left	13 holes	Aleación de titanio
031962004	Distal Lateral Fibular Locking Plate, right	4 holes	Aleación de titanio
031962005	Distal Lateral Fibular Locking Plate, right	5 holes	Aleación de titanio
031962006	Distal Lateral Fibular Locking Plate, right	6 holes	Aleación de titanio
031962007	Distal Lateral Fibular Locking Plate, right	7 holes	Aleación de titanio
031962009	Distal Lateral Fibular Locking Plate, right	9 holes	Aleación de titanio
031962011	Distal Lateral Fibular Locking Plate, right	11 holes	Aleación de titanio
031962013	Distal Lateral Fibular Locking Plate, right	13 holes	Aleación de titanio
031970003	MD Distal Tibial Locking Plate	3 holes	Aleación de titanio
031980004	MD Distal Medial Tibial Locking Plate	4 holes	Aleación de titanio
031990004	MD Distal Anterior Tibial Locking Plate	4 holes	Aleación de titanio
031990006	MD Distal Anterior Tibial Locking Plate	6 holes	Aleación de titanio
032001009	MT Navicular Locking Plate, left		Aleación de titanio
032002009	MT Navicular Locking Plate, right		Aleación de titanio
040070010	3.5mm Locking Screw, Self-tapping	10mm	Aleación de titanio
040070012	3.5mm Locking Screw, Self-tapping	12mm	Aleación de titanio

040070014	3.5mm Locking Screw, Self-tapping	14mm	Aleación de titanio
040070016	3.5mm Locking Screw, Self-tapping	16mm	Aleación de titanio
040070018	3.5mm Locking Screw, Self-tapping	18mm	Aleación de titanio
040070020	3.5mm Locking Screw, Self-tapping	20mm	Aleación de titanio
040070022	3.5mm Locking Screw, Self-tapping	22mm	Aleación de titanio
040070024	3.5mm Locking Screw, Self-tapping	24mm	Aleación de titanio
040070026	3.5mm Locking Screw, Self-tapping	26mm	Aleación de titanio
040070028	3.5mm Locking Screw, Self-tapping	28mm	Aleación de titanio
040070030	3.5mm Locking Screw, Self-tapping	30mm	Aleación de titanio
040070032	3.5mm Locking Screw, Self-tapping	32mm	Aleación de titanio
040070034	3.5mm Locking Screw, Self-tapping	34mm	Aleación de titanio
040070035	3.5mm Locking Screw, Self-tapping	35mm	Aleación de titanio
040070036	3.5mm Locking Screw, Self-tapping	36mm	Aleación de titanio
040070038	3.5mm Locking Screw, Self-tapping	38mm	Aleación de titanio
040070040	3.5mm Locking Screw, Self-tapping	40mm	Aleación de titanio
040070042	3.5mm Locking Screw, Self-tapping	42mm	Aleación de titanio
040070044	3.5mm Locking Screw, Self-tapping	44mm	Aleación de titanio
040070045	3.5mm Locking Screw, Self-tapping	45mm	Aleación de titanio
040070046	3.5mm Locking Screw, Self-tapping	46mm	Aleación de titanio
040070048	3.5mm Locking Screw, Self-tapping	48mm	Aleación de titanio
040070050	3.5mm Locking Screw, Self-tapping	50mm	Aleación de titanio
040070052	3.5mm Locking Screw, Self-tapping	52mm	Aleación de titanio
040070054	3.5mm Locking Screw, Self-tapping	54mm	Aleación de titanio
040070055	3.5mm Locking Screw, Self-tapping	55mm	Aleación de titanio
040070056	3.5mm Locking Screw, Self-tapping	56mm	Aleación de titanio
040070058	3.5mm Locking Screw, Self-tapping	58mm	Aleación de titanio
040070060	3.5mm Locking Screw, Self-tapping	60mm	Aleación de titanio
040070065	3.5mm Locking Screw, Self-tapping	65mm	Aleación de titanio
040070070	3.5mm Locking Screw, Self-tapping	70mm	Aleación de titanio
040070075	3.5mm Locking Screw, Self-tapping	75mm	Aleación de titanio
040080012	3.5mm Locking Screw, Self-drilling	12mm	Aleación de titanio
040080014	3.5mm Locking Screw, Self-drilling	14mm	Aleación de titanio
040080016	3.5mm Locking Screw, Self-drilling	16mm	Aleación de titanio
040080018	3.5mm Locking Screw, Self-drilling	18mm	Aleación de titanio
040080020	3.5mm Locking Screw, Self-drilling	20mm	Aleación de titanio
040080022	3.5mm Locking Screw, Self-drilling	22mm	Aleación de titanio
040080024	3.5mm Locking Screw, Self-drilling	24mm	Aleación de titanio
040080026	3.5mm Locking Screw, Self-drilling	26mm	Aleación de titanio
040080028	3.5mm Locking Screw, Self-drilling	28mm	Aleación de titanio
040080030	3.5mm Locking Screw, Self-drilling	30mm	Aleación de titanio
040080032	3.5mm Locking Screw, Self-drilling	32mm	Aleación de titanio
040080034	3.5mm Locking Screw, Self-drilling	34mm	Aleación de titanio
040080035	3.5mm Locking Screw, Self-drilling	35mm	Aleación de titanio
040080036	3.5mm Locking Screw, Self-drilling	36mm	Aleación de titanio
040080038	3.5mm Locking Screw, Self-drilling	38mm	Aleación de titanio
040080040	3.5mm Locking Screw, Self-drilling	40mm	Aleación de titanio
040080042	3.5mm Locking Screw, Self-drilling	42mm	Aleación de titanio
040080044	3.5mm Locking Screw, Self-drilling	44mm	Aleación de titanio
040080045	3.5mm Locking Screw, Self-drilling	45mm	Aleación de titanio
040080046	3.5mm Locking Screw, Self-drilling	46mm	Aleación de titanio

040080048	3.5mm Locking Screw, Self-drilling	48mm	Aleación de titanio
040080050	3.5mm Locking Screw, Self-drilling	50mm	Aleación de titanio
040080052	3.5mm Locking Screw, Self-drilling	52mm	Aleación de titanio
040080054	3.5mm Locking Screw, Self-drilling	54mm	Aleación de titanio
040080055	3.5mm Locking Screw, Self-drilling	55mm	Aleación de titanio
040080056	3.5mm Locking Screw, Self-drilling	56mm	Aleación de titanio
040080058	3.5mm Locking Screw, Self-drilling	58mm	Aleación de titanio
040080060	3.5mm Locking Screw, Self-drilling	60mm	Aleación de titanio
040090016	5.0mm Locking Screw, Self-tapping	16mm	Aleación de titanio
040090018	5.0mm Locking Screw, Self-tapping	18mm	Aleación de titanio
040090020	5.0mm Locking Screw, Self-tapping	20mm	Aleación de titanio
040090022	5.0mm Locking Screw, Self-tapping	22mm	Aleación de titanio
040090024	5.0mm Locking Screw, Self-tapping	24mm	Aleación de titanio
040090026	5.0mm Locking Screw, Self-tapping	26mm	Aleación de titanio
040090028	5.0mm Locking Screw, Self-tapping	28mm	Aleación de titanio
040090030	5.0mm Locking Screw, Self-tapping	30mm	Aleación de titanio
040090032	5.0mm Locking Screw, Self-tapping	32mm	Aleación de titanio
040090034	5.0mm Locking Screw, Self-tapping	34mm	Aleación de titanio
040090036	5.0mm Locking Screw, Self-tapping	36mm	Aleación de titanio
040090038	5.0mm Locking Screw, Self-tapping	38mm	Aleación de titanio
040090040	5.0mm Locking Screw, Self-tapping	40mm	Aleación de titanio
040090042	5.0mm Locking Screw, Self-tapping	42mm	Aleación de titanio
040090044	5.0mm Locking Screw, Self-tapping	44mm	Aleación de titanio
040090045	5.0mm Locking Screw, Self-tapping	45mm	Aleación de titanio
040090046	5.0mm Locking Screw, Self-tapping	46mm	Aleación de titanio
040090048	5.0mm Locking Screw, Self-tapping	48mm	Aleación de titanio
040090050	5.0mm Locking Screw, Self-tapping	50mm	Aleación de titanio
040090055	5.0mm Locking Screw, Self-tapping	55mm	Aleación de titanio
040090060	5.0mm Locking Screw, Self-tapping	60mm	Aleación de titanio
040090065	5.0mm Locking Screw, Self-tapping	65mm	Aleación de titanio
040090070	5.0mm Locking Screw, Self-tapping	70mm	Aleación de titanio
040090075	5.0mm Locking Screw, Self-tapping	75mm	Aleación de titanio
040090080	5.0mm Locking Screw, Self-tapping	80mm	Aleación de titanio
040090085	5.0mm Locking Screw, Self-tapping	85mm	Aleación de titanio
040090090	5.0mm Locking Screw, Self-tapping	90mm	Aleación de titanio
040100020	5.0mm Locking Screw, Self-drilling	20mm	Aleación de titanio
040100022	5.0mm Locking Screw, Self-drilling	22mm	Aleación de titanio
040100024	5.0mm Locking Screw, Self-drilling	24mm	Aleación de titanio
040100026	5.0mm Locking Screw, Self-drilling	26mm	Aleación de titanio
040100028	5.0mm Locking Screw, Self-drilling	28mm	Aleación de titanio
040100030	5.0mm Locking Screw, Self-drilling	30mm	Aleación de titanio
040100032	5.0mm Locking Screw, Self-drilling	32mm	Aleación de titanio
040100034	5.0mm Locking Screw, Self-drilling	34mm	Aleación de titanio
040100036	5.0mm Locking Screw, Self-drilling	36mm	Aleación de titanio
040100038	5.0mm Locking Screw, Self-drilling	38mm	Aleación de titanio
040100040	5.0mm Locking Screw, Self-drilling	40mm	Aleación de titanio
040100042	5.0mm Locking Screw, Self-drilling	42mm	Aleación de titanio
040100044	5.0mm Locking Screw, Self-drilling	44mm	Aleación de titanio
040100046	5.0mm Locking Screw, Self-drilling	46mm	Aleación de titanio
040100048	5.0mm Locking Screw, Self-drilling	48mm	Aleación de titanio



BIOPROT de bioimplant S.R.L.
IMPLANTES ORTOPÉDICOS

PLACAS BLOQUEADAS DE TITANIO

**INSTRUCCIONES DE
 USO
 ANEXO IIIB**

REV.0
 FECHA: 01/2019

040100050	5.0mm Locking Screw, Self-drilling	50mm	Aleación de titanio
040100055	5.0mm Locking Screw, Self-drilling	55mm	Aleación de titanio
040100060	5.0mm Locking Screw, Self-drilling	60mm	Aleación de titanio
040100065	5.0mm Locking Screw, Self-drilling	65mm	Aleación de titanio
040100070	5.0mm Locking Screw, Self-drilling	70mm	Aleación de titanio
040100075	5.0mm Locking Screw, Self-drilling	75mm	Aleación de titanio
040100080	5.0mm Locking Screw, Self-drilling	80mm	Aleación de titanio
040100085	5.0mm Locking Screw, Self-drilling	85mm	Aleación de titanio
040100090	5.0mm Locking Screw, Self-drilling	90mm	Aleación de titanio
040130060	7.3mm Locking Screw	60mm	Aleación de titanio
040130065	7.3mm Locking Screw	65mm	Aleación de titanio
040130070	7.3mm Locking Screw	70mm	Aleación de titanio
040130075	7.3mm Locking Screw	75mm	Aleación de titanio
040130076	7.3mm Locking Screw, 76.5mm	76.5mm	Aleación de titanio
040130080	7.3mm Locking Screw	80mm	Aleación de titanio
040130083	7.3mm Locking Screw, 83.5mm	83.5mm	Aleación de titanio
040130085	7.3mm Locking Screw	85mm	Aleación de titanio
040130090	7.3mm Locking Screw	90mm	Aleación de titanio
040130095	7.3mm Locking Screw	95mm	Aleación de titanio
040130100	7.3mm Locking Screw	100mm	Aleación de titanio
040130105	7.3mm Locking Screw	105mm	Aleación de titanio
040130110	7.3mm Locking Screw	110mm	Aleación de titanio
040130115	7.3mm Locking Screw	115mm	Aleación de titanio
040140060	6.5mm Cannulated Locking Screw, Fully Threaded	60mm	Aleación de titanio
040140065	6.5mm Cannulated Locking Screw, Fully Threaded	65mm	Aleación de titanio
040140070	6.5mm Cannulated Locking Screw, Fully Threaded	70mm	Aleación de titanio
040140075	6.5mm Cannulated Locking Screw, Fully Threaded	75mm	Aleación de titanio
040140080	6.5mm Cannulated Locking Screw, Fully Threaded	80mm	Aleación de titanio
040140085	6.5mm Cannulated Locking Screw, Fully Threaded	85mm	Aleación de titanio
040140090	6.5mm Cannulated Locking Screw, Fully Threaded	90mm	Aleación de titanio
040140095	6.5mm Cannulated Locking Screw, Fully Threaded	95mm	Aleación de titanio
040140100	6.5mm Cannulated Locking Screw, Fully Threaded	100mm	Aleación de titanio
040140105	6.5mm Cannulated Locking Screw, Fully Threaded	105mm	Aleación de titanio
040140110	6.5mm Cannulated Locking Screw, Fully Threaded	110mm	Aleación de titanio
040140115	6.5mm Cannulated Locking Screw, Fully Threaded	115mm	Aleación de titanio
040140120	6.5mm Cannulated Locking Screw, Fully Threaded	120mm	Aleación de titanio
040240006	2.4mm Locking Screw, Self-tapping, Stardrive	6mm	Aleación de titanio
040240007	2.4mm Locking Screw, Self-tapping, Stardrive	7mm	Aleación de titanio
040240008	2.4mm Locking Screw, Self-tapping, Stardrive	8mm	Aleación de titanio
040240009	2.4mm Locking Screw, Self-tapping, Stardrive	9mm	Aleación de titanio
040240010	2.4mm Locking Screw, Self-tapping, Stardrive	10mm	Aleación de titanio
040240011	2.4mm Locking Screw, Self-tapping, Stardrive	11mm	Aleación de titanio
040240012	2.4mm Locking Screw, Self-tapping, Stardrive	12mm	Aleación de titanio
040240013	2.4mm Locking Screw, Self-tapping, Stardrive	13mm	Aleación de titanio
040240014	2.4mm Locking Screw, Self-tapping, Stardrive	14mm	Aleación de titanio
040240015	2.4mm Locking Screw, Self-tapping, Stardrive	15mm	Aleación de titanio
040240016	2.4mm Locking Screw, Self-tapping, Stardrive	16mm	Aleación de titanio
040240017	2.4mm Locking Screw, Self-tapping, Stardrive	17mm	Aleación de titanio
040240018	2.4mm Locking Screw, Self-tapping, Stardrive	18mm	Aleación de titanio
040240019	2.4mm Locking Screw, Self-tapping, Stardrive	19mm	Aleación de titanio

BIOIMPLANT S.R.L.
 Román Turoldo
 Social Gerente

Ma. VICTORIA COZZI
 FARMACEUTICA
 Mat. 2194

040240020	2.4mm Locking Screw, Self-tapping, Stardrive	20mm	Aleación de titanio
040240022	2.4mm Locking Screw, Self-tapping, Stardrive	22mm	Aleación de titanio
040240024	2.4mm Locking Screw, Self-tapping, Stardrive	24mm	Aleación de titanio
040240026	2.4mm Locking Screw, Self-tapping, Stardrive	26mm	Aleación de titanio
040240028	2.4mm Locking Screw, Self-tapping, Stardrive	28mm	Aleación de titanio
040240030	2.4mm Locking Screw, Self-tapping, Stardrive	30mm	Aleación de titanio
040260070	7.3mm Cannulated Locking Screw, Thread Length 16mm	70mm	Aleación de titanio
040260075	7.3mm Cannulated Locking Screw, Thread Length 16mm	75mm	Aleación de titanio
040260080	7.3mm Cannulated Locking Screw, Thread Length 16mm	80mm	Aleación de titanio
040260085	7.3mm Cannulated Locking Screw, Thread Length 16mm	85mm	Aleación de titanio
040260090	7.3mm Cannulated Locking Screw, Thread Length 16mm	90mm	Aleación de titanio
040260095	7.3mm Cannulated Locking Screw, Thread Length 16mm	95mm	Aleación de titanio
040260100	7.3mm Cannulated Locking Screw, Thread Length 16mm	100mm	Aleación de titanio
040260105	7.3mm Cannulated Locking Screw, Thread Length 16mm	105mm	Aleación de titanio
040260110	7.3mm Cannulated Locking Screw, Thread Length 16mm	110mm	Aleación de titanio
040270070	7.3mm Cannulated Locking Screw, Thread Length 32mm	70mm	Aleación de titanio
040270075	7.3mm Cannulated Locking Screw, Thread Length 32mm	75mm	Aleación de titanio
040270080	7.3mm Cannulated Locking Screw, Thread Length 32mm	80mm	Aleación de titanio
040270085	7.3mm Cannulated Locking Screw, Thread Length 32mm	85mm	Aleación de titanio
040270090	7.3mm Cannulated Locking Screw, Thread Length 32mm	90mm	Aleación de titanio
040270095	7.3mm Cannulated Locking Screw, Thread Length 32mm	95mm	Aleación de titanio
040270100	7.3mm Cannulated Locking Screw, Thread Length 32mm	100mm	Aleación de titanio
040270105	7.3mm Cannulated Locking Screw, Thread Length 32mm	105mm	Aleación de titanio
040270110	7.3mm Cannulated Locking Screw, Thread Length 32mm	110mm	Aleación de titanio
040270115	7.3mm Cannulated Locking Screw, Thread Length 32mm	115mm	Aleación de titanio
040300010	5.0mm Locking Screw, Self-tapping, Flat Head	10mm	Aleación de titanio
040300014	5.0mm Locking Screw, Self-tapping, Flat Head	14mm	Aleación de titanio
040300018	5.0mm Locking Screw, Self-tapping, Flat Head	18mm	Aleación de titanio
040760008	2.4mm Locking Screw II , variable angle, stardrive, self-tapping	8mm	Aleación de titanio
040760010	2.4mm Locking Screw II , variable angle, stardrive, self-tapping	10mm	Aleación de titanio
040760012	2.4mm Locking Screw II , variable angle, stardrive, self-tapping	12mm	Aleación de titanio
040760014	2.4mm Locking Screw II , variable angle, stardrive, self-tapping	14mm	Aleación de titanio
040760016	2.4mm Locking Screw II , variable angle, stardrive, self-tapping	16mm	Aleación de titanio
040760018	2.4mm Locking Screw II , variable angle, stardrive, self-tapping	18mm	Aleación de titanio
040760020	2.4mm Locking Screw II , variable angle, stardrive, self-tapping	20mm	Aleación de titanio
040760022	2.4mm Locking Screw II , variable angle, stardrive, self-tapping	22mm	Aleación de titanio
040760024	2.4mm Locking Screw II , variable angle, stardrive, self-tapping	24mm	Aleación de titanio
040760026	2.4mm Locking Screw II , variable angle, stardrive, self-tapping	26mm	Aleación de titanio
040760028	2.4mm Locking Screw II , variable angle, stardrive, self-tapping	28mm	Aleación de titanio
040760030	2.4mm Locking Screw II , variable angle, stardrive, self-tapping	30mm	Aleación de titanio
040910060	5.5mm Locking Screw, Self-tapping	60mm	Aleación de titanio
040910065	5.5mm Locking Screw, Self-tapping	65mm	Aleación de titanio
040910070	5.5mm Locking Screw, Self-tapping	70mm	Aleación de titanio
040910075	5.5mm Locking Screw, Self-tapping	75mm	Aleación de titanio
040910080	5.5mm Locking Screw, Self-tapping	80mm	Aleación de titanio
040910085	5.5mm Locking Screw, Self-tapping	85mm	Aleación de titanio
040910090	5.5mm Locking Screw, Self-tapping	90mm	Aleación de titanio
040910095	5.5mm Locking Screw, Self-tapping	95mm	Aleación de titanio
040910100	5.5mm Locking Screw, Self-tapping	100mm	Aleación de titanio
040910105	5.5mm Locking Screw, Self-tapping	105mm	Aleación de titanio

040910110	5.5mm Locking Screw, Self-tapping	110mm	Aleación de titanio
040920060	6.5mm Cannulated Locking Screw, Self-tapping	60mm	Aleación de titanio
040920065	6.5mm Cannulated Locking Screw, Self-tapping	65mm	Aleación de titanio
040920070	6.5mm Cannulated Locking Screw, Self-tapping	70mm	Aleación de titanio
040920075	6.5mm Cannulated Locking Screw, Self-tapping	75mm	Aleación de titanio
040920080	6.5mm Cannulated Locking Screw, Self-tapping	80mm	Aleación de titanio
040920085	6.5mm Cannulated Locking Screw, Self-tapping	85mm	Aleación de titanio
040920090	6.5mm Cannulated Locking Screw, Self-tapping	90mm	Aleación de titanio
040920095	6.5mm Cannulated Locking Screw, Self-tapping	95mm	Aleación de titanio
040920100	6.5mm Cannulated Locking Screw, Self-tapping	100mm	Aleación de titanio
040920105	6.5mm Cannulated Locking Screw, Self-tapping	105mm	Aleación de titanio
040920110	6.5mm Cannulated Locking Screw, Self-tapping	110mm	Aleación de titanio
040940014	4,0 mm Locking Screw, self-tapping	14 mm	Aleación de titanio
040940016	4,0 mm Locking Screw, self-tapping	16 mm	Aleación de titanio
040940018	4,0 mm Locking Screw, self-tapping	18 mm	Aleación de titanio
040940020	4,0 mm Locking Screw, self-tapping	20 mm	Aleación de titanio
040940022	4,0 mm Locking Screw, self-tapping	22 mm	Aleación de titanio
040940024	4,0 mm Locking Screw, self-tapping	24 mm	Aleación de titanio
040940026	4,0 mm Locking Screw, self-tapping	26 mm	Aleación de titanio
040940028	4,0 mm Locking Screw, self-tapping	28 mm	Aleación de titanio
040940030	4,0 mm Locking Screw, self-tapping	30 mm	Aleación de titanio
040940032	4,0 mm Locking Screw, self-tapping	32 mm	Aleación de titanio
040940034	4,0 mm Locking Screw, self-tapping	34 mm	Aleación de titanio
040940036	4,0 mm Locking Screw, self-tapping	36 mm	Aleación de titanio
040940038	4,0 mm Locking Screw, self-tapping	38 mm	Aleación de titanio
040940040	4,0 mm Locking Screw, self-tapping	40 mm	Aleación de titanio
040940045	4,0 mm Locking Screw, self-tapping	45 mm	Aleación de titanio
040940050	4,0 mm Locking Screw, self-tapping	50 mm	Aleación de titanio
040940055	4,0 mm Locking Screw, self-tapping	55 mm	Aleación de titanio
040940060	4,0 mm Locking Screw, self-tapping	60 mm	Aleación de titanio
040940065	4,0 mm Locking Screw, self-tapping	65 mm	Aleación de titanio
040940070	4,0 mm Locking Screw, self-tapping	70 mm	Aleación de titanio
040940075	4,0 mm Locking Screw, self-tapping	75 mm	Aleación de titanio
040940080	4,0 mm Locking Screw, self-tapping	80 mm	Aleación de titanio
041050014	1,8 mm Locking Buttress Screw, variable angle, stardrive	14 mm	Aleación de titanio
041050016	1,8 mm Locking Buttress Screw, variable angle, stardrive	16 mm	Aleación de titanio
041050018	1,8 mm Locking Buttress Screw, variable angle, stardrive	18 mm	Aleación de titanio
041050020	1,8 mm Locking Buttress Screw, variable angle, stardrive	20 mm	Aleación de titanio
041050022	1,8 mm Locking Buttress Screw, variable angle, stardrive	22 mm	Aleación de titanio
041050024	1,8 mm Locking Buttress Screw, variable angle, stardrive	24 mm	Aleación de titanio
041050026	1,8 mm Locking Buttress Screw, variable angle, stardrive	26 mm	Aleación de titanio
041170010	MT 2.7mm Locking Screw , stardrive, self-tapping	10mm	Aleación de titanio
041170012	MT 2.7mm Locking Screw , stardrive, self-tapping	12mm	Aleación de titanio
041170014	MT 2.7mm Locking Screw , stardrive, self-tapping	14mm	Aleación de titanio
041170016	MT 2.7mm Locking Screw , stardrive, self-tapping	16mm	Aleación de titanio
041170018	MT 2.7mm Locking Screw , stardrive, self-tapping	18mm	Aleación de titanio
041170020	MT 2.7mm Locking Screw , stardrive, self-tapping	20mm	Aleación de titanio
041170022	MT 2.7mm Locking Screw , stardrive, self-tapping	22mm	Aleación de titanio
041170024	MT 2.7mm Locking Screw , stardrive, self-tapping	24mm	Aleación de titanio
041170026	MT 2.7mm Locking Screw , stardrive, self-tapping	26mm	Aleación de titanio

041170028	MT 2.7mm Locking Screw , stardrive, self-tapping	28mm	Aleación de titanio
041170030	MT 2.7mm Locking Screw , stardrive, self-tapping	30mm	Aleación de titanio
041170032	MT 2.7mm Locking Screw , stardrive, self-tapping	32mm	Aleación de titanio
041170034	MT 2.7mm Locking Screw , stardrive, self-tapping	34mm	Aleación de titanio
041170036	MT 2.7mm Locking Screw , stardrive, self-tapping	36mm	Aleación de titanio
041170038	MT 2.7mm Locking Screw , stardrive, self-tapping	38mm	Aleación de titanio
041170040	MT 2.7mm Locking Screw , stardrive, self-tapping	40mm	Aleación de titanio
041170042	MT 2.7mm Locking Screw , stardrive, self-tapping	42mm	Aleación de titanio
041170044	MT 2.7mm Locking Screw , stardrive, self-tapping	44mm	Aleación de titanio
041170046	MT 2.7mm Locking Screw , stardrive, self-tapping	46mm	Aleación de titanio
041170048	MT 2.7mm Locking Screw , stardrive, self-tapping	48mm	Aleación de titanio
041170050	MT 2.7mm Locking Screw , stardrive, self-tapping	50mm	Aleación de titanio
041180010	MT 3.0mm Locking Screw, stardrive, self-tapping	10mm	Aleación de titanio
041180012	MT 3.0mm Locking Screw, stardrive, self-tapping	12mm	Aleación de titanio
041180014	MT 3.0mm Locking Screw, stardrive, self-tapping	14mm	Aleación de titanio
041180016	MT 3.0mm Locking Screw, stardrive, self-tapping	16mm	Aleación de titanio
041180018	MT 3.0mm Locking Screw, stardrive, self-tapping	18mm	Aleación de titanio
041180020	MT 3.0mm Locking Screw, stardrive, self-tapping	20mm	Aleación de titanio
041180022	MT 3.0mm Locking Screw, stardrive, self-tapping	22mm	Aleación de titanio
041180024	MT 3.0mm Locking Screw, stardrive, self-tapping	24mm	Aleación de titanio
041180026	MT 3.0mm Locking Screw, stardrive, self-tapping	26mm	Aleación de titanio
041180028	MT 3.0mm Locking Screw, stardrive, self-tapping	28mm	Aleación de titanio
041180030	MT 3.0mm Locking Screw, stardrive, self-tapping	30mm	Aleación de titanio
041180032	MT 3.0mm Locking Screw, stardrive, self-tapping	32mm	Aleación de titanio
041180034	MT 3.0mm Locking Screw, stardrive, self-tapping	34mm	Aleación de titanio
041180036	MT 3.0mm Locking Screw, stardrive, self-tapping	36mm	Aleación de titanio
041180038	MT 3.0mm Locking Screw, stardrive, self-tapping	38mm	Aleación de titanio
041180040	MT 3.0mm Locking Screw, stardrive, self-tapping	40mm	Aleación de titanio
041180042	MT 3.0mm Locking Screw, stardrive, self-tapping	42mm	Aleación de titanio
041180044	MT 3.0mm Locking Screw, stardrive, self-tapping	44mm	Aleación de titanio
041180046	MT 3.0mm Locking Screw, stardrive, self-tapping	46mm	Aleación de titanio
041180048	MT 3.0mm Locking Screw, stardrive, self-tapping	48mm	Aleación de titanio
041180050	MT 3.0mm Locking Screw, stardrive, self-tapping	50mm	Aleación de titanio
041190010	MT 3.5mm Locking Screw, stardrive, self-tapping	10mm	Aleación de titanio
041190012	MT 3.5mm Locking Screw, stardrive, self-tapping	12mm	Aleación de titanio
041190014	MT 3.5mm Locking Screw, stardrive, self-tapping	14mm	Aleación de titanio
041190016	MT 3.5mm Locking Screw, stardrive, self-tapping	16mm	Aleación de titanio
041190018	MT 3.5mm Locking Screw, stardrive, self-tapping	18mm	Aleación de titanio
041190020	MT 3.5mm Locking Screw, stardrive, self-tapping	20mm	Aleación de titanio
041190022	MT 3.5mm Locking Screw, stardrive, self-tapping	22mm	Aleación de titanio
041190024	MT 3.5mm Locking Screw, stardrive, self-tapping	24mm	Aleación de titanio
041190026	MT 3.5mm Locking Screw, stardrive, self-tapping	26mm	Aleación de titanio
041190028	MT 3.5mm Locking Screw, stardrive, self-tapping	28mm	Aleación de titanio
041190030	MT 3.5mm Locking Screw, stardrive, self-tapping	30mm	Aleación de titanio
041190032	MT 3.5mm Locking Screw, stardrive, self-tapping	32mm	Aleación de titanio
041190034	MT 3.5mm Locking Screw, stardrive, self-tapping	34mm	Aleación de titanio
041190036	MT 3.5mm Locking Screw, stardrive, self-tapping	36mm	Aleación de titanio
041190038	MT 3.5mm Locking Screw, stardrive, self-tapping	38mm	Aleación de titanio
041190040	MT 3.5mm Locking Screw, stardrive, self-tapping	40mm	Aleación de titanio
041190042	MT 3.5mm Locking Screw, stardrive, self-tapping	42mm	Aleación de titanio



PLACAS BLOQUEADAS DE TITANIO

INSTRUCCIONES DE
USO
ANEXO IIIB

REV.0
FECHA: 01/2019

041190044	MT 3.5mm Locking Screw, stardrive, self-tapping	44mm	Aleación de titanio
041190046	MT 3.5mm Locking Screw, stardrive, self-tapping	46mm	Aleación de titanio
041190048	MT 3.5mm Locking Screw, stardrive, self-tapping	48mm	Aleación de titanio
041190050	MT 3.5mm Locking Screw, stardrive, self-tapping	50mm	Aleación de titanio

INSTRUMENTAL ASOCIADO:

111290000	Locking System 2.4 Instrument Set
111180100	Drill Bit, ϕ 1.8, Length 110mm
111180200	Screwdriver, Stardrive, ϕ 2.4
111180300	Holding Sleeve for ϕ 2.4 Locking Screws
111180400	Drill Sleeve, ϕ 1.8
111180500	Screwdriver Shaft, Stardrive, ϕ 2.4
111180600	Handle, for torque indicating screwdriver, 0.8N.m
111180700	Depth Gauge for ϕ 2.4 Locking Screws
110040100	Guide Wire, ϕ 1.5, Length 150mm
111180900	Drill Sleeve for Guide Wire, ϕ 1.5
111181000	Double Drill Guide, ϕ 2.0/ ϕ 2.7
111181100	Drill Bit, ϕ 2.0, Length 102mm
110222600	Tap for Cortex Screws, ϕ 2.7
111181300	Depth Gauge for ϕ 2.7 Cortex Screws
111181400	Holding Sleeve for ϕ 2.7 Cortex Screws
110224000	T-Handle With Quick Coupling
111181600	Drill Bit, ϕ 2.7, Length 102mm
110223700	Bone Lever, width 6.5mm
110223900	Bone Lever, width 15.5mm
111181900	Periosteal Elevator, Round Edge
111182000	Periosteal Elevator, Straight
111182100	Reduction Forceps With Points, Small
111182200	Reduction Forceps With Points, Large
111182300	Reduction forceps, Curved
111182400	Bending Iron
111290100	Locking System 2.4 Instrument Case
111189900	Variable Angle Drill Sleeve, ϕ 2.4, for ϕ 1.8 Drill Bits

111150000	Proximal Humeral Locking Plate II Instrument Set
111070100	Aiming Arm
111070200	Insertion Handle
111070300	Connecting Screw
111070500	Drill Sleeve, for Load Position
111070600	Outer Sleeve, ϕ 8.0/ ϕ 4.9
111070700	Outer Sleeve, ϕ 8.0/ ϕ 6.0
111070800	Drill Sleeve, ϕ 4.9/ ϕ 2.9
111070900	Drill Sleeve, ϕ 6.0/ ϕ 2.6
111071000	Drill Sleeve for Guide Wire, ϕ 4.9/ ϕ 1.6
111071100	Trocar, ϕ 4.9
111071200	Trocar, ϕ 6.0
111071300	Guide Wire, ϕ 1.5, Length 200mm

BIOMPLANT S.R.L.
Romina M. Voldo
Socia Gerente

Ma. VICTORIA COZZI
FARMACEUTICA
Mat. 2194

Página 25 de 37

140



BIOPROT de bioimplant SRL
 IMPLANTES QUIRÚRGICOS

PLACAS BLOQUEADAS DE TITANIO

INSTRUCCIONES DE
 USO
 ANEXO IIIB

REV.0
 FECHA: 01/2019

111071400	Depth Gauge	
111071500	Direct Measuring Device	
111071600	Screwdriver, Hexagonal, ϕ 2.5	
111071700	Torque Indicating Screwdriver, ϕ 2.5, 1.5Nm	
111071800	Drill Bit, ϕ 2.5, Length 200mm	
111071900	Drill Bit, ϕ 2.8, Length 200mm	
111072000	Fixation Sleeve for Drill Bits, ϕ 3.0	
111072100	Tap for Cortex Screws, ϕ 3.5, Length 200mm	
111072200	Screwdriver Shaft, Hexagonal, ϕ 2.5	
111072500	Drill Sleeve, for Neutral Position	
111072600	L-Wrench, Hexagonal, ϕ 2.5	
111072700	Connecting Screw for Sleeve	
111152300	Proximal Humeral Locking Plate	II Instrument Case

110290000 Low Profile Distal Medial Tibial Locking Plate Instrument Set		
111582600	Aiming Arm, left	
111582500	Aiming Arm, right	
111580600	Outer Sleeve, ϕ 7.5/ ϕ 6.0	
111580700	Drill Sleeve, for Guide Wire, ϕ 6.0/ ϕ 2.0	
111580800	Drill Sleeve, ϕ 6.0/ ϕ 2.8	
111580900	Drill Sleeve, for Drill Bit, ϕ 6.0/ ϕ 2.5	
111581000	Trocar, ϕ 6.0	
111581100	Guide Wire, ϕ 2.0, length 250mm	
111581200	Depth Gauge	
111581300	Depth Gauge	
111581400	Screwdriver SW2.5	
111581500	Torque Indicating Screwdriver, 1.5N.m	
111581600	Drill Bit, ϕ 2.5, length 220mm	
111581700	Drill Bit, ϕ 2.8, length 220mm	
111581800	Fixation Sleeve, for Drill Bits	
111581900	Tap, ϕ 6.5	
111582000	Screwdriver Shaft SW2.5	
111582100	Wrench, for Fixation Sleeve	
111582200	Stopper, for Insertion Guide	
111582300	Pulling Device	
110355800	Wrench, for Fixation Sleeve	
110290100	Low Profile Distal Medial Tibial Locking Plate Instrument Case	

110580000 Small Locking Instrument Set		
110040100	Guide Wire, ϕ 1.5, Length 150mm	
110040200	Drill Bit ϕ 2.8, Length 200mm	
110040300	Drill Sleeve for Drill Bits, ϕ 2.8	
110040400	Depth Gauge for Locking Screws	
110022200	Drill Bit, ϕ 2.5, Length 112mm	
110011100	Bending Iron, Left	
110010700	Tap for Cortex Screws ϕ 3.5	
110040800	Torque Indicating Screwdriver, 1.5Nm	

BIOIMPLANT S.R.L.
 Román Alfredo
 Social Gerente

Ma. VICTORIA COZZI
 FARMACEUTICA
 Mat. 2194

141

110040900	Screwdriver, Hexagonal, ϕ 2.5
110041000	Holding Sleeve
110010100	Double Drill Guide, ϕ 2.5/ ϕ 3.5
110011800	Extraction Screw for Screws, ϕ 2.5
110011900	Broken Screw Extractor
110010500	T-Handle With Quick Coupling
110011200	Screw Forceps, Self-holding
110041600	Drill Sleeve, for Guide Wire ϕ 1.5, Short
110011300	Screwdriver Shaft, Hexagonal, ϕ 2.5
110012600	Bending Iron, Right
110011600	Reduction Forceps With Points
110011700	Bone Holding Forceps, Self-centering
110012200	Periosteal Elevator, Round Edge
110012300	Periosteal Elevator, Straight
110012400	Bone Lever, Small
110012500	Bone Lever, Large
110042500	Holding Sleeve for Torque Indicating Screwdriver
110580001	Small Locking Instrument Case

110550000	Small Fragment Instrument Set
110010100	Double Drill Guide, ϕ 2.5/ ϕ 3.5
110010200	Double Drill Guide, ϕ 2.5/ ϕ 3.0
110010300	DCP Drill Guide ϕ 2.5, for Neutral And Load Position
110010400	Depth Gauge
110010500	T-Handle With Quick Coupling
110010600	Countersink ϕ 6.0, Length 100mm, for Quick Coupling
110010700	Tap for Cortex Screws ϕ 3.5
110010800	Tap for Cancellous Screws ϕ 4.0
110010900	Bending Template for DCP 3.5 and LC-DCP 3.5
110011000	Bending Template for Straight Reconstruction Plates 3.5
110011100	Bending Iron, left
110011200	Screw Forceps, self-holding
110011300	Screwdriver Shaft, Hexagonal, ϕ 2.5
110011400	Holding Sleeve
110011500	Screwdriver, Hexagonal, ϕ 2.5
110011600	Reduction Forceps With Points
110011700	Bone Holding Forceps, Self-centering
110011800	Extraction Screw for Screws, ϕ 2.5
110011900	Broken Screw Extractor
110013400	Drill Bit ϕ 2.5, Length 112mm
110013500	Drill Bit ϕ 3.0, Length 112mm
110012200	Periosteal Elevator, Round Edge
110012300	Periosteal Elevator, Straight
110012400	Bone Lever, Small
110012500	Bone Lever, Large
110012600	Bending Iron, right
110550001	Small Fragment Instrument Case



BIOPROT de bioimplant S.R.L.
(IMPLANTES QUIRÚRGICOS)

PLACAS BLOQUEADAS DE TITANIO

**INSTRUCCIONES DE
 USO
 ANEXO IIIB**

REV.0
 FECHA: 01/2019

110600000	Large Fragment Instrument Set
110060100	Double Drill Guide, $\phi 3.2/\phi 4.5$
110060200	Double Drill Guide, $\phi 3.2/\phi 6.5$
110060300	DCP Drill Guide, $\phi 3.2$
110060400	Drill Guide, Pointed
110060500	Drill Sleeve Insert With Stop for Screws $\phi 4.5$
110060600	Drill Sleeve for Drill Bits, $\phi 3.2$
110060700	Depth Gauge
110060800	T-Handle With Quick Coupling
110060900	Countersink, Large, $\phi 3.0$, Length 112mm, for Cortex Screws $\phi 4.5$
110061000	Tap for Cortex Screws $\phi 4.5$
110061100	Tap for Cancellous Screws $\phi 6.5$
110061200	Bending Template, long
110061300	Bending Template, short
110061400	Bending Iron, left
110061500	Screw Forceps, Self-holding
110061600	Screwdriver Shaft, Hexagonal, Large, $\phi 3.5$
110061700	Screwdriver, hexagonal, large, $\phi 3.5$
110061800	Holding Sleeve
110061900	Reduction Forceps With Points
110062000	Bone Holding forceps, Self-centering, Large
110062100	Extraction Screw for Screws $\phi 3.5$
110062200	Broken Screw Extractor
110063700	Drill Bit $\phi 3.2$, Length 147mm
110063800	Drill Bit $\phi 4.5$, Length 147mm
110062500	Periosteal Elevator, Round Edge
110062600	Periosteal Elevator, Straight
110062700	Bone Lever, Small
110062800	Bone Lever, Large
110062900	Bending Iron, right
110600001	Large Fragment Instrument Case

110610000	Large Locking Instrument Set
110070100	Drill Bit $\phi 4.3$, Length 300mm
110070200	Drill Sleeve for Drill Bits, $\phi 4.3$
110070300	Depth Gauge
110063700	Drill Bit $\phi 3.2$, Length 147mm
110061400	Bending Iron, left
110061000	Tap for Cortex Screws $\phi 4.5$
110070700	Torque Indicating Screwdriver, 4.0N.m
110070800	Screwdriver, Hexagonal, $\phi 3.5$
110070900	Holding Sleeve
110060100	Double Drill Guide, $\phi 3.2/\phi 4.5$
110062100	Extraction Screw for Screws, $\phi 3.5$
110071200	Broken Screw Extractor

BIOIMPLANTS S.R.L.
 Romina F. Boldo
 Socia Gerente

Ma. VICTORIA COZZI
 FARMACEUTICA
 Mat. 2194



BIOPROT de bioimplant SRL
 IMPLANTES QUIRURGICOS

PLACAS BLOQUEADAS DE TITANIO

**INSTRUCCIONES DE
 USO
 ANEXO IIIB**

REV.0
 FECHA: 01/2019

110010500	T-Handle With Quick Coupling
110061500	Screw forceps, Self-holding
110071500	Drill Sleeve for Guide Wire, ϕ 2.5, short
110061600	Screwdriver Shaft, Hexagonal, ϕ 3.5
110071700	Guide Wire, ϕ 2.5, Length 250mm
110062900	Bending Iron, right
110061900	Reduction Forceps With Points
110062000	Bone Holding forceps, Self-centering
110062500	Periosteal Elevator, Round Edge
110062600	Periosteal Elevator, Straight
110062700	Bone Lever, Small
110062800	Bone Lever, Large
110072500	Holding Sleeve for Torque Indicating Screwdriver
110610001	Large Locking Instrument Case

111540000	Distal Femoral Locking Plate Instrument Set
111480100	Insertion Guide for Distal Femur, left
111480200	Centering Sleeve for Guide Wire, ϕ 2.5
111480300	Drill Sleeve for Insertion Guide
111480400	Fixation Bolt for Insertion Guide
111480500	Drill Sleeve, ϕ 4.3
111480600	Trocar
111480700	Depth Gauge
111480800	Insertion Guide for Distal Femur, right
110071700	Guide wire, ϕ 2.5, Length 250mm
110201800	Wrench for Fixation Sleeve
111481100	Fixation Sleeve, for Drill Bit, ϕ 4.3
111481200	Pulling Device
111481300	Wrench for Fixation Bolt
111481500	Stopper for Insertion Guide
111481600	Screwdriver Shaft SW3.5
111481700	Screwdriver SW3.5
110070700	Torque Indicating Screwdriver, 4.0Nm
111481900	Drill Bit ϕ 4.3, Length 300mm
111540100	Distal Femoral Locking Plate Instrument Case

111550000	Proximal Tibial Locking Plate Instrument Set
111500100	Drill Sleeve for Insertion Guide
111500200	Drill Sleeve, ϕ 4.3
111500300	Fixation Bolt for Insertion Guide
111500400	Trocar
111500500	Depth Gauge
111500600	Centering Sleeve for Guide Wire, ϕ 2.5
111500700	Insertion Guide for Proximal Tibia, right
111500800	Insertion Guide for Proximal Tibia, left
110071700	Guide wire, ϕ 2.5, Length 250mm
110201800	Wrench for Fixation Sleeve

BIOIMPLANT S.R.L.
 Romina Volpe
 Socia Gerente

Ma. VICTORIA COZZI
 FARMACEUTICA
 Mat. 2194



BIOPROT
de bioplant SRL
IMPLANTES QUIRURGICOS

PLACAS BLOQUEADAS DE TITANIO

INSTRUCCIONES DE
USO
ANEXO IIIB

REV.0
FECHA: 01/2019

111481100	Fixation Sleeve for Drill Bits
111481200	Pulling Device
111481300	Wrench for Fixation Bolt
111481500	Stopper for Insertion Guide
111501600	Screwdriver Shaft SW3.5
111501700	Screwdriver SW3.5
110070700	Torque Indicating Screwdriver, 4.0Nm
111501900	Drill Bit $\phi 4.3$, Length 300mm
111550100	Proximal Tibial Locking Plate Instrument Case

111900000	Proximal Femoral Locking Plate III Instrument Set
111640100	Insertion Guide for Proximal Femur, left
111640200	Insertion Guide for Proximal Femur, right
111640300	Trocar, $\phi 6.6$
111640400	Drill Sleeve, $\phi 6.6/\phi 4.3$
111640500	Outer Sleeve, $\phi 10.0/\phi 6.6$
111640600	Drill Sleeve for Guide Wire, $\phi 6.6/\phi 2.5$
111640700	Trocar, $\phi 8$
111640800	Drill Sleeve, $\phi 8.0/\phi 3.2$
111640900	Outer Sleeve, $\phi 10.0/\phi 8.0$
111641100	Trocar, $\phi 8.8$
111641200	Drill Sleeve, $\phi 8.8/\phi 5.0$
111641300	Fixation Bolt for Insertion Guide
111641400	Outer Sleeve, $\phi 10.0/\phi 8.8$
111641500	Drill Sleeve for Guide Wire, $\phi 8.8/\phi 2.5$
111641600	Depth Gauge
111196300	Guide Wire, $\phi 2.5$, 320mm
110110700	Guide Wire, $\phi 2.5$, length 320mm, with threaded tip
111481900	Drill Bit $\phi 4.3$, length 300mm
111642700	Drill Bit, $\phi 3.2$, length 300mm
111641700	Drill Bit $\phi 5.0$, length 320mm, cannulated
111641800	Drill Sleeve, $\phi 8.8/\phi 4.5$
111641900	Drill Bit, $\phi 4.5$, length 320mm
110111100	Fixation Sleeve for Drill Bits
111481200	Pulling Device
111481300	Wrench for locking shaft
111481600	Screwdriver Shaft, hexagonal, $\phi 3.5$
111481700	Screwdriver, hexagonal, $\phi 3.5$
110070700	Torque Indicating Screwdriver, 4.0N.m
111480700	Depth Gauge
111642000	Screwdriver, hexagonal, $\phi 4.0$
111642100	Screwdriver, hexagonal, cannulated, $\phi 4.0$
111642200	Tap for Cannulated Screws, $\phi 6.5$, cannulated
111642300	Tap for Cortex Screws, $\phi 4.5$
111642400	Compression Sleeve
110010500	T-Handle With Quick Coupling
111642600	Stopper for Insertion Guide

BIOPRO
de bioplant S.R.L.
Ramiro Turiso
Social Gerente

Ma. VICTORIA COZZI
FARMACEUTICA
Mat. 2194



BIOPROT de bioimplant S.R.L.
 IMPLANTES QUIRURGICOS

PLACAS BLOQUEADAS DE TITANIO

INSTRUCCIONES DE
 USO
 ANEXO IIIB

REV.0
 FECHA: 01/2019

110030700	Wrench, for Fixation Sleeve
111900100	Proximal Femoral Locking Plate III Instrument Case

110650000	Proximal Femoral Locking Instrument Set
110110100	Drill Sleeve for Guide Wire, ϕ 2.5
110110200	Drill Sleeve for Drill Bits, ϕ 4.3
110110300	Drill Sleeve for Drill Bits, Large
110110400	Inner Sleeve, for Large Drill Sleeve
110110500	Drill Bit ϕ 5.0, Cannulated, Length 250mm
110030400	Depth Gauge for Guide Wire
110110700	Guide Wire, ϕ 2.5, Length 320mm, With Threaded tip
110110800	Guide Wire, ϕ 2.5, Length 250mm
110030900	Screwdriver, Hexagonal, Cannulated, ϕ 4.0
110031000	Screwdriver Shaft, Hexagonal, ϕ 4.0
110111100	Fixation Sleeve for Drill Bits
110030700	Wrench for Fixation Sleeve
110111300	Tap for Cannulated Screws ϕ 7.3, Cannulated
110111400	Screwdriver, Hexagonal, ϕ 3.5
110111500	Holding Sleeve
110111600	Tap for Cancellous Bone Screws ϕ 6.5
110111700	Drill Bit, ϕ 4.5, Length 150mm
110060800	T-Handle With Quick Coupling, Cannulated
110650001	Proximal Femoral Locking Instrument Case

110570000	Proximal Femoral Buttress Plate Instrument Set
110030100	Drill Sleeve for Guide Wire
110030200	Guide Wire ϕ 2.0, Length 250mm, With Threaded Tip
110030300	Guide Wire ϕ 2.0, Length 250mm
110030400	Depth Gauge, for Guide Wire
110030500	Drill Bit ϕ 4.5, Cannulated, Length 250mm
110030600	Fixation Sleeve for Drill Bits
110030700	Wrench, for Fixation Sleeve
110030800	Tap for Cannulated Screws, ϕ 6.5, Cannulated
110030900	Screwdriver, Hexagonal, Cannulated, ϕ 4.0
110031000	Screwdriver Shaft, Hexagonal, Cannulated, ϕ 4.0
110570001	Proximal Femoral Buttress Plate Instrument Case

BIOIMPLANT S.R.L.
 Romina Trujillo
 Socia Gerente

[Signature]
 Ma. VICTORIA COZZI
 FARMACEUTICA
 Mat. 2194

- **Parámetros de esterilización (ESTERILIZAR ANTES DE USAR)**

Esterilización: Autoclave por vapor

Tipo de Ciclo: Pre-vacío pulsado

Tiempo de Esterilización: 4 minutos

Temperatura: 132°C

Tiempo de Secado: 30 minutos

Repeticiones Mínimas de Ciclos: Tres veces

- Esterilizador Recomendado: Tuttnauer EHS Series Table-Top Autoclave; (O esterilizador similar con la misma configuración)

3.2. Las prestaciones contempladas en el ítem 3 del Anexo de la Resolución GMC N° 72/98 que dispone sobre los Requisitos Esenciales de Seguridad y Eficacia de los Productos Médicos y los posibles efectos secundarios no deseados;

- **Uso indicado:**

Las Placas Bloqueada Recta se presentan no-estériles, combinadas con tornillos óseos para tratar varias fracturas, incluyendo las de los huesos grandes (la clavícula, la escápula, la pelvis, el húmero, el cúbito, el radio, el fémur, la tibia y el peroné) y de los huesos pequeños (el metacarpo, el metatarso y las falanges).

Las placas óseas son de aleación de titanio Ti6Al4V o aleación de titanio de Ti6Al7Nb

El tornillo: aleación de titanio Ti6Al4V o aleación de titanio Ti6Al7Nb

- **Contraindicaciones**

- La infección general o local, la osteomielitis.
- Los pacientes no tolerantes para la operación o la anestesia debidos al estado de salud.
- La reducción abierta, la fijación interna pueden conducir una infección local o una necrosis de la piel debido al pobre tejido blando local.
- La fijación imposible por osteoporosis severa.


- **Eventos Adversos**

En muchos casos, los adversos estarán más relacionados con la práctica clínica que con el mismo implante.

- La osteoporosis, la inhibición de revascularización, la reabsorción ósea y la formación débil ósea causaría debilidad de fijación, flexión, agrietamiento o fractura del dispositivo o la pérdida prematura de fijación con el hueso, lo que provocaría no unión del hueso.
- La unión demorada, mala unión o no unión de la fractura debido a la alineación impropia.
- Los tejidos fibrosos aumentados alrededor de la fractura debida a la fractura conminuta inestable.
- La infección temprana o tarde, profunda y/o superficial.
- El daño de nervios se puede resultar por trauma quirúrgica.
- Las reacciones de sensibilidad metal por los implantes quirúrgicos se han reportado raramente, y sus significados requieren aun más evaluaciones clínicas.
- El cargo precoz de peso impropio conduciría una fractura del implante por fatiga metal.


- **Advertencias y Precauciones**

- La responsabilidad de formación adecuada, la experiencia en la selección y la colocación de los implantes y la decisión posoperatoria de mantener o retirar los implantes queda a cargo de los cirujanos.
- El cirujano debe discutir con el paciente antes de la operación sobre las expectativas de la cirugía inherente de la utilización del producto. Atención particular debe prestarse a una discusión posoperatoria y el seguimiento médico periódico necesario.
- Se requiere elegir un implante adecuado según la fractura específica y la indicación del implante de acuerdo con los rayos X preoperatorios. Seleccione tornillos e instrumentales de acuerdo con el implante elegido. La cirugía debe terminarse bajo la guía del Manual de Fijación Interna: Técnicas Recomendadas por el Grupo AO-ASIF y otra información relevante.
- La selección correcta del implante es sumamente importante. El producto se debe utilizar en el local anatómico correcto, de acuerdo con las normas aceptadas por la fijación interna. Si no se utiliza un producto adecuado puede resultar en un fracaso clínico precoz. Si no se utiliza un implante adecuado para mantener el suministro de sangre y proporcionar una fijación rígida puede resultar debilidad de fijación, flexión, o fractura del producto y / o hueso.
- Es un requisito realizar una planificación preoperatoria cuidadosamente. Los cirujanos deben elegir los implantes adecuados según las imágenes preoperatorias y mediante el uso de instrumentos de medición, tales como medidores de profundidad y ensayos. Instrumentos destinados a los implantes deben ser aplicados con el fin de tener una combinación segura.
- Con el fin de evitar o reducir al máximo los riesgos específicos asociados a la implantación, los cirujanos deben haber sido entrenados, y la operación debe ser realizada estrictamente de acuerdo con la guía técnica quirúrgica. Los instrumentos deben ser mantenidos y probados de forma regular de acuerdo con las instrucciones de cuidado y mantenimiento.
- Se requiere la manipulación y el almacenaje del producto con cuidado. La abrasión o daños de los componentes puede reducir significativamente la intensidad y la resistencia contra fatiga del producto.
- En general, no hay ninguna limitación entre el implante y las sustancias utilizadas en el ámbito clínico.
- Se tiene que evitar doblar la placa en ángulos extremos con radios pequeños por razón de que este podría provocar la fractura postoperatoria.
- Se tiene que evitar doblar la placa en ángulos extremos con radios pequeños por razón de que este podría provocar la fractura postoperatoria. El uso agresivo de exceso de los instrumentos de doblado puede causar daño macroscópicamente reconocible para el implante (agujeros alargados, bloqueados, etc.). En tales casos, el implante debe ser cambiado por uno nuevo, más cuidadosamente doblado.
- A menos que se indique en la técnica quirúrgica, la modificación del tamaño, forma o estado de superficie está prohibida después del despacho.

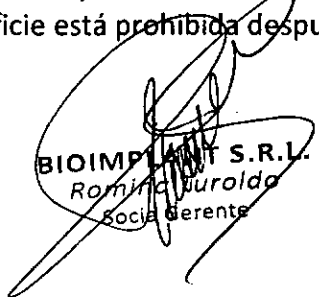

BIOIMPLANT S.R.L.
Romina Jurola
Socia Gerente

Ma. VICTORIA COZZI
FARMACEUTICA
Mat. 2194




Página 33 de 37

- Una vez aplicado, el producto no debe ser reutilizado. Aunque este pueda parecer en buen estado, las tensiones anteriores pueden haber producido imperfecciones que podrían reducir la vida útil del producto.
- Los dispositivos son incompatibles con el ambiente de MRI, a menos que se indique. Los riesgos posibles incluyen pero no limita en:
 - Los dispositivos se calentan o emigran
 - Afecta la imagen de MRI
- Se recomienda consultar al médico antes de entrar en ambientes que potencialmente podrían afectar el funcionamiento del implante.
- Los cirujanos deben verificar si los implantes están colocados correctamente usando la tecnología de imagen adecuada.
- Se prohíbe la actividad física o carga prematura después de la cirugía.
- El uso a largo plazo puede provocar aflojamiento o fractura del implante; Después de la consolidación del hueso, se sugiere sacar el implante, de lo contrario pueda provocar un accidente médico.
- Se debe aconsejar al paciente a informar cualquier cambio inusual en el sitio de operación suyo al cirujano. Los pacientes deben ser estrechamente monitorizados si un cambio en el sitio de fijación ha sido detectado. El cirujano debe evaluar la posibilidad del fracaso clínico posterior y discutir con el paciente sobre la necesidad de las medidas consideradas para ayudar a la curación.
- Las consideraciones para la extracción del implante después de la curación. Si el dispositivo no se quita después de la finalización de su uso previsto, cualquiera de las siguientes complicaciones se puede producir: (1) la corrosión, con reacción tisular o dolor; (2) la migración de la posición del implante dando lesiones, (3) el riesgo de lesiones adicionales de un trauma postoperatorio; (4) el doblado, aflojamiento, y / o la rotura, por los cuales podrían resultar dificultad al retirarse; (5) el dolor, la incomodidad o sensaciones anormales debido a la presencia del dispositivo; (6) la posibilidad de aumentar del riesgo de infección; (7) la osteopenia. El cirujano debe evaluar cuidadosamente los riesgos y beneficios al momento de decidir si se quita o no el implante. La extracción del implante debe ser seguido por el tratamiento postoperatorio adecuado para evitar una nueva fractura. Si el paciente es mayor de edad y tiene un bajo nivel de actividad, el cirujano puede decidir no retirar el implante así evitando los riesgos de una segunda cirugía.
- Los implantes sacados del cuerpo deben ser tratados de una manera anti-contaminación según estipulaciones de cada país con el fin de prevenir la infección cruzada.



BIOIMPLANT S.R.L.
Romina Vialto
Socia Gerente



Ma. VICTORIA COZZI
FARMACEUTICA
Mat. 2194

Nota: Las técnicas quirúrgicas para cada modelo se entregan junto con las placas correspondientes.

La Extracción del Implante

Exponga la placa y los tornillos de los tejidos blandos en su superficie mediante la incisión original utilizando el atornillador hexágono. Retire la placa al remover todos los tornillos.

La Eliminación del Implante

Los implantes retirados deben ser manejado de una manera anti-contaminación de acuerdo con las estipulaciones de cada país con el fin de evitar la infección cruzada.

3.4. Todas las informaciones que permitan comprobar si el producto médico está bien instalado y pueda funcionar correctamente y con plena seguridad, así como los datos relativos a la naturaleza y frecuencia de las operaciones de mantenimiento y calibrado que haya que efectuar para garantizar permanentemente el buen funcionamiento y la seguridad de los productos médicos: (Ver 3.3)

3.5. La información útil para evitar ciertos riesgos relacionados con la implantación del producto médico: (Ver 3.3)

3.6. La información relativa a los riesgos de interferencia recíproca relacionados con la presencia del producto médico en investigaciones o tratamientos específicos: (No aplica)

3.7. Las instrucciones necesarias en caso de rotura del envase protector de la esterilidad y si corresponde la indicación de los métodos adecuados de reesterilización:

Las placas bloqueadas son productos médicos no estériles que precisan ser esterilizados antes de su uso.


3.8. Si un producto médico está destinado a reutilizarse, los datos sobre los procedimientos apropiados para la reutilización, incluida la limpieza, desinfección, el acondicionamiento y, en su caso, el método de esterilización si el producto debe ser reesterilizado, así como cualquier limitación respecto al número posible de reutilizaciones. En caso de que los productos médicos deban ser esterilizados antes de su uso, las instrucciones de limpieza y esterilización deberán estar formuladas de modo que si se siguen correctamente, el producto siga cumpliendo los requisitos previstos en la Sección I (Requisitos Generales) del Anexo de la Resolución GMC N° 72/98 que dispone sobre los Requisitos Esenciales de Seguridad y Eficacia de los productos médicos:

Las placas bloqueadas son productos médicos de un solo uso que deben ser esterilizadas antes de usar.

- **Esterilización:** Autoclave al vacío por un mínimo de 4 minutos a 132°C. Tres ciclos mínimos.



BIOIMPLANT S.R.L.
Román Turoldo
Socia Gerente

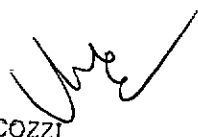

Ma. VICTORIA COZZI
FARMACEUTICA
Mat. 2194

- **Almacenamiento:**

1. El dispositivo debe ser almacenado en un lugar fresco, bien ventilado y seco a una temperatura por debajo de los 40°C y a una humedad relativa de entre el 30% y el 80%.
2. Proteja el dispositivo de la luz solar directa.

3.9. Información sobre cualquier tratamiento o procedimiento adicional que deba realizarse antes de utilizar el producto médico (por ejemplo, esterilización, montaje final, entre otros): Esterilizar antes de usar

- **Esterilización**

Los implantes e instrumentales deben ser esterilizados antes de que sean utilizados para la cirugía. Los métodos de esterilización para ellos son iguales. Todos los instrumentos enumerados en los manuales pueden ser esterilizados directamente sin ser descompuestos. Pueden ser re-esterilizados hasta 25 veces.


- **Parámetros de Esterilización:**
Esterilización: autoclave a vapor
Tipo de Ciclo: Pre-vacío pulsado
Tiempo de Esterilización: 4 minutos
Temperatura: 132°C
Tiempo de Secado: 30 minutos
Repeticiones Mínimas de Ciclos: Tres veces
 - **Esterilizador Recomendado:** Tuttnauer EHS Series Table-Top Autoclave; (O esterilizador similar con la misma configuración)
 - **Método de Envoltura:** la Bolsa
- Adicionalmente, por favor note los puntos siguientes:
- Si se lleva a cabo un método de esterilización que no sea el que se menciona por arriba, el usuario debe validarlo antes de la esterilización de acuerdo con la norma ANSI / AAMI ST 79: 2010 "Guía completa para esterilización grupal y el aseguramiento de esterilidad en centros de atención de la salud";
 - No esterilice el dispositivo con su empaque original
 - No esterilice el dispositivo contaminado por los fluidos corporales

LOS ITEMS 3.10, 3.11, 3.12, 3.13, 3.14, 3.15, 3.16 DE LA DISPOSICION: 2318/02 T.O. 2004, ANEXO IIIB

No aplican



BIOIMPLANT S.R.L.
Romina Wroldo
Socia Gerente



Ma. VICTORIA COZZI
FARMACEUTICA
Mat. 2194



República Argentina - Poder Ejecutivo Nacional
2019 - Año de la Exportación

Hoja Adicional de Firmas
Anexo

Número:

Referencia: Rot de Inst de Uso Bio Implant SRL

El documento fue importado por el sistema GEDO con un total de 38 pagina/s.

Digitally signed by GESTION DOCUMENTAL ELECTRONICA - GDE
DN: cn=GESTION DOCUMENTAL ELECTRONICA - GDE, c=AR, o=SECRETARIA DE GOBIERNO DE MODERNIZACION,
ou=SECRETARIA DE MODERNIZACION ADMINISTRATIVA, serialNumber=CUIT 30715117564
Date: 2019.07.19 11:01:31 -0300'

Digitally signed by GESTION DOCUMENTAL ELECTRONICA - GDE
DN: cn=GESTION DOCUMENTAL ELECTRONICA - GDE, c=AR, o=SECRETARIA DE GOBIERNO DE MODERNIZACION,
ou=SECRETARIA DE MODERNIZACION ADMINISTRATIVA, serialNumber=CUIT 30715117564
Date: 2019.07.19 11:01:33 -0300'



Secretaría de
Gobierno de Salud



Ministerio de Salud y Desarrollo Social
Presidencia de la Nación

CERTIFICADO DE AUTORIZACIÓN E INSCRIPCIÓN

Expediente N°: 1-47-3110-1051-19-0

La Administración Nacional de Medicamentos, Alimentos y Tecnología Médica (ANMAT) certifica que de acuerdo con lo solicitado por BIOIMPLANT S.R.L., se autoriza la inscripción en el Registro Nacional de Productores y Productos de Tecnología Médica (RPPTM), de un nuevo producto con los siguientes datos identificatorios característicos:

Nombre descriptivo: PLACAS BLOQUEADAS DE TITANIO

Código de identificación y nombre técnico UMDNS: 13-050- Placas, para Huesos

Marca(s) de (los) producto(s) médico(s): DOUBLE MEDICAL

Clase de Riesgo: III

Indicación/es autorizada/s: para tratar fracturas de huesos grandes (clavícula, escapula, pelvis, húmero, cúbito, radio, fémur, tibia y peroné) y pequeños (metacarpo, metatarso y falanges).

Modelo/s:

Nº Referencia	Descripción	Perforac.
020140003	T-Locking Plate, Combi-holes	3 Holes
020140004	T-Locking Plate, Combi-holes	4 Holes
020140005	T-Locking Plate, Combi-holes	5 Holes
020140006	T-Locking Plate, Combi-holes	6 Holes

Sedes y Delegaciones

Tel. (+54-11) 4340-0800 - <http://www.anmat.gov.ar> - República Argentina

Sede Central
Av. de Mayo 869, CABA

Sede Alsina
Alsina 665/671, CABA

Sede INAME
Av. Caseros 2161, CABA

Sede INAL
Estados Unidos 25, CABA

Sede Prod. Médicos
Av. Belgrano 1480, CABA

Deleg. Mendoza
Remedios de Escalada de
San Martín 1909, Mendoza
Prov. de Mendoza

Deleg. Córdoba
Obispo Trejo 635,
Córdoba,
Prov. de Córdoba

Deleg. Paso de los Libres
Ruta Nacional 117, km.10,
CO.TE.CAR., Paso de los
Libres, Prov. de Corrientes

Deleg. Posadas
Roque González 1137,
Posadas, Prov. de
Misiones

Deleg. Santa Fé
Eva Perón 2456,
Santa Fé,
Prov. de Santa Fé

020140007 T-Locking Plate, Combi-holes 7 Holes
020140008 T-Locking Plate, Combi-holes 8 Holes
020150004 T-Buttress Locking Plate, Combi-holes 4 Holes
020150005 T-Buttress Locking Plate, Combi-holes 5 Holes
020150006 T-Buttress Locking Plate, Combi-holes 6 Holes
020150003 T-Buttress Locking Plate, Combi-holes 3 Holes
020150007 T-Buttress Locking Plate, Combi-holes 7 Holes
020150008 T-Buttress Locking Plate, Combi-holes 8 Holes
020161003 L-Buttress Locking Plate left 3 Holes
020161004 L-Buttress Locking Plate left 4 Holes
020161005 L-Buttress Locking Plate left 5 Holes
020161006 L-Buttress Locking Plate left 6 Holes
020161007 L-Buttress Locking Plate left 7 Holes
020161008 L-Buttress Locking Plate left 8 Holes
020162003 L-Buttress Locking Plate, right 3 Holes
020162004 L-Buttress Locking Plate, right 4 Holes
020162005 L-Buttress Locking Plate, right 5 Holes
020162006 L-Buttress Locking Plate, right 6 Holes
020162007 L-Buttress Locking Plate, right 7 Holes
020162008 L-Buttress Locking Plate, right 8 Holes
020171004 Distal Femoral Sub-Condylar Locking Plate, Left 4 Holes

MB



Secretaría de
Gobierno de Salud



Ministerio de Salud y Desarrollo Social
Presidencia de la Nación

020171005	Distal Femoral Sub-Condylar Locking Plate, Left	5 Holes
020171006	Distal Femoral Sub-Condylar Locking Plate, Left	6 Holes
020171007	Distal Femoral Sub-Condylar Locking Plate, Left	7 Holes
020171008	Distal Femoral Sub-Condylar Locking Plate, Left	8 Holes
020171009	Distal Femoral Sub-Condylar Locking Plate, Left	9 Holes
020171010	Distal Femoral Sub-Condylar Locking Plate, Left	10 Holes
020171011	Distal Femoral Sub-Condylar Locking Plate, Left	11 Holes
020171012	Distal Femoral Sub-Condylar Locking Plate, Left	12 Holes
020172004	Distal Femoral Sub-Condylar Locking Plate, Right	4 Holes
020172005	Distal Femoral Sub-Condylar Locking Plate, Right	5 Holes
020172006	Distal Femoral Sub-Condylar Locking Plate, Right	6 Holes
020172007	Distal Femoral Sub-Condylar Locking Plate, Right	7 Holes
020172008	Distal Femoral Sub-Condylar Locking Plate, Right	8 Holes
020172009	Distal Femoral Sub-Condylar Locking Plate, Right	9 Holes
020172010	Distal Femoral Sub-Condylar Locking Plate, Right	10 Holes
020172011	Distal Femoral Sub-Condylar Locking Plate, Right	11 Holes
020172012	Distal Femoral Sub-Condylar Locking Plate, Right	12 Holes
020211060	Calcaneal Locking Plate I , Left	8 Holes
020211075	Calcaneal Locking Plate I , Left	10 Holes
020212060	Calcaneal Locking Plate I , Right	8 Holes

Sedes y Delegaciones

Tel. (+54-11) 4340-0800 - <http://www.anmat.gov.ar> - República Argentina

Sede Central
Av. de Mayo 869, CABA

Sede Alsina
Alsina 665/671, CABA

Sede INAME
Av. Caseros 2161, CABA

Sede INAL
Estados Unidos 25, CABA

Sede Prod. Médicos
Av. Belgrano 1480, CABA

Deleg. Mendoza
Remedios de Escalada de
San Martín 1909, Mendoza
Prov. de Mendoza

Deleg. Córdoba
Obispo Trejo 635,
Córdoba,
Prov. de Córdoba

Deleg. Paso de los Libres
Ruta Nacional 117, km.10,
CO.TE.CAR., Paso de los
Libres, Prov. de Corrientes

Deleg. Posadas
Roque González 1137,
Posadas, Prov. de
Misiones

Deleg. Santa Fé
Eva Perón 2456,
Santa Fé,
Prov. de Santa Fé

020212075 Calcaneal Locking Plate I , Right 10 Holes
020221015 Calcaneal Locking Plate II , Left 15 Holes
020222015 Calcaneal Locking Plate II , Right 15 Holes
020231015 Calcaneal Locking Plate II , Left 15 Holes
020232015 Calcaneal Locking Plate II , Right 15 Holes
020251003 Proximal Femoral Locking Plate, Multi-holes, Left 3 Holes
020251004 Proximal Femoral Locking Plate, Multi-holes, Left 4 Holes
020251005 Proximal Femoral Locking Plate, Multi-holes, Left 5 Holes
020251006 Proximal Femoral Locking Plate, Multi-holes, Left 6 Holes
020251007 Proximal Femoral Locking Plate, Multi-holes, Left 7 Holes
020251008 Proximal Femoral Locking Plate, Multi-holes, Left 8 Holes
020251009 Proximal Femoral Locking Plate, Multi-holes, Left 9 Holes
020251010 Proximal Femoral Locking Plate, Multi-holes, Left 10 Holes
020251011 Proximal Femoral Locking Plate, Multi-holes, Left 11 Holes
020251012 Proximal Femoral Locking Plate, Multi-holes, Left 12 Holes
020252003 Proximal Femoral Locking Plate, Multi-holes, Right 3 Holes
020252004 Proximal Femoral Locking Plate, Multi-holes, Right 4 Holes
020252005 Proximal Femoral Locking Plate, Multi-holes, Right 5 Holes
020252006 Proximal Femoral Locking Plate, Multi-holes, Right 6 Holes
020252007 Proximal Femoral Locking Plate, Multi-holes, Right 7 Holes

UMP



Secretaría de
Gobierno de Salud



Ministerio de Salud y Desarrollo Social
Presidencia de la Nación

- 020252008 Proximal Femoral Locking Plate, Multi-holes, Right 8 Holes
- 020252009 Proximal Femoral Locking Plate, Multi-holes, Right 9 Holes
- 020252010 Proximal Femoral Locking Plate, Multi-holes, Right 10 Holes
- 020252011 Proximal Femoral Locking Plate, Multi-holes, Right 11 Holes
- 020252012 Proximal Femoral Locking Plate, Multi-holes, Right 12 Holes
- 020280004 Pelvic Reconstruction Locking Plate 4 Holes
- 020280005 Pelvic Reconstruction Locking Plate 5 Holes
- 020280006 Pelvic Reconstruction Locking Plate 6 Holes
- 020280007 Pelvic Reconstruction Locking Plate 7 Holes
- 020280008 Pelvic Reconstruction Locking Plate 8 Holes
- 020280009 Pelvic Reconstruction Locking Plate 9 Holes
- 020280010 Pelvic Reconstruction Locking Plate 10 Holes
- 020280011 Pelvic Reconstruction Locking Plate 11 Holes
- 020280012 Pelvic Reconstruction Locking Plate 12 Holes
- 020280013 Pelvic Reconstruction Locking Plate 13 Holes
- 020280014 Pelvic Reconstruction Locking Plate 14 Holes
- 020320000 3.5mm Locking Spacer
- 020330000 5.0mm Locking Spacer
- 020360000 3.5mm Cable Spacer II
- 020370000 5.0mm Cable Spacer II

WAS

Sedes y Delegaciones

Tel. (+54-11) 4340-0800 - <http://www.anmat.gov.ar> - República Argentina

Sede Central
Av. de Mayo 869, CABA

Sede Alsina
Alsina 665/671, CABA

Sede INAME
Av. Caseros 2161, CABA

Sede INAL
Estados Unidos 25, CABA

Sede Prod. Médicos
Av. Belgrano 1480, CABA

Deleg. Mendoza
Remedios de Escalada de
San Martín 1909, Mendoza
Prov. de Mendoza

Deleg. Córdoba
Obispo Trejo 635,
Córdoba,
Prov. de Córdoba

Deleg. Paso de los Libres
Ruta Nacional 117, km.10,
CO.TE.CAR., Paso de los
Libres, Prov. de Corrientes

Deleg. Posadas
Roque González 1137,
Posadas, Prov. de
Misiones

Deleg. Santa Fé
Eva Perón 2456,
Santa Fé,
Prov. de Santa Fé

020380003	Proximal Humeral Locking Plate II, Combi-holes	3 Holes
020380004	Proximal Humeral Locking Plate II, Combi-holes	4 Holes
020380005	Proximal Humeral Locking Plate II, Combi-holes	5 Holes
020380006	Proximal Humeral Locking Plate II, Combi-holes	6 Holes
020380007	Proximal Humeral Locking Plate II, Combi-holes	7 Holes
020380008	Proximal Humeral Locking Plate II, Combi-holes	8 Holes
020380009	Proximal Humeral Locking Plate II, Combi-holes	9 Holes
020380010	Proximal Humeral Locking Plate II, Combi-holes	10 Holes
020380011	Proximal Humeral Locking Plate II, Combi-holes	11 Holes
020380012	Proximal Humeral Locking Plate II, Combi-holes	12 Holes
020380013	Proximal Humeral Locking Plate II, Combi-holes	13 Holes
020400006	Broad Limited Contact Locking Plate, Combi-holes	6 Holes
020400007	Broad Limited Contact Locking Plate, Combi-holes	7 Holes
020400008	Broad Limited Contact Locking Plate, Combi-holes	8 Holes
020400009	Broad Limited Contact Locking Plate, Combi-holes	9 Holes
020400010	Broad Limited Contact Locking Plate, Combi-holes	10 Holes
020400011	Broad Limited Contact Locking Plate, Combi-holes	11 Holes
020400012	Broad Limited Contact Locking Plate, Combi-holes	12 Holes
020400013	Broad Limited Contact Locking Plate, Combi-holes	13 Holes
020400014	Broad Limited Contact Locking Plate, Combi-holes	14 Holes

WFB



Secretaría de
Gobierno de Salud



Ministerio de Salud y Desarrollo Social
Presidencia de la Nación

- 020400015 Broad Limited Contact Locking Plate, Combi-holes 15 Holes
- 020400016 Broad Limited Contact Locking Plate, Combi-holes 16 Holes
- 020410004 Narrow Limited Contact Locking Plate, Combi-holes 4 Holes
- 020410005 Narrow Limited Contact Locking Plate, Combi-holes 5 Holes
- 020410006 Narrow Limited Contact Locking Plate, Combi-holes 6 Holes
- 020410007 Narrow Limited Contact Locking Plate, Combi-holes 7 Holes
- 020410008 Narrow Limited Contact Locking Plate, Combi-holes 8 Holes
- 020410009 Narrow Limited Contact Locking Plate, Combi-holes 9 Holes
- 020410010 Narrow Limited Contact Locking Plate, Combi-holes 10 Holes
- 020410011 Narrow Limited Contact Locking Plate, Combi-holes 11 Holes
- 020410012 Narrow Limited Contact Locking Plate, Combi-holes 12 Holes
- 020410013 Narrow Limited Contact Locking Plate, Combi-holes 13 Holes
- 020410014 Narrow Limited Contact Locking Plate, Combi-holes 14 Holes
- 020410015 Narrow Limited Contact Locking Plate, Combi-holes 15 Holes
- 020410016 Narrow Limited Contact Locking Plate, Combi-holes 16 Holes
- 020421003 Proximal Femoral Locking Plate, Combi-holes, Left 3 Holes
- 020421004 Proximal Femoral Locking Plate, Combi-holes, Left 4 Holes
- 020421005 Proximal Femoral Locking Plate, Combi-holes, Left 5 Holes
- 020421006 Proximal Femoral Locking Plate, Combi-holes, Left 6 Holes
- 020421007 Proximal Femoral Locking Plate, Combi-holes, Left 7 Holes
- 020421008 Proximal Femoral Locking Plate, Combi-holes, Left 8 Holes

Sedes y Delegaciones

Tel. (+54-11) 4340-0800 - <http://www.anmat.gov.ar> - República Argentina

Sede Central

Av. de Mayo 869, CABA

Sede Alsina

Alsina 665/671, CABA

Sede INAME

Av. Caseros 2161, CABA

Sede INAL

Estados Unidos 25, CABA

Sede Prod. Médicos

Av. Belgrano 1480, CABA

Deleg. Mendoza

Remedios de Escalada de
San Martín 1909, Mendoza
Prov. de Mendoza

Deleg. Córdoba

Obispo Trejo 635,
Córdoba,
Prov. de Córdoba

Deleg. Paso de los Libres

Ruta Nacional 117, km.10,
CO.TE.CAR., Paso de los
Libres, Prov. de Corrientes

Deleg. Posadas

Roque González 1137,
Posadas, Prov. de
Misiones

Deleg. Santa Fé

Eva Perón 2456,
Santa Fé,
Prov. de Santa Fé

020421009	Proximal Femoral Locking Plate, Combi-holes, Left	9 Holes
020421010	Proximal Femoral Locking Plate, Combi-holes, Left	10 Holes
020421011	Proximal Femoral Locking Plate, Combi-holes, Left	11 Holes
020421012	Proximal Femoral Locking Plate, Combi-holes, Left	12 Holes
020421013	Proximal Femoral Locking Plate, Combi-holes, Left	13 Holes
020422003	Proximal Femoral Locking Plate, Combi-holes, Right	3 Holes
020422004	Proximal Femoral Locking Plate, Combi-holes, Right	4 Holes
020422005	Proximal Femoral Locking Plate, Combi-holes, Right	5 Holes
020422006	Proximal Femoral Locking Plate, Combi-holes, Right	6 Holes
020422007	Proximal Femoral Locking Plate, Combi-holes, Right	7 Holes
020422008	Proximal Femoral Locking Plate, Combi-holes, Right	8 Holes
020422009	Proximal Femoral Locking Plate, Combi-holes, Right	9 Holes
020422010	Proximal Femoral Locking Plate, Combi-holes, Right	10 Holes
020422011	Proximal Femoral Locking Plate, Combi-holes, Right	11 Holes
020422012	Proximal Femoral Locking Plate, Combi-holes, Right	12 Holes
020422013	Proximal Femoral Locking Plate, Combi-holes, Right	13 Holes
020451003	Distal Medial Tibial Locking Plate, Combi-holes, Left	3 Holes
020451004	Distal Medial Tibial Locking Plate, Combi-holes, Left	4 Holes
020451005	Distal Medial Tibial Locking Plate, Combi-holes, Left	5 Holes
020451006	Distal Medial Tibial Locking Plate, Combi-holes, Left	6 Holes
020451007	Distal Medial Tibial Locking Plate, Combi-holes, Left	7 Holes

WJB



Secretaría de
Gobierno de Salud



Ministerio de Salud y Desarrollo Social
Presidencia de la Nación

020451008	Distal Medial Tibial Locking Plate, Combi-holes, Left	8 Holes
020451009	Distal Medial Tibial Locking Plate, Combi-holes, Left	9 Holes
020451010	Distal Medial Tibial Locking Plate, Combi-holes, Left	10 Holes
020451011	Distal Medial Tibial Locking Plate, Combi-holes, Left	11 Holes
020451012	Distal Medial Tibial Locking Plate, Combi-holes, Left	12 Holes
020451013	Distal Medial Tibial Locking Plate, Combi-holes, Left	13 Holes
020452003	Distal Medial Tibial Locking Plate, Combi-holes, Right	3 Holes
020452004	Distal Medial Tibial Locking Plate, Combi-holes, Right	4 Holes
020452005	Distal Medial Tibial Locking Plate, Combi-holes, Right	5 Holes
020452006	Distal Medial Tibial Locking Plate, Combi-holes, Right	6 Holes
020452007	Distal Medial Tibial Locking Plate, Combi-holes, Right	7 Holes
020452008	Distal Medial Tibial Locking Plate, Combi-holes, Right	8 Holes
020452009	Distal Medial Tibial Locking Plate, Combi-holes, Right	9 Holes
020452010	Distal Medial Tibial Locking Plate, Combi-holes, Right	10 Holes
020452011	Distal Medial Tibial Locking Plate, Combi-holes, Right	11 Holes
020452012	Distal Medial Tibial Locking Plate, Combi-holes, Right	12 Holes
020452013	Distal Medial Tibial Locking Plate, Combi-holes, Right	13 Holes
020470004	Small Limited Contact Locking Plate, Combi-holes	4 Holes
020470005	Small Limited Contact Locking Plate, Combi-holes	5 Holes
020470006	Small Limited Contact Locking Plate, Combi-holes	6 Holes
020470007	Small Limited Contact Locking Plate, Combi-holes	7 Holes

WMB

Sedes y Delegaciones

Tel. (+54-11) 4340-0800 - <http://www.anmat.gov.ar> - República Argentina

Sede Central

Av. de Mayo 869, CABA

Sede Alsina

Alsina 665/671, CABA

Sede INAME

Av. Caseros 2161, CABA

Sede INAL

Estados Unidos 25, CABA

Sede Prod. Médicos

Av. Belgrano 1480, CABA

Deleg. Mendoza

Remedios de Escalada de
San Martín 1909, Mendoza
Prov. de Mendoza

Deleg. Córdoba

Obispo Trejo 635,
Córdoba,
Prov. de Córdoba

Deleg. Paso de los Libres

Ruta Nacional 117, km.10,
CO.TE.CAR., Paso de los
Libres, Prov. de Corrientes

Deleg. Posadas

Roque González 1137,
Posadas, Prov. de
Misiones

Deleg. Santa Fé

Eva Perón 2456,
Santa Fé,
Prov. de Santa Fé

020470008 Small Limited Contact Locking Plate, Combi-holes 8 Holes
020470009 Small Limited Contact Locking Plate, Combi-holes 9 Holes
020470010 Small Limited Contact Locking Plate, Combi-holes 10 Holes
020470011 Small Limited Contact Locking Plate, Combi-holes 11 Holes
020470012 Small Limited Contact Locking Plate, Combi-holes 12 Holes
020470013 Small Limited Contact Locking Plate, Combi-holes 13 Holes
020470014 Small Limited Contact Locking Plate, Combi-holes 14 Holes
020480004 Reconstruction Locking Plate, combi-holes 4 Holes
020480005 Reconstruction Locking Plate, combi-holes 5 Holes
020480006 Reconstruction Locking Plate, combi-holes 6 Holes
020480007 Reconstruction Locking Plate, combi-holes 7 Holes
020480008 Reconstruction Locking Plate, combi-holes 8 Holes
020480009 Reconstruction Locking Plate, combi-holes 9 Holes
020480010 Reconstruction Locking Plate, combi-holes 10 Holes
020480011 Reconstruction Locking Plate, combi-holes 11 Holes
020480012 Reconstruction Locking Plate, combi-holes 12 Holes
020480013 Reconstruction Locking Plate, combi-holes 13 Holes
020480014 Reconstruction Locking Plate, combi-holes 14 Holes
020490004 One-third Tubular Locking Plate, Combi-holes 4 Holes
020490005 One-third Tubular Locking Plate, Combi-holes 5 Holes
020490006 One-third Tubular Locking Plate, Combi-holes 6 Holes

WFB



Secretaría de
Gobierno de Salud



Ministerio de Salud y Desarrollo Social
Presidencia de la Nación

020490007	One-third Tubular Locking Plate, Combi-holes	7 Holes
020490008	One-third Tubular Locking Plate, Combi-holes	8 Holes
020490009	One-third Tubular Locking Plate, Combi-holes	9 Holes
020490010	One-third Tubular Locking Plate, Combi-holes	10 Holes
020490011	One-third Tubular Locking Plate, Combi-holes	11 Holes
020490012	One-third Tubular Locking Plate, Combi-holes	12 Holes
020501004	Clavicular Locking Plate, Combi-holes, Left	4 Holes
020501005	Clavicular Locking Plate, Combi-holes, Left	5 Holes
020501006	Clavicular Locking Plate, Combi-holes, Left	6 Holes
020501007	Clavicular Locking Plate, Combi-holes, Left	7 Holes
020501008	Clavicular Locking Plate, Combi-holes, Left	8 Holes
020501009	Clavicular Locking Plate, Combi-holes, Left	9 Holes
020501010	Clavicular Locking Plate, Combi-holes, Left	10 Holes
020501011	Clavicular Locking Plate, Combi-holes, Left	11 Holes
020501012	Clavicular Locking Plate, Combi-holes, Left	12 Holes
020502004	Clavicular Locking Plate, Combi-holes, Right	4 Holes
020502005	Clavicular Locking Plate, Combi-holes, Right	5 Holes
020502006	Clavicular Locking Plate, Combi-holes, Right	6 Holes
020502007	Clavicular Locking Plate, Combi-holes, Right	7 Holes
020502008	Clavicular Locking Plate, Combi-holes, Right	8 Holes
020502009	Clavicular Locking Plate, Combi-holes, Right	9 Holes

WAB
Sedes y Delegaciones

Tel. (+54-11) 4340-0800 - <http://www.anmat.gov.ar> - República Argentina

Sede Central

Av. de Mayo 869, CABA

Sede Alsina

Alsina 665/671, CABA

Sede INAME

Av. Caseros 2161, CABA

Sede INAL

Estados Unidos 25, CABA

Sede Prod. Médicos

Av. Belgrano 1480, CABA

Deleg. Mendoza

Remedios de Escalada de
San Martín 1909, Mendoza
Prov. de Mendoza

Deleg. Córdoba

Obispo Trejo 635,
Córdoba,
Prov. de Córdoba

Deleg. Paso de los Libres

Ruta Nacional 117, km.10,
CO.TE.CAR., Paso de los
Libres, Prov. de Corrientes

Deleg. Posadas

Roque González 1137,
Posadas, Prov. de
Misiones

Deleg. Santa Fé

Eva Perón 2456,
Santa Fé,
Prov. de Santa Fé

020502010	Clavicular Locking Plate, Combi-holes, Right	10 Holes
020502011	Clavicular Locking Plate, Combi-holes, Right	11 Holes
020502012	Clavicular Locking Plate, Combi-holes, Right	12 Holes
020511002	Distal Fibular Locking Plate, Combi-holes, Left	2 Holes
020511003	Distal Fibular Locking Plate, Combi-holes, Left	3 Holes
020511004	Distal Fibular Locking Plate, Combi-holes, Left	4 Holes
020511005	Distal Fibular Locking Plate, Combi-holes, Left	5 Holes
020511006	Distal Fibular Locking Plate, Combi-holes, Left	6 Holes
020511007	Distal Fibular Locking Plate, Combi-holes, Left	7 Holes
020511008	Distal Fibular Locking Plate, Combi-holes, Left	8 Holes
020511009	Distal Fibular Locking Plate, Combi-holes, Left	9 Holes
020512002	Distal Fibular Locking Plate, Combi-holes, Right	2 Holes
020512003	Distal Fibular Locking Plate, Combi-holes, Right	3 Holes
020512004	Distal Fibular Locking Plate, Combi-holes, Right	4 Holes
020512005	Distal Fibular Locking Plate, Combi-holes, Right	5 Holes
020512006	Distal Fibular Locking Plate, Combi-holes, Right	6 Holes
020512007	Distal Fibular Locking Plate, Combi-holes, Right	7 Holes
020512008	Distal Fibular Locking Plate, Combi-holes, Right	8 Holes
020512009	Distal Fibular Locking Plate, Combi-holes, Right	9 Holes
020520003	Olecranon Locking Plate III, Combi-holes	3 Holes
020520004	Olecranon Locking Plate III, Combi-holes	4 Holes

WAB



Secretaría de
Gobierno de Salud



Ministerio de Salud y Desarrollo Social
Presidencia de la Nación

020520005	Olecranon Locking Plate III, Combi-holes	5 Holes
020520006	Olecranon Locking Plate III, Combi-holes	6 Holes
020520007	Olecranon Locking Plate III, Combi-holes	7 Holes
020520008	Olecranon Locking Plate III, Combi-holes	8 Holes
020541000	Proximal Femoral Locking Plate for 7.3mm Cannulated Locking Screw, Left 51mm	
020542000	Proximal Femoral Locking Plate for 7.3mm Cannulated Locking Screw, Right 51mm	
020621004	Distal Humeral Sub-Condylar Locking Plate, Left	4 Holes
020621007	Distal Humeral Sub-Condylar Locking Plate, Left	7 Holes
020621010	Distal Humeral Sub-Condylar Locking Plate, Left	10 Holes
020622004	Distal Humeral Sub-Condylar Locking Plate, Right	4 Holes
020622007	Distal Humeral Sub-Condylar Locking Plate, Right	7 Holes
020622010	Distal Humeral Sub-Condylar Locking Plate, Right	10 Holes
020631002	Distal Radius Locking Plate I, Small, Combi-holes, Left	2 Holes
020631003	Distal Radius Locking Plate I, Small, Combi-holes, Left	3 Holes
020631004	Distal Radius Locking Plate I, Small, Combi-holes, Left	4 Holes
020631005	Distal Radius Locking Plate I, Small, Combi-holes, Left	5 Holes
020632002	Distal Radius Locking Plate I, Small, Combi-holes, Right	2 Holes
020632003	Distal Radius Locking Plate I, Small, Combi-holes, Right	3 Holes

Handwritten signature/initials

Sedes y Delegaciones

Tel. (+54-11) 4340-0800 - <http://www.anmat.gov.ar> - República Argentina

Sede Central
Av. de Mayo 869, CABA

Sede Alsina
Alsina 665/671, CABA

Sede INAME
Av. Caseros 2161, CABA

Sede INAL
Estados Unidos 25, CABA

Sede Prod. Médicos
Av. Belgrano 1480, CABA

Deleg. Mendoza
Remedios de Escalada de
San Martín 1909, Mendoza
Prov. de Mendoza

Deleg. Córdoba
Obispo Trejo 635,
Córdoba,
Prov. de Córdoba

Deleg. Paso de los Libres
Ruta Nacional 117, km.10,
CO.TE.CAR., Paso de los
Libres, Prov. de Corrientes

Deleg. Posadas
Roque González 1137,
Posadas, Prov. de
Misiones

Deleg. Santa Fé
Eva Perón 2456,
Santa Fé,
Prov. de Santa Fé

020632004 Distal Radius Locking Plate I , Small, Combi-holes, Right 4 Holes
020632005 Distal Radius Locking Plate I , Small, Combi-holes, Right 5 Holes
020641002 Distal Radius Locking Plate I , Large, Combi-holes, Left 2 Holes
020641003 Distal Radius Locking Plate I , Large, Combi-holes, Left 3 Holes
020641004 Distal Radius Locking Plate I , Large, Combi-holes, Left 4 Holes
020641005 Distal Radius Locking Plate I , Large, Combi-holes, Left 5 Holes
020642002 Distal Radius Locking Plate I , Large, Combi-holes, Right 2 Holes
020642003 Distal Radius Locking Plate I , Large, Combi-holes, Right 3 Holes
020642004 Distal Radius Locking Plate I , Large, Combi-holes, Right 4 Holes
020642005 Distal Radius Locking Plate I , Large, Combi-holes, Right 5 Holes
020651002 Proximal Radius Locking Plate, left 2 Holes
020651003 Proximal Radius Locking Plate, left 3 Holes
020651004 Proximal Radius Locking Plate, left 4 Holes
020652002 Proximal Radius Locking Plate, right 2 Holes
020652003 Proximal Radius Locking Plate, right 3 Holes
020652004 Proximal Radius Locking Plate, right 4 Holes
020661003 Distal Clavicular Locking Plate, Combi-holes, Left 3 Holes
020661004 Distal Clavicular Locking Plate, Combi-holes, Left 4 Holes
020661005 Distal Clavicular Locking Plate, Combi-holes, Left 5 Holes
020661006 Distal Clavicular Locking Plate, Combi-holes, Left 6 Holes

WAB



Secretaría de
Gobierno de Salud



Ministerio de Salud y Desarrollo Social
Presidencia de la Nación

- 020661007 Distal Clavicular Locking Plate, Combi-holes, Left 7 Holes
- 020661008 Distal Clavicular Locking Plate, Combi-holes, Left 8 Holes
- 020662003 Distal Clavicular Locking Plate, Combi-holes, Right 3 Holes
- 020662004 Distal Clavicular Locking Plate, Combi-holes, Right 4 Holes
- 020662005 Distal Clavicular Locking Plate, Combi-holes, Right 5 Holes
- 020662006 Distal Clavicular Locking Plate, Combi-holes, Right 6 Holes
- 020662007 Distal Clavicular Locking Plate, Combi-holes, Right 7 Holes
- 020662008 Distal Clavicular Locking Plate, Combi-holes, Right 8 Holes
- 020811003 Clavicular Hook Locking Plate III, Height 15mm, Combi-holes, Left 3 Holes
- 020811004 Clavicular Hook Locking Plate III, Height 15mm, Combi-holes, Left 4 Holes
- 020811005 Clavicular Hook Locking Plate III, Height 15mm, Combi-holes, Left 5 Holes
- 020811006 Clavicular Hook Locking Plate III, Height 15mm, Combi-holes, Left 6 Holes
- 020811007 Clavicular Hook Locking Plate III, Height 15mm, Combi-holes, Left 7 Holes
- 020811008 Clavicular Hook Locking Plate III, Height 15mm, Combi-holes, Left 8 Holes

WIB

Sedes y Delegaciones

Tel. (+54-11) 4340-0800 - <http://www.anmat.gov.ar> - República Argentina

Sede Central
Av. de Mayo 869, CABA

Sede Alsina
Alsina 665/671, CABA

Sede INAME
Av. Caseros 2161, CABA

Sede INAL
Estados Unidos 25, CABA

Sede Prod. Médicos
Av. Belgrano 1480, CABA

Deleg. Mendoza
Remedios de Escalada de
San Martín 1909, Mendoza
Prov. de Mendoza

Deleg. Córdoba
Obispo Trejo 635,
Córdoba,
Prov. de Córdoba

Deleg. Paso de los Libres
Ruta Nacional 117, km.10,
CO.TE.CAR., Paso de los
Libres, Prov. de Corrientes

Deleg. Posadas
Roque González 1137,
Posadas, Prov. de
Misiones

Deleg. Santa Fé
Eva Perón 2456,
Santa Fé,
Prov. de Santa Fé

020812003 Clavicular Hook Locking Plate III, Height 15mm, Combi-holes, Right 3
Holes

020812004 Clavicular Hook Locking Plate III, Height 15mm, Combi-holes, Right 4
Holes

020812005 Clavicular Hook Locking Plate III, Height 15mm, Combi-holes, Right 5
Holes

020812006 Clavicular Hook Locking Plate III, Height 15mm, Combi-holes, Right 6
Holes

020812007 Clavicular Hook Locking Plate III, Height 15mm, Combi-holes, Right 7
Holes

020812008 Clavicular Hook Locking Plate III, Height 15mm, Combi-holes, Right 8
Holes

020821003 Clavicular Hook Locking Plate III, Height 18mm, Combi-holes, Left 3
Holes

020821004 Clavicular Hook Locking Plate III, Height 18mm, Combi-holes, Left 4
Holes

020821005 Clavicular Hook Locking Plate III, Height 18mm, Combi-holes, Left 5
Holes

020821006 Clavicular Hook Locking Plate III, Height 18mm, Combi-holes, Left 6
Holes

WAB



Secretaría de
Gobierno de Salud



Ministerio de Salud y Desarrollo Social
Presidencia de la Nación

020821007	Clavicular Hook Locking Plate III, Height 18mm, Combi-holes, Left	7
	Holes	
020821008	Clavicular Hook Locking Plate III, Height 18mm, Combi-holes, Left	8
	Holes	
020822003	Clavicular Hook Locking Plate III, Height 18mm, Combi-holes, Right	3
	Holes	
020822004	Clavicular Hook Locking Plate III, Height 18mm, Combi-holes, Right	4
	Holes	
020822005	Clavicular Hook Locking Plate III, Height 18mm, Combi-holes, Right	5
	Holes	
020822006	Clavicular Hook Locking Plate III, Height 18mm, Combi-holes, Right	6
	Holes	
020822007	Clavicular Hook Locking Plate III, Height 18mm, Combi-holes, Right	7
	Holes	
020822008	Clavicular Hook Locking Plate III, Height 18mm, Combi-holes, Right	8
	Holes	
020961005	Proximal Humeral Greater Tuberosity Locking Compression Plate, left	5
	Holes	
020962005	Proximal Humeral Greater Tuberosity Locking Compression Plate, right	5
	Holes	

WJB

Sedes y Delegaciones

Tel. (+54-11) 4340-0800 - <http://www.anmat.gov.ar> - República Argentina

Sede Central
Av. de Mayo 869, CABA

Sede Alsina
Alsina 665/671, CABA

Sede INAME
Av. Caseros 2161, CABA

Sede INAL
Estados Unidos 25, CABA

Sede Prod. Médicos
Av. Belgrano 1480, CABA

Deleg. Mendoza
Remedios de Escalada de
San Martín 1909, Mendoza
Prov. de Mendoza

Deleg. Córdoba
Obispo Trejo 635,
Córdoba,
Prov. de Córdoba

Deleg. Paso de los Libres
Ruta Nacional 117, km.10,
CO.TE.CAR., Paso de los
Libres, Prov. de Corrientes

Deleg. Posadas
Roque González 1137,
Posadas, Prov. de
Misiones

Deleg. Santa Fé
Eva Perón 2456,
Santa Fé,
Prov. de Santa Fé

020971000 Proximal Femoral Locking Plate, for 7.3mm Cannulated Locking Screw,
small, left 43mm

020972000 Proximal Femoral Locking Plate, for 7.3mm Cannulated Locking Screw,
small, right 43mm

021301003 Distal Lateral Femoral Condylar Locking Plate, left 3 Holes

021302003 Distal Lateral Femoral Condylar Locking Plate, right 3 Holes

021311003 Distal Medial Femoral Condylar Locking Plate, left 3 Holes

021312003 Distal Medial Femoral Condylar Locking Plate, right 3 Holes

021321004 Distal Medial Femoral Sub-condylar Locking Plate, large, left 4 Holes

021321007 Distal Medial Femoral Sub-condylar Locking Plate, large, left 7 Holes

021322004 Distal Medial Femoral Sub-condylar Locking Plate, large, right 4 Holes

021322007 Distal Medial Femoral Sub-condylar Locking Plate, large, right 7 Holes

021331004 Distal Medial Femoral Sub-condylar Locking Plate, small, left 4 Holes

021331007 Distal Medial Femoral Sub-condylar Locking Plate, small, left 7 Holes

021332004 Distal Medial Femoral Sub-condylar Locking Plate, small, right 4 Holes

021332007 Distal Medial Femoral Sub-condylar Locking Plate, small, right 7 Holes

021371003 Distal Lateral Humeral Locking Plate, Combi-holes, Left 3 Holes

021371004 Distal Lateral Humeral Locking Plate, Combi-holes, Left 4 Holes

021371005 Distal Lateral Humeral Locking Plate, Combi-holes, Left 5 Holes

021371006 Distal Lateral Humeral Locking Plate, Combi-holes, Left 6 Holes

021371007 Distal Lateral Humeral Locking Plate, Combi-holes, Left 7 Holes

WMB



Secretaría de
Gobierno de Salud



Ministerio de Salud y Desarrollo Social
Presidencia de la Nación

- 021371008 Distal Lateral Humeral Locking Plate, Combi-holes, Left 8 Holes
- 021371009 Distal Lateral Humeral Locking Plate, Combi-holes, Left 9 Holes
- 021372003 Distal Lateral Humeral Locking Plate, Combi-holes, Right 3 Holes
- 021372004 Distal Lateral Humeral Locking Plate, Combi-holes, Right 4 Holes
- 021372005 Distal Lateral Humeral Locking Plate, Combi-holes, Right 5 Holes
- 021372006 Distal Lateral Humeral Locking Plate, Combi-holes, Right 6 Holes
- 021372007 Distal Lateral Humeral Locking Plate, Combi-holes, Right 7 Holes
- 021372008 Distal Lateral Humeral Locking Plate, Combi-holes, Right 8 Holes
- 021372009 Distal Lateral Humeral Locking Plate, Combi-holes, Right 9 Holes
- 021381003 Distal Medial Humeral Locking Plate, Combi-holes, Left 3 Holes
- 021381004 Distal Medial Humeral Locking Plate, Combi-holes, Left 4 Holes
- 021381005 Distal Medial Humeral Locking Plate, Combi-holes, Left 5 Holes
- 021381006 Distal Medial Humeral Locking Plate, Combi-holes, Left 6 Holes
- 021381007 Distal Medial Humeral Locking Plate, Combi-holes, Left 7 Holes
- 021381008 Distal Medial Humeral Locking Plate, Combi-holes, Left 8 Holes
- 021381009 Distal Medial Humeral Locking Plate, Combi-holes, Left 9 Holes
- 021382003 Distal Medial Humeral Locking Plate, Combi-holes, Right 3 Holes
- 021382004 Distal Medial Humeral Locking Plate, Combi-holes, Right 4 Holes
- 021382005 Distal Medial Humeral Locking Plate, Combi-holes, Right 5 Holes
- 021382006 Distal Medial Humeral Locking Plate, Combi-holes, Right 6 Holes
- 021382007 Distal Medial Humeral Locking Plate, Combi-holes, Right 7 Holes

Sedes y Delegaciones

Tel. (+54-11) 4340-0800 - <http://www.anmat.gov.ar> - República Argentina

Sede Central
Av. de Mayo 869, CABA

Sede Alsina
Alsina 665/671, CABA

Sede INAME
Av. Caseros 2161, CABA

Sede INAL
Estados Unidos 25, CABA

Sede Prod. Médicos
Av. Belgrano 1480, CABA

Deleg. Mendoza
Remedios de Escalada de
San Martín 1909, Mendoza
Prov. de Mendoza

Deleg. Córdoba
Obispo Trejo 635,
Córdoba,
Prov. de Córdoba

Deleg. Paso de los Libres
Ruta Nacional 117, km.10,
CO.TE.CAR., Paso de los
Libres, Prov. de Corrientes

Deleg. Posadas
Roque González 1137,
Posadas, Prov. de
Misiones

Deleg. Santa Fé
Eva Perón 2456,
Santa Fé,
Prov. de Santa Fé

W/B

021382008	Distal Medial Humeral Locking Plate, Combi-holes, Right	8 Holes
021382009	Distal Medial Humeral Locking Plate, Combi-holes, Right	9 Holes
021390003	T-Locking Plate, Right Angled, Head 3 Holes, Combi-holes	3 Holes
021400004	T-Locking Plate, Right Angled, Head 4 Holes, Combi-holes	4 Holes
021390005	T-Locking Plate, Right Angled, Head 3 Holes, Combi-holes	5 Holes
021400006	T-Locking Plate, Right Angled, Head 4 Holes, Combi-holes	6 Holes
021411003	T-Locking Plate, Oblique Angled, Combi-holes, Left	3 Holes
021411004	T-Locking Plate, Oblique Angled, Combi-holes, Left	4 Holes
021411005	T-Locking Plate, Oblique Angled, Combi-holes, Left	5 Holes
021411006	T-Locking Plate, Oblique Angled, Combi-holes, Left	6 Holes
021412003	T-Locking Plate, Oblique Angled, Combi-holes, Right	3 Holes
021412004	T-Locking Plate, Oblique Angled, Combi-holes, Right	4 Holes
021412005	T-Locking Plate, Oblique Angled, Combi-holes, Right	5 Holes
021412006	T-Locking Plate, Oblique Angled, Combi-holes, Right	6 Holes
021421003	Olecranon Locking Plate II, Combi-holes, Left	3 Holes
021421004	Olecranon Locking Plate II, Combi-holes, Left	4 Holes
021421005	Olecranon Locking Plate II, Combi-holes, Left	5 Holes
021421006	Olecranon Locking Plate II, Combi-holes, Left	6 Holes
021421007	Olecranon Locking Plate II, Combi-holes, Left	7 Holes
021421008	Olecranon Locking Plate II, Combi-holes, Left	8 Holes

WJB



Secretaría de
Gobierno de Salud



Ministerio de Salud y Desarrollo Social
Presidencia de la Nación

021422003	Olecranon Locking Plate II, Combi-holes, Right	3 Holes
021422004	Olecranon Locking Plate II, Combi-holes, Right	4 Holes
021422005	Olecranon Locking Plate II, Combi-holes, Right	5 Holes
021422006	Olecranon Locking Plate II, Combi-holes, Right	6 Holes
021422007	Olecranon Locking Plate II, Combi-holes, Right	7 Holes
021422008	Olecranon Locking Plate II, Combi-holes, Right	8 Holes
021431005	Distal Femoral Locking Plate II, Combi-holes, Left	5 Holes
021431006	Distal Femoral Locking Plate II, Combi-holes, Left	6 Holes
021431007	Distal Femoral Locking Plate II, Combi-holes, Left	7 Holes
021431008	Distal Femoral Locking Plate II, Combi-holes, Left	8 Holes
021431009	Distal Femoral Locking Plate II, Combi-holes, Left	9 Holes
021431010	Distal Femoral Locking Plate II, Combi-holes, Left	10 Holes
021431011	Distal Femoral Locking Plate II, Combi-holes, Left	11 Holes
021431012	Distal Femoral Locking Plate II, Combi-holes, Left	12 Holes
021431013	Distal Femoral Locking Plate II, Combi-holes, Left	13 Holes
021432005	Distal Femoral Locking Plate II, Combi-holes, Right	5 Holes
021432006	Distal Femoral Locking Plate II, Combi-holes, Right	6 Holes
021432007	Distal Femoral Locking Plate II, Combi-holes, Right	7 Holes
021432008	Distal Femoral Locking Plate II, Combi-holes, Right	8 Holes

WMB

Sedes y Delegaciones

Tel. (+54-11) 4340-0800 - <http://www.anmat.gov.ar> - República Argentina

Sede Central

Av. de Mayo 869, CABA

Sede Alsina

Alsina 665/671, CABA

Sede INAME

Av. Caseros 2161, CABA

Sede INAL

Estados Unidos 25, CABA

Sede Prod. Médicos

Av. Belgrano 1480, CABA

Deleg. Mendoza

Remedios de Escalada de
San Martín 1909, Mendoza
Prov. de Mendoza

Deleg. Córdoba

Obispo Trejo 635,
Córdoba,
Prov. de Córdoba

Deleg. Paso de los Libres

Ruta Nacional 117, km.10,
CO.TE.CAR., Paso de los
Libres, Prov. de Corrientes

Deleg. Posadas

Roque González 1137,
Posadas, Prov. de
Misiones

Deleg. Santa Fé

Eva Perón 2456,
Santa Fé,
Prov. de Santa Fé

021432009	Distal Femoral Locking Plate II, Combi-holes, Right	9 Holes
021432010	Distal Femoral Locking Plate II, Combi-holes, Right	10 Holes
021432011	Distal Femoral Locking Plate II, Combi-holes, Right	11 Holes
021432012	Distal Femoral Locking Plate II, Combi-holes, Right	12 Holes
021432013	Distal Femoral Locking Plate II, Combi-holes, Right	13 Holes
021441003	Lateral Tibial Head Buttress Locking Plate II, Combi-holes, Left	3 Holes
021441004	Lateral Tibial Head Buttress Locking Plate II, Combi-holes, Left	4 Holes
021441005	Lateral Tibial Head Buttress Locking Plate II, Combi-holes, Left	5 Holes
021441006	Lateral Tibial Head Buttress Locking Plate II, Combi-holes, Left	6 Holes
021441007	Lateral Tibial Head Buttress Locking Plate II, Combi-holes, Left	7 Holes
021441008	Lateral Tibial Head Buttress Locking Plate II, Combi-holes, Left	8 Holes
021441009	Lateral Tibial Head Buttress Locking Plate II, Combi-holes, Left	9 Holes
021441010	Lateral Tibial Head Buttress Locking Plate II, Combi-holes, Left	10 Holes
021441011	Lateral Tibial Head Buttress Locking Plate II, Combi-holes, Left	11 Holes
021441012	Lateral Tibial Head Buttress Locking Plate II, Combi-holes, Left	12 Holes
021441013	Lateral Tibial Head Buttress Locking Plate II, Combi-holes, Left	13 Holes

WZ



Secretaría de
Gobierno de Salud



Ministerio de Salud y Desarrollo Social
Presidencia de la Nación

021442003	Lateral Tibial Head Buttress Locking Plate II, Combi-holes, Right Holes	3
021442004	Lateral Tibial Head Buttress Locking Plate II, Combi-holes, Right Holes	4
021442005	Lateral Tibial Head Buttress Locking Plate II, Combi-holes, Right Holes	5
021442006	Lateral Tibial Head Buttress Locking Plate II, Combi-holes, Right Holes	6
021442007	Lateral Tibial Head Buttress Locking Plate II, Combi-holes, Right Holes	7
021442008	Lateral Tibial Head Buttress Locking Plate II, Combi-holes, Right Holes	8
021442009	Lateral Tibial Head Buttress Locking Plate II, Combi-holes, Right Holes	9
021442010	Lateral Tibial Head Buttress Locking Plate II, Combi-holes, Right 10 Holes	
021442011	Lateral Tibial Head Buttress Locking Plate II, Combi-holes, Right 11 Holes	
021442012	Lateral Tibial Head Buttress Locking Plate II, Combi-holes, Right 12 Holes	

WMB

Sedes y Delegaciones

Tel. (+54-11) 4340-0800 - <http://www.anmat.gov.ar> - República Argentina

Sede Central
Av. de Mayo 869, CABA

Sede Alsina
Alsina 665/671, CABA

Sede INAME
Av. Caseros 2161, CABA

Sede INAL
Estados Unidos 25, CABA

Sede Prod. Médicos
Av. Belgrano 1480, CABA

Deleg. Mendoza
Remedios de Escalada de
San Martín 1909, Mendoza
Prov. de Mendoza

Deleg. Córdoba
Obispo Trejo 635,
Córdoba,
Prov. de Córdoba

Deleg. Paso de los Libres
Ruta Nacional 117, km.10,
CO.TE.CAR., Paso de los
Libres, Prov. de Corrientes

Deleg. Posadas
Roque González 1137,
Posadas, Prov. de
Misiones

Deleg. Santa Fé
Eva Perón 2456,
Santa Fé,
Prov. de Santa Fé

021442013 Lateral Tibial Head Buttress Locking Plate II, Combi-holes, Right

13 Holes

021450006 1.5S-Mini Wed Locking Plate 34.4mm

021460008 2.0S-Mini Wed Locking Plate 47.7mm

021470010 2.4S-Mini Wed Locking Plate 77.2mm

021501004 Low Profile Distal medial Tibial Locking Plate, Minimally Invasive

Technique left 4 Holes

021501006 Low Profile Distal medial Tibial Locking Plate, Minimally Invasive

Technique left 6 Holes

021501008 Low Profile Distal medial Tibial Locking Plate, Minimally Invasive

Technique left 8 Holes

021501010 Low Profile Distal medial Tibial Locking Plate, Minimally Invasive

Technique left 10 Holes

021501012 Low Profile Distal medial Tibial Locking Plate, Minimally Invasive

Technique left 12 Holes

021501014 Low Profile Distal medial Tibial Locking Plate, Minimally Invasive

Technique left 14 Holes

021502004 Low Profile Distal medial Tibial Locking Plate, Minimally Invasive

Technique right 4 Holes

021502006 Low Profile Distal medial Tibial Locking Plate, Minimally Invasive

Technique right 6 Holes

UMB



Secretaría de
Gobierno de Salud



Ministerio de Salud y Desarrollo Social
Presidencia de la Nación

- 021502008 Low Profile Distal medial Tibial Locking Plate, Minimally Invasive
Technique right 8 Holes
- 021502010 Low Profile Distal medial Tibial Locking Plate, Minimally Invasive
Technique right 10 Holes
- 021502012 Low Profile Distal medial Tibial Locking Plate, Minimally Invasive
Technique right 12 Holes
- 021502014 Low Profile Distal medial Tibial Locking Plate, Minimally Invasive
Technique right 14 Holes
- 021510008 Curved Locking Compression Plate, broad 8 Holes
- 021510009 Curved Locking Compression Plate, broad 9 Holes
- 021510010 Curved Locking Compression Plate, broad 10 Holes
- 021510011 Curved Locking Compression Plate, broad 11 Holes
- 021510012 Curved Locking Compression Plate, broad 12 Holes
- 021510013 Curved Locking Compression Plate, broad 13 Holes
- 021510014 Curved Locking Compression Plate, broad 14 Holes
- 021510015 Curved Locking Compression Plate, broad 15 Holes
- 021510016 Curved Locking Compression Plate, broad 16 Holes
- 021510017 Curved Locking Compression Plate, broad 17 Holes
- 021510018 Curved Locking Compression Plate, broad 18 Holes
- 021520002 Posterior Medial Proximal Tibial Plate, combi-holes 2 Holes
- 021520004 Posterior Medial Proximal Tibial Plate, combi-holes 4 Holes

UMB

Sedes y Delegaciones

Tel. (+54-11) 4340-0800 - <http://www.anmat.gov.ar> - República Argentina

Sede Central
Av. de Mayo 869, CABA

Sede Alsina
Alsina 665/671, CABA

Sede INAME
Av. Caseros 2161, CABA

Sede INAL
Estados Unidos 25, CABA

Sede Prod. Médicos
Av. Belgrano 1480, CABA

Deleg. Mendoza
Remedios de Escalada de
San Martín 1909, Mendoza
Prov. de Mendoza

Deleg. Córdoba
Obispo Trejo 635,
Córdoba,
Prov. de Córdoba

Deleg. Paso de los Libres
Ruta Nacional 117, km.10,
CO.TE.CAR., Paso de los
Libres, Prov. de Corrientes

Deleg. Posadas
Roque González 1137,
Posadas, Prov. de
Misiones

Deleg. Santa Fé
Eva Perón 2456,
Santa Fé,
Prov. de Santa Fé

021520006	Posterior Medial Proximal Tibial Plate, combi-holes	6 Holes
021520008	Posterior Medial Proximal Tibial Plate, combi-holes	8 Holes
021531002	Distal Radius Locking Plate I , variable angle, small, left	2 Holes
021531003	Distal Radius Locking Plate I , variable angle, small, left	3 Holes
021531004	Distal Radius Locking Plate I , variable angle, small, left	4 Holes
021531005	Distal Radius Locking Plate I , variable angle, small, left	5 Holes
021532002	Distal Radius Locking Plate I , variable angle, small, right	2 Holes
021532003	Distal Radius Locking Plate I , variable angle, small, right	3 Holes
021532004	Distal Radius Locking Plate I , variable angle, small, right	4 Holes
021532005	Distal Radius Locking Plate I , variable angle, small, right	5 Holes
021541002	Distal Radius Locking Plate I , variable angle, large, left	2 Holes
021541003	Distal Radius Locking Plate I , variable angle, large, left	3 Holes
021541004	Distal Radius Locking Plate I , variable angle, large, left	4 Holes
021541005	Distal Radius Locking Plate I , variable angle, large, left	5 Holes
021542002	Distal Radius Locking Plate I , variable angle, large, right	2 Holes
021542003	Distal Radius Locking Plate I , variable angle, large, right	3 Holes

WAB



Secretaría de
Gobierno de Salud



Ministerio de Salud y Desarrollo Social
Presidencia de la Nación

- 021542004 Distal Radius Locking Plate I , variable angle, large, right 4 Holes
- 021542005 Distal Radius Locking Plate I , variable angle, large, right 5 Holes
- 021550006 Cable Locking Plate 6 Holes
- 021550008 Cable Locking Plate 8 Holes
- 021550010 Cable Locking Plate 10 Holes
- 021561002 Cable Grip Locking Plate, Left 2 Holes
- 021561004 Cable Grip Locking Plate, Left 4 Holes
- 021561006 Cable Grip Locking Plate, Left 6 Holes
- 021561008 Cable Grip Locking Plate, Left 8 Holes
- 021561010 Cable Grip Locking Plate, Left 10 Holes
- 021561012 Cable Grip Locking Plate, Left 12 Holes
- 021562002 Cable Grip Locking Plate, Right 2 Holes
- 021562004 Cable Grip Locking Plate, Right 4 Holes
- 021562006 Cable Grip Locking Plate, Right 6 Holes
- 021562008 Cable Grip Locking Plate, Right 8 Holes
- 021562010 Cable Grip Locking Plate, Right 10 Holes
- 021562012 Cable Grip Locking Plate, Right 12 Holes
- 021601003 Superior Anterior Distal Clavicle Locking Plate, with lateral extension, combi-holes, left 3 Holes

WAS

Sedes y Delegaciones

Tel. (+54-11) 4340-0800 - <http://www.anmat.gov.ar> - República Argentina

Sede Central

Av. de Mayo 869, CABA

Sede Alsina

Alsina 665/671, CABA

Sede INAME

Av. Caseros 2161, CABA

Sede INAL

Estados Unidos 25, CABA

Sede Prod. Médicos

Av. Belgrano 1480, CABA

Deleg. Mendoza

Remedios de Escalada de
San Martín 1909, Mendoza
Prov. de Mendoza

Deleg. Córdoba

Obispo Trejo 635,
Córdoba,
Prov. de Córdoba

Deleg. Paso de los Libres

Ruta Nacional 117, km.10,
CO.TE.CAR., Paso de los
Libres, Prov. de Corrientes

Deleg. Posadas

Roque González 1137,
Posadas, Prov. de
Misiones

Deleg. Santa Fé

Eva Perón 2456,
Santa Fé,
Prov. de Santa Fé

021601004 Superior Anterior Distal Clavicle Locking Plate, with lateral extension,
combi-holes, left 4 Holes

021601005 Superior Anterior Distal Clavicle Locking Plate, with lateral extension,
combi-holes, left 5 Holes

021601006 Superior Anterior Distal Clavicle Locking Plate, with lateral extension,
combi-holes, left 6 Holes

021602003 Superior Anterior Distal Clavicle Locking Plate, with lateral extension,
combi-holes, right 3 Holes

021602004 Superior Anterior Distal Clavicle Locking Plate, with lateral extension,
combi-holes, right 4 Holes

021602005 Superior Anterior Distal Clavicle Locking Plate, with lateral extension,
combi-holes, right 5 Holes

021602006 Superior Anterior Distal Clavicle Locking Plate, with lateral extension,
combi-holes, right 6 Holes

021611004 Distal Medial Radius Locking Plate, left 4 Holes

021612004 Distal Medial Radius Locking Plate, right 4 Holes

021620005 Distal Lateral Radius Locking Plate 5 Holes

021620006 Distal Lateral Radius Locking Plate 6 Holes

021631003 Proximal Femoral Locking Plate III, left 3 Holes

021631004 Proximal Femoral Locking Plate III, Left 4 Holes

LMB



Secretaría de
Gobierno de Salud



Ministerio de Salud y Desarrollo Social
Presidencia de la Nación

021631005	Proximal Femoral Locking Plate III, Left	5 Holes
021631006	Proximal Femoral Locking Plate III, Left	6 Holes
021631007	Proximal Femoral Locking Plate III, Left	7 Holes
021631008	Proximal Femoral Locking Plate III, Left	8 Holes
021631009	Proximal Femoral Locking Plate III, Left	9 Holes
021631010	Proximal Femoral Locking Plate III, Left	10 Holes
021631011	Proximal Femoral Locking Plate III, Left	11 Holes
021631012	Proximal Femoral Locking Plate III, Left	12 Holes
021631013	Proximal Femoral Locking Plate III, Left	13 Holes
021632003	Proximal Femoral Locking Plate III, Right	3 Holes
021632004	Proximal Femoral Locking Plate III, Right	4 Holes
021632005	Proximal Femoral Locking Plate III, Right	5 Holes
021632006	Proximal Femoral Locking Plate III, Right	6 Holes
021632007	Proximal Femoral Locking Plate III, Right	7 Holes
021632008	Proximal Femoral Locking Plate III, Right	8 Holes
021632009	Proximal Femoral Locking Plate III, Right	9 Holes
021632010	Proximal Femoral Locking Plate III, Right	10 Holes
021632011	Proximal Femoral Locking Plate III, Right	11 Holes
021632012	Proximal Femoral Locking Plate III, Right	12 Holes

W/B

Sedes y Delegaciones

Tel. (+54-11) 4340-0800 - <http://www.anmat.gov.ar> - República Argentina

Sede Central
Av. de Mayo 869, CABA

Sede Alsina
Alsina 665/671, CABA

Sede INAME
Av. Caseros 2161, CABA

Sede INAL
Estados Unidos 25, CABA

Sede Prod. Médicos
Av. Belgrano 1480, CABA

Deleg. Mendoza
Remedios de Escalada de
San Martín 1909, Mendoza
Prov. de Mendoza

Deleg. Córdoba
Obispo Trejo 635,
Córdoba,
Prov. de Córdoba

Deleg. Paso de los Libres
Ruta Nacional 117, km.10,
CO.TE.CAR., Paso de los
Libres, Prov. de Corrientes

Deleg. Posadas
Roque González 1137,
Posadas, Prov. de
Misiones

Deleg. Santa Fé
Eva Perón 2456,
Santa Fé,
Prov. de Santa Fé

021632013 Proximal Femoral Locking Plate III, Right 13 Holes
021641004 Ankle Fusion Locking Plate, left 4 Holes
021642004 Ankle Fusion Locking Plate, right 4 Holes
021671005 Proximal Lateral Tibial Locking Plate II left 5 Holes
021671009 Proximal Lateral Tibial Locking Plate II left 9 Holes
021671013 Proximal Lateral Tibial Locking Plate II left 13 Holes
021671017 Proximal Lateral Tibial Locking Plate II left 17 Holes
021672005 Proximal Lateral Tibial Locking Plate II right 5 Holes
021672009 Proximal Lateral Tibial Locking Plate II right 9 Holes
021672013 Proximal Lateral Tibial Locking Plate II right 13 Holes
021672017 Proximal Lateral Tibial Locking Plate II right 17 Holes
021681004 Proximal Medial Tibial Locking Plate II left 4 Holes
021681006 Proximal Medial Tibial Locking Plate II left 6 Holes
021681008 Proximal Medial Tibial Locking Plate II left 8 Holes
021682004 Proximal Medial Tibial Locking Plate II right 4 Holes
021682006 Proximal Medial Tibial Locking Plate II right 6 Holes
021682008 Proximal Medial Tibial Locking Plate II right 8 Holes
021691003 Proximal Posterior Tibial Locking Plate left, II 3 Holes
021691005 Proximal Posterior Tibial Locking Plate left, II 5 Holes
021691007 Proximal Posterior Tibial Locking Plate left, II 7 Holes
021692003 Proximal Posterior Tibial Locking Plate right, II 3 Holes

LMB



021692005	Proximal Posterior Tibial Locking Plate right, II	5 Holes
021692007	Proximal Posterior Tibial Locking Plate right, II	7 Holes
021800004	Humeral Shaft Straight Locking Plate	4 Holes
021800005	Humeral Shaft Straight Locking Plate	5 Holes
021800006	Humeral Shaft Straight Locking Plate	6 Holes
021800007	Humeral Shaft Straight Locking Plate	7 Holes
021800008	Humeral Shaft Straight Locking Plate	8 Holes
021800009	Humeral Shaft Straight Locking Plate	9 Holes
021800010	Humeral Shaft Straight Locking Plate	10 Holes
021800011	Humeral Shaft Straight Locking Plate	11 Holes
021800012	Humeral Shaft Straight Locking Plate	12 Holes
021800013	Humeral Shaft Straight Locking Plate	13 Holes
021800014	Humeral Shaft Straight Locking Plate	14 Holes
021800015	Humeral Shaft Straight Locking Plate	15 Holes
021800016	Humeral Shaft Straight Locking Plate	16 Holes
021810127	Sacroiliac Joint Plate II	127mm
021810137	Sacroiliac Joint Plate II	137mm
021810147	Sacroiliac Joint Plate II	147mm
021821003	Proximal Lateral Tibial Locking Plate, left	3 Holes
021821004	Proximal Lateral Tibial Locking Plate, left	4 Holes

UMB

Sedes y Delegaciones

Tel. (+54-11) 4340-0800 - <http://www.anmat.gov.ar> - República Argentina

Sede Central
 Av. de Mayo 869, CABA

Sede Alsina
 Alsina 665/671, CABA

Sede INAME
 Av. Caseros 2161, CABA

Sede INAL
 Estados Unidos 25, CABA

Sede Prod. Médicos
 Av. Belgrano 1480, CABA

Deleg. Mendoza
 Remedios de Escalada de
 San Martín 1909, Mendoza
 Prov. de Mendoza

Deleg. Córdoba
 Obispo Trejo 635,
 Córdoba,
 Prov. de Córdoba

Deleg. Paso de los Libres
 Ruta Nacional 117, km.10,
 CO.TE.CAR., Paso de los
 Libres, Prov. de Corrientes

Deleg. Posadas
 Roque González 1137,
 Posadas, Prov. de
 Misiones

Deleg. Santa Fé
 Eva Perón 2456,
 Santa Fé,
 Prov. de Santa Fé

021821005	Proximal Lateral Tibial Locking Plate, left	5 Holes
021821006	Proximal Lateral Tibial Locking Plate, left	6 Holes
021821007	Proximal Lateral Tibial Locking Plate, left	7 Holes
021821008	Proximal Lateral Tibial Locking Plate, left	8 Holes
021821009	Proximal Lateral Tibial Locking Plate, left	9 Holes
021821010	Proximal Lateral Tibial Locking Plate, left	10 Holes
021821011	Proximal Lateral Tibial Locking Plate, left	11 Holes
021821012	Proximal Lateral Tibial Locking Plate, left	12 Holes
021821013	Proximal Lateral Tibial Locking Plate, left	13 Holes
021821014	Proximal Lateral Tibial Locking Plate, left	14 Holes
021821015	Proximal Lateral Tibial Locking Plate, left	15 Holes
021822003	Proximal Lateral Tibial Locking Plate, right	3 Holes
021822004	Proximal Lateral Tibial Locking Plate, right	4 Holes
021822005	Proximal Lateral Tibial Locking Plate, right	5 Holes
021822006	Proximal Lateral Tibial Locking Plate, right	6 Holes
021822007	Proximal Lateral Tibial Locking Plate, right	7 Holes
021822008	Proximal Lateral Tibial Locking Plate, right	8 Holes
021822009	Proximal Lateral Tibial Locking Plate, right	9 Holes
021822010	Proximal Lateral Tibial Locking Plate, right	10 Holes
021822011	Proximal Lateral Tibial Locking Plate, right	11 Holes
021822012	Proximal Lateral Tibial Locking Plate, right	12 Holes

WFB



Secretaría de
Gobierno de Salud



Ministerio de Salud y Desarrollo Social
Presidencia de la Nación

- 021822013 Proximal Lateral Tibial Locking Plate, right 13 Holes
- 021822014 Proximal Lateral Tibial Locking Plate, right 14 Holes
- 021822015 Proximal Lateral Tibial Locking Plate, right 15 Holes
- 021881003 Posterior Lateral Fibular Locking Plate combi holes, left 3 Holes
- 021881004 Posterior Lateral Fibular Locking Plate combi holes, left 4 Holes
- 021881005 Posterior Lateral Fibular Locking Plate combi holes, left 5 Holes
- 021881006 Posterior Lateral Fibular Locking Plate combi holes, left 6 Holes
- 021881007 Posterior Lateral Fibular Locking Plate combi holes, left 7 Holes
- 021881008 Posterior Lateral Fibular Locking Plate combi holes, left 8 Holes
- 021882003 Posterior Lateral Fibular Locking Plate combi holes, right 3 Holes
- 021882004 Posterior Lateral Fibular Locking Plate combi holes, right 4 Holes
- 021882005 Posterior Lateral Fibular Locking Plate combi holes, right 5 Holes
- 021882006 Posterior Lateral Fibular Locking Plate combi holes, right 6 Holes
- 021882007 Posterior Lateral Fibular Locking Plate combi holes, right 7 Holes
- 021882008 Posterior Lateral Fibular Locking Plate combi holes, right 8 Holes
- 021891003 Distal Lateral Humeral Locking Plate, head 2.7 hole, combi-holes, Left 3 Holes
- 021891004 Distal Lateral Humeral Locking Plate, head 2.7 hole, combi-holes, Left 4 Holes
- WMS* 021891005 Distal Lateral Humeral Locking Plate, head 2.7 hole, combi-holes, Left 5 Holes

Sedes y Delegaciones

Tel. (+54-11) 4340-0800 - <http://www.anmat.gov.ar> - República Argentina

Sede Central
Av. de Mayo 869, CABA

Sede Alsina
Alsina 665/671, CABA

Sede INAME
Av. Caseros 2161, CABA

Sede INAL
Estados Unidos 25, CABA

Sede Prod. Médicos
Av. Belgrano 1480, CABA

Deleg. Mendoza
Remedios de Escalada de
San Martín 1909, Mendoza
Prov. de Mendoza

Deleg. Córdoba
Obispo Trejo 635,
Córdoba,
Prov. de Córdoba

Deleg. Paso de los Libres
Ruta Nacional 117, km.10,
CO.TE.CAR., Paso de los
Libres, Prov. de Corrientes

Deleg. Posadas
Roque González 1137,
Posadas, Prov. de
Misiones

Deleg. Santa Fé
Eva Perón 2456,
Santa Fé,
Prov. de Santa Fé

021891006 Distal Lateral Humeral Locking Plate, head 2.7 hole, combi-holes, Left 6
Holes

021891007 Distal Lateral Humeral Locking Plate, head 2.7 hole, combi-holes, Left 7
Holes

021891008 Distal Lateral Humeral Locking Plate, head 2.7 hole, combi-holes, Left 8
Holes

021891009 Distal Lateral Humeral Locking Plate, head 2.7 hole, combi-holes, Left 9
Holes

021892003 Distal Lateral Humeral Locking Plate, head 2.7 hole, combi-holes, Right
3 Holes

021892004 Distal Lateral Humeral Locking Plate, head 2.7 hole, combi-holes, Right
4 Holes

021892005 Distal Lateral Humeral Locking Plate, head 2.7 hole, combi-holes, Right
5 Holes

021892006 Distal Lateral Humeral Locking Plate, head 2.7 hole, combi-holes, Right
6 Holes

021892007 Distal Lateral Humeral Locking Plate, head 2.7 hole, combi-holes, Right
7 Holes

021892008 Distal Lateral Humeral Locking Plate, head 2.7 hole, combi-holes, Right
8 Holes

WLB



021892009 Distal Lateral Humeral Locking Plate, head 2.7 hole, combi-holes, Right
9 Holes

021901003 Distal Medial Humeral Locking Plate, head 2.7 hole, combi-holes, Left 3
Holes

021901004 Distal Medial Humeral Locking Plate, head 2.7 hole, combi-holes, Left 4
Holes

021901005 Distal Medial Humeral Locking Plate, head 2.7 hole, combi-holes, Left 5
Holes

021901006 Distal Medial Humeral Locking Plate, head 2.7 hole, combi-holes, Left 6
Holes

021901007 Distal Medial Humeral Locking Plate, head 2.7 hole, combi-holes, Left 7
Holes

021901008 Distal Medial Humeral Locking Plate, head 2.7 hole, combi-holes, Left 8
Holes

021901009 Distal Medial Humeral Locking Plate, head 2.7 hole, combi-holes, Left 9
Holes

021902003 Distal Medial Humeral Locking Plate, head 2.7 hole, combi-holes, Right
3 Holes

021902004 Distal Medial Humeral Locking Plate, head 2.7 hole, combi-holes, Right
4 Holes

ANMAT

Sedes y Delegaciones

Tel. (+54-11) 4340-0800 - <http://www.anmat.gov.ar> - República Argentina

Sede Central
Av. de Mayo 869, CABA

Sede Alsina
Alsina 665/671, CABA

Sede INAME
Av. Caseros 2161, CABA

Sede INAL
Estados Unidos 25, CABA

Sede Prod. Médicos
Av. Belgrano 1480, CABA

Deleg. Mendoza
Remedios de Escalada de
San Martín 1909, Mendoza
Prov. de Mendoza

Deleg. Córdoba
Obispo Trejo 635,
Córdoba,
Prov. de Córdoba

Deleg. Paso de los Libres
Ruta Nacional 117, km.10,
CO.TE.CAR., Paso de los
Libres, Prov. de Corrientes

Deleg. Posadas
Roque González 1137,
Posadas, Prov. de
Misiones

Deleg. Santa Fé
Eva Perón 2456,
Santa Fé,
Prov. de Santa Fé

021902005 Distal Medial Humeral Locking Plate, head 2.7 hole, combi-holes, Right
5 Holes

021902006 Distal Medial Humeral Locking Plate, head 2.7 hole, combi-holes, Right
6 Holes

021902007 Distal Medial Humeral Locking Plate, head 2.7 hole, combi-holes, Right
7 Holes

021902008 Distal Medial Humeral Locking Plate, head 2.7 hole, combi-holes, Right
8 Holes

021902009 Distal Medial Humeral Locking Plate, head 2.7 hole, combi-holes, Right
9 Holes

022100008 Wrist Fusion Locking Plate, length 117mm, standard, combi-holes 8
Holes

022110008 Wrist Fusion Locking Plate, length 117mm, short, combi-holes 8 Holes

022120009 Wrist Fusion Locking Plate, length 112mm, straight, combi-holes 9
Holes

022181002 Distal Radius Locking Plate, head 7 holes variable angle, left 2 Holes

022181003 Distal Radius Locking Plate, head 7 holes variable angle, left 3 Holes

022181004 Distal Radius Locking Plate, head 7 holes variable angle, left 4 Holes

022181005 Distal Radius Locking Plate, head 7 holes variable angle, left 5 Holes

022182002 Distal Radius Locking Plate, head 7 holes variable angle, right 2 Holes

022182003 Distal Radius Locking Plate, head 7 holes variable angle, right 3 Holes

W3



Secretaría de
Gobierno de Salud



Ministerio de Salud y Desarrollo Social
Presidencia de la Nación

- 022182004 Distal Radius Locking Plate, head 7 holes variable angle, right 4 Holes
- 022182005 Distal Radius Locking Plate, head 7 holes variable angle, right 5 Holes
- 022360004 Proximal Medial Tibial Osteotomy Locking Plate, combi-holes 4 Holes
- 022981003 Distal Lateral Femoral Osteotomy Locking Plate, combi-holes, left 3 Holes
- 022981004 Distal Lateral Femoral Osteotomy Locking Plate, combi-holes, left 4 Holes
- 022981005 Distal Lateral Femoral Osteotomy Locking Plate, combi-holes, left 5 Holes
- 022981006 Distal Lateral Femoral Osteotomy Locking Plate, combi-holes, left 6 Holes
- 022981007 Distal Lateral Femoral Osteotomy Locking Plate, combi-holes, left 7 Holes
- 022981008 Distal Lateral Femoral Osteotomy Locking Plate, combi-holes, left 8 Holes
- 022981009 Distal Lateral Femoral Osteotomy Locking Plate, combi-holes, left 9 Holes
- 022981010 Distal Lateral Femoral Osteotomy Locking Plate, combi-holes, left 10 Holes
- 022981011 Distal Lateral Femoral Osteotomy Locking Plate, combi-holes, left 11 Holes

Sedes y Delegaciones

Tel. (+54-11) 4340-0800 - <http://www.anmat.gov.ar> - República Argentina

Sede Central
Av. de Mayo 869, CABA

Sede Alsina
Alsina 665/671, CABA

Sede INAME
Av. Caseros 2161, CABA

Sede INAL
Estados Unidos 25, CABA

Sede Prod. Médicos
Av. Belgrano 1480, CABA

Deleg. Mendoza
Remedios de Escalada de
San Martín 1909, Mendoza
Prov. de Mendoza

Deleg. Córdoba
Obispo Trejo 635,
Córdoba,
Prov. de Córdoba

Deleg. Paso de los Libres
Ruta Nacional 117, km.10,
CO.TE.CAR., Paso de los
Libres, Prov. de Corrientes

Deleg. Posadas
Roque González 1137,
Posadas, Prov. de
Misiones

Deleg. Santa Fé
Eva Perón 2456,
Santa Fé,
Prov. de Santa Fé

022981012 Distal Lateral Femoral Osteotomy Locking Plate, combi-holes,left
12 Holes

022981013 Distal Lateral Femoral Osteotomy Locking Plate, combi-holes,left
13 Holes

022982003 Distal Lateral Femoral Osteotomy Locking Plate, combi-holes,right 3
Holes

022982004 Distal Lateral Femoral Osteotomy Locking Plate, combi-holes,right 4
Holes

022982005 Distal Lateral Femoral Osteotomy Locking Plate, combi-holes,right 5
Holes

022982006 Distal Lateral Femoral Osteotomy Locking Plate, combi-holes,right 6
Holes

022982007 Distal Lateral Femoral Osteotomy Locking Plate, combi-holes,right 7
Holes

022982008 Distal Lateral Femoral Osteotomy Locking Plate, combi-holes,right 8
Holes

022982009 Distal Lateral Femoral Osteotomy Locking Plate, combi-holes,right 9
Holes

022982010 Distal Lateral Femoral Osteotomy Locking Plate, combi-holes,right
10 Holes

WLB



Secretaría de
Gobierno de Salud



Ministerio de Salud y Desarrollo Social
Presidencia de la Nación

- 022982011 Distal Lateral Femoral Osteotomy Locking Plate, combi-holes,right
11 Holes
- 022982012 Distal Lateral Femoral Osteotomy Locking Plate, combi-holes,right
12 Holes
- 022982013 Distal Lateral Femoral Osteotomy Locking Plate, combi-holes,right
13 Holes
- 022991003 Proximal Lateral Tibial Osteotomy Locking Plate,left 3 Holes
- 022991005 Proximal Lateral Tibial Osteotomy Locking Plate,left 5 Holes
- 022991007 Proximal Lateral Tibial Osteotomy Locking Plate,left 7 Holes
- 022991009 Proximal Lateral Tibial Osteotomy Locking Plate,left 9 Holes
- 022991011 Proximal Lateral Tibial Osteotomy Locking Plate,left 11 Holes
- 022991013 Proximal Lateral Tibial Osteotomy Locking Plate,left 13 Holes
- 022992003 Proximal Lateral Tibial Osteotomy Locking Plate,right 3 Holes
- 022992005 Proximal Lateral Tibial Osteotomy Locking Plate,right 5 Holes
- 022992007 Proximal Lateral Tibial Osteotomy Locking Plate,right 7 Holes
- 022992009 Proximal Lateral Tibial Osteotomy Locking Plate,right 9 Holes
- 022992011 Proximal Lateral Tibial Osteotomy Locking Plate,right 11 Holes
- 022992013 Proximal Lateral Tibial Osteotomy Locking Plate,right 13 Holes
- 025331000 Coronoid Locking Plate, Left left
- 025332000 Coronoid Locking Plate, Right right
- 025341003 Distal Lateral Tibial Locking Plate II, combi-holes,left 3 Holes

WLS

Sedes y Delegaciones

Tel. (+54-11) 4340-0800 - <http://www.anmat.gov.ar> - República Argentina

Sede Central
Av. de Mayo 869, CABA

Sede Alsina
Alsina 665/671, CABA

Sede INAME
Av. Caseros 2161, CABA

Sede INAL
Estados Unidos 25, CABA

Sede Prod. Médicos
Av. Belgrano 1480, CABA

Deleg. Mendoza
Remedios de Escalada de
San Martín 1909, Mendoza
Prov. de Mendoza

Deleg. Córdoba
Obispo Trejo 635,
Córdoba,
Prov. de Córdoba

Deleg. Paso de los Libres
Ruta Nacional 117, km.10,
CO.TE.CAR., Paso de los
Libres, Prov. de Corrientes

Deleg. Posadas
Roque González 1137,
Posadas, Prov. de
Misiones

Deleg. Santa Fé
Eva Perón 2456,
Santa Fé,
Prov. de Santa Fé

025341005	Distal Lateral Tibial Locking Plate II , combi-holes,Left	5Holes
025341007	Distal Lateral Tibial Locking Plate II , combi-holes,left	7 Holes
025341009	Distal Lateral Tibial Locking Plate II , combi-holes,Left	9Holes
025341011	Distal Lateral Tibial Locking Plate II , combi-holes,left	11 Holes
025341013	Distal Lateral Tibial Locking Plate II , combi-holes,Left	13Holes
025341015	Distal Lateral Tibial Locking Plate II , combi-holes,left	15 Holes
025341017	Distal Lateral Tibial Locking Plate II , combi-holes,Left	17Holes
025342003	Distal Lateral Tibial Locking Plate II , combi-holes,right	3 Holes
025342005	Distal Lateral Tibial Locking Plate II , combi-holes,Right	5Holes
025342007	Distal Lateral Tibial Locking Plate II , combi-holes,right	7 Holes
025342009	Distal Lateral Tibial Locking Plate II , combi-holes,Right	9Holes
025342011	Distal Lateral Tibial Locking Plate II , combi-holes,right	11 Holes
025342013	Distal Lateral Tibial Locking Plate II , combi-holes,Right	13Holes
025342015	Distal Lateral Tibial Locking Plate II , combi-holes,right	15 Holes
025342017	Distal Lateral Tibial Locking Plate II , combi-holes,Right	17Holes
025351008	Distal Diaphyseal Metacarpal Radius Locking Plate, left	8 Holes
025351010	Distal Diaphyseal Metacarpal Radius Locking Plate, left	10 Holes
025351012	Distal Diaphyseal Metacarpal Radius Locking Plate, left	12 Holes
025352008	Distal Diaphyseal Metacarpal Radius Locking Plate, right	8 Holes
025352010	Distal Diaphyseal Metacarpal Radius Locking Plate, right	10 Holes

WJB



- 025352012 Distal Diaphyseal Metacarpal Radius Locking Plate, right 12 Holes
- 025551005 Proximal Humeral Greater Tuberosity Locking Compression Plate,
suture hole, left 5 Holes
- 025552005 Proximal Humeral Greater Tuberosity Locking Compression Plate,
suture hole, right 5 Holes
- 025681006 Superior Anterior Clavicle Locking Plate II, combi-holes, left 6 Holes
- 025681007 Superior Anterior Clavicle Locking Plate II, combi-holes, left 7 Holes
- 025681008 Superior Anterior Clavicle Locking Plate II, combi-holes, left 8 Holes
- 025681009 Superior Anterior Clavicle Locking Plate II, combi-holes, left 9 Holes
- 025681010 Superior Anterior Clavicle Locking Plate II, combi-holes, left
Holes 10
- 025681011 Superior Anterior Clavicle Locking Plate II, combi-holes, left
Holes 11
- 025681012 Superior Anterior Clavicle Locking Plate II, combi-holes, left
Holes 12
- 025682006 Superior Anterior Clavicle Locking Plate II, combi-holes, right 6 Holes
- 025682007 Superior Anterior Clavicle Locking Plate II, combi-holes, right 7 Holes
- 025682008 Superior Anterior Clavicle Locking Plate II, combi-holes, right 8 Holes
- 025682009 Superior Anterior Clavicle Locking Plate II, combi-holes, right 9 Holes

LMZ

Sedes y Delegaciones

Tel. (+54-11) 4340-0800 - <http://www.anmat.gov.ar> - República Argentina

Sede Central
Av. de Mayo 869, CABA

Sede Alsina
Alsina 665/671, CABA

Sede INAME
Av. Caseros 2161, CABA

Sede INAL
Estados Unidos 25, CABA

Sede Prod. Médicos
Av. Belgrano 1480, CABA

Deleg. Mendoza
Remedios de Escalada de
San Martín 1909, Mendoza
Prov. de Mendoza

Deleg. Córdoba
Obispo Trejo 635,
Córdoba,
Prov. de Córdoba

Deleg. Paso de los Libres
Ruta Nacional 117, km.10,
CO.TE.CAR., Paso de los
Libres, Prov. de Corrientes

Deleg. Posadas
Roque González 1137,
Posadas, Prov. de
Misiones

Deleg. Santa Fé
Eva Perón 2456,
Santa Fé,
Prov. de Santa Fé

025682010 Superior Anterior Clavicle Locking Plate II, combi-holes, right 10

Holes

025682011 Superior Anterior Clavicle Locking Plate II, combi-holes, right 11

Holes

025682012 Superior Anterior Clavicle Locking Plate II, combi-holes, right 12

Holes

025850006 Metaphyseal Locking Plate, 3,5 6 Holes

025850007 Metaphyseal Locking Plate, 3,5 7 Holes

025850008 Metaphyseal Locking Plate, 3,5 8 Holes

025850009 Metaphyseal Locking Plate, 3,5 9 Holes

025850010 Metaphyseal Locking Plate, 3,5 10 Holes

025850011 Metaphyseal Locking Plate, 3,5 11 Holes

025850012 Metaphyseal Locking Plate, 3,5 12 Holes

025850014 Metaphyseal Locking Plate, 3,5 14 Holes

025870004 Distal Medial Humeral Reconstruction Locking Plate 4 Holes

025870006 Distal Medial Humeral Reconstruction Locking Plate 6 Holes

025870008 Distal Medial Humeral Reconstruction Locking Plate 8 Holes

025870010 Distal Medial Humeral Reconstruction Locking Plate 10 Holes

025870012 Distal Medial Humeral Reconstruction Locking Plate 12 Holes

026012008 Olecranon Locking Plate, curved, riht 8 holes

WRB



Secretaría de
Gobierno de Salud



Ministerio de Salud y Desarrollo Social
Presidencia de la Nación

- 026020007 Olecranon Locking Plate, straight 7 holes
- 026110004 MD Distal Tibial T-Locking Plate 4 holes
- 026121004 MD Distal Tibial L-Locking Plate, left 4 holes
- 026122004 MD Distal Tibial L-Locking Plate, right 4 holes
- 026131006 Distal Diaphyseal Radius Locking Plate II, left 6 Holes
- 026131008 Distal Diaphyseal Radius Locking Plate II, left 8 Holes
- 026131010 Distal Diaphyseal Radius Locking Plate II, left 10 Holes
- 026131012 Distal Diaphyseal Radius Locking Plate II, left 12 Holes
- 026131014 Distal Diaphyseal Radius Locking Plate II, left 14 Holes
- 026132006 Distal Diaphyseal Radius Locking Plate II, right 6 Holes
- 026132008 Distal Diaphyseal Radius Locking Plate II, right 8 Holes
- 026132010 Distal Diaphyseal Radius Locking Plate II, right 10 Holes
- 026132012 Distal Diaphyseal Radius Locking Plate II, right 12 Holes
- 026132014 Distal Diaphyseal Radius Locking Plate II, right 14 Holes
- 026141000 Acromion Locking Plate, left
- 026142000 Acromion Locking Plate, right
- 026151000 Scapula Medial Locking Plate, left
- 026152000 Scapula Medial Locking Plate, right
- 026161000 Scapula Lateral Locking Plate, left
- 026162000 Scapula Lateral Locking Plate, right
- 026331000 Glenoid Locking Plate, left

WLB

Sedes y Delegaciones

Tel. (+54-11) 4340-0800 - <http://www.anmat.gov.ar> - República Argentina

Sede Central
Av. de Mayo 869, CABA

Sede Alsina
Alsina 665/671, CABA

Sede INAME
Av. Caseros 2161, CABA

Sede INAL
Estados Unidos 25, CABA

Sede Prod. Médicos
Av. Belgrano 1480, CABA

Deleg. Mendoza
Remedios de Escalada de
San Martín 1909, Mendoza
Prov. de Mendoza

Deleg. Córdoba
Obispo Trejo 635,
Córdoba,
Prov. de Córdoba

Deleg. Paso de los Libres
Ruta Nacional 117, km.10,
CO.TE.CAR., Paso de los
Libres, Prov. de Corrientes

Deleg. Posadas
Roque González 1137,
Posadas, Prov. de
Misiones

Deleg. Santa Fé
Eva Perón 2456,
Santa Fé,
Prov. de Santa Fé

026332000 Glenoid Locking Plate, right
026611012 Olecranon U-locking Plate, left 12 Holes
026612012 Olecranon U-locking Plate, right 12 Holes
026620006 MT Calcaneal Step Locking Plate
026620008 MT Calcaneal Step Locking Plate
026620010 MT Calcaneal Step Locking Plate
026651000 MT First Metatarsal Osteotomy Locking Plate, left
026651003 MT First Metatarsal Osteotomy Locking Plate, left 3 Holes
026651004 MT First Metatarsal Osteotomy Locking Plate, left 4 Holes
026651005 MT First Metatarsal Osteotomy Locking Plate, left 5 Holes
026651006 MT First Metatarsal Osteotomy Locking Plate, left 6 Holes
026651007 MT First Metatarsal Osteotomy Locking Plate, left 7 Holes
026652000 MT First Metatarsal Osteotomy Locking Plate, right
026652003 MT First Metatarsal Osteotomy Locking Plate, right 3 Holes
026652004 MT First Metatarsal Osteotomy Locking Plate, right 4 Holes
026652005 MT First Metatarsal Osteotomy Locking Plate, right 5 Holes
026652006 MT First Metatarsal Osteotomy Locking Plate, right 6 Holes
026652007 MT First Metatarsal Osteotomy Locking Plate, right 7 Holes
026660000 MT Lateral Column Lengthening Locking Plate
026660008 MT Lateral Column Lengthening Locking Plate
026660010 MT Lateral Column Lengthening Locking Plate

WMS



Secretaría de
Gobierno de Salud



Ministerio de Salud y Desarrollo Social
Presidencia de la Nación

026660012	MT Lateral Column Lengthening Locking Plate	
030300006	1,5 mm Locking Screw, stardrive, self tapping	6 mm
030300007	1,5 mm Locking Screw, stardrive, self tapping	7 mm
030300008	1,5 mm Locking Screw, stardrive, self tapping	8 mm
030300009	1,5 mm Locking Screw, stardrive, self tapping	9 mm
030300010	1,5 mm Locking Screw, stardrive, self tapping	10 mm
030300011	1,5 mm Locking Screw, stardrive, self tapping	11 mm
030300012	1,5 mm Locking Screw, stardrive, self tapping	12 mm
030300013	1,5 mm Locking Screw, stardrive, self tapping	13 mm
030300014	1,5 mm Locking Screw, stardrive, self tapping	14 mm
030300015	1,5 mm Locking Screw, stardrive, self tapping	15 mm
030300016	1,5 mm Locking Screw, stardrive, self tapping	16 mm
030300017	1,5 mm Locking Screw, stardrive, self tapping	17 mm
030300018	1,5 mm Locking Screw, stardrive, self tapping	18 mm
030300019	1,5 mm Locking Screw, stardrive, self tapping	19 mm
030300020	1,5 mm Locking Screw, stardrive, self tapping	20 mm
030320006	2,0mm Locking Screw, stardrive, self tapping	6 mm
030320007	2,0mm Locking Screw, stardrive, self tapping	7 mm
030320008	2,0mm Locking Screw, stardrive, self tapping	8 mm
030320009	2,0mm Locking Screw, stardrive, self tapping	9 mm
030320010	2,0mm Locking Screw, stardrive, self tapping	10 mm

W/B

Sedes y Delegaciones

Tel. (+54-11) 4340-0800 - <http://www.anmat.gov.ar> - República Argentina

Sede Central
Av. de Mayo 869, CABA

Sede Alsina
Alsina 665/671, CABA

Sede INAME
Av. Caseros 2161, CABA

Sede INAL
Estados Unidos 25, CABA

Sede Prod. Médicos
Av. Belgrano 1480, CABA

Deleg. Mendoza
Remedios de Escalada de
San Martín 1909, Mendoza
Prov. de Mendoza

Deleg. Córdoba
Obispo Trejo 635,
Córdoba,
Prov. de Córdoba

Deleg. Paso de los Libres
Ruta Nacional 117, km.10,
CO.TE.CAR., Paso de los
Libres, Prov. de Corrientes

Deleg. Posadas
Roque González 1137,
Posadas, Prov. de
Misiones

Deleg. Santa Fé
Eva Perón 2456,
Santa Fé,
Prov. de Santa Fé

030320011	2,0mm Locking Screw, stardrive, self tapping	11 mm
030320012	2,0mm Locking Screw, stardrive, self tapping	12 mm
030320013	2,0mm Locking Screw, stardrive, self tapping	13 mm
030320014	2,0mm Locking Screw, stardrive, self tapping	14 mm
030320015	2,0mm Locking Screw, stardrive, self tapping	15 mm
030320016	2,0mm Locking Screw, stardrive, self tapping	16 mm
030320017	2,0mm Locking Screw, stardrive, self tapping	17 mm
030320018	2,0mm Locking Screw, stardrive, self tapping	18 mm
030320019	2,0mm Locking Screw, stardrive, self tapping	19 mm
030320020	2,0mm Locking Screw, stardrive, self tapping	20 mm
030320022	2,0mm Locking Screw, stardrive, self tapping	22 mm
030320024	2,0mm Locking Screw, stardrive, self tapping	24 mm
030320026	2,0mm Locking Screw, stardrive, self tapping	26 mm
030320028	2,0mm Locking Screw, stardrive, self tapping	28 mm
030320030	2,0mm Locking Screw, stardrive, self tapping	30 mm
030340006	2,7 mm Locking Screw, stardrive, selft-tapping	6 mm
030340007	2,7 mm Locking Screw, stardrive, selft-tapping	7 mm
030340008	2,7 mm Locking Screw, stardrive, selft-tapping	8 mm
030340009	2,7 mm Locking Screw, stardrive, selft-tapping	9 mm
030340010	2,7 mm Locking Screw, stardrive, selft-tapping	10 mm
030340011	2,7 mm Locking Screw, stardrive, selft-tapping	11 mm

WIB



Secretaría de
Gobierno de Salud



Ministerio de Salud y Desarrollo Social
Presidencia de la Nación

030340012	2,7 mm Locking Screw, stardrive, self-tapping	12 mm
030340013	2,7 mm Locking Screw, stardrive, self-tapping	13 mm
030340014	2,7 mm Locking Screw, stardrive, self-tapping	14 mm
030340015	2,7 mm Locking Screw, stardrive, self-tapping	15 mm
030340016	2,7 mm Locking Screw, stardrive, self-tapping	16 mm
030340017	2,7 mm Locking Screw, stardrive, self-tapping	17 mm
030340018	2,7 mm Locking Screw, stardrive, self-tapping	18 mm
030340019	2,7 mm Locking Screw, stardrive, self-tapping	19 mm
030340020	2,7 mm Locking Screw, stardrive, self-tapping	20 mm
030340022	2,7 mm Locking Screw, stardrive, self-tapping	22 mm
030340024	2,7 mm Locking Screw, stardrive, self-tapping	24 mm
030340026	2,7 mm Locking Screw, stardrive, self-tapping	26 mm
030340028	2,7 mm Locking Screw, stardrive, self-tapping	28 mm
030340030	2,7 mm Locking Screw, stardrive, self-tapping	30 mm
030340035	2.7mm Locking Screw, Stardrive, Self-tapping	35mm
030340040	2.7mm Locking Screw, Stardrive, Self-tapping	40mm
030340045	2.7mm Locking Screw, Stardrive, Self-tapping	45mm
030340050	2.7mm Locking Screw, Stardrive, Self-tapping	50mm
031271205	MD Anterior Process Calcaneal Locking Plate, small, left	5 Holes
031271705	MD Anterior Process Calcaneal Locking Plate, large, left	5 Holes
031272205	MD Anterior Process Calcaneal Locking Plate, small, right	5 Holes

WAB

Sedes y Delegaciones

Tel. (+54-11) 4340-0800 - <http://www.anmat.gov.ar> - República Argentina

Sede Central
Av. de Mayo 869, CABA

Sede Alsina
Alsina 665/671, CABA

Sede INAME
Av. Caseros 2161, CABA

Sede INAL
Estados Unidos 25, CABA

Sede Prod. Médicos
Av. Belgrano 1480, CABA

Deleg. Mendoza
Remedios de Escalada de
San Martín 1909, Mendoza
Prov. de Mendoza

Deleg. Córdoba
Obispo Trejo 635,
Córdoba,
Prov. de Córdoba

Deleg. Paso de los Libres
Ruta Nacional 117, km.10,
CO.TE.CAR., Paso de los
Libres, Prov. de Corrientes

Deleg. Posadas
Roque González 1137,
Posadas, Prov. de
Misiones

Deleg. Santa Fé
Eva Perón 2456,
Santa Fé,
Prov. de Santa Fé

031272705 MD Anterior Process Calcaneal Locking Plate, large, right 5 Holes
031281005 MD Posterior Tuberosity Calcaneal Locking Plate, small, left 5 Holes
031281006 MD Posterior Tuberosity Calcaneal Locking Plate, large, left 6 Holes
031282005 MD Posterior Tuberosity Calcaneal Locking Plate, small, right 5 Holes
031282006 MD Posterior Tuberosity Calcaneal Locking Plate, large, right 6 Holes
031291008 MD AP and PT Calcaneal Locking Plate, medium, left 8 Holes
031291009 MD AP and PT Calcaneal Locking Plate, large, left 9 Holes
031291308 MD AP and PT Calcaneal Locking Plate, small, left 8 Holes
031292008 MD AP and PT Calcaneal Locking Plate, medium, right 8 Holes
031292009 MD AP and PT Calcaneal Locking Plate, large, right 9 Holes
031292308 MD AP and PT Calcaneal Locking Plate, small, right 8 Holes
031301011 MD Calcaneal Locking Plate, small, left 11Holes
031302011 MD Calcaneal Locking Plate, small, right 11Holes
031311011 MD Calcaneal Locking Plate, medium, left 11Holes
031312011 MD Calcaneal Locking Plate, medium, right 11Holes
031321011 MD Calcaneal Locking Plate, large, left 11Holes
031322011 MD Calcaneal Locking Plate, large, right 11Holes
031360002 MT T-Locking Plate for First Metatarsa 2 Holes
031851006 MD Fifth Metatarsal Locking Plate
031852006 MD Fifth Metatarsal Locking Plate
031861009 MT Cuboid Locking Plate, left

WFB



Secretaría de
Gobierno de Salud



Ministerio de Salud y Desarrollo Social
Presidencia de la Nación

- 031862009 MT Cuboid Locking Plate, right
- 031870005 MT Talus Locking Plate
- 031941005 MT First MTC Plantar Fusion Locking Plate, large, left
- 031942005 MT First MTC Plantar Fusion Locking Plate, large, right
- 031951005 MT First MTC Plantar Fusion Locking Plate, small, left
- 031952005 MT First MTC Plantar Fusion Locking Plate, small, right
- 031961004 Distal Lateral Fibular Locking Plate, left 4 holes
- 031961005 Distal Lateral Fibular Locking Plate, left 5 holes
- 031961006 Distal Lateral Fibular Locking Plate, left 6 holes
- 031961007 Distal Lateral Fibular Locking Plate, left 7 holes
- 031961009 Distal Lateral Fibular Locking Plate, left 9 holes
- 031961011 Distal Lateral Fibular Locking Plate, left 11 holes
- 031961013 Distal Lateral Fibular Locking Plate, left 13 holes
- 031962004 Distal Lateral Fibular Locking Plate, right 4 holes
- 031962005 Distal Lateral Fibular Locking Plate, right 5 holes
- 031962006 Distal Lateral Fibular Locking Plate, right 6 holes
- 031962007 Distal Lateral Fibular Locking Plate, right 7 holes
- 031962009 Distal Lateral Fibular Locking Plate, right 9 holes
- 031962011 Distal Lateral Fibular Locking Plate, right 11 holes
- 031962013 Distal Lateral Fibular Locking Plate, right 13 holes
- 031970003 MD Distal Tibial Locking Plate 3 holes

UAB

Sedes y Delegaciones

Tel. (+54-11) 4340-0800 - <http://www.anmat.gov.ar> - República Argentina

Sede Central
Av. de Mayo 869, CABA

Sede Alsina
Alsina 665/671, CABA

Sede INAME
Av. Caseros 2161, CABA

Sede INAL
Estados Unidos 25, CABA

Sede Prod. Médicos
Av. Belgrano 1480, CABA

Deleg. Mendoza
Remedios de Escalada de
San Martín 1909, Mendoza
Prov. de Mendoza

Deleg. Córdoba
Obispo Trejo 635,
Córdoba,
Prov. de Córdoba

Deleg. Paso de los Libres
Ruta Nacional 117, km.10,
CO.TE.CAR., Paso de los
Libres, Prov. de Corrientes

Deleg. Posadas
Roque González 1137,
Posadas, Prov. de
Misiones

Deleg. Santa Fé
Eva Perón 2456,
Santa Fé,
Prov. de Santa Fé

031980004	MD Distal Medial Tibial Locking Plate	4 holes
031990004	MD Distal Anterior Tibial Locking Plate	4 holes
031990006	MD Distal Anterior Tibial Locking Plate	6 holes
032001009	MT Navicular Locking Plate,	left
032002009	MT Navicular Locking Plate,	right
040070010	3.5mm Locking Screw, Self-tapping	10mm
040070012	3.5mm Locking Screw, Self-tapping	12mm
040070014	3.5mm Locking Screw, Self-tapping	14mm
040070016	3.5mm Locking Screw, Self-tapping	16mm
040070018	3.5mm Locking Screw, Self-tapping	18mm
040070020	3.5mm Locking Screw, Self-tapping	20mm
040070022	3.5mm Locking Screw, Self-tapping	22mm
040070024	3.5mm Locking Screw, Self-tapping	24mm
040070026	3.5mm Locking Screw, Self-tapping	26mm
040070028	3.5mm Locking Screw, Self-tapping	28mm
040070030	3.5mm Locking Screw, Self-tapping	30mm
040070032	3.5mm Locking Screw, Self-tapping	32mm
040070034	3.5mm Locking Screw, Self-tapping	34mm
040070035	3.5mm Locking Screw, Self-tapping	35mm
040070036	3.5mm Locking Screw, Self-tapping	36mm
040070038	3.5mm Locking Screw, Self-tapping	38mm

WJB



Secretaría de
Gobierno de Salud



Ministerio de Salud y Desarrollo Social
Presidencia de la Nación

040070040	3.5mm Locking Screw, Self-tapping	40mm
040070042	3.5mm Locking Screw, Self-tapping	42mm
040070044	3.5mm Locking Screw, Self-tapping	44mm
040070045	3.5mm Locking Screw, Self-tapping	45mm
040070046	3.5mm Locking Screw, Self-tapping	46mm
040070048	3.5mm Locking Screw, Self-tapping	48mm
040070050	3.5mm Locking Screw, Self-tapping	50mm
040070052	3.5mm Locking Screw, Self-tapping	52mm
040070054	3.5mm Locking Screw, Self-tapping	54mm
040070055	3.5mm Locking Screw, Self-tapping	55mm
040070056	3.5mm Locking Screw, Self-tapping	56mm
040070058	3.5mm Locking Screw, Self-tapping	58mm
040070060	3.5mm Locking Screw, Self-tapping	60mm
040070065	3.5mm Locking Screw, Self-tapping	65mm
040070070	3.5mm Locking Screw, Self-tapping	70mm
040070075	3.5mm Locking Screw, Self-tapping	75mm
040080012	3.5mm Locking Screw, Self-drilling	12mm
040080014	3.5mm Locking Screw, Self-drilling	14mm
040080016	3.5mm Locking Screw, Self-drilling	16mm
040080018	3.5mm Locking Screw, Self-drilling	18mm
040080020	3.5mm Locking Screw, Self-drilling	20mm

WLB

Sedes y Delegaciones

Tel. (+54-11) 4340-0800 - <http://www.anmat.gov.ar> - República Argentina

Sede Central
Av. de Mayo 869, CABA

Sede Alsina
Alsina 665/671, CABA

Sede INAME
Av. Caseros 2161, CABA

Sede INAL
Estados Unidos 25, CABA

Sede Prod. Médicos
Av. Belgrano 1480, CABA

Deleg. Mendoza
Remedios de Escalada de
San Martín 1909, Mendoza
Prov. de Mendoza

Deleg. Córdoba
Obispo Trejo 635,
Córdoba,
Prov. de Córdoba

Deleg. Paso de los Libres
Ruta Nacional 117, km.10,
CO.TE.CAR., Paso de los
Libres, Prov. de Corrientes

Deleg. Posadas
Roque González 1137,
Posadas, Prov. de
Misiones

Deleg. Santa Fé
Eva Perón 2456,
Santa Fé,
Prov. de Santa Fé

040080022	3.5mm Locking Screw, Self-drilling	22mm
040080024	3.5mm Locking Screw, Self-drilling	24mm
040080026	3.5mm Locking Screw, Self-drilling	26mm
040080028	3.5mm Locking Screw, Self-drilling	28mm
040080030	3.5mm Locking Screw, Self-drilling	30mm
040080032	3.5mm Locking Screw, Self-drilling	32mm
040080034	3.5mm Locking Screw, Self-drilling	34mm
040080035	3.5mm Locking Screw, Self-drilling	35mm
040080036	3.5mm Locking Screw, Self-drilling	36mm
040080038	3.5mm Locking Screw, Self-drilling	38mm
040080040	3.5mm Locking Screw, Self-drilling	40mm
040080042	3.5mm Locking Screw, Self-drilling	42mm
040080044	3.5mm Locking Screw, Self-drilling	44mm
040080045	3.5mm Locking Screw, Self-drilling	45mm
040080046	3.5mm Locking Screw, Self-drilling	46mm
040080048	3.5mm Locking Screw, Self-drilling	48mm
040080050	3.5mm Locking Screw, Self-drilling	50mm
040080052	3.5mm Locking Screw, Self-drilling	52mm
040080054	3.5mm Locking Screw, Self-drilling	54mm
040080055	3.5mm Locking Screw, Self-drilling	55mm
040080056	3.5mm Locking Screw, Self-drilling	56mm

WKB



Secretaría de
Gobierno de Salud



Ministerio de Salud y Desarrollo Social
Presidencia de la Nación

040080058	3.5mm Locking Screw, Self-drilling	58mm
040080060	3.5mm Locking Screw, Self-drilling	60mm
040090016	5.0mm Locking Screw, Self-tapping	16mm
040090018	5.0mm Locking Screw, Self-tapping	18mm
040090020	5.0mm Locking Screw, Self-tapping	20mm
040090022	5.0mm Locking Screw, Self-tapping	22mm
040090024	5.0mm Locking Screw, Self-tapping	24mm
040090026	5.0mm Locking Screw, Self-tapping	26mm
040090028	5.0mm Locking Screw, Self-tapping	28mm
040090030	5.0mm Locking Screw, Self-tapping	30mm
040090032	5.0mm Locking Screw, Self-tapping	32mm
040090034	5.0mm Locking Screw, Self-tapping	34mm
040090036	5.0mm Locking Screw, Self-tapping	36mm
040090038	5.0mm Locking Screw, Self-tapping	38mm
040090040	5.0mm Locking Screw, Self-tapping	40mm
040090042	5.0mm Locking Screw, Self-tapping	42mm
040090044	5.0mm Locking Screw, Self-tapping	44mm
040090045	5.0mm Locking Screw, Self-tapping	45mm
040090046	5.0mm Locking Screw, Self-tapping	46mm
040090048	5.0mm Locking Screw, Self-tapping	48mm
040090050	5.0mm Locking Screw, Self-tapping	50mm

Handwritten signature

Sedes y Delegaciones

Tel. (+54-11) 4340-0800 - <http://www.anmat.gov.ar> - República Argentina

Sede Central
Av. de Mayo 869, CABA

Sede Alsina
Alsina 665/671, CABA

Sede INAME
Av. Caseros 2161, CABA

Sede INAL
Estados Unidos 25, CABA

Sede Prod. Médicos
Av. Belgrano 1480, CABA

Deleg. Mendoza
Remedios de Escalada de
San Martín 1909, Mendoza
Prov. de Mendoza

Deleg. Córdoba
Obispo Trejo 635,
Córdoba,
Prov. de Córdoba

Deleg. Paso de los Libres
Ruta Nacional 117, km.10,
CO.TE.CAR., Paso de los
Libres, Prov. de Corrientes

Deleg. Posadas
Roque González 1137,
Posadas, Prov. de
Misiones

Deleg. Santa Fé
Eva Perón 2456,
Santa Fé,
Prov. de Santa Fé

040090055	5.0mm Locking Screw, Self-tapping	55mm
040090060	5.0mm Locking Screw, Self-tapping	60mm
040090065	5.0mm Locking Screw, Self-tapping	65mm
040090070	5.0mm Locking Screw, Self-tapping	70mm
040090075	5.0mm Locking Screw, Self-tapping	75mm
040090080	5.0mm Locking Screw, Self-tapping	80mm
040090085	5.0mm Locking Screw, Self-tapping	85mm
040090090	5.0mm Locking Screw, Self-tapping	90mm
040100020	5.0mm Locking Screw, Self-drilling	20mm
040100022	5.0mm Locking Screw, Self-drilling	22mm
040100024	5.0mm Locking Screw, Self-drilling	24mm
040100026	5.0mm Locking Screw, Self-drilling	26mm
040100028	5.0mm Locking Screw, Self-drilling	28mm
040100030	5.0mm Locking Screw, Self-drilling	30mm
040100032	5.0mm Locking Screw, Self-drilling	32mm
040100034	5.0mm Locking Screw, Self-drilling	34mm
040100036	5.0mm Locking Screw, Self-drilling	36mm
040100038	5.0mm Locking Screw, Self-drilling	38mm
040100040	5.0mm Locking Screw, Self-drilling	40mm
040100042	5.0mm Locking Screw, Self-drilling	42mm
040100044	5.0mm Locking Screw, Self-drilling	44mm

WMB



Secretaría de
Gobierno de Salud



Ministerio de Salud y Desarrollo Social
Presidencia de la Nación

040100046	5.0mm Locking Screw, Self-drilling	46mm
040100048	5.0mm Locking Screw, Self-drilling	48mm
040100050	5.0mm Locking Screw, Self-drilling	50mm
040100055	5.0mm Locking Screw, Self-drilling	55mm
040100060	5.0mm Locking Screw, Self-drilling	60mm
040100065	5.0mm Locking Screw, Self-drilling	65mm
040100070	5.0mm Locking Screw, Self-drilling	70mm
040100075	5.0mm Locking Screw, Self-drilling	75mm
040100080	5.0mm Locking Screw, Self-drilling	80mm
040100085	5.0mm Locking Screw, Self-drilling	85mm
040100090	5.0mm Locking Screw, Self-drilling	90mm
040130060	7.3mm Locking Screw	60mm
040130065	7.3mm Locking Screw	65mm
040130070	7.3mm Locking Screw	70mm
040130075	7.3mm Locking Screw	75mm
040130076	7.3mm Locking Screw, 76.5mm	76.5mm
040130080	7.3mm Locking Screw	80mm
040130083	7.3mm Locking Screw, 83.5mm	83.5mm
040130085	7.3mm Locking Screw	85mm
040130090	7.3mm Locking Screw	90mm
040130095	7.3mm Locking Screw	95mm

W/B

Sedes y Delegaciones

Tel. (+54-11) 4340-0800 - <http://www.anmat.gov.ar> - República Argentina

Sede Central
Av. de Mayo 869, CABA

Sede Alsina
Alsina 665/671, CABA

Sede INAME
Av. Caseros 2161, CABA

Sede INAL
Estados Unidos 25, CABA

Sede Prod. Médicos
Av. Belgrano 1480, CABA

Deleg. Mendoza
Remedios de Escalada de
San Martín 1909, Mendoza
Prov. de Mendoza

Deleg. Córdoba
Obispo Trejo 635,
Córdoba,
Prov. de Córdoba

Deleg. Paso de los Libres
Ruta Nacional 117, km.10,
CO.TE.CAR., Paso de los
Libres, Prov. de Corrientes

Deleg. Posadas
Roque González 1137,
Posadas, Prov. de
Misiones

Deleg. Santa Fé
Eva Perón 2456,
Santa Fé,
Prov. de Santa Fé

040130100	7.3mm Locking Screw	100mm
040130105	7.3mm Locking Screw	105mm
040130110	7.3mm Locking Screw	110mm
040130115	7.3mm Locking Screw	115mm
040140060	6.5mm Cannulated Locking Screw, Fully Threaded	60mm
040140065	6.5mm Cannulated Locking Screw, Fully Threaded	65mm
040140070	6.5mm Cannulated Locking Screw, Fully Threaded	70mm
040140075	6.5mm Cannulated Locking Screw, Fully Threaded	75mm
040140080	6.5mm Cannulated Locking Screw, Fully Threaded	80mm
040140085	6.5mm Cannulated Locking Screw, Fully Threaded	85mm
040140090	6.5mm Cannulated Locking Screw, Fully Threaded	90mm
040140095	6.5mm Cannulated Locking Screw, Fully Threaded	95mm
040140100	6.5mm Cannulated Locking Screw, Fully Threaded	100mm
040140105	6.5mm Cannulated Locking Screw, Fully Threaded	105mm
040140110	6.5mm Cannulated Locking Screw, Fully Threaded	110mm
040140115	6.5mm Cannulated Locking Screw, Fully Threaded	115mm
040140120	6.5mm Cannulated Locking Screw, Fully Threaded	120mm
040240006	2.4mm Locking Screw, Self-tapping, Stardrive	6mm
040240007	2.4mm Locking Screw, Self-tapping, Stardrive	7mm
040240008	2.4mm Locking Screw, Self-tapping, Stardrive	8mm
040240009	2.4mm Locking Screw, Self-tapping, Stardrive	9mm

WKS



Secretaría de
Gobierno de Salud



Ministerio de Salud y Desarrollo Social
Presidencia de la Nación

040240010	2.4mm Locking Screw, Self-tapping, Stardrive	10mm
040240011	2.4mm Locking Screw, Self-tapping, Stardrive	11mm
040240012	2.4mm Locking Screw, Self-tapping, Stardrive	12mm
040240013	2.4mm Locking Screw, Self-tapping, Stardrive	13mm
040240014	2.4mm Locking Screw, Self-tapping, Stardrive	14mm
040240015	2.4mm Locking Screw, Self-tapping, Stardrive	15mm
040240016	2.4mm Locking Screw, Self-tapping, Stardrive	16mm
040240017	2.4mm Locking Screw, Self-tapping, Stardrive	17mm
040240018	2.4mm Locking Screw, Self-tapping, Stardrive	18mm
040240019	2.4mm Locking Screw, Self-tapping, Stardrive	19mm
040240020	2.4mm Locking Screw, Self-tapping, Stardrive	20mm
040240022	2.4mm Locking Screw, Self-tapping, Stardrive	22mm
040240024	2.4mm Locking Screw, Self-tapping, Stardrive	24mm
040240026	2.4mm Locking Screw, Self-tapping, Stardrive	26mm
040240028	2.4mm Locking Screw, Self-tapping, Stardrive	28mm
040240030	2.4mm Locking Screw, Self-tapping, Stardrive	30mm
040260070	7.3mm Cannulated Locking Screw, Thread Length 16mm	70mm
040260075	7.3mm Cannulated Locking Screw, Thread Length 16mm	75mm
040260080	7.3mm Cannulated Locking Screw, Thread Length 16mm	80mm
040260085	7.3mm Cannulated Locking Screw, Thread Length 16mm	85mm
040260090	7.3mm Cannulated Locking Screw, Thread Length 16mm	90mm

WMB

Sedes y Delegaciones

Tel. (+54-11) 4340-0800 - <http://www.anmat.gov.ar> - República Argentina

Sede Central
Av. de Mayo 869, CABA

Sede Alsina
Alsina 665/671, CABA

Sede INAME
Av. Caseros 2161, CABA

Sede INAL
Estados Unidos 25, CABA

Sede Prod. Médicos
Av. Belgrano 1480, CABA

Deleg. Mendoza
Remedios de Escalada de
San Martín 1909, Mendoza
Prov. de Mendoza

Deleg. Córdoba
Obispo Trejo 635,
Córdoba,
Prov. de Córdoba

Deleg. Paso de los Libres
Ruta Nacional 117, km.10,
CO.TE.CAR., Paso de los
Libres, Prov. de Corrientes

Deleg. Posadas
Roque González 1137,
Posadas, Prov. de
Misiones

Deleg. Santa Fé
Eva Perón 2456,
Santa Fé,
Prov. de Santa Fé

040260095 7.3mm Cannulated Locking Screw, Thread Length 16mm 95mm

040260100 7.3mm Cannulated Locking Screw, Thread Length 16mm

100mm

040260105 7.3mm Cannulated Locking Screw, Thread Length 16mm

105mm

040260110 7.3mm Cannulated Locking Screw, Thread Length 16mm

110mm

040270070 7.3mm Cannulated Locking Screw, Thread Length 32mm 70mm

040270075 7.3mm Cannulated Locking Screw, Thread Length 32mm 75mm

040270080 7.3mm Cannulated Locking Screw, Thread Length 32mm 80mm

040270085 7.3mm Cannulated Locking Screw, Thread Length 32mm 85mm

040270090 7.3mm Cannulated Locking Screw, Thread Length 32mm 90mm

040270095 7.3mm Cannulated Locking Screw, Thread Length 32mm 95mm

040270100 7.3mm Cannulated Locking Screw, Thread Length 32mm

100mm

040270105 7.3mm Cannulated Locking Screw, Thread Length 32mm

105mm

040270110 7.3mm Cannulated Locking Screw, Thread Length 32mm

110mm

040270115 7.3mm Cannulated Locking Screw, Thread Length 32mm

115mm

WJ2



Secretaría de
Gobierno de Salud



Ministerio de Salud y Desarrollo Social
Presidencia de la Nación

040300010	5.0mm Locking Screw, Self-tapping, Flat Head	10mm
040300014	5.0mm Locking Screw, Self-tapping, Flat Head	14mm
040300018	5.0mm Locking Screw, Self-tapping, Flat Head	18mm
040760008	2.4mm Locking Screw II, variable angle, stardrive, self-tapping	8mm
040760010	2.4mm Locking Screw II, variable angle, stardrive, self-tapping	10mm
040760012	2.4mm Locking Screw II, variable angle, stardrive, self-tapping	12mm
040760014	2.4mm Locking Screw II, variable angle, stardrive, self-tapping	14mm
040760016	2.4mm Locking Screw II, variable angle, stardrive, self-tapping	16mm
040760018	2.4mm Locking Screw II, variable angle, stardrive, self-tapping	18mm
040760020	2.4mm Locking Screw II, variable angle, stardrive, self-tapping	20mm
040760022	2.4mm Locking Screw II, variable angle, stardrive, self-tapping	22mm
040760024	2.4mm Locking Screw II, variable angle, stardrive, self-tapping	24mm
040760026	2.4mm Locking Screw II, variable angle, stardrive, self-tapping	26mm
040760028	2.4mm Locking Screw II, variable angle, stardrive, self-tapping	28mm
040760030	2.4mm Locking Screw II, variable angle, stardrive, self-tapping	30mm
040910060	5.5mm Locking Screw, Self-tapping	60mm
040910065	5.5mm Locking Screw, Self-tapping	65mm
040910070	5.5mm Locking Screw, Self-tapping	70mm
040910075	5.5mm Locking Screw, Self-tapping	75mm
040910080	5.5mm Locking Screw, Self-tapping	80mm

WAB

Sedes y Delegaciones

Tel. (+54-11) 4340-0800 - <http://www.anmat.gov.ar> - República Argentina

Sede Central
Av. de Mayo 869, CABA

Sede Alsina
Alsina 665/671, CABA

Sede INAME
Av. Caseros 2161, CABA

Sede INAL
Estados Unidos 25, CABA

Sede Prod. Médicos
Av. Belgrano 1480, CABA

Deleg. Mendoza
Remedios de Escalada de
San Martín 1909, Mendoza
Prov. de Mendoza

Deleg. Córdoba
Obispo Trejo 635,
Córdoba,
Prov. de Córdoba

Deleg. Paso de los Libres
Ruta Nacional 117, km.10,
CO.TE.CAR., Paso de los
Libres, Prov. de Corrientes

Deleg. Posadas
Roque González 1137,
Posadas, Prov. de
Misiones

Deleg. Santa Fé
Eva Perón 2456,
Santa Fé,
Prov. de Santa Fé

040910085	5.5mm Locking Screw, Self-tapping	85mm
040910090	5.5mm Locking Screw, Self-tapping	90mm
040910095	5.5mm Locking Screw, Self-tapping	95mm
040910100	5.5mm Locking Screw, Self-tapping	100mm
040910105	5.5mm Locking Screw, Self-tapping	105mm
040910110	5.5mm Locking Screw, Self-tapping	110mm
040920060	6.5mm Cannulated Locking Screw, Self-tapping	60mm
040920065	6.5mm Cannulated Locking Screw, Self-tapping	65mm
040920070	6.5mm Cannulated Locking Screw, Self-tapping	70mm
040920075	6.5mm Cannulated Locking Screw, Self-tapping	75mm
040920080	6.5mm Cannulated Locking Screw, Self-tapping	80mm
040920085	6.5mm Cannulated Locking Screw, Self-tapping	85mm
040920090	6.5mm Cannulated Locking Screw, Self-tapping	90mm
040920095	6.5mm Cannulated Locking Screw, Self-tapping	95mm
040920100	6.5mm Cannulated Locking Screw, Self-tapping	100mm
040920105	6.5mm Cannulated Locking Screw, Self-tapping	105mm
040920110	6.5mm Cannulated Locking Screw, Self-tapping	110mm
040940014	4,0 mm Locking Screw, self-tapping	14 mm
040940016	4,0 mm Locking Screw, self-tapping	16 mm
040940018	4,0 mm Locking Screw, self-tapping	18 mm
040940020	4,0 mm Locking Screw, self-tapping	20 mm

LM3



Secretaría de
Gobierno de Salud



Ministerio de Salud y Desarrollo Social
Presidencia de la Nación

040940022	4,0 mm Locking Screw, self-tapping	22 mm
040940024	4,0 mm Locking Screw, self-tapping	24 mm
040940026	4,0 mm Locking Screw, self-tapping	26 mm
040940028	4,0 mm Locking Screw, self-tapping	28 mm
040940030	4,0 mm Locking Screw, self-tapping	30 mm
040940032	4,0 mm Locking Screw, self-tapping	32 mm
040940034	4,0 mm Locking Screw, self-tapping	34 mm
040940036	4,0 mm Locking Screw, self-tapping	36 mm
040940038	4,0 mm Locking Screw, self-tapping	38 mm
040940040	4,0 mm Locking Screw, self-tapping	40 mm
040940045	4,0 mm Locking Screw, self-tapping	45 mm
040940050	4,0 mm Locking Screw, self-tapping	50 mm
040940055	4,0 mm Locking Screw, self-tapping	55 mm
040940060	4,0 mm Locking Screw, self-tapping	60 mm
040940065	4,0 mm Locking Screw, self-tapping	65 mm
040940070	4,0 mm Locking Screw, self-tapping	70 mm
040940075	4,0 mm Locking Screw, self-tapping	75 mm
040940080	4,0 mm Locking Screw, self-tapping	80 mm
041050014	1,8 mm Locking Buttress Screw, variable angle, stardrive	14 mm
041050016	1,8 mm Locking Buttress Screw, variable angle, stardrive	16 mm
041050018	1,8 mm Locking Buttress Screw, variable angle, stardrive	18 mm

WAB

Sedes y Delegaciones

Tel. (+54-11) 4340-0800 - <http://www.anmat.gov.ar> - República Argentina

Sede Central
Av. de Mayo 869, CABA

Sede Alsina
Alsina 665/671, CABA

Sede INAME
Av. Caseros 2161, CABA

Sede INAL
Estados Unidos 25, CABA

Sede Prod. Médicos
Av. Belgrano 1480, CABA

Deleg. Mendoza
Remedios de Escalada de
San Martín 1909, Mendoza
Prov. de Mendoza

Deleg. Córdoba
Obispo Trejo 635,
Córdoba,
Prov. de Córdoba

Deleg. Paso de los Libres
Ruta Nacional 117, km.10,
CO.TE.CAR., Paso de los
Libres, Prov. de Corrientes

Deleg. Posadas
Roque González 1137,
Posadas, Prov. de
Misiones

Deleg. Santa Fé
Eva Perón 2456,
Santa Fé,
Prov. de Santa Fé

041050020 1,8 mm Locking Buttress Screw, variable angle, stardrive20 mm
041050022 1,8 mm Locking Buttress Screw, variable angle, stardrive22 mm
041050024 1,8 mm Locking Buttress Screw, variable angle, stardrive24 mm
041050026 1,8 mm Locking Buttress Screw, variable angle, stardrive26 mm
041170010 MT 2.7mm Locking Screw , stardrive, self-tapping 10mm
041170012 MT 2.7mm Locking Screw , stardrive, self-tapping 12mm
041170014 MT 2.7mm Locking Screw , stardrive, self-tapping 14mm
041170016 MT 2.7mm Locking Screw , stardrive, self-tapping 16mm
041170018 MT 2.7mm Locking Screw , stardrive, self-tapping 18mm
041170020 MT 2.7mm Locking Screw , stardrive, self-tapping 20mm
041170022 MT 2.7mm Locking Screw , stardrive, self-tapping 22mm
041170024 MT 2.7mm Locking Screw , stardrive, self-tapping 24mm
041170026 MT 2.7mm Locking Screw , stardrive, self-tapping 26mm
041170028 MT 2.7mm Locking Screw , stardrive, self-tapping 28mm
041170030 MT 2.7mm Locking Screw , stardrive, self-tapping 30mm
041170032 MT 2.7mm Locking Screw , stardrive, self-tapping 32mm
041170034 MT 2.7mm Locking Screw , stardrive, self-tapping 34mm
041170036 MT 2.7mm Locking Screw , stardrive, self-tapping 36mm
041170038 MT 2.7mm Locking Screw , stardrive, self-tapping 38mm
041170040 MT 2.7mm Locking Screw , stardrive, self-tapping 40mm
041170042 MT 2.7mm Locking Screw , stardrive, self-tapping 42mm

W/B



Secretaría de
Gobierno de Salud



Ministerio de Salud y Desarrollo Social
Presidencia de la Nación

041170044 MT 2.7mm Locking Screw , stardrive, self-tapping 44mm
041170046 MT 2.7mm Locking Screw , stardrive, self-tapping 46mm
041170048 MT 2.7mm Locking Screw , stardrive, self-tapping 48mm
041170050 MT 2.7mm Locking Screw , stardrive, self-tapping 50mm
041180010 MT 3.0mm Locking Screw, stardrive, self-tapping 10mm
041180012 MT 3.0mm Locking Screw, stardrive, self-tapping 12mm
041180014 MT 3.0mm Locking Screw, stardrive, self-tapping 14mm
041180016 MT 3.0mm Locking Screw, stardrive, self-tapping 16mm
041180018 MT 3.0mm Locking Screw, stardrive, self-tapping 18mm
041180020 MT 3.0mm Locking Screw, stardrive, self-tapping 20mm
041180022 MT 3.0mm Locking Screw, stardrive, self-tapping 22mm
041180024 MT 3.0mm Locking Screw, stardrive, self-tapping 24mm
041180026 MT 3.0mm Locking Screw, stardrive, self-tapping 26mm
041180028 MT 3.0mm Locking Screw, stardrive, self-tapping 28mm
041180030 MT 3.0mm Locking Screw, stardrive, self-tapping 30mm
041180032 MT 3.0mm Locking Screw, stardrive, self-tapping 32mm
041180034 MT 3.0mm Locking Screw, stardrive, self-tapping 34mm
041180036 MT 3.0mm Locking Screw, stardrive, self-tapping 36mm
041180038 MT 3.0mm Locking Screw, stardrive, self-tapping 38mm
041180040 MT 3.0mm Locking Screw, stardrive, self-tapping 40mm
041180042 MT 3.0mm Locking Screw, stardrive, self-tapping 42mm

LMB

Sedes y Delegaciones

Tel. (+54-11) 4340-0800 - <http://www.anmat.gov.ar> - República Argentina

Sede Central
Av. de Mayo 869, CABA

Sede Alsina
Alsina 665/671, CABA

Sede INAME
Av. Caseros 2161, CABA

Sede INAL
Estados Unidos 25, CABA

Sede Prod. Médicos
Av. Belgrano 1480, CABA

Deleg. Mendoza
Remedios de Escalada de
San Martín 1909, Mendoza
Prov. de Mendoza

Deleg. Córdoba
Obispo Trejo 635,
Córdoba,
Prov. de Córdoba

Deleg. Paso de los Libres
Ruta Nacional 117, km.10,
CO.TE.CAR., Paso de los
Libres, Prov. de Corrientes

Deleg. Posadas
Roque González 1137,
Posadas, Prov. de
Misiones

Deleg. Santa Fé
Eva Perón 2456,
Santa Fé,
Prov. de Santa Fé

041180044	MT 3.0mm Locking Screw, stardrive, self-tapping	44mm
041180046	MT 3.0mm Locking Screw, stardrive, self-tapping	46mm
041180048	MT 3.0mm Locking Screw, stardrive, self-tapping	48mm
041180050	MT 3.0mm Locking Screw, stardrive, self-tapping	50mm
041190010	MT 3.5mm Locking Screw, stardrive, self-tapping	10mm
041190012	MT 3.5mm Locking Screw, stardrive, self-tapping	12mm
041190014	MT 3.5mm Locking Screw, stardrive, self-tapping	14mm
041190016	MT 3.5mm Locking Screw, stardrive, self-tapping	16mm
041190018	MT 3.5mm Locking Screw, stardrive, self-tapping	18mm
041190020	MT 3.5mm Locking Screw, stardrive, self-tapping	20mm
041190022	MT 3.5mm Locking Screw, stardrive, self-tapping	22mm
041190024	MT 3.5mm Locking Screw, stardrive, self-tapping	24mm
041190026	MT 3.5mm Locking Screw, stardrive, self-tapping	26mm
041190028	MT 3.5mm Locking Screw, stardrive, self-tapping	28mm
041190030	MT 3.5mm Locking Screw, stardrive, self-tapping	30mm
041190032	MT 3.5mm Locking Screw, stardrive, self-tapping	32mm
041190034	MT 3.5mm Locking Screw, stardrive, self-tapping	34mm
041190036	MT 3.5mm Locking Screw, stardrive, self-tapping	36mm
041190038	MT 3.5mm Locking Screw, stardrive, self-tapping	38mm
041190040	MT 3.5mm Locking Screw, stardrive, self-tapping	40mm
041190042	MT 3.5mm Locking Screw, stardrive, self-tapping	42mm

WB



041190044 MT 3.5mm Locking Screw, stardrive, self-tapping 44mm

041190046 MT 3.5mm Locking Screw, stardrive, self-tapping 46mm

041190048 MT 3.5mm Locking Screw, stardrive, self-tapping 48mm

041190050 MT 3.5mm Locking Screw, stardrive, self-tapping 50mm

Instrumental Asociado

Período de vida útil: NC

Forma de presentación: no estéril, por unidad

Condición de uso: Venta exclusiva a profesionales e instituciones sanitarias

Fabricante:

DOUBLE MEDICAL TECHNOLOGY INC.

Lugar de Elaboración:

No. 18, Shanbianhong East Road, Haicang District, Xiamen, 361026, P.R. China.

Se extiende el presente Certificado de Autorización e Inscripción del PM-2342-01, con una vigencia de cinco (5) años a partir de la fecha de la Disposición autorizante.

Expediente N° 1-47-3110-1051-19-0

Disposición N°

6789 | 22 AGO. 2019

Dr. Waldo Belloso
Subadministrador Nacional
ANMAT

Sedes y Delegaciones

Tel. (+54-11) 4340-0800 - <http://www.anmat.gov.ar> - República Argentina

Sede Central

Av. de Mayo 869, CABA

Sede Alsina

Alsina 665/671, CABA

Sede INAME

Av. Caseros 2161, CABA

Sede INAL

Estados Unidos 25, CABA

Sede Prod. Médicos

Av. Belgrano 1480, CABA

Deleg. Mendoza

Remedios de Escalada de
San Martín 1909, Mendoza
Prov. de Mendoza

Deleg. Córdoba

Obispo Trejo 635,
Córdoba,
Prov. de Córdoba

Deleg. Paso de los Libres

Ruta Nacional 117, km.10,
CO.TE.CAR., Paso de los
Libres, Prov. de Corrientes

Deleg. Posadas

Roque González 1137,
Posadas, Prov. de
Misiones

Deleg. Santa Fé

Eva Perón 2456,
Santa Fé,
Prov. de Santa Fé